Capio™ SLIM

OPEN ACCESS

Suture Capturing Device

Directions for Use	2
Instrucciones de uso	8
Mode d'emploi	14
Gebrauchsanweisung	20
Istruzioni per l'uso	26
Gebruiksaanwijzing	32
Instruções de Utilização	38

50635526-01

2019-12

TABLE OF CONTENTS

WARNING 3

PRIOR TO USE..... 3

DEVICE DESCRIPTION..... 3

INTENDED USE/INDICATIONS FOR USE 3

WARNINGS..... 4

PRECAUTIONS..... 4

CONTRAINDICATIONS..... 4

ADVERSE EVENTS 4

DIRECTIONS FOR USE..... 4

TECHNIQUE TIPS 6

HOW SUPPLIED..... 6

WARRANTY 7

Capio™ SLIM

OPEN ACCESS

Suture Capturing Device

Rx ONLY

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PRIOR TO USE

Carefully inspect the unit to verify that neither the contents nor the sterile package has been damaged in shipment. If damage is found contact your Boston Scientific representative.

DEVICE DESCRIPTION

The Capio SLIM Suture Capturing Device was designed to be used by the physician to facilitate the consistent placement of sutures in difficult to access locations during open surgical procedures. The Capio SLIM Device was designed to throw, catch and retrieve sutures in one step. The device consists of a one-hand-activated plunger at the proximal end, a tubular shaft and a curved needle driver at the distal end. Depressing the plunger actuates the carrier and drives the suture/dart through the tissue. The needle is automatically caught by the device's needle catch mechanism for easy suture tying, or the device may be re-loaded for additional suture throws.

Sutures sold separately.

Contents

(1) Capio SLIM Suture Capturing Device

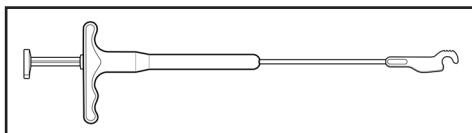
Figure 1 Capio SLIM Suture Capturing Device

INTENDED USE/INDICATIONS FOR USE

The Capio SLIM Suture Capturing Device is intended for use in general suturing applications during open surgery to assist in the placement of suture material in tissues at the operative site. The Capio SLIM Suture Capturing Device is to be used with sensory and/or direct visual control.

The Capio SLIM Suture Capturing Device is also intended for use as an aid in the insertion, placement, fixation or anchoring of surgical mesh during urogynecologic procedures.

WARNINGS

- Punctures or lacerations of vessels, nerves, bladder, urethra or bowel may occur during placement and may require surgical repair.
- This product is intended for use only by clinicians with adequate training and experience in open techniques and procedures.
- Once the suturing action has begun, do not turn, rotate or torque the Capio™ SLIM Device. To do so, may result in injury.

PRECAUTIONS

- Use only size 0 Capio Sutures with the Capio SLIM Suture Capturing Device.
- The physician is advised to consult the medical literature regarding techniques, complications and hazards associated with the intended procedures.
- The Capio SLIM Suture Capturing Device is not intended to be used to manipulate tissue, organ or bone.
- Do not bend or alter the device as this may result in mechanical damage to the device.

CONTRAINDICATIONS

The use of the Capio SLIM Suture Capturing Device is contraindicated as follows:

- For placing sutures into or through bone.

ADVERSE EVENTS

Potential complications associated with the use of the Capio SLIM Device include, but are not limited to:

- allergic reaction
- bleeding
- dyspareunia
- hematoma
- infection
- inflammation
- injury to internal organs/tissue including perforation and occlusion
- injury to vessels and nerves
- pain
- unretrieved device fragment

DIRECTIONS FOR USE

Directions for Open Surgical Procedure

Refer to Technique Tips section for added assistance.

1. Place the suture dart into the carrier at the distal end of the Capio SLIM Suture Capturing Device. Engage the suture dart into the carrier by gently pulling down on the suture, ensuring that the dart is pulled completely into the carrier. (Refer to Figure 2.)

Note: When loading a dart into the device, verify that the dart is properly positioned in the carrier. The tip of the dart should not protrude from the Capio Device tip.


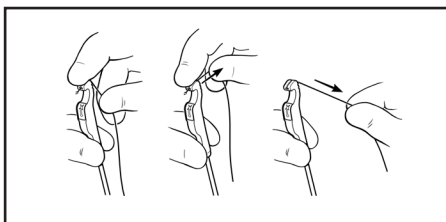
Figure 2.

2. Position the device with the distal tip in the desired location against the tissue to be sutured as shown in Figure 3 and hold firmly in place. Use slight counter-traction on the suture while positioning the Capio™ SLIM Suture Capturing Device tip in order to maintain the position of the dart in the carrier. Avoid excessive counter-traction on the suture as this has the potential to damage the suture during deployment.

Note: The technique of maintaining slight counter-traction on the suture is best described as placing just enough tension on the suture to keep the dart in the carrier. Excessive counter-traction is best described as preventing free movement of the suture during deployment.

3. Fully depress the plunger in one continuous motion, as shown in Figures 3 and 4, to pass the suture through the tissue.


Note: Once the suturing action has been initiated, completely depress the plunger. Do not stop in mid-stitch and allow the dart to reverse direction. It is important to maintain pressure against the tissue with the device to ensure that the dart penetrates the tissue and latches into the catch.

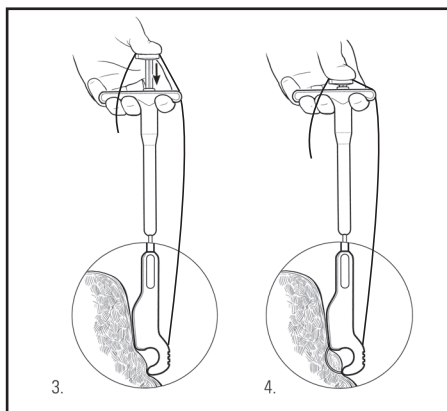Avoid excessive counter-traction that would prevent free movement of the suture during deployment

Figures 3 and 4.

4. Unplunge. Release plunger and verify it is back in its original position. This allows the carrier to retract into the housing, leaving the dart in the catch and the suture through the tissue as shown in Figure 5 and 6.

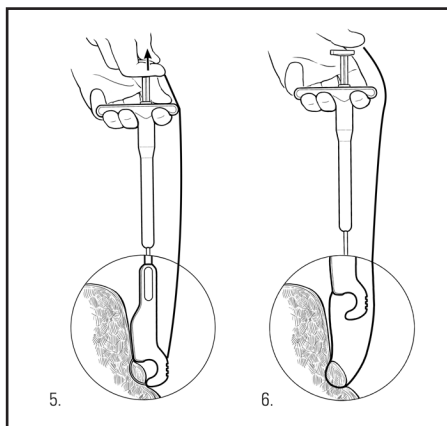Figures 5 and 6.

5. Remove the Capio™ SLIM Suture Capturing Device (Figure 6) from the incision.
6. Remove the dart from the catch by holding onto the suture and sliding the dart up toward the handle to the key hole slot. If additional suturing is required, the surgeon may repeat Steps 1 through 6.
7. Complete surgical procedure as desired. Ensure darts are properly discarded.

TECHNIQUE TIPS

Read and understand the following technique tips for the Capio SLIM Suture Capturing Device.

- Use slight counter-traction on the suture while positioning the Capio SLIM Suture Capturing Device tip, in order to maintain the position of the dart in the carrier. Avoid excessive counter-traction on the suture as this has the potential to damage the suture during deployment. The technique of maintaining slight counter-traction on the suture is best described as placing just enough tension on the suture to keep the dart in the carrier. Excessive counter-traction is best described as preventing free movement of the suture during deployment.

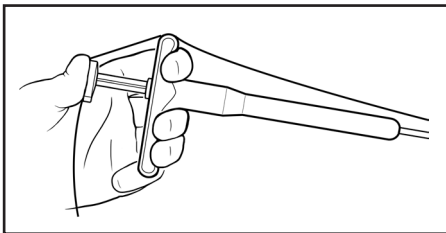
Figure 7.

- Between each throw, rinse/soak the head of device in sterile saline while plunging and unplunging the driver button to help clear the distal tip of tissue and blood.

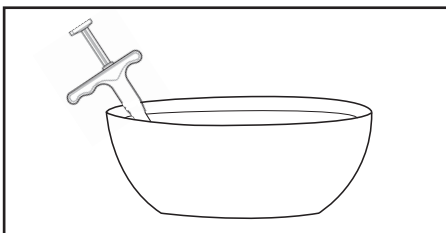
Figure 8.

Steps to Remove a Device Fragment

1. In the event the device breaks and a fragment is in the patient, assess to determine if it can be located and removed without damage to the surrounding anatomy. Use caution not to disturb the operative site until it is located to avoid pushing it further from view.
2. If the fragment is in direct view it should be grasped with a fine needle driver or forceps and returned to the manufacturer along with the device for investigation.
3. If the fragment is not in direct view, use of x-ray studies may help to localize its location. Radiologic assistance should be utilized to guide dissection and retrieve the fragment, minimizing injury to surrounding tissue.

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WARRANTY

Boston Scientific Corporation (BSC) warrants that reasonable care has been used in the design and manufacture of this instrument. **This warranty is in lieu of and excludes all other warranties not expressly set forth herein, whether express or implied by operation of law or otherwise, including, but not limited to, any implied warranties of merchantability or fitness for a particular purpose.** Handling, storage, cleaning and sterilization of this instrument as well as other factors relating to the patient, diagnosis, treatment, surgical procedures and other matters beyond BSC's control directly affect the instrument and the results obtained from its use. BSC's obligation under this warranty is limited to the repair or replacement of this instrument and BSC shall not be liable for any incidental or consequential loss, damage or expense directly or indirectly arising from the use of this instrument. BSC neither assumes, nor authorizes any other person to assume for it, any other or additional liability or responsibility in connection with this instrument. **BSC assumes no liability with respect to instruments reused, reprocessed or resterilized and makes no warranties, express or implied, including but not limited to merchantability or fitness for a particular purpose, with respect to such instruments.**

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Urological Guidewires

Rx ONLY

Caution: Federal Law (USA) restricts this device to sale by or on the order of a physician.

WARNING

Contents supplied STERILE using an ethylene oxide (EO) process. Do not use if sterile barrier is damaged. If damage is found, call your Boston Scientific representative.

For single use only. Do not reuse, reprocess or resterilize. Reuse, reprocessing or resterilization may compromise the structural integrity of the device and/or lead to device failure which, in turn, may result in patient injury, illness or death. Reuse, reprocessing or resterilization may also create a risk of contamination of the device and/or cause patient infection or cross-infection, including, but not limited to, the transmission of infectious disease(s) from one patient to another. Contamination of the device may lead to injury, illness or death of the patient.

After use, dispose of product and packaging in accordance with hospital, administrative and/or local government policy.

DEVICE DESCRIPTION

Urological Guidewires are used to provide a foundation for diagnostic and therapeutic catheters. These guidewires are offered in a variety of sizes (diameters and lengths), tip designs, coatings and stiffness.

Sensor™ Nitinol Wires with Hydrophilic Tip are designed with a 5 cm hydrophilically coated urethane tip to facilitate access and PTFE-coated body to aid in handling and manipulation of the guidewire.

Zebra™ Guidewires are designed with a radiopaque distal spring tip encapsulated in a flexible outer jacket and are striped for enhanced visualization through the scope. A Torque Vise is provided within the package for Zebra Guidewires for ease of handling.

PTFE Guidewires are available for use in standard urological procedures. PTFE coating facilitates trackability of devices as well as advancement through the urinary system

Amplatz Super Stiff™ Guidewires are constructed of a flat wire which produces a larger inner core, increased stiffness and added strength for easier access and better stability during catheter placement and exchange and are ideal for use in severely obstructed tracts.

INTENDED USE/INDICATIONS FOR USE

Urological Guidewires are intended to facilitate the placement of endourological instruments during diagnostic or interventional procedures. These guidewires are not intended for coronary artery, vascular or neurological use.

CONTRAINDICATIONS

None known.

WARNINGS

- A thorough understanding of the technical principles, clinical applications, and risks associated with the use of guidewires is necessary before using this product. Use of this device should be restricted to use by or under the supervision of physicians trained in urologic endoscopic procedures. Care should be exercised to prevent perforation or trauma of the linings and associated tissue, channels or ducts.
- Failure to abide by the following warnings might result in damage to the channel or duct, abrasion of the hydrophilic coating, release of plastic fragments from the guidewire, damage to or breakage/separation of the guidewire, that may necessitate intervention.

- Use caution if using with a metal needle or cannula. If the guidewire is being used with a metal cannula or needle and the guidewire needs to be withdrawn, remove the guidewire and metal cannula/needle as a unit, to reduce potential damage to the wire. If a needle is used for initial placement, a plastic entry needle is recommended when using the guidewire. Extreme caution should be observed when used with a one-wall puncture style needle.
- Use extreme caution when using a laser, making sure to avoid contact with the guidewire. Direct contact may cause damage to the wire and/or sever the wire.
- Do not reshape the guidewire by any means. Attempting to reshape the wire may cause damage, resulting in the release of wire fragments to the urinary system.
- When exchanging or withdrawing a catheter over the guidewire, secure and maintain the guidewire in place under fluoroscopy to avoid unexpected guidewire advancement. Otherwise damage to the urinary channel by the wire's tip may occur.
- Manipulate the guidewire slowly and carefully in the urinary system while confirming the behavior and location of the wire's tip under fluoroscopy. Excessive manipulation of the guidewire without fluoroscopic confirmation may result in perforation or trauma of the linings or associated tissues, channels or ducts. If any resistance is felt or if the tip's behavior and/or location seems improper, **STOP** manipulating the guidewire and/or the catheter and determine the cause by fluoroscopy. Failure to exercise proper caution may result in bending, kinking, separation of the guidewire's tip, damage to the catheter, or damage to the urinary system. If necessary, remove the guidewire and ancillary device or scope as a complete unit to avoid complications.
- Do not attempt to use the guidewire if it has been bent, kinked or damaged. Use of a damaged wire may result in damage to the linings and associated tissue, channels or ducts or release of wire fragments into the urinary system.
- A retrieving device, such as a gripper or basket forceps, should only be used after the guidewire has been removed from the patient's channel or duct. Using a retrieving device while the guidewire is in place may cause the guidewire to break.

PRECAUTIONS

- Do not use this product without reading and understanding the complete instructions enclosed herein.
- The entire operation should be carried out in a sterile field.
- Product is sterile in an unopened and undamaged unit package. Do not use if the unit package or the guidewire is broken, damaged or soiled. Return any defective product to Boston Scientific. All guidewires should be disposed of safely and properly after use, following local regulations for medical waste management.
- When using a drug or a device concurrently with a wire, the operator should have a full understanding of the properties/ characteristics of the drug or device so as to avoid damage to the wire.
- For hydrophilic wires only, the surface of the guidewire is not lubricious unless it is wet. Before taking it out of its holder and inserting it through a catheter, fill the holder and the catheter with physiological saline solution.
- The guidewire should be advanced through the scope using short, deliberate 2-3 cm movements to prevent inadvertent damage to the device or patient.
- When reinserting the guidewire back into the holder, take care not to damage the wire's coating with the edge of the holder.

- Do not use a metal torque device with the guidewire. Use of a metal torque device may result in damage to the wire. Also do not slip a tightened up torque device over the wire, as this may result in damage to the wire.
- Due to variations of certain catheter tip inner diameters, abrasion of the hydrophilic coating may occur during manipulation. If any resistance is felt during introduction of the catheter, it is advisable to stop using such catheters.
- After removal from the patient's urinary system, and prior to reinserting it into the same patient during the same catheterization, hydrophilic guidewires should be rinsed in a bowl full of physiological saline solution. Use of alcohol, antiseptic solutions or other solvents must be avoided because they may adversely affect the surface of the guidewire.

ADVERSE EVENTS

Complications which can result from the use of guidewires in urological applications include:

- Perforation of the urinary tract
- Acute Bleeding
- Hemorrhage
- Tissue Trauma
- Edema
- Foreign Object in Body
- Infection
- Hemoglobinuria
- Peritonitis
- Ureter Avulsion

DIRECTIONS FOR USE

Prior to Use

Carefully examine the unit to verify that neither the contents nor the sterile package has been damaged in shipment. Guidewires are delicate instruments and should be handled carefully. Prior to use and when possible during the procedure, inspect the guidewire carefully for bends or kinks which may have occurred.

Hydrophilic Coated Guidewires only

Before removing the hydrophilic guidewire from the protective hoop (by its distal end), inject physiological saline solution into the hub end of the holder in order to completely wet the surface of the guidewire. Use of alcohol, antiseptic solutions or other solvents must be avoided, as they may adversely affect the surface of the guidewire. Do not wipe guidewire with dry gauze.

Before using, inspect the guidewire for separation of the wound coil spring.

1. Remove the guidewire from protective hoop. Save the hoop to store the guidewire if it will be used again during the same procedure on the same patient.
2. Before using, inspect the guidewire for the following:
Wire rotation (see Step 4 in Recommended Technique)
Roughness or abrasion at the tip
Kinking along the length of the guidewire

RECOMMENDED TECHNIQUE

Introduction

1. The guidewire may be introduced into the patient in any of the following manners:

Ureteroscopic

- a. Place the guidewire via a scope into the ureter first, to gain initial access before placing a catheter over the guidewire.

- b. Preload a catheter over the guidewire and place as a complete unit into the ureter.
- c. Back-load the guidewire through a preplaced catheter.

Percutaneous

- a. Once the percutaneous access needle has been placed into the desired portion of the collecting system, advance the guidewire through the needle to the desired location.

Warning: Do not withdraw the guidewire through a metal cannula or needle. If the guidewire needs to be withdrawn, remove the guidewire and metal cannula or needle as a unit. Failure to do so may result in destruction and/or separation of the outer polymer jacket requiring retrieval. If a needle is used for initial placement, a plastic entry needle is recommended when using the guidewire. Extreme caution should be observed when used with a one-wall puncture style needle.

Note: Keep at least 5 cm of the wire extended out of the proximal end of the scope or catheter during the introduction at all times.

Steering

2. If steerability is desired, thread the proximal end of the guidewire through the tapered end of a “torque vise” (such as the Boston Scientific Torque Vise). Tighten the “torque vise” on the guidewire until it firmly grasps the guidewire.

Note: For the optimal distal guidewire response, the “torque vise” should be positioned as close to the catheter hub as possible.

Advancement

3. To tighten, hold the proximal end of the “torque vise” stationary. Slide the distal, tapered end of the “torque vise” back, pushing the two ends together, tightening the grip of the “torque vise” on the guidewire.
4. To steer the guidewire, slowly rotate the “torque vise” in either direction. Continue to rotate until the distal tip is oriented to the selected direction and advance.
5. To release the “torque vise”, reverse steps 2 & 3 above.
6. Carefully advance guidewire to the desired position.
7. Fluoroscopically confirm the position of the distal radiopaque tip of the guidewire.

Withdrawal

8. Carefully and slowly withdraw the guidewire from the patient to avoid any kinking.

HOW SUPPLIED

Handling and Storage

Store in a cool, dry, dark place. Do not use if package is opened or damaged. Use the device prior to the “Use by” date noted on the product label. Do not use if labeling is incomplete or illegible.

WARRANTY

Boston Scientific Corporation (BSC) warrants that reasonable care has been used in the design and manufacture of this instrument. **This warranty is in lieu of and excludes all other warranties not expressly set forth herein, whether express or implied by operation of law or otherwise, including, but not limited to, any implied warranties of merchantability or fitness for a particular purpose.** Handling, storage, cleaning and sterilization of this instrument as well as other factors relating to the patient, diagnosis, treatment, surgical procedures and other matters beyond BSC’s control directly affect the instrument and the results obtained from its use. BSC’s obligation under this warranty is limited to the repair or replacement of this instrument and BSC shall not be liable for any incidental or consequential loss, damage or expense directly or indirectly arising from the use of this instrument. BSC neither assumes, nor authorizes any other person to assume for it, any other or additional liability or responsibility in connection with this instrument. **BSC assumes no liability with respect to instruments reused, reprocessed or resterilized and makes no warranties, express or implied, including but not limited to merchantability or fitness for a particular purpose, with respect to such instruments.**

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All rights reserved.