

SHENZHEN BONA PHARMA TECHNOLOGY CO., LTD
BONA SMART PLANT, NO. 6, GAOXINYUAN GUANSHENG 1ST ROAD, LUHU COMMUNITY,
GUANHU SUBDISTRICT, LONGHUA DISTRICT, SHENZHEN CITY, GUANGDONG PROVINCE, P.R.
CHINA
Tel: +86 755 2699 4208 Fax: +86 755 2652 7737
Website: www.bona-cn.com

INSTRUCTIONS FOR USING

1. Object:


Spray pumps 18Din G#, 18/415 G#, 18/415 K#, 18Din L#

2. How to use: (See below picture)

Step 1 Take off the dust cap and spacing ring (if have).

Step 2 Make sure the orifice of actuator facing toward the affected part and container upright as well, press actuator quick and deeply, medicines in bottle will be spray out to affected part.

Step 3 Clean the actuator and cover the dust cap.

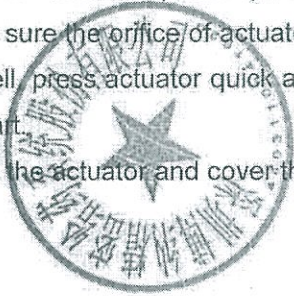

«ШЭНЬЧЖЭНЬ БОНА ФАРМА ТЕХНОЛОДЖИ КО., ЛТД.»

BONA SMART PLANT, №6, Gaoxinyuan Guansheng 1st road, Luhua Community, субрайон
Гуанху, район Лонгхуа, Шэньчжэнь, провинция Гуандун, Китайская народная республика

Тел.: +86 755 2699 4208 Факс: +86 755 2652 7737

Веб-сайт: www.bona-cn.com

ИНСТРУКЦИИ ПО ИСПОЛЬЗОВАНИЮ

1. Объект:

Распыляющие помпы 18Din G#, 18/415 G#, 18/415 K#, 18Din L#

2. Как использовать:

Этап 1 Снять пылезащитный колпачок и распорное кольцо (если имеется)

Этап 2 Удостовериться, что отверстие активатора направлено на поврежденный участок, а контейнер расположен вертикально, нажать на активатор быстро и сильно, лекарственные средства во флаконе распылятся на поврежденный участок тела.

Этап 3 Очищают активатор и закрывают пылезащитным колпачком.

(печать)

