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
[illegible][illegible]

HARMONIC SYNERGY®

Dissecting Hook, Curved Blade and Torque Wrench

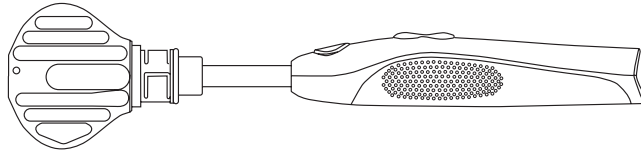Crochet dissecteur, lame courbe et clé dynamométrique
Dissektionshaken, gebogener scharfer Spatel
und Klingenschlüssel
Uncino per dissezione, lama curva e chiave torsionometrica
Gancho Dissector, Lâmina Curva e Chave de Torção
Disector de gancho, lâmina curva y llave de apriete
Dissectiehaak, gebogen blad en torsiesleutel
Dissektionshook, buet blad og momentnøgle
Dissekaatiokoukku, kaareva terä ja momenttiavain
Αγκιστρο Αποκόλλησης, κυρτή λεπίδα
και κλειδί ροπής στρέψης
Dissektionskrok, böjt blad och momentnyckel
Haczyk do cięcia, zakrzywione ostrze
i klucz dynamometryczny
Bontóhorog, hajlított penge és nyomatékkulcs
Řezací hák, zakřivená čepel a momentový klíč
Prerezávací háčik, zakřivená čepel a momentový klúč
Dissekeringskrok, buet blad og momentnøkkel

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





Important: This package insert is designed to provide Instructions For Use of the HARMONIC SYNERGY® instruments and Gray Torque Wrench. It is not a reference to surgical techniques.

HARMONIC, HARMONIC SYNERGY, and ULTRACISION are trademarks of Ethicon Endo-Surgery.

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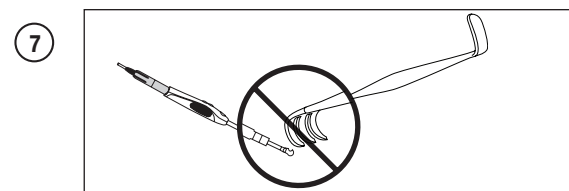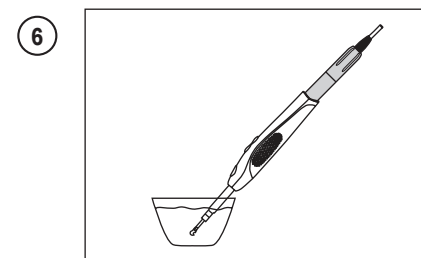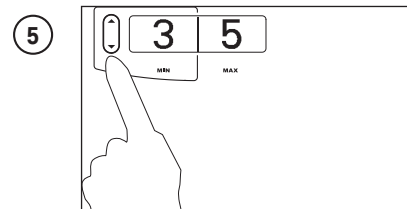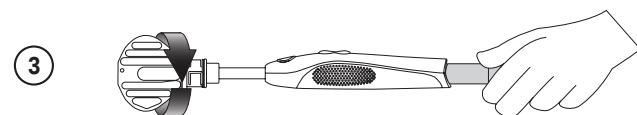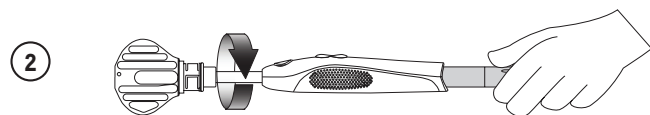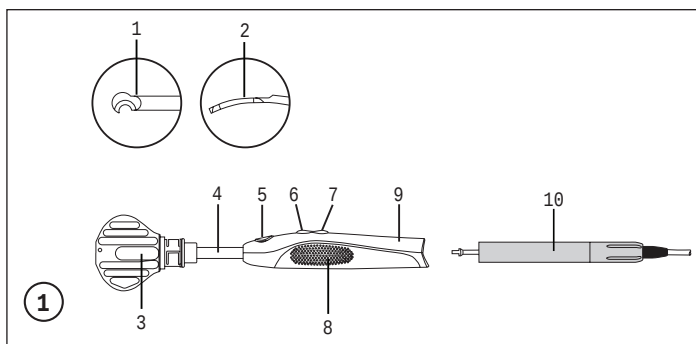[illegible]

HARMONIC, HARMONIC SYNERGY and ULTRACISION are trademarks of Ethicon Endo-Surgery.

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ETHICON
PART OF THE *Johnson & Johnson* FAMILY OF COMPANIES

Ethicon Endo-Surgery, LLC
475 Calle C
Guaynabo, PR 00969 USA


© Ethicon Endo-Surgery, Inc. 2018

Rev. 2018-09-31

P40575P09

ETHICON

HARMONIC FOCUS® Shears + Adaptive Tissue Technology
Ciseaux HARMONIC FOCUS® + Technologie Tissulaire Adaptive
HARMONIC FOCUS® Koagulationsschere + Adaptive Tissue Technology
Forbici HARMONIC FOCUS® + tecnologia adattativa al tessuto
Tesoura HARMONIC FOCUS® + Tecnologia Tecidual Adaptável
Cizallas HARMONIC FOCUS® + Tecnología adaptativa tisular
HARMONIC FOCUS®-schaar + adaptieve weefseltechnologie
HARMONIC FOCUS® saks + adaptiv vævsteknologi
HARMONIC FOCUS®-sakset + mukautuva kudostekniikka
Ψαλίδι HARMONIC FOCUS® + Προσαρμοζόμενη τεχνολογία ιστού
HARMONIC FOCUS® sax + adaptiv vävnadsteknologi
Nożyce HARMONIC FOCUS® + adaptacyjna technologia tkankowa
HARMONIC FOCUS® metsző + adaptív szövetechnológia
Nůžky HARMONIC FOCUS®+ s adaptivní tkáňovou technologií
Nožnice HARMONIC FOCUS®+ s adaptivnou tkanivovou technológiou
HARMONIC FOCUS® skjærere + adaptiv vevs-teknologi
HARMONIC FOCUS® Makas + Adaptif Doku Teknolojisi
Ножницы с технологией адаптации к тканям HARMONIC FOCUS®
Pensă foarfecă HARMONIC FOCUS® + tehnologia de adaptare la țesut
Gunting HARMONIC FOCUS® + Teknologi Jaringan Adaptif
Kéo HARMONIC FOCUS® + Công nghệ Mô thích ứng
Kohanduva koetehnoloogiaga käärid HARMONIC FOCUS®+
HARMONIC FOCUS® škéres + adaptivā audu tehnoloģija
HARMONIC FOCUS® + žirkles su prisitaikancio audinio technologija
HARMONIC FOCUS® ножици + адаптивна към тъкани технология
Škare HARMONIC FOCUS®+ s prilagodljivom tehnologijom tkiva
Škarje HARMONIC FOCUS®+ z adaptivno tehnologijo rezanja tkiva
HARMONIC FOCUS® ножици + технологија за приспособување на ткивото
HARMONIC FOCUS® маказе + адаптивна технологија за ткиво
HARMONIC FOCUS®+ маказе са адаптивном технологијом ткива
HARMONIC FOCUS® +




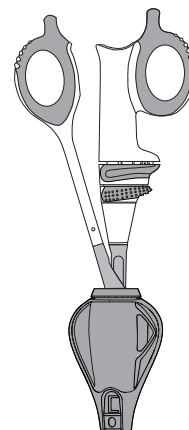
Please read all information carefully.

Failure to properly follow the instructions may lead to serious surgical consequences, such as failure to ligate.


Important: This package insert is designed to provide instructions for use of the HARMONIC FOCUS® Shears + Adaptive Tissue Technology. It is not a reference to ligation techniques.

HARMONIC, HARMONIC FOCUS and ULTRACISION are trademarks of Ethicon Endo-Surgery.


Instructions, Instrucciones, Gebrauchsanweisung, Istruzioni, Instruções, Instrucciones, Instruções, Brugervejledning, Ohjeet, Oðnysic, Anvisningar, Instrukcije, Utasítások, Návod k použití, Návod na použitie, Instruksjoner, Talimatlar, Инструкции, Instrucțiuni, Instruksi, Hướng dẫn, Kasutusjuhend, Instruccijas, Instrukcija, Инструкции, Upute, Navodila, Упатства, Упутства, Упутства,

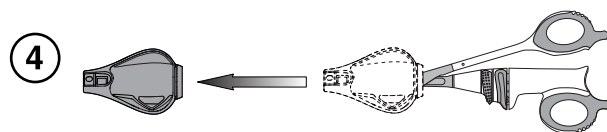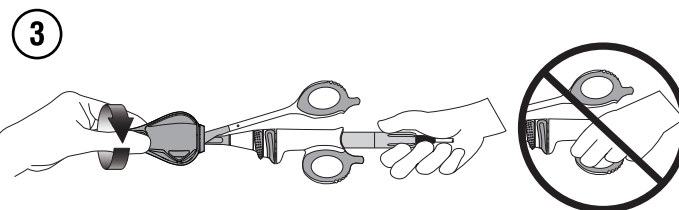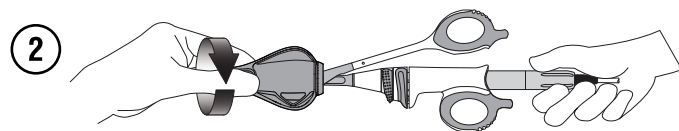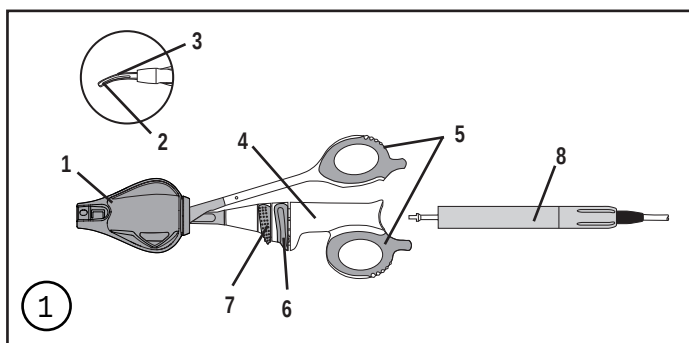
5

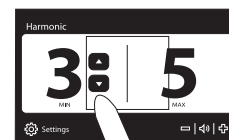
6

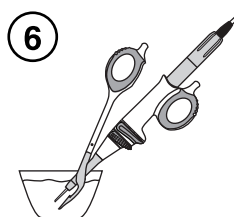7

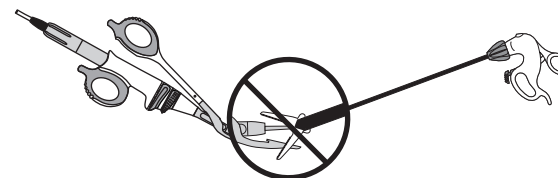8

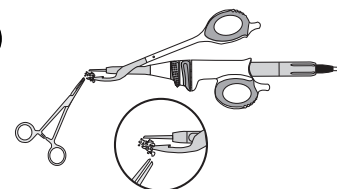[illegible][illegible]

[illegible][illegible]

STERILE	EO
---------	----

[illegible]

	<p>Lot N° de lot Ch.-B. Lotto N° do lote N° de lote Lotnr. Parti Erän koodi J J. "Üj üüÜÜg Batchnummer Numer partii produkcyjnej Tétel Gärge Gärğa Lot</p>	<p>Lot u j k d w Lot Lot L.¹ s[n xuít Partii Partija Partija 1 0 u k d H 0 Serija Serija u j k d 0 1 0 u k d 0 1 u s " h j k d j "</p>
	<p>Use Until Date À utiliser avant Verw. bis Utilizzare entro Validade A utilizar antes de Gebruik vóór Holdbar til angivne dato Käytettävä viimeistään j j d u s e s " e s ü ü Ü e y c j s Ü q z Använd före Koniec okresu przydatności do użytku A feltüntetett dátumig használható fel Použit do data Použiteľnô do</p>	<p>Brukes før-dato Son Kullanma Tarihi u k s s s s s s s h j k d A se utiliza înainte de data Digunakan Hingga Tanggal S d ng L j n ng "y Kasutamise lõpptähtaeg Deriguma termiäg Tinka naudoti iki R j f s d s e s 0 H s H 0 s 0 Datum „Uporabiti do“ Rok uporabe r 0 h j k f s s s j w f H s r 0 k s d s f h j s 0 u s s 0 k s 0 0 " 0 s k s f h j k s d f H s H 0 k s d 0</p>

[illegible][illegible]

[illegible]

P40730P01

REF
HAR9F

Ethicon Endo-Surgery (Europe) GmbH
Hummelsbuetteler Steindamm 71
22851 Norderstedt
GERMANY

ETHICON ENDO-SURGERY, INC.
4545 Creek Road
Cincinnati, OH 45242-2839 USA
1-877-ETHICON

ETHICON


ETHICON ENDO-SURGERY, LLC ©EES, LLC 2014
475 Calle C
Guaynabo, PR 00969 USA

Rev. 2014-02

P40730P01

HARMONIC FOCUS® Long Shears + Adaptive Tissue Technology
 Ciseaux longs HARMONIC FOCUS® + technologie tissulaire adaptative
 HARMONIC FOCUS® Lange Koagulationsschere + Adaptive Gewebetechnologie
 Forbici lunghe HARMONIC FOCUS® + tecnologia adattativa al tessuto
 Tesoura Comprida + Tecnologia Adaptativa aos Tecidos HARMONIC FOCUS®
 Tijeras largas + Tecnología de adaptación tisular HARMONIC FOCUS®
 HARMONIC FOCUS® lange schaar + adaptieve weefseltechnologie
 HARMONIC FOCUS® lang saks + adaptiv vævsteknologi
 HARMONIC FOCUS® pitkä leikkausinstrumentti + mukautuva kudosteknologia
 Μακρύ ψαλίδι HARMONIC FOCUS® + με Προσαρμόζόμενη τεχνολογία ιστού
 HARMONIC FOCUS® lång sax + adaptiv vävnadsteknik
 Długie nożyce HARMONIC FOCUS® + Adaptacyjna technologia tkankowa
 HARMONIC FOCUS® hosszú olló + intelligens szövét-technológia
 Dlouhé nůžky HARMONIC FOCUS® s adaptivní tkáňovou technologií
 Dlhé nožnice HARMONIC FOCUS® s adaptívnou tkanivovou technológiou
 HARMONIC FOCUS® Lange skjærere + Adaptiv vevsteknologi
 HARMONIC FOCUS® Uzun Makas + Adaptif Doku Teknolojisi
 Длинные ножницы с технологией адаптации к тканям HARMONIC FOCUS® +
 Pensă foarfecă lungă HARMONIC FOCUS® + cu tehnologie de adaptare la țesut
 Teknologi Jaringan Adaptif + Gunting Besar Panjang HARMONIC FOCUS®
 Dao cắt HARMONIC FOCUS® Long Shears + Công nghệ Mô thích ứng
 HARMONIC FOCUS® pikad käärid + adaptiivne koetehnoloogia
 HARMONIC FOCUS® garás škéres + adaptiva audu tehnológia
 HARMONIC FOCUS® ilgos žirkklės su audinių pritaikymo technologija
 Дълги ножници HARMONIC FOCUS® + адаптивна технология за тъкани
 HARMONIC FOCUS® duge škare + tehnologija adaptivnog tkiva
 Dolge skarje HARMONIC FOCUS® + prilagodljiva tkivna tehnologija
 HARMONIC FOCUS® долги ножници + технология за приспособуване на ткивото
 HARMONIC FOCUS® дуге маказе + адаптивна технология за ткиво
 HARMONIC FOCUS® дугачке маказе са адаптивном технологијом ткива
 HARMONIC FOCUS® +

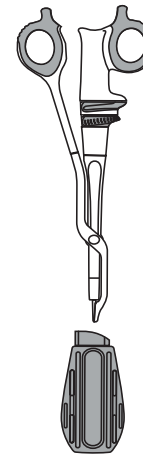Please read all information carefully.

Failure to properly follow the instructions may lead to serious surgical consequences, such as failure to ligate.










Important: This package insert is designed to provide instructions for use of the HARMONIC FOCUS® Long Shears + Adaptive Tissue Technology. It is not a reference to ligation techniques.

HARMONIC, HARMONIC FOCUS and ULTRACISION are trademarks of Ethicon Endo-Surgery.

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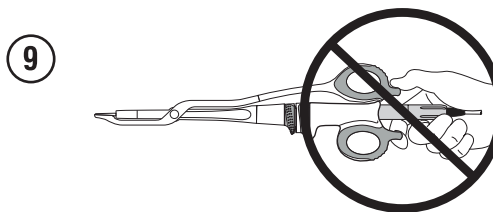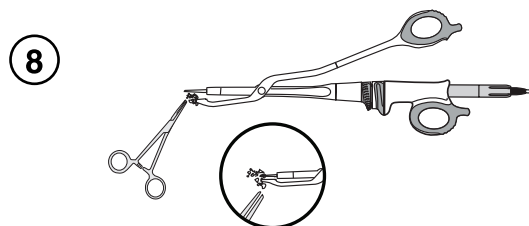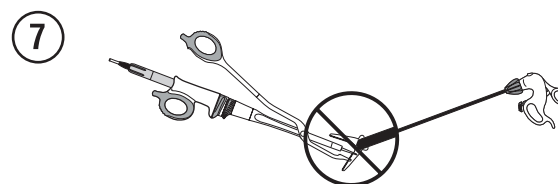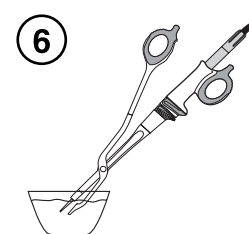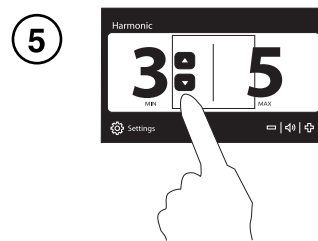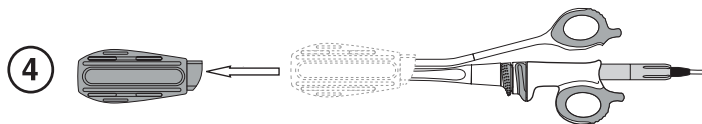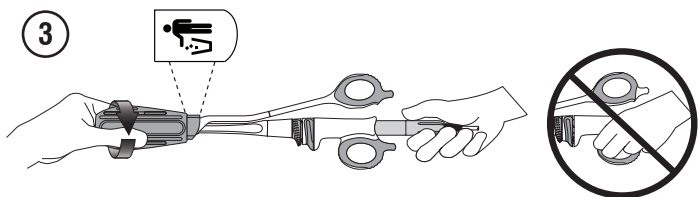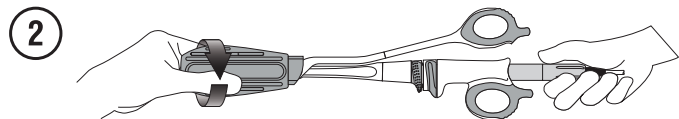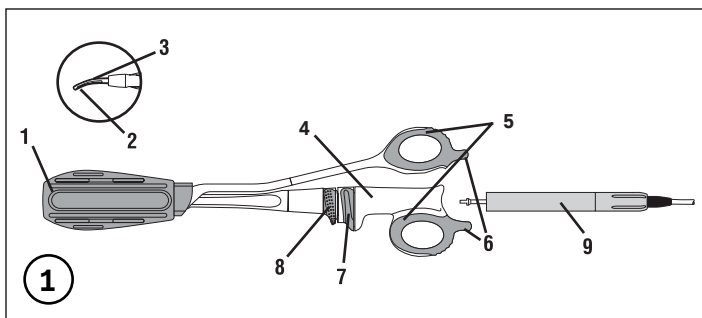
REP

[illegible]








HARMONIC ACE® 5 mm Diameter Shears + Adaptive Tissue Technology

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

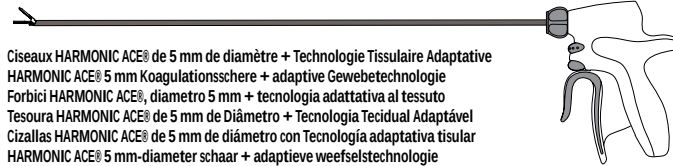
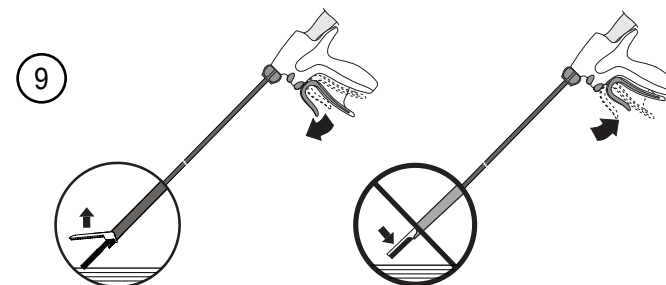
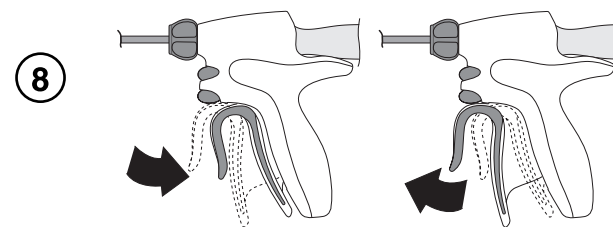
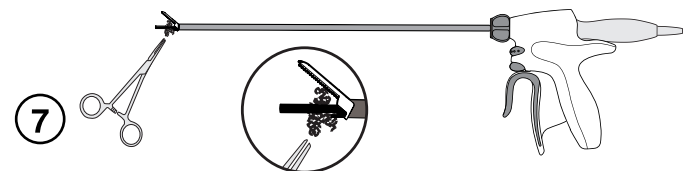
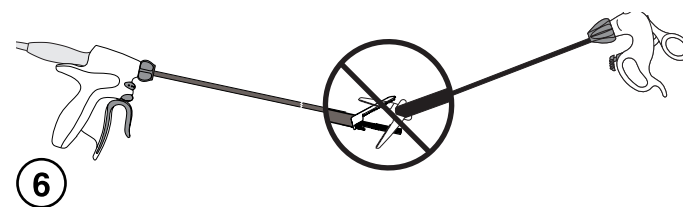
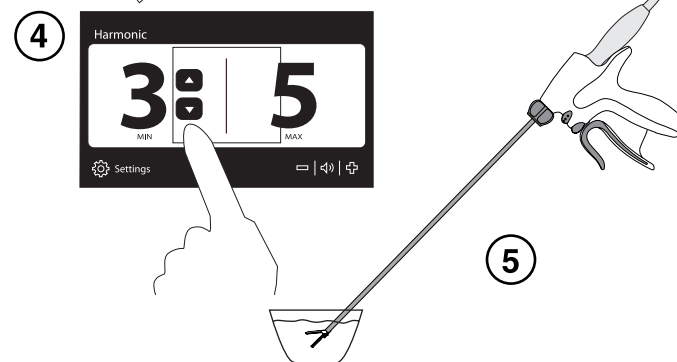
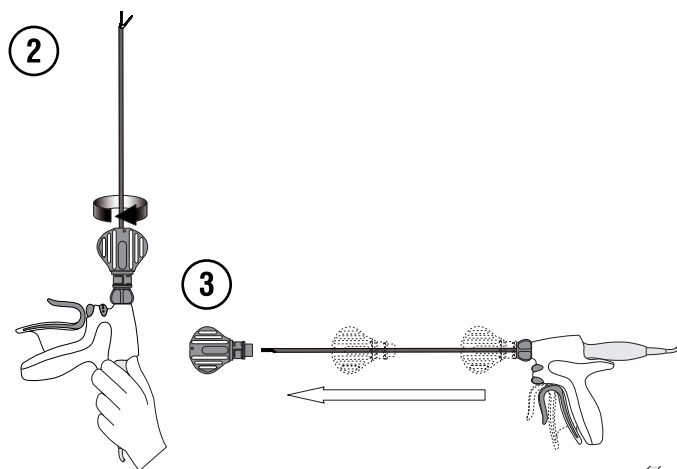
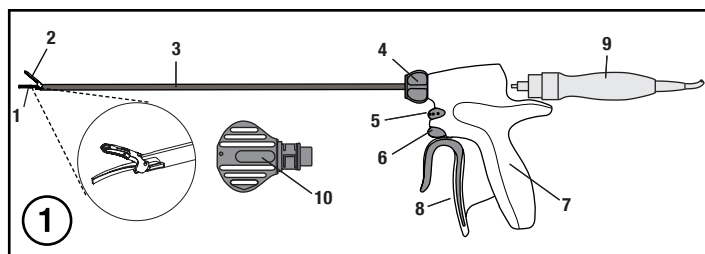
Please read all information carefully.

Failure to properly follow the instructions may lead to serious surgical consequences.

Important: This package insert is designed to provide instructions for use of the HARMONIC ACE®+7, 5 mm Diameter Shears with Advanced Hemostasis. It is not a reference to surgical techniques.

HARMONIC® and HARMONIC ACE® are trademarks of Ethicon Endo-Surgery.

Instructions, Instrucciones, Gebrauchsanweisung, Istruzioni, Instruções, Instrucciones, Instructies, Brugervejledning, Ohjeet, Oðnyisc, Anvisningar, Instrukcje, Utasítások, Návod k použití, Návod na použitie, Instruksjoner, Talimatlar, Инструкции, Instrucţiuni, Instruksi, Hướng dẫn, Kasutusjuhend, Instrukcijas, Instrukcija, Инструкции, Upute, Navodila, Упутства, Упутства, Упутства.

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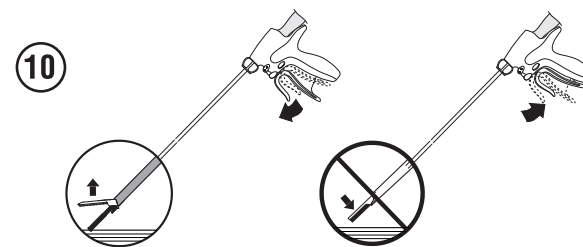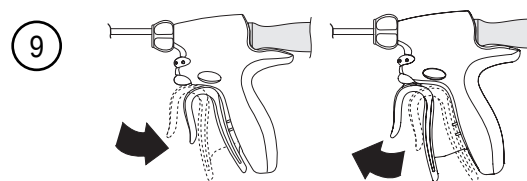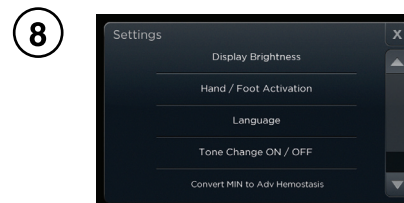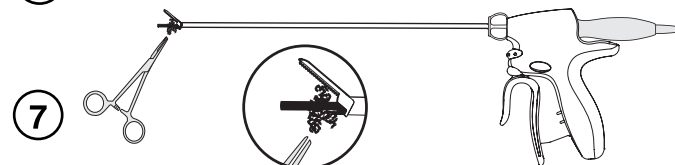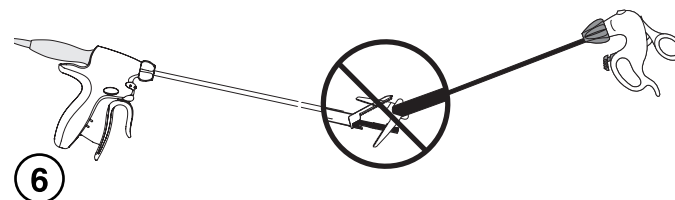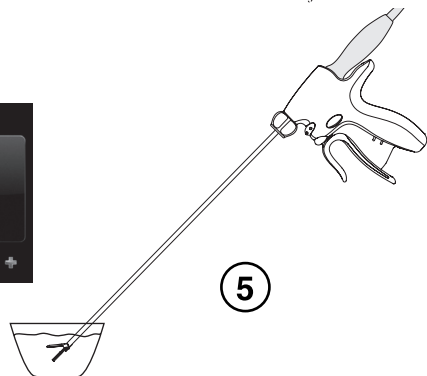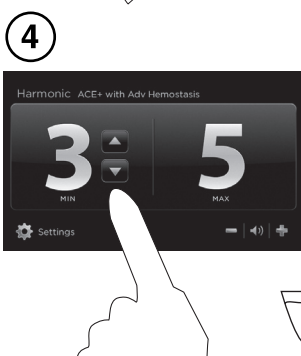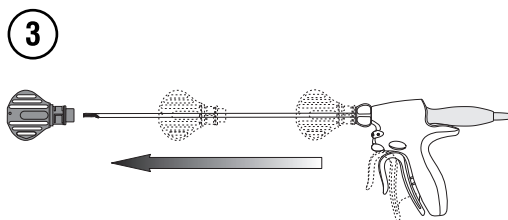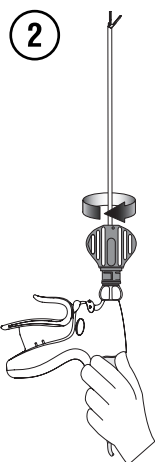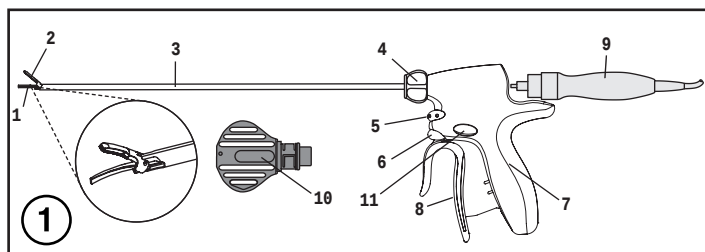Please read all information carefully.








Failure to properly follow the instructions may lead to serious surgical consequences.

Important: This package insert is designed to provide instructions for use of the HARMONIC® HD 1000i Shears. It is not a reference to surgical techniques.

ETHICON

PART OF THE **Johnson & Johnson** FAMILY OF COMPANIES

Instructions, Instrucciones, Gebrauchsanweisung, Istruzioni, Instruções, Instrucciones, Instruções, Brugervejledning, Ohjeet, Oðnýðir, Anvisningar, Instrukcje, Utasítások, Návod k použití, Návod na použitie, Instruksjoner, Talimatlar, Инструкции, Instrucțiuni, Instruksi, Hướng dẫn, Kasutusjuhend, Instrukcijas, Instrukcija, Инструкции, Upute, Navodila, Упатства, Упутства, Упутства,

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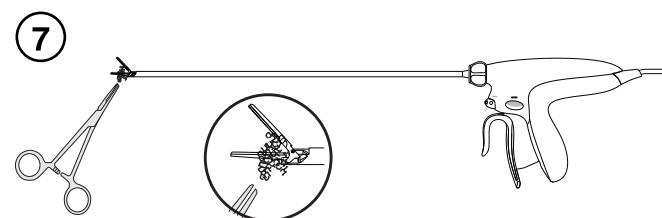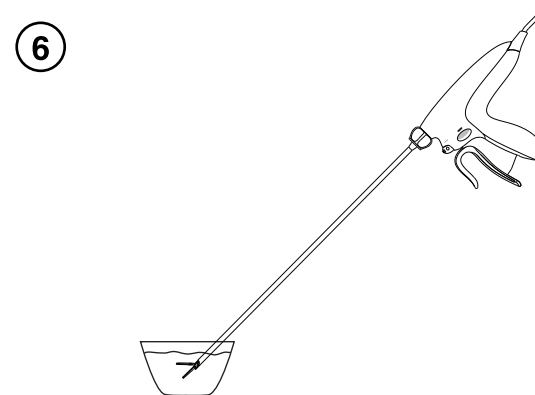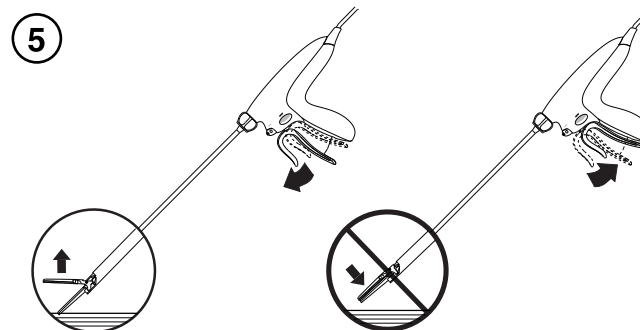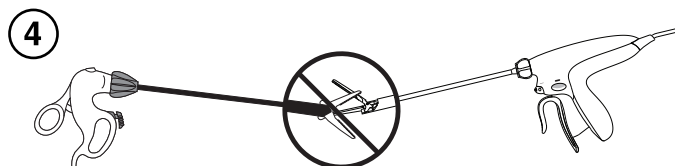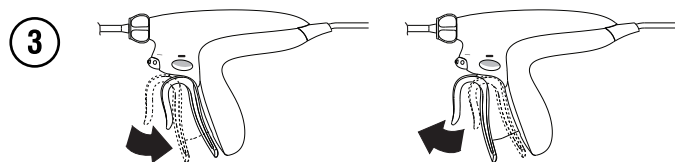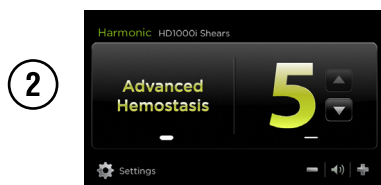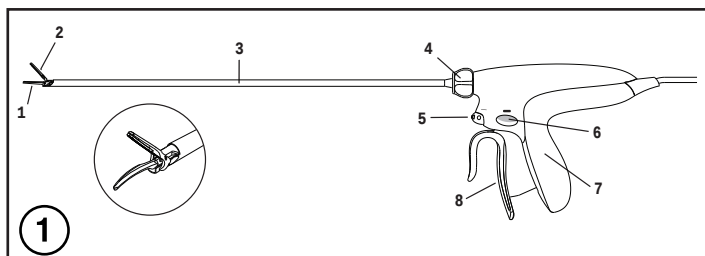ETHICON

PART OF THE *JOHNSON & JOHNSON* FAMILY OF COMPANIES


Ethicon Endo-Surgery, LLC ©Ethicon Endo-Surgery, Inc. 2017
475 Calle C
Guaynabo, PR 00969 USA

Rev. 2017-02-28

P000162P02

ETHICON Enseal®

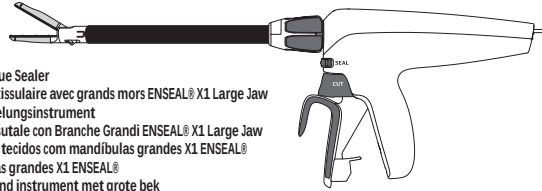ENSEAL® X1 Large Jaw Tissue Sealer
Dispositif de coagulation tissulaire avec grands mors ENSEAL® X1 Large Jaw
ENSEAL® X1 Makroversiegelungsinstrument
Dispositivo di chiusura tissutale con Branche Grandi ENSEAL® X1 Large Jaw
Dispositivo de selagem de tecidos com mandíbulas grandes X1 ENSEAL®
Sellador de tejido de ramas grandes X1 ENSEAL®
ENSEAL® X1 Weefseldichtend instrument met grote bek
ENSEAL® X1 vævsforsegler med store kæber
ENSEAL® X1 Large Jaw -kudosfuusolaite
Στεγανοποιητής ιστών ENSEAL® X1 μεγάλων σιαγόνων
ENSEAL® X1 tillslutningsinstrument för vävnad med stora käftar
Zgrzewarka do tkanek ENSEAL® X1 Large Jaw
ENSEAL® X1 nagy pófájú szövetlezáró
Tkáňový spojovací nástroj ENSEAL® X1 s veľkou čelistí
Nástroj na uzatváranie tkaniv ENSEAL® X1 Large Jaw
ENSEAL® X1 Stor kjevetang for vevsforsegling
ENSEAL® X1 Büyük Çene Doku Kapatıcı
Коагулятор-диссектор тканей с браншами большего размера ENSEAL® X1
Dispozitiv de etanșare a țesuturilor ENSEAL® X1 cu fălci mari
Penutup Jaringan Rahang Besar ENSEAL® X1
Dao hàn mach hong kep lớn ENSEAL® X1
Suure haaratsiga koesulgur ENSEAL® X1
ENSEAL® X1 audu nosiešanas instruments ar lielām spilēm
ENSEAL® audīnių koaguliacijos instrumentas su X1 dydžio žiotimis
ENSEAL® X1 Коагулятор-трансектор на тъкани с голяма челюст
Uređaj za spajanje tkiva s velikim čeljustima ENSEAL® X1
Pripomoček za zatesnitev tkiva ENSEAL® X1 z veliko čeljustjo
Уред за спојување ткива со широка вилица ENSEAL® X1
Инструмент ENSEAL® X1 за затварање ткива са великом челюсти
ENSEAL® X1









Please read all information carefully.

Failure to properly follow the instructions may lead to serious surgical consequences.

Important: This package insert is designed to provide instructions for use of the ENSEAL® X1 Large Jaw Tissue Sealer. It is not a reference to surgical techniques.

ETHICON
PART OF THE **JOHNSON & JOHNSON** FAMILY OF COMPANIES

Instructions, Instrucciones, Gebrauchsanweisung, Istruzioni, Instruções, Instrucciones, Instructies, Brugervejledning, Ohjeet, Oðnysic, Anvisningar, Instrukcje, Utasítások, Návod k použití, Návod na použitie, Instruksjoner, Talimatlar, Инструкции, Instrucțiuni, Instruksi, Hướng dẫn, Kasutusjuhend, Instrukcijas, Instrukcija, Инструкции, Upute, Navodila, Упутства, Упутства,

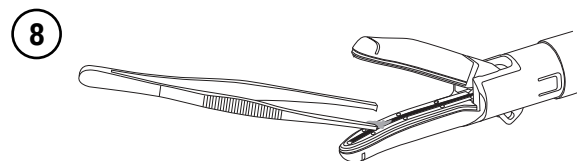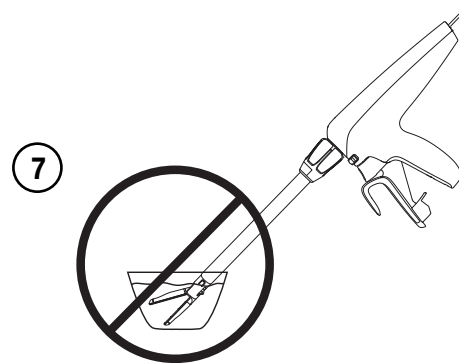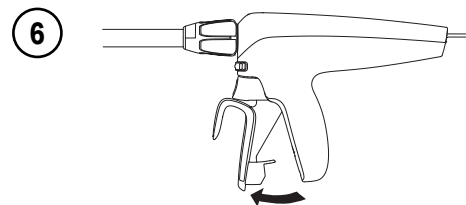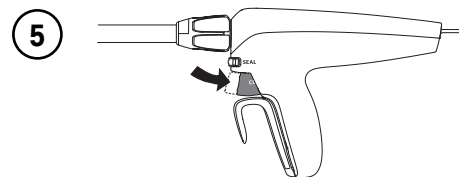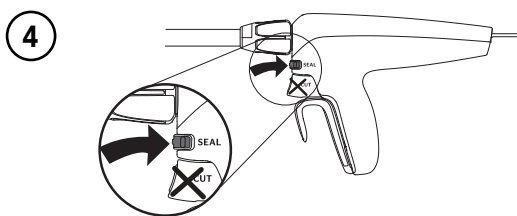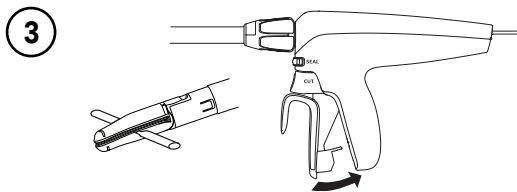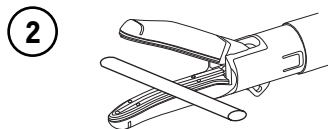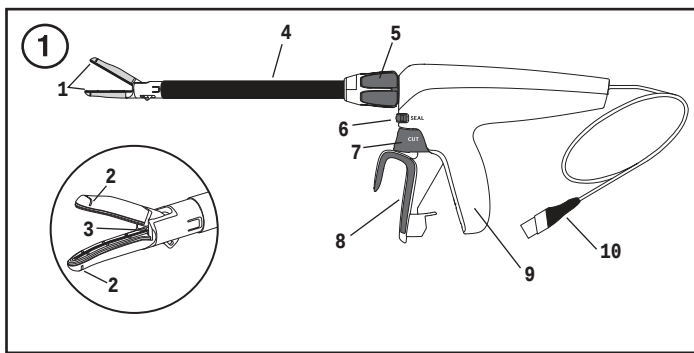	<p>Manufacturer/Date of Manufacture Fabricant/Date de fabrication Hersteller/Herstellungsdatum Fabbicante/Data di fabbricazione Fabricante/Data de Fabrico Fabricante/Fecha de fabricaci6n Fabrikant/Fabrigagedatum Producent/Produktionsdato Valmistaja/Valmistuspäivamäärä TÜÜÜüälgÜüüēg/ ē ēļēēēāü aÜÜÜüälgēg Tiliverkare/Tillverkningsdatum Product/Data produkci Gyártó/Gyártás ideje Výrobce/Datum výroby Výrobcia/Dátum výroby</p>	<p>Produsent/Produktsjonsdato Ūretitci/Ūretim Tarihi ᠮᠤᠰᠣᠳᠤᠨᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ Productor/Data fabricației Producers/Tanggal Produksi Nh~ s'n xuiłt/Ng'y s'n xuiłt Tootja/tootmiskupaapēv Rağotljs/Rağoğanas datums Gamintojas/pagaminimo data ᠮᠤᠰᠣᠳᠤᠨᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ProizvoLaI/datum proizvodnje Proizvajalec/Datum proizvodnje ᠮᠤᠰᠣᠳᠤᠨᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠮᠤᠰᠣᠳᠤ' Ōuᠯ ᠾᠥ᠋᠋᠋᠋᠋ ᠬᠤ᠋᠋᠋᠋᠋ ᠶ᠋ᠢ ᠶ᠋ ᠑ - ᠶ᠋</p>
<div>EC</div> <div>REP</div>	<p>Authorized Representative in the European Community Représentant autorisé dans la Communauté européenne Bevollmächtigter in der Europäischen Gemeinschaft Rappresentante autorizzato per la Comunità Europea Representante autorizado na Comunidade Europeia Representante autorizado en la Comunidad Europea Bevoegd vertegenwoordiger bij de Europese Gemeenschap Autoriseret repræsentant i det europæiske fællesskab Valtuutettu edustaja Euroopan yhteisön alueille Aзғуығасөйлепдес әңгілесі жән х әңгіледә Әлҗ х Үйег Тезәһүдүү Auktoriserad representant i Europeiska gemenskapen Autoryzowany przedstawiciel w Unii Europejskiej Az Eur-pai K“z’ssög meghatalmazott képviselője AutorizovanĹ zřstupce v Evropském společenstvĹ Autorizovaný zástupca EU Autoriseret representant i Det europeiske fellesskap Avrupa Topluluğunda Yetkilili Temsilci ᠮᠤᠰᠣᠳᠤᠨᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ Representant autorizat ‘n Comunitatea Europeană. Perwakilan Resmi di Komunitas Eropa ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ Volitatud esindaja Euroopa Ühenduses Pilnvarotaits pūrstūvis Eiropas Kopienū jgaliotasis atstovas Europos Bendrijoje ᠮᠤᠰᠣᠳᠤᠨᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ ᠶ᠋ᠢᠵᠡᠭᠦᠷᠲᠢᠨ Ovlažteni predstavnik u Europskoj zajednici Pooblašeni zastopnik za Evropsko skupnost fēdōñsjđ fēkjēsñsēdēq fēt cēñsfññsōkō lōjñdāqfō fēdō ñ jđqđ fēkjēsñsēdēq fē kēñsfñññsē kēd qđ</p>	

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ETHICON

PART OF THE Johnson & Johnson FAMILY OF COMPANIES


Ethicon Endo-Surgery, LLC ©Ethicon Endo-Surgery, Inc. 2020
475 Calle C
Guaynabo, PR 00969 USA

Rev. 2020-10-31

CE 0123
500564871

ETHICON Enseal®

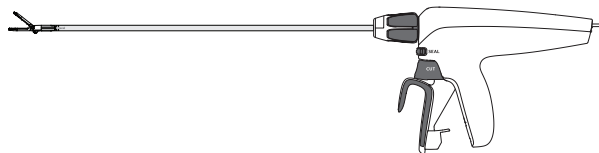en: X1 Tissue Sealer	lv: X1 audu nosiešanas instruments
bg: X1 коагулатор-трансектор на тъкани	mk: Уред за спојување ткива X1
cs: Tkáňový spojovací nástroj X1	nl: X1 Weefseldichtend instrument
da: X1-vævsforsegler	no: X1 vevsforsegler
de: X1 Versiegelungsinstrument	pl: Zgrzewarka do tkanek X1
el: Στεγανοποιητής ιστών X1	pt: Dispositivo de selagem de tecidos X1
es: Sellador de tejido X1	ro: Dispozitiv de sigilare a țesuturilor X1
et: Koesulgur X1	ru: Коагулятор-диссектор тканей X1
fi: X1 -kudosfuusolaite	sk: Nástroj na uzatváranie tkanív X1
fr: Dispositif de coagulation tissulaire X1	sl: Pripomoček za zatesnitev tkiva X1
hr: Uredaj za spajanje tkiva X1	sr: X1 инструмент за затварање ткива
hu: X1 szövetlezáró	sv: X1 tillslutningsinstrument för vävnad
id: Penutup Jaringan X1	tr: X1 Doku Kapatıcı
it: Dispositivo di chiusura tissutale X1	vi: Dao hàn mạch X1
lt: X1 audinio koaguliacijos instrumentas	zh: X1

Please read all information carefully.

Failure to properly follow the instructions may lead to serious surgical consequences.

Important: This package insert is designed to provide instructions for use of the ENSEAL® X1 Tissue Sealer. It is not a reference to surgical techniques.

ETHICON™
PART OF THE **Johnson & Johnson** FAMILY OF COMPANIES

Instructions, Инструкции, Návod k použití, Brugervejledning, Gebrauchsanweisung, Οδηγίες, Instrucciones, Kasutusjuhend, Ohjeet, Instructions, Upute, Utasítások, Instruksi, Istruzioni, Instrukcija, Instrukcijas, Упатства, Instruções, Instruksjoner, Instrukcje, Instruções, Instrucțiuni, Инструкции, Návod na použitie, Navodila, Упутства, Anvisningar, Talimatlar, Hướng dẫn,

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HARMONIC is a trademark of Ethicon Endo-Surgery.