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
Anhui Hongyu Wuzhou Medical manufacturer Co., Ltd

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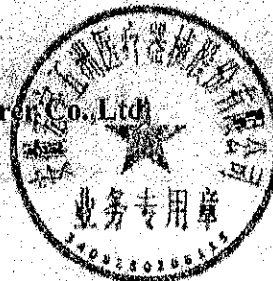ANHUI HONGYU WUZHOU MEDICAL MANUFACTURER CO.,LTD

Add: No.2, Guanyin Road, Economy Development Zone, Taihu, Anhui Province, China

INSTRUCTION FOR USE OF MEDICAL DEVICE

The disposable infusion set is for single use, Sterile, Non Toxic Non pyrogenic.
And is constructed by spike, drip chamber, tubing roller clamp, needle and etc

Direction of use:

1. Prepare the Iv container
2. Remove set from pouch and close roller clamp
3. Remove spike protector from spike
4. Insert spike into container
5. Hang container fill drip chamber by squeezing until approximate half full.
6. Remove protector from adapter
7. Open roller clamp prime set and purge air from tubing
8. Close roller clamp until rollers meets bottom of frame
9. Attach adapter to vein puncture device
10. Adjust fluid flow with roller clamp
11. Checking maintenance of proper flow rate at regular intervals.

Caution:

1. Sterile unless package has been opened or damaged
2. Do not use if tip protectors are not in place
3. Avoid freezing and excessive heat.
4. The product should be used immediately after envelope opened, discard after use.
5. Not for insertion into blood or plasma containers
- 6 For single use only. Do not resterilize or re-use.

Shelf life and storage temperature:

Production from the date of production is valid for three years , must be used within the period of validity. The batch number and expiry date was printing on the seal of individual packing.

This product should be stored in the relatibe humidity of not more than 80%, no corrosive gas , cool , dry and clean indoor.

ANHUI HONGYU WUZHOU MEDICAL MANUFACTURER CO., LTD.

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ANHUI HONGYU WUZHOU MEDICAL MANUFACTURER CO., LTD

Add: No. 2, Quanyin Road, Economy Development Zone, Taihu, Anhui Province, China

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Anhui Hongyu Wuzhou Medical Manufacturer Co., Ltd.

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Add: No.2, Quanyin Road, Economy Development Zone, Taihu, Anhui Province, China

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The disposable transfusion set with needle is for single use, Sterile, Non Toxic Non pyrogenic. And is constructed by ABS spike, flexible drip chamber with filters, Air vent, roller clamp, Y-injection port, slip or lock adapter.

Direction of use:

- 1: Along the packing tear tear small mouth
- 2: Insert the cork puncture outfit bags of puncture
- 3: Open the flow regulator, adjust the drop in the bucket surface in drops of $1/2 \sim 2/3$, to vent air from liquid road, with the regulator to adjust according to the need of blood transfusion blood transfusion speed, can be used normally

Caution:

1. Sterile unless package has been opened or damaged
2. Do not use if tip protectors are not in place
3. Avoid freezing and excessive heat.
4. The product should be used immediately after envelope opened, discard after use.
- 6 For single use only. Do not re-sterilize or re-use.

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Production from the date of production is valid for three years, must be used within the period of validity. The batch number and expiry date was printing on the seal of individual packing.

This product should be stored in the relative humidity of not more than 80%, suggested storage temperature should be -5°C to $+30^{\circ}\text{C}$, no corrosive gas, a cool, dry, well ventilated and clean room. Storage warehouse should be prevention and pest prevention measures.

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