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issued 2021-01

INGREDIENTS: PA Poliamide + medical soft silicone gel

Spycra Contact: sterile, bi-elastic soft silicone adhesive wound contact layer developed with "ReSkin Silicone Skin Technology", a unique patented adhesive technology

Product description

Elastic, one-sided silicone adhesive wound contact layer with self-adhering properties made of soft and flexible polyamide fibers. Repositionable and rounded edges to stay better in place. When being used on exuding wounds, the wide scaled pore structure of the wound contact layer allows an optimal vertical migration of wound fluid towards an absorbent secondary dressing. Does not absorb.

ReSkin Silicone Skin Technology

"ReSkin Silicone Skin Technology" is a patented combination of elastic fibers and a soft, self-adhering silicone layer on one side of the dressing. "ReSkin Silicone Skin Technology" is integrated in Spycra Contact and prevents the dressing from sticking to the wound. The soft and perfect adherence on the skin allows changes without loss of epidermis and dressing changes as pain free as possible. When used on fragile skin, Spycra Contact secures a soft and perfect adherence on the skin allowing dressing changes without pain and loss of epidermis. Thanks to the "Reskin Silicone Skin Technology", Spycra Contact can be used on exuding or bleeding wounds whereby the comfort of the patient is crucial (painless removal of the dressing and optimal support of the wound healing processes).

Mode of action

The wide scaled pore structure of the dressing allows excess of wound fluid to pass vertically towards an outer absorbent dressing. The soft silicone layer seals off the wound edges avoiding maceration and keeping the surrounding skin dry and healthy limiting the risk of maceration. The outer absorbent dressing needs to be changed when saturated limiting the risk of maceration. Spycra Contact may be left in place up to several days with a maximum of 14 days depending on the condition of the wound and the amount of wound fluid. The soft adherence secures no loss of epidermis during dressing changes and makes these changes as pain free as possible.

Intended use

Can be used as wound contact layer to take care of heavy exuding wounds such as skin tears, abrasions, surgical incisions, partial thickness burns, fixation of transplants, diabetic ulcers, venous and arterial ulcers, pressure sores, Negative Pressure Wound Therapy,... Spycra Contact can also be used as a protective layer on non-exuding wounds and on areas with fragile skin. Where clinically indicated, treatments such as topical steroids, topical antimicrobials or hydrogels can be applied over the top or under Spycra Contact.

Instructions for us

- Clean the wound carefully and dry the surrounding skin thoroughly.
- Choose a Spycra Contact that covers the wound and the surrounding skin by at least 2 cm. In case of large wounds, it is preferable to get a larger protection of the surrounding skin. Before removing the protective films, Spycra Contact can be cut to size, if needed. Choose a larger size instead of several – overlapping – small sizes to secure a perfect adherence. (Spycra Contact does not secure a proper adherence when overlapping; it is not recommended to overlap several small sizes).

- Hold Spycra Contact with the protective film and adhesive layer upwards. Before removing the protective film, it is recommended to gently stretch Spycra Contact in order to have a better grip on the cutting line. The protective paper can now be easily be gripped at both sides of the cutting line.
- Hold, at the same time, a top of the undulation at the left side of the cutting line (with the left hand) and one at the right side (with the right hand). Remove the protective film by tearing evenly the left side and the right side and make sure that approximately 1 to 2 cm of the protective film is still left. Moisten gloves to avoid adherence of Spycra Contact onto the gloves.
- Apply Spycra Contact on the wound and gently smooth and ensure a good seal onto the surrounding skin.
- Remove the remaining protective film. Make sure that Spycra Contact is applied in a normal way (not stretched). If necessary, Spycra Contact or a part of Spycra Contact can be removed and repositioned. Where clinically indicated, treatments such as topical steroids, topical antimicrobials or hydrogels can be applied above or under Spycra Contact.
- Apply a secondary absorbing dressing on top Spycra Contact.
- A large fixation dressing or bandage can be used to have an optimal fixation of the secondary dressing. Can be used under compression bandages.

Frequency of change

Spycra Contact can be left in place up to 14 days or as long as the commonly used protocol proposes and depends on the condition of the wound (exudate should pass freely through the dressing and the pores should not be blocked). If saturated, the outer secondary absorbent dressing should be changed.

Precautions

- The wound should be inspected for signs of infection according to clinical practice. Consult a healthcare professional for the appropriate medical treatment.
- Spycra Contact may be used on Epidermolysis Bullosa patients (E.B.) after consulting a qualified health care professional
- Spycra Contact used on partial thickness burns with high risk of rapid granulation: avoid placing unnecessary pressure upon the dressing and lift or reposition the dressing at least every second day.
- When Spycra Contact is used for the fixation of skin grafts and protection of blisters, the dressing should not be changed before the fifth day post application.
- when used on bleeding wounds or wounds with high viscosity exudate, Spycra Contact should be covered with a moist absorbent dressing.
- Moist skin or heavy transpiration will diminish the adherence of the Spycra Contact. Before applying the dressing, make sure the skin is completely dry.

Check the packaging before use. Sterility is guaranteed unless inner package is damaged or opened prior to use. Do not re-sterilize

When re-using the same dressing, it is obvious that the dressing isn't sterile anymore; the dressing can be contaminated. As a consequence, applying the dressing on a wound increases the risk of wound infection.

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issued 2021-01

INGREDIENTS: Polyester + medical soft silicone gel

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EN SPYCRAS Contact TWO

Spycra Contact Two: sterile, soft double sided silicone adhesive wound contact layer

Product description

Elastic, double-sided silicone adhesive wound contact layer with self-adhering properties made of soft and flexible polyamide fibers. When used on exuding wounds, the wide scaled pore structure of the dressing allows excess of wound fluid to pass vertically towards an outer absorbent dressing. Does not absorb.

When used on fragile skin, Spycra Contact Two secures a soft and perfect adherence on the skin allowing dressing changes without pain and loss of epidermis.

Does not stick onto the wound. Repositionable
Spycra Contact Two can be used on exuding or bleeding wounds whereby the comfort of the patient is crucial (painless removal of the dressing and optimal support of the wound healing processes).

Mode of action

The wide scaled pore structure of the dressing allows excess of wound fluid to pass vertically towards an outer absorbent dressing. The soft silicone layer seals off the wound edges avoiding maceration and keeping the surrounding skin dry and healthy limiting the risk of maceration. The outer absorbent dressing needs to be changed when saturated limiting the risk of maceration. Spycra Contact Two may be left in place up to several days with a maximum of 14 days depending on the condition of the wound and the amount of wound fluid. The soft adherence secures no loss of epidermis during dressing changes and makes these changes as pain free as possible.

Intended use

Can be used as wound contact layer to take care of heavy exuding wounds such as skin tears, abrasions, surgical incisions, partial thickness burns, fixation of transplants, diabetic ulcers, venous and arterial ulcers, pressure sores, Negative Pressure Wound Therapy, ...

Spycra Contact Two can also be used as a protective layer on non-exuding wounds and on areas with fragile skin.

Where clinically indicated, treatments such as topical steroids, topical antimicrobials or hydrogels can be applied over the top or under Spycra Contact Two.

Instructions for use

1. Clean the wound carefully and dry the surrounding skin thoroughly.
2. Choose a Spycra Contact Two that covers the wound and the surrounding skin by at least 2 cm. In case of large wounds, it is preferable to get a larger protection of the surrounding skin.
3. When you hold Spycra Contact Two, make sure the large protective foil is upwards; take off the protective foil. Moisten gloves to avoid adherence of Spycra Contact Two onto the gloves.
4. Apply Spycra Contact Two on the wound and gently smooth and ensure a good seal onto the surrounding skin.
5. Remove the remaining protective foil. If necessary, Spycra Contact Two can be lifted and repositioned.
6. Apply a secondary absorbing dressing on top Spycra Contact Two.

7. A large fixation dressing or bandage can be used to have an optimal fixation of the secondary dressing. Can be used under compression bandages.

Frequency of change

Spycra Contact Two can be left in place up to 14 days or as long as the commonly used protocol proposes and depends on the condition of the wound (exudate should pass freely through the dressing and the pores should not be blocked). If saturated, the outer secondary absorbent dressing should be changed.

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- Spycra Contact Two may be used on Epidermolysis Bullosa patients (E.B.) after consulting a qualified health care professional

- Spycra Contact Two used on partial thickness burns with high risk of rapid granulation: avoid placing unnecessary pressure upon the dressing and lift or reposition the dressing at least every second day.

- When Spycra Contact Two is used for the fixation of skin grafts and protection of blisters, the dressing should not be changed before the fifth day post application.

- when used on bleeding wounds or wounds with high viscosity exudate, Spycra Contact Two should be covered with a moist absorbent dressing.

- Moist skin or heavy transpiration will diminish the adherence of Spycra Contact Two. Before applying the dressing, make sure the skin is completely dry.

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

INGREDIENTS: PA polyurethane film + medical soft silicone gel

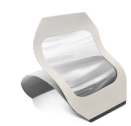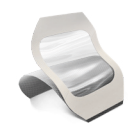Spycra Film: sterile, flexible soft silicone adherent polyurethane film developed with "ReSkin Silicone Skin Technology", a unique patented adhesive technology.

Product description
Spycra Film is a transparent thin polyurethane film with a soft, self-adhering silicone surface designed to protect fragile and sensitive skin. Spycra Film minimizes the risk of skin breakdown. Spycra Film is very conformable and adapts well to surface irregularities and different body contours. It provides a flexible transparent covering, to protect your patient's skin or wound from microbial contamination, fluid strikethrough or other external contamination. Spycra Film has breathable properties and supports the wound healing process based on the principle of moist wound healing. Spycra Film does not absorb. Spycra Film can be used as a stand-alone dressing or in combination with other products.

ReSkin Silicone Skin Technology
"ReSkin Silicone Skin Technology" is a patented combination of elastic material and a soft, self-adhering silicone layer on one side of the dressing. "ReSkin Silicone Skin Technology" is integrated in Spycra Film, an elastic polyurethane film which does not stick to the wound. Causes minimal trauma and pain when changing the dressing. When used on fragile skin, Spycra Film secures a soft and perfect adherence on the skin allowing dressing changes without pain and loss of epidermis. Thanks to the "ReSkin Silicone Skin Technology", Spycra Film can be used on 'skin at risk' to support the fragile skin. Spycra Film can also be used on superficial non-to-low exuding wounds whereby the comfort of the patient is crucial (painless removal of the dressing and optimal support of the wound healing processes).

Mode of action
Spycra Film is designed to protect fragile and sensitive skin. Spycra Film minimizes the risk of skin breakdown. Spycra Film is very conformable and adapts well to surface irregularities and different body contours. It provides a flexible transparent covering, to protect your patient's skin or wound from microbial contamination, fluid strikethrough, or other external contamination. Spycra Film has breathable properties and supports the wound healing process based on the principle of moist wound healing.

Intended use
Spycra Film is designed to protect fragile and sensitive skin. Spycra Film minimizes the risk of skin breakdown. Spycra Film can be used as a secondary dressing fixating e.g. gauzes or non-woven swabs or as a landing zone for aggressive tapes to secure devices to skin. Spycra Film can also be used as a secondary protective layer over IV sites when the IV device is secured through adhesive tape or sutures.

Instructions for use
1. Clean the skin or wound carefully and dry the skin thoroughly.
2. Choose a size large enough to cover the entire skin area to protect. In case of a wound, choose a size which covers the wound and the surrounding skin by at least 2 cm.
3. Hold Spycra Film with the protective transparent layer upwards and remove the protective layer.
4. Apply Spycra Film on the fragile skin and gently smooth Spycra Film in place ensuring a good seal.

Spycra Film is based on 'ReSkin Silicone Skin Technology' developed by Reskin Medical N.V.

 Reskin Medical NV, Poort West-Limburg 1399, Industrieweg 114, B-3980 Tessenderlo, Belgium
T: +32 (0)13 67 13 67, F: +32 (0)13 67 13 68

5. Remove the white upper protective tabs/borders.
Note: Make sure that Spycra Film is applied in a normal way (not stretched). If necessary, the dressing or part of the dressing can be removed and repositioned.
Where clinically indicated, treatments such as topical steroids, topical antimicrobials or hydrogels can be applied under Spycra Film. When applying hydrogels, creams or ointments under Spycra Film, the adherence of the dressing can be compromised.

Frequency of change
Spycra Film can be left in place for several days depending on the wound/skin condition or as indicated by good clinical practice.
Lift a corner of the dressing and gently pull off the film along with the skin.

Precautions
- Spycra Film should not be used directly on full thickness wounds involving muscle, tendon or bone or on third degree burns.
- Moist skin or heavy transpiration will diminish the adherence of Spycra Film. Before applying the dressing, make sure the skin is completely dry.
- When taking a bath, make sure that Spycra Film does not come in contact with water.

Warnings
Spycra Film should not be used as a primary fixation for IV, cannulae, ports or other infusion and/or life-sustaining devices.

Check the packaging before use.
Sterility is guaranteed unless inner package is damaged or opened prior to use. Do not re-sterilize

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EN SPYCRA[®]

PROTECT

Spycra Protect: sterile, bi-elastic soft silicone adhesive dressing developed with "ReSkin Silicone Skin Technology"; a unique patented adhesive technology

Product description

Spycra Protect is a soft, silken-like one-sided silicone adhesive dressing made of bi-elastic lycra fibers with self-adhering properties. Spycra Protect is bi-elastic; elastic in both directions - in the length and in the width -, can be repositioned and has rounded edges to stay better in place.

Spycra Protect has breathable properties and supports the wound healing process based on the principle of moist wound healing.

Spycra Protect does not absorb.

ReSkin Silicone Skin Technology

"ReSkin Silicone Skin Technology" is a patented combination of elastic fibers and a soft, self-adhering silicone layer on one side of the dressing.

"ReSkin Silicone Skin Technology" is integrated in Spycra Protect, a bi-elastic dressing which does not stick to the wound. Causes minimal trauma and pain when changing the dressing. When used on fragile skin, Spycra Protect secures a soft and perfect adherence on the skin allowing dressing changes without pain and loss of epidermis. Thanks to the "ReSkin Silicone Skin Technology", Spycra Protect can be used 'on skin at risk' to support the fragile skin. Spycra Protect can be used on superficial non-to-low exuding wounds whereby the comfort of the patient is crucial (painless removal of the dressing and optimal support of the wound healing processes).

Mode of action

Used as protective layer to support fragile skin (skin at risk):

The soft silicone adherent layer in combination with the bi-elastic top layer secures a perfect support of fragile skin.

The bi-elasticity minimizes the shearing and traction forces on the skin, limiting the risk of secondary skin lesions. The wide scaled pore structure of the dressing allows the dressing to breathe; transpiration can evaporate immediately limiting the risk of maceration. The dressing can easily be re-positioned without loss of adherence. The dressing can stay in place for several days. Dressing changes are painless and a-traumatic.

Used for wound care:

The soft silicone adherent layer seals off the wound edges avoiding wound fluid to come in contact with the surrounding skin limiting the risk of maceration and keeps the surrounding skin dry and healthy. The dressing can stay in place for several days. Dressing changes are painless and a-traumatic.

Intended use

1. Support of fragile skin.

Spycra Protect can be used preventively as a protective layer when applied on 'skin at risk' to support fragile skin. The bi-elastic properties of the dressing protect the fragile skin tissues minimizing the risk of skin tears.

Spycra Protect can also be used to minimize the friction and shearing forces among E.B patients, bedridden patients and wheelchair and prostheses users.

2. Care of superficial non-to-low exuding wounds.

Spycra Protect can be used on different superficial non-to-low exuding wounds such as: skin tears, scratch injuries, surgical incisions, partial thickness burns, blisters, shearing injuries, abrasions, intertrigo, eczema,...

Instructions for use

1. Clean the skin or wound carefully and dry the skin thoroughly.

2. Choose a size that covers the wound and the surrounding skin by at least 2 cm. If needed, Spycra Protect can, before removing the protective film, be cut to size. It is preferable to use a larger size instead of several - overlapping - small sizes to secure a perfect adherence. (Spycra Protect does not secure a proper adherence when overlapping; it is not recommended to overlap several small sizes).

3. Hold Spycra Protect with the protective - film and adhesive layer upwards. Before removing the protective - film, it is recommended to gently stretch the dressing in order to have a better grip on the cutting line. Fold the dressing with the opened cutting line upwards and rub softly the already removed parts of the protective layer.

The rubbing facilitates the removal of the protective layer from the dressing. The protective paper can now be easily be gripped at both sides of the cutting line.

Spycra Protect is based on 'ReSkin Silicone Skin Technology' developed by Reskin Medical N.V.

4. Hold, at the same time, a top of the undulation at the left side of the cutting line (with the left hand) and one at the right side (with the right hand). Remove the protective film by tearing evenly the left side and the right side and make sure that approximately 1 to 2 cm of the protective film is still left.

Moisten gloves to avoid adherence to Spycra Protect.

5. Apply Spycra Protect on the wound (or on the fragile skin which needs protection) and gently smooth Spycra Protect in place onto the surrounding skin, ensuring a good seal.

6. Remove the remaining protective film. Make sure that Spycra Protect is applied in a normal way (not stretched). If necessary, the dressing or part of the dressing can be removed and repositioned. Where clinically indicated, treatments such as topical steroids, topical antimicrobials or hydrogels can be applied under Spycra Protect. Applying hydrogels, creams or ointments under Spycra Protect compromises the adherence of the dressing; it is recommended to put an extra outer -fixation dressing on top of Spycra Protect.

7. If used on a dry skin surface or if used on a moist non-to-low exuding wound, no extra secondary dressing is needed.

The soft, silken-like top layer does not stick to the bed linen.

The soft silken-like top layer and the rounded edges ensure a perfect stay on the wound or on the fragile skin; the dressing does not roll.

8. The bi-elasticity of Spycra Protect ensures a perfect fit if applied on jointed or contoured areas (e.g. under arm, under breast, inner elbow, groin,...). It is recommended to use a dressing which has a larger size to have a perfect seal. The bi-elasticity of Spycra Protect limits the risk of traction on the skin avoiding secondary skin lesions.

Frequency of change

Spycra Protect can be left in place for several days.

If used as protective layer on a skin at risk or to limit the consequences of shearing forces, it is recommended to change the dressing once a week.

Precautions

- The wound should be inspected for signs of infection according to clinical practice. Consult a healthcare professional for the appropriate medical treatment.

- Spycra Protect may be used on Epidermolysis Bullosa patients (E.B.) after consulting a qualified health care professional

- Spycra Protect used on partial thickness burns with high risk of rapid granulation: avoid placing unnecessary pressure upon the dressing and lift or reposition the dressing at least every second day.

- When Spycra Protect is used for the fixation of skin grafts and protection of blisters, the dressing should not be changed before the fifth day post application.

- Moist skin or heavy transpiration will diminish the adherence of Spycra Protect. Before applying the dressing, make sure the skin is completely dry.

Check the packaging before use.

Sterility is guaranteed unless inner package is damaged or opened prior to use. Do not re-sterilize

When re-using the same dressing, it is obvious that the dressing isn't sterile anymore; the dressing can be contaminated.

As a consequence, applying the dressing on a wound increases the risk of wound infection.

 Reskin Medical NV, Poort West-Limburg 1399, Industrieweg 114, B-3980 Tessenderlo, Belgium
T: +32 (0)13 67 13 67, F: +32 (0)13 67 13 68

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Transparent film on roll with silicone adhesive

Spycra Roll: transparent polyurethane film with an silicone adhesive developed with "ReSkin Silicone Skin Technology"; a unique patented adhesive technology.

Product description

Spycra Roll is a transparent thin polyurethane film with a soft, self-adhering silicone surface designed to protect fragile and sensitive skin.

Spycra Roll minimizes the risk of skin breakdown. Spycra Roll is very conformable and adapts well to surface irregularities and different body contours. It provides a flexible transparent covering and to protects your patient's skin or wound from microbial contamination, fluid strikethrough or other external contamination.

Spycra Roll has breathable properties and supports the wound healing process based on the principle of moist wound healing.

Spycra Roll does not absorb. Spycra Roll can be used as a stand-alone dressing or in combination with other products. Spycra Roll can be cut to size and is therefore a perfect solution to protect small or larger areas.

ReSkin Silicone Skin Technology

"ReSkin Silicone Skin Technology" is a patented combination of elastic fibers and a soft, self-adhering silicone layer on one side of the dressing. "ReSkin Silicone Skin Technology" is integrated in Spycra Roll, an elastic polyurethane film which does not stick to the wound and causes minimal trauma and pain when changing the dressing. When used on (fragile) skin, Spycra Roll secures a soft and perfect adherence on the skin allowing dressing changes without pain and loss of epidermis. Thanks to the "ReSkin Silicone Skin Technology", Spycra Roll can be used on 'skin at risk' to support the fragile skin. Spycra Roll can be used on superficial non-to-low exuding wounds whereby the comfort of the patient is crucial (painless removal of the dressing and optimal support of the wound healing processes).

Mode of action

Spycra Roll is designed to protect fragile and sensitive skin. Spycra Roll minimizes the risk of skin breakdown.

Spycra Roll is very conformable and adapts well to surface irregularities and different body contours. It provides a flexible transparent covering, to protect your patient's skin or wound from microbial contamination, fluid strikethrough, or other external contamination.

Spycra Roll has breathable properties and supports the wound healing process based on the principle of moist wound healing.

Intended use

- as a primary dressing on fragile and sensitive skin to minimize skin breakdown.
- as a protective film to protect skin or wounds from microbial or other contamination and from fluid strikethrough.
- as a secondary dressing fixating gauzes and non-woven swabs.

Instructions for use

1. Wash your hands before starting the care
2. Clean the skin or wound carefully and dry the skin thoroughly.
3. Apply an absorbent patch if needed
4. Choose a size large enough to cover the entire skin area to protect. In case of a wound, choose a size which covers the wound and the surrounding skin by at least 2 cm.
5. Hold Spycra Roll with the protective transparent layer upwards and remove the protective layer.
6. Hold the transparent Film on Roll by the two green border protections and apply Spycra Roll on the fragile skin; gently smooth Spycra Roll in place ensuring a good seal.
7. Remove (tear off) the green border protections and gently smooth Spycra Roll once more to secure a good seal.

Note: Make sure that Spycra Roll is applied in a normal way (not stretched).

Where clinically indicated, treatments such as topical steroids, topical antimicrobials or hydrogels can be applied under Spycra Roll. When applying hydrogels, creams or ointments under Spycra Roll, the adherence of Spycra Roll can be compromised.

Removal of Spycra Roll

Lift a corner of the Film and gently pull off the film along with the skin.

Precautions

- Moist skin or heavy transpiration will diminish the adherence of Spycra Roll. Before applying the dressing, make sure the skin is completely dry.
- Spycra Roll may not be used on open wounds as a primary dressing.
- When taking a bath, make sure that Spycra Roll does not come in contact with water.

Warnings

Spycra Roll should not be used as a primary fixation for IV, cannulae, ports or other infusion and/or life-sustaining devices.

 **ReSkin Medical NV, Poort West Limburg 1399, Industrieweg 114, B-3980 Tessenderlo, Belgium**
T: +32 (0)13 67 13 67 F: +32 (0)13 67 13 68

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issued 2020-03

INGREDIENTS: polyurethane film + medical soft silicone gel

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