# STERRAD™ Booster STERRAD™ Adaptor

en

## INSTRUCTIONS FOR USE

### INDICATIONS FOR USE

The STERRAD Booster (REF 15400-020) is designed to deliver hydrogen peroxide into the interior of a lumen when used with the STERRAD™ Sterilization System (STERRAD™ 100S). The STERRAD Booster is inserted into a STERRAD Adaptor. The STERRAD Adaptor, (REF 15401, 15402, 15403, & 15407) (Adaptors I, II, III, IV) are to be attached on medical device ports containing lumens. Please read these instructions carefully before proceeding.

### WARNING

DO NOT ATTEMPT TO STERILIZE ITEMS OR MATERIALS THAT DO NOT COMPLY WITH THESE INSTRUCTIONS. IN ADDITION, YOU SHOULD READ THE MEDICAL DEVICE MANUFACTURER'S REPROCESSING INSTRUCTIONS OR CALL YOUR ASP REPRESENTATIVE TO DETERMINE WHETHER AN ITEM CAN BE STERILIZED BY THE STERILIZER.

There are three general steps to follow:

- Step 1 – Determine the dimensions (Inside Diameter and Total Channel Length) of each lumen in your medical device using information from the Medical Device Manufacturer. Note that "Working Length" is typically shorter than the "Total Lumen Length" for most medical devices.
- Step 2 – Determine if a STERRAD Booster Adaptor is required for processing in the STERRAD Sterilization System by referring to Section A. Note that more than one STERRAD Booster /Adaptor may be required for some devices. A STERRAD Booster/Adaptor is required for each proximal port on the device if one of the lumens of the device require a Booster/Adaptor.
- Step 3 – If a device requires a STERRAD Booster/Adaptor, follow steps and illustrations in section B for selection and attachment of the Booster/Adaptor.

Use Personal Protective Equipment (e.g. Chemical Resistant Gloves) at all times when handling.

### WARNING

FAILING TO USE A STERRAD BOOSTER / ADAPTOR WHEN REQUIRED, COULD RESULT IN A LUMEN WHICH IS NOT STERILE.

### Section A:

Determine if lumen dimensions of the medical device being processed need a STERRAD Booster/Adaptor for processing in the STERRAD Sterilization System. There are four different medical device lumen configurations listed in this section:

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**ADVANCED STERILIZATION PRODUCTS, INC.**

33 Technology Drive, Irvine, CA 92618, USA

+1(888)-STERRAD

ASP, The Netherlands BV  
BIC 1, 5657 BX, Eindhoven, The Netherlands

## 1) Medical devices that have a single Stainless Steel lumen

### **Use only Adaptors I, II, and III**

For STERRAD 100S:

Cycle type: Select a Short Cycle in sterilizer, if applicable

Use a STERRAD Booster/Adaptor for each lumen with the following:

Lumen size:

An inside diameter of  $\geq 1$  mm but  $< 2$  mm and a length  $> 125$  mm and  $\leq 500$  mm,

An inside diameter of  $\geq 2$  mm but  $< 3$  mm and a length  $> 250$  mm and  $\leq 500$  mm,

An inside diameter of  $\geq 3$  mm and a length of  $> 400$  mm and  $\leq 500$  mm.

## 2) Medical Tubing – single polyethylene or Teflon (PTFE) lumen

### **Use only Adaptors I, II, and III**

For STERRAD 100S:

Cycle type: Select the Long Cycle in sterilizer

Use a STERRAD Booster/Adaptor for each lumen with the following:

An inside diameter of 1 mm or greater ( $\geq 1$  mm) and a length of greater than 1000 mm ( $> 1000$  mm) but less than or equal to 2000 mm ( $\leq 2000$  mm)

## 3) Flexible Endoscopes with lumens made from polyethylene or Teflon (PTFE)

**Adaptors I, II, and III:** Use on endoscope ports that are not recessed.

**Adaptor IV:** Only use on recessed ports of the endoscope.

For STERRAD 100S:

**Cycle type:** Select the Long Cycle in sterilizer

Use Adaptors I, II, and III with single-channel flexible endoscopes having:

An inside diameter of 1 mm or larger ( $\geq 1$  mm) and a length of greater than 500 mm ( $> 500$  mm) and less than or equal to 2000 mm ( $\leq 2000$  mm)

Use Adaptors I, II, and III with multi-channel flexible endoscopes having:

An inside diameter of 1 mm or larger ( $\geq 1$  mm) and a length of less than or equal to 2000 mm ( $\leq 2000$  mm)

Use Adaptor IV with single-channel flexible endoscopes having:

An inside diameter of 2 mm or larger ( $\geq 2$  mm) and a length of greater than 500 mm ( $> 500$  mm) and less than or equal to 959 mm ( $\leq 959$  mm),

or,

An inside diameter of 1.2 mm or larger ( $\geq 1.2$  mm) and a length of greater than 500 mm ( $> 500$  mm) and less than or equal to 850 mm ( $\leq 850$  mm).

## 4) Medical devices with a lumen made of copper containing alloys, such as Monel

### **Use only Adaptors I, II, and III**

For STERRAD 100S:

Cycle type: Select a Short Cycle in sterilizer, if applicable

Use a STERRAD Booster/Adaptor for each lumen with the following:

Lumen size: Inside diameter of 3 mm or larger ( $\geq 3$  mm), and a length of equal to or less than 500 mm ( $\leq 500$  mm).

Please refer to the STERRAD™ User's Guides for more information on the lumen requirements for Boosters.

### **Cautions**

STERRAD Booster contains a minimum of 145  $\mu$ l of 58% hydrogen peroxide solution. Confirm that the indicator strip in the STERRAD Booster bag is yellow. If the indicator strip is red, a booster may be damaged. DO NOT open the package if the indicator strip is red. Discard entire contents of the bag. *Boosters should be used within 24 hours of activation.*

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### **WARNING**

Thorough cleaning of medical devices, especially at the interface area, is essential for assurance of sterilization by all sterilization methods, including the STERRAD™ Sterilizer as there is a risk of occlusion of organisms at the interface between the Adaptor and the medical device when the Booster/Adaptor is in place. Position the Booster(s) so that the hydrogen peroxide vapor can diffuse down the entire lumen without escaping prior to reaching all internal surfaces. Some endoscopes are so complex that the STERRAD Booster cannot attach properly. If the attachment cannot be securely attached, do not process the device in the STERRAD Sterilizer. If endoscope requires an ethylene oxide /vent cap for sterilization, attach vent cap prior to processing in the STERRAD Sterilizer. Should a cancellation occur when Booster and Adaptor are in the chamber, they should be discarded and the devices should be repackaged using new STERRAD Booster and Adaptor when restarting the cycle.

**Booster and Adaptor: Single Use** 

**DO NOT RE-USE THE BOOSTER AND ADAPTOR. RE-USE OF BOOSTER AND ADAPTOR MAY RESULT IN A NON-STERILE DEVICE.**

### Storage

Store below 25°C (77°F).

**Note:** If refrigerated, allow Booster to reach room temperature prior to use.

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
## Section B:

### 1 Select a Proper Adaptor

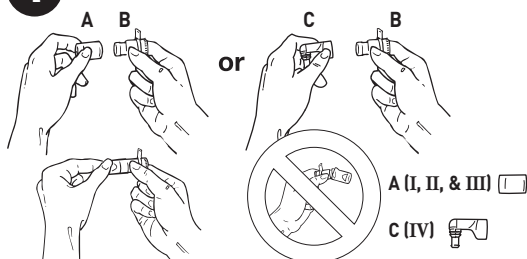STERRAD Adaptors I, II, & III are provided in three color-coded sizes. STERRAD Adaptor IV is clear in color and shape/size is distinctly different from Adaptors I, II, & III. For Adaptors I, II, & III, select the largest size that, when connected to a Booster, will not fall off the device. Connect one STERRAD Adaptor (A/C) to a Booster (B) by placing the open end of the Booster into the larger opening of the Adaptor. Be sure the Booster is fully inserted into the Adaptor. For devices with complex/multiple ports, see suggestions in Step 2.

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

## 2

# Attach STERRAD Adaptor-Booster Complex to Device

## Flexible Endoscope
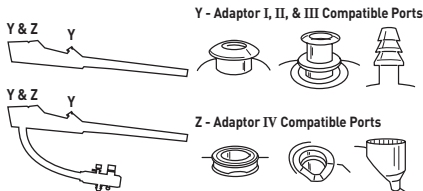
## Rigid Device

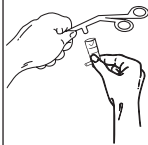For devices requiring STERRAD Booster and STERRAD Adaptor I, II, & III, attach the STERRAD Adaptor/Booster complex to the device by pushing the end of the STERRAD Adaptor over the end of the lumen or lumen port of the device. For endoscopes with a recessed port, attach the STERRAD Adaptor IV/Booster complex by inserting the Adaptor into the endoscope port as far as the Adaptor will allow.

Suggestions: More than one STERRAD Booster may be required for some devices. Some devices have multiple lumens or lumen ports that sometimes converge into one lumen. One STERRAD Adaptor/Booster complex should be attached to each lumen port opening to ensure adequate sterilant exposure to each channel.

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
### 3 Assure Proper Fit

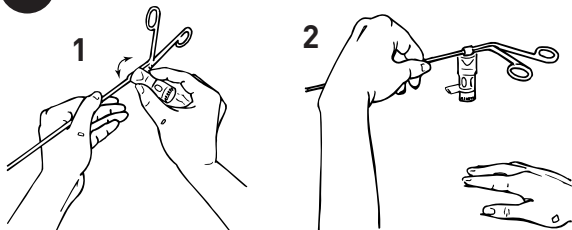1. To ensure a good fit for STERRAD Adaptors I, II, & III, grasp the STERRAD Adaptor and rotate it one-quarter turn in each direction. For STERRAD Adaptor IV, firmly push the insertion end into the recessed port as far as the Adaptor will allow. Do not use excessive force.
2. Turn the device over to test the Adaptor's grip on the device. A properly fitted STERRAD Adaptor will grip the device securely, without falling off.

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
## 6 Reattach STERRAD Adaptor-Booster Complex onto Device

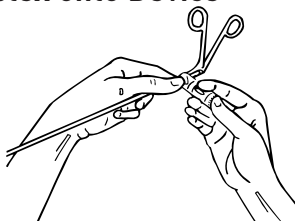Replace the STERRAD Adaptor/Booster complex onto the device (see steps **1** to **3**). Ensure that the attachment of the STERRAD Adaptor to the device is not compromised when placed in the sterilization packaging. After the sterilization is successfully completed, prior to using the device, remove the STERRAD Adaptor/Booster complex and discard it according to healthcare facility's procedures or following local waste regulations.

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**Danger** Contains: Hydrogen peroxide solution 58%

**Hazard statement** May intensify fire; oxidizer. Harmful if swallowed. Causes severe skin burns and eye damage. Harmful if inhaled. May cause respiratory irritation. Toxic to aquatic life.

**Precautionary statement**

**Prevention** Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. Keep/Store away from clothing and other combustible materials. Take any precaution to avoid mixing with combustibles. Do not breathe mist/vapors. Wash thoroughly after handling. Do not eat, drink or smoke when using this product. Use only outdoors or in a well-ventilated area. Avoid release to the environment. Wear protective gloves/protective clothing/eye protection/face protection.

**Response** If swallowed: Rinse mouth. Do NOT induce vomiting. If on skin (or hair): Take off immediately all contaminated clothing. Rinse skin with water/shower. Wash contaminated clothing before reuse. If inhaled: Remove person to fresh air and keep comfortable for breathing. If in eyes: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a poison center/doctor. In case of fire: Use appropriate media to extinguish.

**Storage** Store in a well-ventilated place. Keep container tightly closed. Store locked up.

**Disposal** Dispose of contents/container in accordance with local/regional/national/international regulations.

**\*\*For country specific GHS and SDS information, refer to country specific SDS located at [asp.com/service-and-support/sds-datasheets](http://asp.com/service-and-support/sds-datasheets).**

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[www.nevacert.ru](http://www.nevacert.ru) info@nevacert.ru © **Isuuntoisia vaiheita:**