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Instruction for Use

2. WCT Brackets Description

The World Class Technology Corporation (WCT) metal brackets included in this report are the WCT Mini-Twin Conventional Brackets (Figure 3), and WCT H4 Self-Ligating Brackets (Figure 4), both marketed under WCT's Ortho Classic brand. These brackets are provided in a variety of sizes and prescriptions; for the complete list of metal brackets encompassed by this report, refer to Table 8 and Table 9 in Section 13.1.

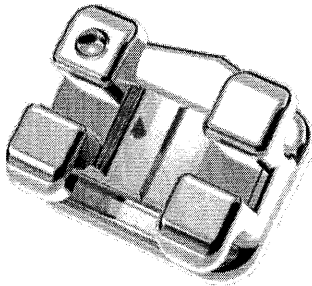
*Figure 3: Mini-Twin Brackets
Ligating Bracket*

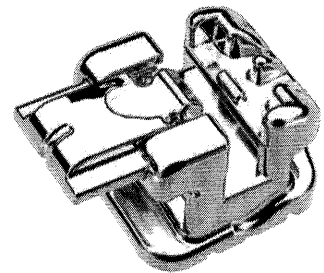
Figure 4: H4 Self-Ligating Bracket

As with other bracket systems, the WCT brackets are bonded to the teeth with commercially available materials, and linked together by archwire that applies continuous pressure to achieve the desired tooth movement. Also as with other systems, the brackets are invasive products designed to be used in the oral cavity, and are also intended for long-term use (continuous use over 30 days); the WCT brackets may be adhered to patients' teeth for several years.

The H4, as a self-ligating bracket, consists of two primary components: a base, which bonds (via adhesives) to the tooth, and the door, which slides open to allow capture and entrapment of the archwire during bracket placement. Both the Mini-Twin and the H4 also have specially engineered base footprints to optimize bonding to the tooth via adhesives. Both bracket types are produced using WCT's proprietary Base-Lok™ technology, which is designed to provide a very robust adhesive bond via a simulated mesh base. The Mini-Twin Brackets are available in simulated 100-mesh (Figure 5) and simulated 80-mesh base varieties, while the H4 Brackets are available in the simulated 80-mesh base.

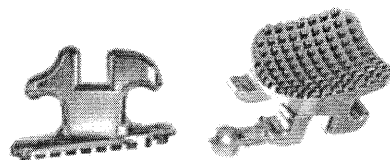
Figure 5: Base-Lok 100-Mesh Base

WCT Mini-Twin and H4 Brackets are both manufactured from 17-4 PH stainless steel, a standard material for orthodontic products. 17-4 PH stainless steel has a long history biocompatibility and use in a wide variety of invasive applications. These brackets are also both manufactured using a metal injection molding (MiM) process, also a widely accepted standard

Revision History

105-7300-10

Metal Buccal Tubes Instructions for Use

Rev	Description of Changes	Release Date
A	CO#20030: Initial Release	2020-05-22
B	CO#20084: Updated section 11 and 12	2020-08-20

1. Product Family & Part Numbers

Product Family	Product Part Numbers
BioMIM Single, Double, Triple (Convertible & Non-convertible), H4 and Pitts21 Buccal Tubes	*907.XXXX, 908.XXXX and 947.XXXX

*XXXX = Multiple part numbers in product family

2. Intended Use

WCT Metal Buccal Tubes are medical devices for *Professional Use Only* by dental professionals in the recommended indications. Devices are supplied:

- Non-sterile
- Designed for single-use
- For use by dentists and orthodontists only.

WCT devices are intended to be used in conjunction with other orthodontic devices.

WCT orthodontic appliances are used to aid in the movement of natural teeth in patients with malocclusion during orthodontic treatment.

3. Contraindications

- Patient's inability or unwillingness to cooperate/or follow the treatment plan.
- Known allergies to any of the components or materials in the system (see Table1).
- Any illness and/or underlying conditions that preclude orthodontic treatment.
- Any existing tooth root resorption.
- Any existing decalcification of the tooth enamel.
- Patients with deficient oral hygiene.

4. Materials

All part numbers referenced on sec. 1 are manufactured using material per Table 1

Table 1

Stainless Steel 17-4PH per ASTM B883-19

7. General Information for the Dentist/Orthodontist

- As part of developing a treatment plan, and before appliance placement, assess the need for interdisciplinary coordination with other professionals, such as speech pathologists, otolaryngologists, physicians, dentists and/or orthodontists.
- Follow the manufacturer's instructions for any orthodontic bonding agents, instruments or other materials used in the orthodontic treatment.
- Orthodontic training in standard procedures will determine appropriate instruments for use during appliance placement and removal.
- Do not touch bonding surfaces with bare fingers since skin oils may diminish adhesion.
- Oral hygiene is of particular importance for immunocompromised patients. Closely monitor oral hygiene on immunocompromised patients.
- Assess whether further orthodontic treatment is advisable in the presence of root resorption.


8. Disposal (if applicable)

- Disposal of all orthodontic appliances must follow regional and national regulations.

9. Storage and Handling for medical devices (if applicable)

- The device should be stored in a dry environment under ambient conditions.

10. Name and address of labeler

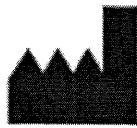World Class Technology

1300 NE Alpha Dr.
McMinnville. OR 97128

This document contains information considered PROPRIETARY and CONFIDENTIAL to World Class Technology.
Documents viewed from the server are Controlled. Printed documents are for Reference Only unless stamped Controlled.

Item Number: 105-7300-10	Revision: B	CO#20084	2020-08-20	Page 2 of 4
--------------------------	-------------	----------	------------	-------------

1300 NE Alpha Drive | McMinnville, Oregon 97128 USA

1-866-WCT-CORP | Fax: 503-435-2432 | Email: info@worldclasstech.com | www.worldclasstech.com

Instructions for Use		
Document Title	Document Description	Revision
105-7300-07	Metal Archwire and Straight Wire Instructions for Use	C

Revision History

Rev	Description of Changes	Release Date
A	CO#20030: Initial Release	2020-05-22
B	CO#20084: Updated section 12 and 13	2020-08-20
C	CO#21054: Add product family 300.XXXX; Updated section 6 (nitinol warning); Updated header/footer to Greenlight Guru template.	2021-07-27

1. Product Family & Part Numbers

Product Family	Product Part Numbers
Beta Titanium Archwire Forms	646.XXXX 61.40.170.XXXXX 61.40.270.XXXXX
Stainless Steel Archwire Forms	651.XXXX 61.40.180.XXXXX 61.40.280.XXXXX
Super Elastic Archwire Forms	625.XXXX 627.XXXX 41.40.120.XXXXX 41.40.220.XXXXX 41.40.320.XXXXX 41.40.420.XXXXX 41.40.520.XXXXX 61.40.120.XXXXX 61.40.125.XXXXX 61.40.130.XXXXX 61.40.135.XXXXX 61.40.220.XXXXX 61.40.225.XXXXX 61.40.230.XXXXX 61.40.235.XXXXX 61.40.420.XXXXX 61.40.425.XXXXX 61.40.520.XXXXX 61.40.525.XXXXX 61.40.150.11622 61.40.156.01622 61.40.250.11622 61.40.256.01622 61.40.150.11725 61.40.156.01725 61.40.250.11725 61.40.256.01725

WCT CONFIDENTIAL

Once printed, this document becomes UNCONTROLLED
Please verify latest version before each use

Instructions for Use		
Document Title	Document Description	Revision
105-7300-07	Metal Archwire and Straight Wire Instructions for Use	C

5. Materials

All part numbers referenced on sec. 1 are manufactured using materials per Table 1

Table 1

Materials List	
Stainless Steel Archwires and Straight Wires	304V per ASTM A313
Beta Titanium Archwires and Straight Wires	Ti Beta 3 (Ti 11.5Mo-6Zr-4.5Sn)
Super Elastic Archwires	Nickel Titanium #1 per ASTM F2063
Thermal Activated Archwires	Nickel Titanium #1 per ASTM F2063
Thermal Activated Ultra-Soft Archwires	Nickel Titanium #8 per ASTM F2063

6. Warnings and Precautionary Measures

- When gripping and/or holding archwire or straight wires, only use instruments without sharp edges or serrated surface. Scratches and scuffing of the wire caused by such instruments can lead to the wire breakage during treatment..
- Excessive force and/or repeated bending of the archwire and straight wire can lead to breakage.
- Do not use on patients with known allergies to any of the materials in the system (see Table 1).
- Super Elastic, Thermal Activated and Thermal Activated Ultra-Soft Archwires contain nitinol, an alloy of nickel and titanium and should not be used for individuals with known allergic sensitivity to these metals. Prior to use, patients should be counseled on the materials contained in the device, as well as the potential for allergy/hypersensitivity to these materials. (see Table 1).
- Immediately remove Orthodontic Appliance(s) in the event of an allergic reaction.
- Follow all regional and national standards regarding the use of orthodontic appliances.
- Do not use any products which are damaged, or do not comply with the labeling specifications.

7. Patient Information

- There is no information available which would preclude the use of commonly available oral healthcare products.
- Chewing hard foods can cause appliances to be damaged.
- Some sports may cause damage to orthodontic appliances, and their presence may increase the risk of harm in the event of certain sports-related injuries.
- When participating in sports, always wear appropriate mouth and/or bracket guards as recommended by the dental professional treating the patient.
- Always inform the MRI or radiology staff that braces are in place before any procedure to ensure appropriate measures are taken for the procedure.

WCT CONFIDENTIAL

Once printed, this document becomes UNCONTROLLED
Please verify latest version before each use

Page 3 of 4

1300 NE Alpha Drive | McMinnville, Oregon 97128 USA
1-866-WCT-CORP | Fax: 503-435-2432 | Email: info@worldclasstech.com | www.worldclasstech.com

Rev	Description of Changes	Release Date
A	CO#20030: Initial Release	2020-05-22
B	CO#20048: Updated product part numbers sequence in Sec. 1	2020-06-15
C	CO#20084: Updated section 12 and 13.	2020-08-20

1. Product Family & Part Numbers

Product Family	Product Part Numbers
Duraband	33.30.XXX.XXXXX
Fierro Bands	34.30.XXX.XXXXX

*XXXX = Multiple part numbers in product family

2. Intended Use

WCT Orthodontic bands wrap around the molar and bicuspid (premolar) teeth of the maxillary or mandibular dental arch and function in conjunction with a variety of orthodontic appliances supporting orthodontic treatment. Bands are pressed and bonded to bicuspid and molar teeth to serve as a sturdy anchorage point for supporting auxiliary devices in orthodontic treatment.

3. Instructions for use

Bands are pressed and bonded to bicuspid and molar teeth to serve as a sturdy anchorage point for supporting auxiliary devices. Weld or solder the required attachments in accordance with your treatment plan. Use conventional orthodontic band cement following manufacturer's instructions. Place it using the instrument of your choice. Remove excess cement from around the band and from the occlusal surface of the tooth.

4. Contraindications

- Patient's inability or unwillingness to cooperate/or follow the treatment plan.
- Known allergies to any of the components or materials in the system (see Table1).
- Any illness and/or underlying conditions that preclude orthodontic treatment.
- Any existing tooth root resorption.
- Any existing decalcification of the tooth enamel.
- Patients with deficient oral hygiene.

5. Materials

All part numbers referenced on sec. 1 are manufactured using materials per Table 1

Table 1

305 Stainless Steel per AMS5514H		
Element Name	Min. (% by weight)	Max. (% by weight)
Carbon (C)	--	0.12
Manganese (Mn)	--	2.00
Silicon (Si)	--	1.00
Phosphorus (P)	--	0.040
Sulfur (S)	--	0.030
Chromium (Cr)	17.00	19.00
Nickel (Ni)	10.00	13.00
Molybdenum	--	0.75
Copper (Cu)	--	0.75

This document contains information considered PROPRIETARY and CONFIDENTIAL to World Class Technology. Documents viewed from the server are Controlled. Printed documents are for Reference Only unless stamped Controlled.

Item Number: 105-7300-08 Revision: C CO#20084 2020-08-20 Page 1 of 4

1300 NE Alpha Drive | McMinnville, Oregon 97128 USA

1-866-WCT-CORP | Fax: 503-435-2432 | Email: info@worldclasstech.com | www.worldclasstech.com