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SHOFU All Ceramic Preparation Kit

**For Preparation of
All Ceramic Restorations**

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SHOFU All Ceramic Preparation Kit

The demand for ceramics in dentistry is at an all time high. Newer ceramic materials available today produce superb ceramic restorations. This is challenging the dental professional to consider newer materials, esthetics, function, and clinical expertise when preparing an all ceramic restoration.

Under the guidance of Dr. Takao Maruyama of Osaka University, Dr. Yoshihiko Mutoke and Shofu have developed an all ceramic preparation system utilizing diamond burs to prepare an all ceramic restoration.

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(CAUTIONS)

●Never exceed the maximum permissible rotation speeds described in Table A. ●Insert the shank completely into the chuck of the handpiece. ●Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation and dirt and exposure of shank due to repeated use. ●Perform test run to check safety and smooth operation without run-out prior to use. ●Use intermittently with gentle touch and water spray. ●Avoid any use at acute angles or with excessive pressure to thin and long abrasives. ●Use local dust extractor and protective equipment while using. ●This product is intended for use by dental professionals only. ●Do not use this product for any purposes other than the intended use.

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(STERILIZATION)

●Perform autoclave under the condition described in Table B. ●When using intraorally, sterilize before each use. ●For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. ●Do not use chlorine disinfectant liquid. ●Follow the instructions for use of washer disinfectant and sterilizer, respectively. ●For the details about REPROCESSING procedures, visit www.shofu.de.

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0833-1, 833V-1

Head length 9.3 mm, diameter 1.1 mm (at the top of working part) with tapered 4° angle. Use to slice cut an interproximal area or box preparation at the interproximal area of class II inlay. Use 1/2 of the diameter, except when preparing the lingual of the cervical area of a lower anterior tooth then use 2/3 of the diameter.

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0836-1, 836V-1

Head length 9.0 mm, diameter 1.3 mm (at the top of working part) with tapered 4° angle. Use to prep anterior tooth, buccal side of upper and lower premolar or cervical area on lingual side. Use 1/2 -2/3 of the diameter.

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0839-1, 839V-1


Head length 9.0 mm, diameter 1.5 mm (at the top of working part) with tapered 6° angle. Use to prep upper and lower molar on buccal and lingual side. For cervical area, use 2/3 of the diameter.

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0825-1, 825V-1.....

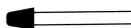
Head length 3.5 mm, diameter 1.5 mm (at the top of working part) with tapered 20° angle. Use to prep the occlusal box-form of an upper and lower premolar and molar.

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0826-1, 826V-1.....

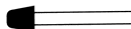
Head length 3.5 mm, diameter 2.0 mm (at the top of working part) with tapered 20°. Use for occlusal box-form of upper and lower molar.

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0893-1, 893V-1.....

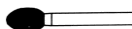Head length 4.9 mm. For cutting upper and lower anterior tooth on lingual side (to maintain occlusal space) and lower premolar on occlusal area.

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
0899-1, 899V-1.....

Tapered 20° angle (at inverted conical shape). Use to reduce an upper premolar and the occlusal area of an upper and lower molar to maintain the occlusal space. Keep the upper axis of the diamond bur parallel to the occlusal surface to prepare a 140° angle, the average occlusal cusp angle of a molar. By altering the inclination of the diamond bur, adjust the cusp top angle to accommodate the opposing tooth, the upheaval shape of the diamond bur produces a concave prep making the necessary space for the points of a crown. It also smoothes the finish of the overall edge of the occlusal and axial angle of the molar.

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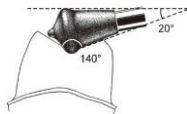808V-1

Head length 8.0 mm. Use this diamond bur after the prep is finished to eliminate the free enamel in the cervical area. User a lower speed to avoid chipping the enamel.

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0896-1

Use to prep a 0.5 mm depth groove in a laminate veneer preparation.

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
1. Slice cut the interproximal area of the anterior tooth (Fig. 1-1). Use the long axis of the bur (Fig. 1-2). After slice cut, retract tissue with retraction cord.

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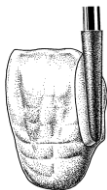1-1

1-2

2. Reduction of incisal edge (Fig. 2-1). Reduce the incisal edge 1.5 – 2.0 mm. Angle the diamond toward the lingual to stay parallel with the incisal line.

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3-1


3-2

4. Reduce the labial surface to the depth the lateral guide grooves (Fig 4-1 and 4-2)

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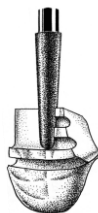4-1


4-2

5. Prepare the labial and lingual cervical area. Cut from the lowest point of the cervical scallop mesially and distally on the labial and lingual surfaces. Make sure to reduce enough in the center of the labial and lingual surfaces. Reduce enough in the center of the labial surface and not too much in the interproximal subgingival cervical areas (Fig. 5-1 and 5-2).

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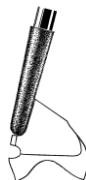
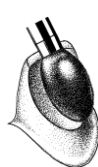
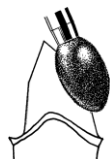
7

8. Prepare the shoulder margin. Smooth and finish the enamel at the cervical area (Fig. 8).

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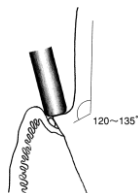
8

9. Smooth and finish the shoulder prep by using the outside 1/3 of the bur cutting edge that is away from the long axis of the prep. Move the bur in the direction of the rotation (Fig. 9).

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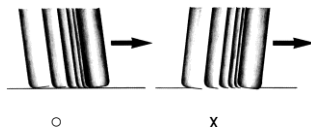
9

10. Place the bur vertically on the shoulder margin surface, not parallel to the long axis of the tooth preparation (Fig. 10).

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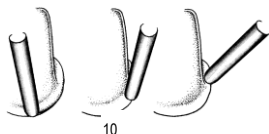
10

11. Completed preparation (Fig. 11). Imagine the final shape while cutting the preparation (Fig. 12 and 13)

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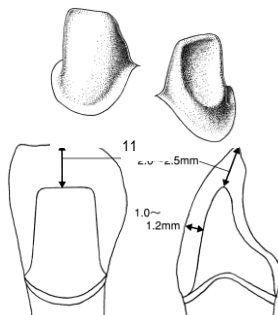
1. Preparation of lower premolar (1st bicuspid). Interproximal view (Fig. 14). Buccal view (Fig. 15).

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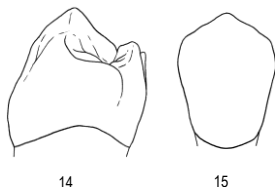2. Reduce the four surfaces of the premolar. A slight undercut will occur while following the lingual inclination of the premolar. The bur, when placed correctly against the occlusal plane, will tip toward the buccal surface (Fig 16).

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3. Reduce the apex of the buccal cusp following the inclined plane angles (Fig. 17).

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
4. Utilizing the premolar anatomy. The occlusal surface is prepared similar to an anterior preparation (Fig. 18).

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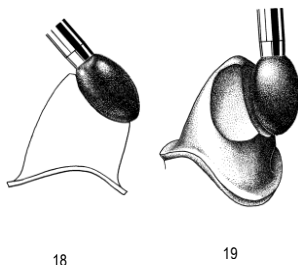
5. To make broad mesial and distal pits in the occlusal surface, reduce toward each interproximal direction (Fig. 19).

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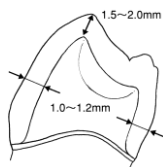
20

6. Completed preparation. Interproximal view (Fig. 20).

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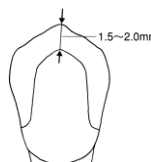21

7. Completed preparation. Buccal view (Fig. 21).

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22


23

2. Reduce occlusal surface to create occlusal space (Fig. 24-1). Use shape of the bur top to cut the central groove and to shape the side areas for the pits of the crown (Fig. 24-2).

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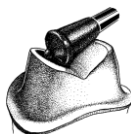24-1


24-2


3. Use the nipple part of the bur top to smooth and finish the overall edge of the occlusal and axial angle (Fig. 25-1 and 25-2).

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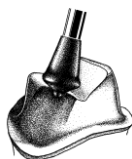25-1

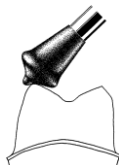
25-2


4. Completed preparation. Interproximal view (Fig. 26).

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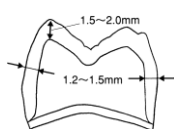26

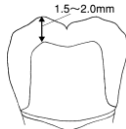27

5. Completed preparation. Buccal view (Fig. 27).

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
1. Cut the step for the box preparation in the interproximal area (Fig.28).

To make a clear margin on the buccal and lingual surface, cut with the bur counter-clockwise toward the apex of the tooth. Move from the buccal toward the lingual (check the occlusal contact to determine the form of the preparation prior to cutting the prep).

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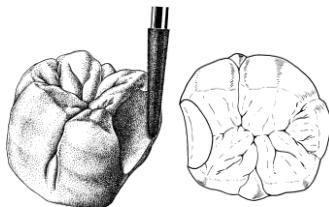
30

2. Form the box preparation on the occlusal surface (Fig. 31). Prepare a gently curved preparation, nota dove-tail preparation (Fig.32).

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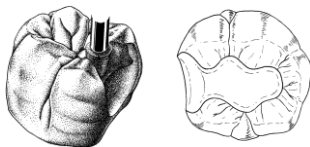The occlusal portion of the preparation must be at least 1.5 mm from the fossa to the floor of the preparation (Fig. 33). The box depth of the prep is cut 2.5 - 3.0 mm considering the length of the bur (3.5 mm) (Fig. 34).

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33


34


3. Round the corner of pulpal and axial walls. (Fig. 35-1 and 35-2). Rounding of the box prep will strengthen and produce an excellent inlay (Fig. 36).

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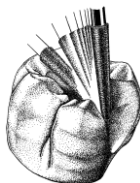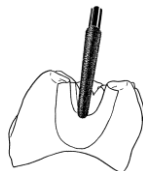
1. Retract tissue prior to preparation with retraction cord. Cut a lateral guide groove 0.5 mm in depth, then outline the groove with a marking pen (Fig. 37).

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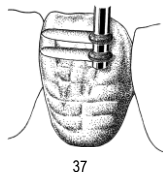
37

2. Following the drawn outline, prep the surface (Fig. 38).

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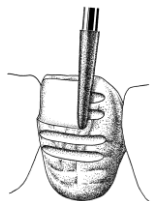
38

3. In order to maintain proper thickness of the veneer, keep the bur vertical to the cervical tissue (Fig. 39).

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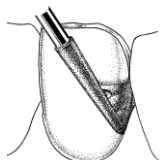
39

4. Prepare the embrasure at the lower interproximal area (Fig. 40). Prep the interproximals to the contact point. Prep from the facial in order not to expose the margin on the mesial or distal sides (Fig. 41). Prep at the correct angle to avoid undercuts mesially and distally.

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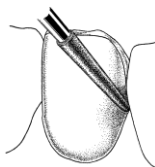5. Prepare the embrasure in the upper interproximal area (Fig. 42).

Prep the upper interproximal areas the same as the lower interproximals to avoid enamel tipping, move the bur in the direction of rotation.

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

This assortment is recommended for systematic finishing and polishing of ceramic restorations. It includes three Dura-Green Stones, three Dura-White Stones, three CeraMaster Coarse Polishers (blue), and three CeraMaster Polishers (white and blue bands). First, use Dura-Green or Dura-White Stones for grinding and contouring the ceramic restoration. Then use CeraMaster Coarse for finishing and polishing and CeraMaster for super polishing. The CeraMaster, which are filled with a high density of pure diamond particles, are durable and designed for quick finishing and polishing of ceramic restorations to provide a smooth surface.

PRECAUTIONS AND WARNINGS

● Never exceed the maximum permissible rotation speeds described in Table A. ● Insert the shank completely into the chuck of the handpiece. ● Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use. ● Perform test run to check safety and smooth operation without run-out prior to use. ● Avoid any use at acute angles or with excessive pressure. ● Use local dust extractor and protective equipment while using. ● This product is intended for use by dental professionals only. ● Do not use this product for any purposes other than the intended use. ● Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

● Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination. ● Perform autoclave under the condition described in Table B. ● When using intraorally, sterilize before each use. ● For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. ● Follow the instructions for use of washer disinfectant and sterilizer, respectively. ● Do not use chlorine disinfectant liquid. ● Dura-Green Stones and Dura-White Stones: Do not perform plasma sterilization. ● For the details about REPROCESSING procedures, visit www.shofu.de.

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CeraMaster Finishing & Polishing Kit

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CompoMaster Finishing & Polishing Kit

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

This assortment is recommended for systematic finishing and polishing of composite restorations. It includes four T&F Hybrid Points, three CompoMaster Coarse Polishers (yellow), and three CompoMaster Polishers (white and yellow bands). First use T&F Hybrid Points for grinding and contouring the composite restoration. Then use CompoMaster Coarse for finishing and polishing and CompoMaster for superpolishing. The CompoMaster, which are filled with a high density of pure diamond particles, are durable and designed for quick finishing and polishing of composite restorations to provide a smooth surface. All polishers are autoclavable for reuse. When additional superpolishing is required, use DirectDia Paste or SuperBuff polishing discs.

PRECAUTIONS AND WARNINGS

● Never exceed the maximum permissible rotation speeds described in Table A. ● Insert the shank completely into the chuck of the handpiece. ● Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use. ● Perform test run to check safety and smooth operation without run-out prior to use. ● Avoid any use at acute angles or with excessive pressure. ● Use local dust extractor and protective equipment while using. ● Use T&F Hybrid Points with water spray. ● This product is intended for use by dental professionals only. ● Do not use this product for any purposes other than the intended use. ● Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

● Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination. ● Perform autoclave under the condition described in Table B. ● When using intraorally, sterilize before each use. ● For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. ● Follow the instructions for use of washer disinfectant and sterilizer, respectively. ● Do not use chlorine disinfectant liquid. ● For the details about REPROCESSING procedures, visit www.shofu.de.

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nl

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PRODUCT DESCRIPTION

This kit is designed for finishing and polishing all types of composite restorations. It includes two Dura-Green and two Dura-White Stones for finishing, six CompoSite for prepolishing and polishing and two CompoSite Fine (white band) Polishers for superpolishing. Dura-Green and Dura-White Stones, CompoSite and CompoSite Fine Polishers are autoclavable for reuse.

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STORAGE AND STERILIZATION

- Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination.
- Perform autoclave under the condition described in Table B.
- When using intraorally, sterilize before each use.
- For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize.
- Follow the instructions for use of washer disinfectant and sterilizer, respectively.
- Do not use chlorine disinfectant liquid.
- Dura-Green Stones and Dura-White Stones: Do not perform plasma sterilization.
- For the details about REPROCESSING procedures, visit www.shofu.de.

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nl

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PRODUCT DESCRIPTION

This kit is designed for finishing and polishing composite restorations. It has two Dura-White Stones, eight CompoSite and eight CompoSite Fine Polishers. Dura-White Stones will finish composite restorations smoothly. CompoSite Polishers will prepolish and polish and CompoSite Fine (with white banded shank) Polishers will superpolish composite restorations. Dura-White Stones, CompoSite and CompoSite Fine Polishers are autoclavable for reuse.

PRECAUTIONS AND WARNINGS

- Never exceed the maximum permissible rotation speeds described in Table A.
- Insert the shank completely into the chuck of the handpiece.
- Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use.
- Perform test run to check safety and smooth operation without run-out prior to use.
- Avoid any use at acute angles or with excessive pressure.
- Use local dust extractor and protective equipment while using.
- This product is intended for use by dental professionals only.
- Do not use this product for any purposes other than the intended use.
- Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

- Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination.
- Perform autoclave under the condition described in Table B.
- When using intraorally, sterilize before each use.
- For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize.
- Follow the instructions for use of washer disinfectant and sterilizer, respectively.
- Do not use chlorine disinfectant liquid.
- Dura-White Stones: Do not perform plasma sterilization.
- For the details about REPROCESSING procedures, visit www.shofu.de.

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

CompoSite and CompoSite Fine polishers will polish and final polish composite restorations without scratching surrounding tooth structure. Use CompoSite for finishing and pre-polishing and CompoSite Fine for polishing.

PRECAUTIONS AND WARNINGS

- Never exceed the maximum permissible rotation speeds described in Table A.
- Use FG type with lower speed airmotor handpiece. High speed air turbine units can break cups or points.
- Insert the shank completely into the chuck of the handpiece.
- Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use.
- Perform test run to check safety and smooth operation without run-out prior to use.
- Avoid any use at acute angles or with excessive pressure.
- Use local dust extractor and protective equipment while using.
- This product is intended for use by dental professionals only.
- Do not use this product for any purposes other than the intended use.
- Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

- Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination.
- Perform autoclave under the condition described in Table B.
- When using intraorally, sterilize before each use.
- For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize.
- Follow the instructions for use of washer disinfectant and sterilizer, respectively.
- Do not use chlorine disinfectant liquid.
- For the details about REPROCESSING procedures, visit www.shofu.de.

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

This kit is designed for reducing and finishing composite fillings, also interproximal and subgingival. It has six Dura-Green and six Dura-White Stones in CN1 and FL2 shapes. Dura-Green Stones will fast finish composite restorations. Dura-White Stones will fine finish composite restorations. Dura-Green and Dura-White Stones are autoclavable for reuse.

PRECAUTIONS AND WARNINGS

- Never exceed the maximum permissible rotation speeds described in Table A. ●Insert the shank completely into the chuck of the handpiece. ●Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use. ●Perform test run to check safety and smooth operation without run-out prior to use. ●Avoid any use at acute angles or with excessive pressure. ●Use local dust extractor and protective equipment while using. ●This product is intended for use by dental professionals only. ●Do not use this product for any purposes other than the intended use. ●Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

- Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination. ●Perform autoclave under the condition described in Table B. ●When using intraorally, sterilize before each use. ●For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. ●Follow the instructions for use of washer disinfectant and sterilizer, respectively. ●Do not use chlorine disinfectant liquid. ●Do not perform plasma sterilization. ●For the details about REPROCESSING procedures, visit www.shofu.de.

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Location of the symbols glossary:
https://www.shofu.com/en/symbols_glossary/

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

SHOFU GingiCurettage Diamonds are designed to establish a chamfered finish to the tooth surface in crown and bridge preparations.

PRECAUTIONS AND WARNINGS

●Never exceed the maximum permissible rotation speeds described in Table A. ●Insert the shank completely into the chuck of the handpiece. ●Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use. ●Perform test run to check safety and smooth operation without run-out prior to use. ●Use intermittently with gentle touch and water spray. ●Avoid any use at acute angles or with excessive pressure. ●Use local dust extractor and protective equipment while using. ●This product is intended for use by dental professionals only. ●Do not use this product for any purposes other than the intended use. ●Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

●Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination. ●Perform autoclave under the condition described in Table B. ●When using intraorally, sterilize before each use. ●For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. ●Follow the instructions for use of washer disinfectant and sterilizer, respectively. ●Do not use chlorine disinfectant liquid. ●For the details about REPROCESSING procedures, visit www.shofu.de.

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nl

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PRODUCT DESCRIPTION

SHOFU Hybrid Points, diamond-coated abrasives are used for cavity, crown and bridge preparation.

PRECAUTIONS AND WARNINGS

- Never exceed the maximum permissible rotation speeds described in Table A. ● Insert the shank completely into the chuck of the handpiece. ● Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use. ● Perform test run to check safety and smooth operation without run-out prior to use. ● Use intermittently with gentle touch and water spray. ● Avoid any use at acute angles or with excessive pressure. ● Use local dust extractor and protective equipment while using. ● This product is intended for use by dental professionals only. ● Do not use this product for any purposes other than the intended use. ● Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

- Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination. ● Perform autoclave under the condition described in Table B. ● When using intraorally, sterilize before each use. ● For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. ● Follow the instructions for use of washer disinfectant and sterilizer, respectively. ● Do not use chlorine disinfectant liquid. ● For the details about REPROCESSING procedures, visit www.shofu.de.

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

This stainless steel mandrels, Straight Hand Piece will securely hold unmounted wheels for optimal control.

PRECAUTIONS AND WARNINGS

●Never exceed the maximum permissible rotation speeds described in Table A. ●When maximum permissible rotation speed of unmounted wheel is lower than 50,000 min⁻¹, follow the speed for your safety. ●Follow the instructions for use of the unmounted wheel to be used. ●Mount an unmounted wheel onto the mandrel securely. ●Insert the shank completely into the chuck of the handpiece. ●Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt. ●Perform test run to check safety and smooth operation without run-out prior to use. ●Always use an unmounted wheel in clockwise motion. ●Use local dust extractor and protective equipment while using. ●This product is intended for use by dental professionals only. ●Do not use this product for any purposes other than the intended use. ●Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

●Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination. ●When necessary, perform autoclave under the condition described in Table B. ●When necessary, sterilize before use. ●When necessary, remove any contaminants with washer disinfectant, then sterilize. ●Follow the instructions for use of washer disinfectant and sterilizer, respectively. ●Do not use chlorine disinfectant liquid. ●For the details about REPROCESSING procedures, visit www.shofu.de.

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

This kit is designed for making porcelain veneer placement. Kit has one Regular-grit and two Fine-grit FG Robot diamond points, three T&F Hybrid points, two Dura-White stones in FL2 and RE1 shapes, eight Ceramisté Standard (for prepolishing), Ultra (for polishing) and Ultra II (for superpolishing) in knife, wheel and cup shapes. Robot Points, T&F Hybrid Points and Dura-White Stones will finish porcelain smoothly. Ceramisté Polishers will polish and final polish porcelain veneers. Robot Points, T&F Hybrid Points, Dura-White Stones and Ceramisté Polishers are autoclavable for reuse.

PRECAUTIONS AND WARNINGS

• Never exceed the maximum permissible rotation speeds described in Table A. • Insert the shank completely into the chuck of the handpiece. • Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use. • Perform test run to check safety and smooth operation without run-out prior to use. • Avoid any use at acute angles or with excessive pressure. • Use local dust extractor and protective equipment while using. • This product is intended for use by dental professionals only. • Do not use this product for any purposes other than the intended use. • Dispose in accordance with relevant local and national regulation.

STORAGE AND STERILIZATION

• Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination. • Perform autoclave under the condition described in Table B. • When using intraorally, sterilize before each use. • For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. • Follow the instructions for use of washer disinfectant and sterilizer, respectively. • Do not use chlorine disinfectant liquid. • Dura-White Stones: Do not perform plasma sterilization. • For the details about REPROCESSING procedures, visit www.shofu.de.

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READ CAREFULLY BEFORE USE

en PREPARATION and REPROCESSING INSTRUCTIONS according to EN ISO 17664 for users in the European Union

Scope

These instructions are designed according to EN ISO 17664 for users in the European Union and applicable to SHOFU reusable rotary dental instruments for intraoral use. All SHOFU abrasives and polishers should be cleaned and sterilized according to these instructions before first use and each subsequent reuse.

Warnings

- Chlorine disinfectant and peracetic acid must not be used to avoid corrosion of the instruments.
- Do not expose the instruments to non-approved chemical solutions (i.e. solutions not listed in these instructions).

Limitations on reprocessing

Instruments do not have an indefinite functional life. All reusable instruments are subjected to repeated stresses related to routine use, cleaning, disinfection and/or sterilization processes. The product life is determined by wear and damage caused during use. Any damaged or defective instruments must not be used (i.g. instruments with corrosion, surface flaws, deformation, dirt and exposed shank, etc.).

Instructions

1. Initial treatment after use

Perform the following steps within one (1) hour of the patient procedure to prevent contamination from drying on the instruments:

- Detach single-use components from the reusable components and discard the single-use components.
- Prepare an enzymatic bath using 8 mL of Advanced Sterilization Products Enzol®/Cidezyme® Enzymatic Detergent per 1 L of tap water.
- Immerse instruments for a minimum of three (3) minutes.

2. Cleaning and disinfection

Method A Mechanical cleaning and disinfection (recommended)

A-1. Preparation

- The following will be required for cleaning and disinfection of the instruments:
 - Soft-bristled toothbrush
 - Critical water (water that is extensively treated, usually by a multistep treatment process that could include a carbon bed, softening, DI and RO or distillation)
 - Non-linting wipes or cloths for drying
 - neodisher® MediClean forte
 - Suitable bur block
 - ISO 15883 compliant washer disinfectant

A-2 Manual pre-cleaning

- Using a soft-bristled toothbrush, brush the instruments in enzymatic bath for a minimum of 30 seconds and until no visible contamination is present.
- Remove the instruments from enzymatic bath.
- Rinse the instruments under running tap water, for a minimum of 10 seconds.

A-3. Mechanical cleaning and disinfection

- Place the instruments into a bur block. Leave the block lid open and place the block on the shelf of an ISO 15883 compliant washer disinfectant.

- Perform automated cleaning and disinfection using the following parameters:

Phase	Minimum time	Temperature	Type of detergent/water
Pre-Wash	2 minutes	Cold	Tap water
Draining			
Wash	5 minutes	Heated (55 °C)	Tap water and 0.5% (5 mL/L) neodisher® MediClean forte
Draining			
Rinse 1	3 minutes	Cold	Critical water
Draining			
Rinse 2	2 minutes	Cold	Critical water
Draining			
Thermal Disinfection	5 minutes	Heated (≥ 90 °C)	Critical water
Dry	15 minutes	Heated	

Notes: Cleaning and disinfection validation has been performed using Miele Professional G7836 CD washer disinfector with vario-TD program. Parameters may vary depending on the washer disinfectors.

- If needed, thoroughly dry all surfaces of the instruments using non-linting wipes or cloths, changing wipes/cloths when necessary to ensure that the instruments are completely dry.

A-4. Inspection

- Visually inspect each instrument for the absence or presence of remaining contamination in a well-lit area. If contamination is present, repeat the mechanical cleaning and disinfection until all visible contamination is removed.

Note: Discard the instruments if any damage or defects (e.g. corrosion, surface flaws and/or deformations), which would prevent proper operation, are found.

Method B] Manual cleaning and disinfection

B-1. Manual cleaning

B-1-1. Preparation

- The following will be required for manual cleaning of the instruments:
- Advanced Sterilization Products Enzol®/Cidezyme® Enzymatic Detergent
 - Soft-bristled toothbrush
 - Critical water (water that is extensively treated, usually by a multistep treatment process that could include a carbon bed, softening, DI and RO or distillation)
 - Non-linting wipes or cloths for drying
 - Sonicator
 - Glass beaker

B-1-2. Manual cleaning

- Remove instruments from enzymatic bath.
- Using a soft-bristled toothbrush, brush the instruments for a minimum of 30 seconds and until no visible contamination is present.
- Rinse the instruments under running tap water, for a minimum of 10 seconds.
- Prepare enzymatic detergent solution using Advanced Sterilization Products Enzol® or Cidezyme® Enzymatic Detergent according to the detergent manufacturer's instructions using 8 to 16 mL per liter of tap water.
- Pour the enzymatic detergent solution in a glass beaker and place the beaker in a sonicator. Degas the sonicator.
- Immerse the instruments in the solution of the beaker and sonicate the instruments for 15 to 20 minutes at a temperature range of 35 °C to 45 °C and a frequency of 44±6 kHz.
- Rinse the instruments in a critical water bath for a minimum of 30 seconds.

- Thoroughly dry all surfaces of the instruments using non-linting wipes or cloths, changing wipes/ cloths when necessary to ensure that the instruments are completely dry.

B-1-3. Inspection

- Visually inspect each instrument for the absence or presence of remaining contamination in a well-lit area. If contamination is present, repeat the manual cleaning until all visible contamination is removed.

Note: Discard the instruments if any damage or defects (e.g. corrosion, surface flaws and/or deformations), which would prevent proper operation, are found.

B-2. Manual disinfection

B-2-1. Preparation

- The following will be required for manual disinfection of the instruments:
 - Trays or basins for immersion of instruments during disinfection and rinsing
 - At least limited virucidal instrument disinfectant on VAH List (e.g. CIDEX® OPA High level disinfection solution with a concentration of 0.55% of ortho-phthalaldehyde)
 - Critical water (water that is extensively treated, usually by a multistep treatment process that could include a carbon bed, softening, DI and RO or distillation) for final rinsing.
 - Sterile, non-linting wipes or cloths for drying

B-2-2. Manual disinfection

- Prepare the disinfectant solution according to the manufacturer's instructions.
- Completely immerse the instruments in the disinfectant solution.
- Remove the instruments from the disinfectant solution.
- Rinse the instruments thoroughly with critical water.
- Thoroughly dry all surfaces of the instruments using sterile, non-linting wipes or cloths, changing wipes/cloths when necessary to ensure that the instruments are completely dry.

Notes: Perform manual disinfection following the disinfectant manufacturer's instruction.

3. Storage after disinfection

Instruments should be protected from contamination until sterilization. Disinfected and dried instruments should be handled and stored in a manner that protects them from recontamination.

4. Sterilization

4-1. Preparation

The following will be required for sterilization:

- EN 13060 compliant autoclave
- Suitable bur block
- ISO 11607-1 compliant pouches for steam sterilization, as applicable

4-2. Packaging

Handle instrument(s) as follows:

- Place a single instrument in a sterilization pouch in accordance with local procedures (e.g. AAMI ST79).
- Place multiple instruments in a bur block and then place the block in a sterilization pouch in accordance with local procedures (e.g. AAMI ST79).

4-3. Sterilization

Sterilize the instruments using the appropriate parameters listed below:

Procedure	Dynamic-air-removal Steam Sterilization Cycle
Exposure time	≥ 3 minutes
Temperature	134 °C
Minimum drying time	20 minutes

5. Storage after sterilization

- Keep instruments in sterilization packaging in a dry and clean environment.
- Sterility cannot be guaranteed if packaging becomes open, damaged, or wet.
- Check the packaging and the instruments before use (packaging integrity, no excessive humidity and validity period).

Customer service

- The cleaning, disinfection and sterilization information is provided in accordance with EN ISO 17664, AAMI TIR12, and AAMI TIR30.
- These instructions have been validated by SHOFU INC. as being capable of preparing the rotary dental instruments for reuse. It remains the responsibility of the processor to ensure that the processing is actually performed using equipment, materials and personnel in the reprocessing facility to achieve the desired result. This requires validation and routine monitoring of the process. Likewise, any deviation by the processor from the recommended process in these instructions should be properly evaluated for effectiveness with potential adverse consequences.

 SHOFU INC.
11 Kamitakamatsu-cho, Fukuine, Higashiyama-ku,
Kyoto 605-0983, Japan
www.shofu.com

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READ CAREFULLY BEFORE USE

INDICATIONS FOR USE

Finishing natural teeth and restorations

PRECAUTIONS

1. This product is intended for used by dental professionals only.
2. Do not use this product for any purpose other than specifically outlined in these instructions for use.

DIRECTIONS FOR USE

1. Sterilize diamond burs before each use.
 - Sterilization method: autoclave this product for 3 min at 134 °C (273 °F).

Cycle Type	Holding time	Temperature	Dry time
Pre-vacuum (3 times)	3'00"	134°C/ 273°F	10'00"

- Remove any contaminants, using washer disinfectant.
 - Do not immerse diamond burs into chlorine disinfectants.
 - Follow the instructions for use of washer disinfectant and sterilizer, respectively.
2. Mount a diamond bur on a dental handpiece (dental air turbine, micro motor etc.).
 - Insert the shank completely into the chuck of the handpiece, following the manufacturer's instructions.
 - Perform test run to check safety and smooth operation without run-out prior to use.
 - Make sure that head of a diamond bur is in good condition before use. Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt, exposure of shank due to repeated use.
 3. Use following the standard manner.
 - Never exceed the maximum permissible rotation speed.
 - Use local dust extractor, protective mask, etc while finishing to avoid adverse effect of dust to the human body.
 - Use protective equipment such as glasses while using.
 - Use intermittently with gentle touch and water spray.
 - Avoid any use at acute angles or with excessive pressure to thin and long abrasives. They are fragile and may bend or fracture under excessive loads.
 4. For reuse after using them intraorally, remove any contaminants immediately using washer disinfectant and sterilize.

PRODUCTS

- Robot Points
- Hybrid Points
- T&F Hybrid Points
- TPE Diamonds
- GingiCurettage

Manufacturer

SHOFU INC.

11 Kamitakamatsu-cho, Fukuine,
Higashiyama-ku, Kyoto 605-0983, Japan
www.shofu.com

EC REP Authorized Representative

SHOFU DENTAL GmbH

Am Brüll 17, 40878 Ratingen, Germany

SHOFU DENTAL CORPORATION

1225 Stone Drive, San Marcos, California
92078-4059, USA

**For the details about REPROCESSING
procedures, visit www.shofu.de**

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nl

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PRODUCT DESCRIPTION

For composite, porcelain or gold restorations. Super-Snap Buff Disks have no metal centers or exposed mandrels that can discolor or gouge composite materials. It features a synthetic velvet polishing cloth outer surface which retains polishing pastes with minimum splatter or waste.

PRECAUTIONS AND WARNINGS

- Never exceed the maximum permissible rotation speeds described in Table A. • Do not re-use Super-Snap Buff Disk.
- Mount a disk onto the metal mandrel provided securely. • Insert the shank completely into the chuck of the handpiece. • Do not use damaged or defective disk or mandrel i.e. deformation, dirt corrosion or surface flaw. • Perform test run to check safety and smooth operation without run-out prior to use. • Use local dust extractor and protective equipment while using. • This product is intended for use by dental professionals only. • Do not use this product for any purposes other than the intended use. • Dispose in accordance with relevant local and national regulation.

STORAGE

Avoid exposure to moisture, corrosive chemical and its vapor. Keep away from external pressure (physical load) and avoid contamination.

STERILIZATION

Metal mandrel • Perform autoclave under the condition described in Table B. • When using intraorally, sterilize before each use. • For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize. • Follow the instructions for use of washer disinfectant and sterilizer, respectively. • Do not use chlorine disinfectant liquid. • For the details about REPROCESSING procedures, visit www.shofu.de.

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PRODUCT DESCRIPTION

Super-Snap is an easy to use, color co-ordinated, 4-step disposable disk system with metal-free center, designed for contouring, finishing and polishing of all types of direct aesthetic resins (microfilled, micro-hybrid, nano and nano-hybrid). Dedicated Super-Snap metal and plastic mandrels are available.

INDICATIONS FOR USE

Contouring, finishing and polishing of all types of direct aesthetic resin materials

PRECAUTIONS AND WARNINGS

- Do not use damaged, defective or worn metal mandrels, i.e. corrosion, surface degradation due to repeated use.
- Use protective mask and eyewear when operating handpieces.
- Avoid contacting soft tissues to avoid injury during use.
- Super-Snap Disks and plastic mandrels are designed for single use only. Discard the Super-Snap Disks and plastic mandrel after use in accordance with the disposal protocol of your healthcare facility.
- Do not use this product for any purposes other than specifically outlined in these instructions for use.
- The product is intended for use by dental professionals only.

DIRECTIONS FOR USE

- Select proper abrasive for each case, Super-Snap Black, Super-Snap Violet to remove excess restorative material and contour to the desired shape. When using Super-Snap Disks, mount the disk securely on the dedicated Super-Snap mandrel.
- Insert the shank completely into the chuck of the handpiece (micro motor), following the manufacturer's instruction. Perform test run to check safety and smooth operation without run-out prior to use.
- Apply uniform, steady pressure for finishing and feather light pressure for polishing. Unidirectional movement of the disks prevents formation of scratches and a visible step at the restoration-tooth junction.
 - Follow the recommended finishing and polishing sequence with the color coded disks for the best results. Reducing or skipping the recommended steps may compromise the quality of the final polish.
 - Use saliva ejector or suction to eliminate moisture and debris during the procedure. Use protective mask and eyewear while using the product.
 - Avoid exceeding maximum permissible speed.

Product name	Maximum permissible rotation speed	Water spray
Super-Snap Disks with metal mandrel	15,000 min ⁻¹	No
Super-Snap Disks with plastic mandrel	5,000 min ⁻¹	No

- Clean the debris from restoration surface after each finishing/ polishing step.
- After using dedicated metal mandrel, remove any contaminants immediately, disinfect and sterilize. Dispose the Super-Snap Disks and plastic mandrel after each use.

STORAGE AND STERILIZATION

Super-Snap Disks and mandrels should be stored away from direct sunlight, moisture and corrosive chemicals/vapor. Sterilize Super-Snap metal mandrels before use. For reuse after using, remove any contaminants immediately with washer disinfectant, then sterilize.

(Autoclave)

Product name	Cycle Type	Temperature	Holding time	Dry time
Super-Snap metal mandrel	Pre-vacuum 3 times	134 °C (273 °F)	3'00"	10'00"

Follw manufacturer's instructions for use of washer disinfectant and sterilizer, respectively. Chloride disinfectant and dry heat should be avoided

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For the details about REPROCESSING procedures, visit **www.shofu.de**.

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

Super-Snap Plastic Mandrel is specifically designed for use with the Super-Snap Disks.

PRECAUTIONS AND WARNINGS

- Do not re-use Super-Snap Plastic Mandrels.
- Never exceed the maximum permissible rotation speeds.
- Mount Super-Snap Disk onto the mandrel securely.
- Insert the shank completely into the chuck of the handpiece.
- Do not use damaged or defective instruments i.e. deformation, dirt.
- Perform test run to check safety and smooth operation without run-out prior to use.
- Always use in clockwise motion.
- Use local dust extractor and protective equipment while using.
- This product is intended for use by dental professionals only.
- Do not use this product for any purposes other than the intended use.
- Dispose in accordance with relevant local and national regulation.

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PRODUCT DESCRIPTION

Super-Snap is an easy to use, color co-ordinated, 4-step disposable disk system with metal-free center, designed for hassle-free, safe and efficient contouring, finishing and polishing of all types of direct aesthetic resins (microfilled, micro-hybrid, nano and nano-hybrid) to a natural gloss. The metal-free center eliminate the risk of restoration gouging or discoloration. Available in 2 sizes (diameter) 12 mm disk with larger surface area and 8 mm mini disk for safer application on cervical (gingival) area of the restoration. Super-Snap X-TREME Green and Red disks are thicker and more resilient for a better tactile sensation and optimal, faster polishing results. Super-Snap X-TREME Red disks have a 3D surface structure to minimize friction, heat and debris build-up for a smoother, glossier restoration surface. Super-Snap Polystrips are double-ended interproximal strips with the same color code as the Super-Snap Disks and incorporate two abrasive grits per strip, for finishing and polishing of the interproximal restoration surface with similar results as the Super-Snap Disks. These polystrips come with an abrasive-free plastic center for easy insertion and removal. Dura-White Stones CA, CompoSite CA and CompoSite Fine CA polishers complement the kit by facilitating the finishing and polishing of hard to reach labial/lingual surfaces as well as occlusal surfaces.

- en

READ CAREFULLY BEFORE USE

INDICATIONS FOR USE

Finishing and polishing of all types of direct aesthetic resin materials

PRECAUTIONS AND WARNINGS

- Do not use damaged, defective or worn instruments/mandrels, i.e. corrosion, surface degradation due to repeated use.
- Use protective mask and eyewear when operating handpieces.
- Avoid contacting soft tissues to avoid injury during use.
- Super-Snap Disks and Polystrips are designed for single use only. Discard the Super-Snap Disks and polystrips after use in accordance with the disposal protocol of your healthcare facility.
- Do not use this product for any purposes other than specifically outlined in these instructions for use.
- The product is intended for use by dental professionals only.

DIRECTIONS FOR USE

- Select proper abrasive for each case, Super-Snap Black, Super-Snap Violet or Dura-White Stones CA, to remove excess restorative material and contour to the desired shape. When using Super-Snap Disks, mount the disk securely on the dedicated Super-Snap mandrel.
- Insert the shank completely into the chuck of the handpiece (micro motor), following the manufacturer's instruction. Perform test run to check safety and smooth operation without run-out prior to use.
- Apply uniform, steady pressure for finishing and feather light pressure for polishing. Unidirectional movement of the disks prevents formation of secondary scratches and a visible step at the restoration-tooth junction.

➤

Follow the recommended finishing and polishing sequence with the color coded disks for the best results. Use CompoSite followed by CompoSite Fine to finish and polish hard to reach labial/lingual surfaces as well as occlusal surfaces. Reducing or skipping the recommended steps may compromise the quality of the final polish.

➤

Use saliva ejector or suction to eliminate moisture and debris during the procedure. Use protective mask and eyewear while using the product.

➤

Avoid exceeding maximum permissible speed for each abrasive.

Product name	Maximum permissible rotation speed	Water spray
Super-Snap Disks	15.000 min ⁻¹	No
Dura-White Stones	30.000 min ⁻¹	No
CompoSite	20.000 min ⁻¹	Yes
CompoSite Fine	20.000 min ⁻¹	Yes

- Clean the debris from restoration surface after each finishing/polishing step.
- After using dedicated metal mandrel, Dura-White Stone, CompoSite and CompoSite Fine polishers, remove any contaminants immediately, disinfect and sterilize. Dispose the Super-Snap Disks after each use.

Super-Snap Polystrips

Insert the abrasive free center of Super-Snap Polystrip in the interproximal area. Follow the same color coded abrasive grit as Super-Snap Disks. Adapt over the restoration surface to be finished and polished and apply uniform pressure in back and forth motion, starting with the coarse section followed by the finer section of the strip. Eg. Black followed by Violet and then, Green followed by Red.

- Do not reuse Super-Snap Polystrips.
- Exercise caution during use near the cervical area to avoid injuring the gum with the edge of the Polystrip.
- Avoid bending the strip at sharp angles.

STORAGE AND STERILIZATION

Super-Snap Disks and Polystrips should be stored away from direct sunlight, moisture and corrosive chemicals/vapor. Sterilize Dura-White Stones, CompoSite, CompoSite Fine and Super-Snap metal mandrels before use.

Remove any contaminants using suitable cleaning liquid or disinfectant and sterilize.

Product name	Autoclave
Super-Snap metal mandrel	134 °C (273 °F) – 3'00"
Dura-White Stones	or
CompoSite	121 °C (250 °F) – 30'00"
CompoSite Fine	

Follow manufacturer's instructions for use of cleaning liquid, disinfectant or sterilizer. Chloride disinfectant and dry heat should be avoided. Plasma Sterilization should also be avoided for sterilization of Dura-White Stones as it leads to degradation of the bond between the head and shank.

CAUTION: US Federal law restricts this device to sale by or on the order of a dental professional.

Location of the symbols glossary: https://www.shofu.com/en/symbols_glossary/

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nl

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PRODUCT DESCRIPTION

SHOFU T&F Hybrid Points, superfine diamond coated abrasives are used for finishing enamel, porcelain and composites with best results.

PRECAUTIONS AND WARNINGS

- Never exceed the maximum permissible rotation speeds described in Table A.
- Insert the shank completely into the chuck of the handpiece.
- Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use.
- Perform test run to check safety and smooth operation without run-out prior to use.
- Use intermittently with gentle touch and water spray.
- Avoid any use at acute angles or with excessive pressure.
- Use local dust extractor and protective equipment while using.
- This product is intended for use by dental professionals only.
- Do not use this product for any purposes other than the intended use.
- Dispose in accordance with relevant local and national regulation.

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READ CAREFULLY BEFORE USE

PRODUCT DESCRIPTION

SHOFU TPE, Tissue-Protective End-cutting,
Diamonds are recommended for time saving
preparation of subgingival shoulder margins.

PRECAUTIONS AND WARNINGS

- Never exceed the maximum permissible rotation speeds described in Table A.
- Do not use damaged or defective instruments i.e. corrosion, surface flaw, deformation, dirt and exposure of shank due to repeated use.
- Perform test run to check safety and smooth operation without run-out prior to use.
- Use intermittently with gentle touch and water spray.
- Avoid any use at acute angles or with excessive pressure.
- Use local dust extractor and protective equipment while using.
- This product is intended for use by dental professionals only.
- Do not use this product for any purposes other than the intended use.
- Dispose in accordance with relevant local and national regulation.

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