
Implant Restoration | Permanent Cementation

Check fit, occlusion, proximal contacts and complete all finishing procedures.

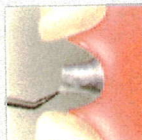Dry and isolate the prepared abutments. Absolute dryness is not necessary, but obvious contamination should be avoided.

Make sure the interior of the restoration is clean and dry. To facilitate future removal of implant restorations, lubricate the surface of the abutment with a thin layer of KY Jelly® or other water-based lubricant.

Bleed the syringe before installation of the catalyst. Depending on the size of the crown, dispense 10-20 drops of the catalyst into the syringe and bleed the syringe again. This ensures complete flow of base and catalyst components. Squeeze the thin layer of the cement directly into the preparation. You have 45-60 seconds working time.

Immediately seat the crown(s) and either hold in place or have the patient gently close on a cotton roll until the initial set which is 2-2.5 minutes. The initial set is characterized by a distinct rubbery phase in which excess cement can be easily removed using an explorer or periodontal probe.

Final cement set is between 4-5 minutes, depending on conditions.

Precautions: Though rare, some individuals may be sensitive to this product. If an allergic reaction occurs, discontinue use, remove from patient and refer patient to a physician. Important Application Information: Only a small bead of Crown Set around the margin is necessary to securely hold your crown. This is especially true with highly retentive restorations. Using less cement saves money and allows for easy retrieval!

MIS Implant Technologies Ltd.
P.O. Box 7, Van Lee Industrial Park,
2015000, Taipei
Website: www.mis-implant.com

CE
0483

DS REP BV
Schipholweg 55, 2718 SG
Zoetermeer, Netherlands
Tel: +31 79 347 1288

Manufacturer

Use by

Date of Manufacture

Batch Code

Catalogue Number

Caution, consult
accompanying documents

Upper Limit of Temperature

© MIS Implant Technologies Ltd. All Rights Reserved

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