EnSite Precision™ Surface Electrode Kit

REF

EN0020-P

Instructions for Use

This warranty is given by Abbott Medical. To make inquiries regarding this warranty, use the contact information for either the manufacturer or Australian Sponsor located on the back cover.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

The benefits to you given by this warranty are in addition to your other rights and remedies under the Australian Consumer Law.

Abbott Medical warrants to the buyer that the goods, for a period equal to the validated shelf life of the goods (the "Warranty Period"), shall meet the product specifications established by the manufacturer when used in accordance with the manufacturer's instructions for use and shall be free from defects in materials and workmanship. Unless otherwise obligated by law or mandatory consumer guarantees, under this warranty Abbott Medical will, at its option, replace or repair the goods at its factory, if the buyer, at the buyer's expense, returns the goods within the Warranty Period to Abbott Medical to the address below and after Abbott Medical confirms that the goods are defective.

Abbott Medical
Roseville Distribution Center
Attn: Postmarket Surveillance
2305 Walnut Street
Roseville, MN 55113
Tel.: (651) 756-2000

To the maximum extent permitted by law, if a mandatory term is implied by law or a mandatory consumer guarantee applies to the goods, and the goods are not of a kind ordinarily acquired for personal, domestic or household use or consumption, Abbott Medical's liability for the breach of the term or guarantee is limited to, at Abbott Medical's option, either replacement or repair of the goods or payment of the costs of replacing or repairing the goods.

EXCEPT FOR MANDATORY TERMS IMPLIED BY LAW, MANDATORY CONSUMER GUARANTEES, OR AS EXPRESSLY PROVIDED IN THIS WARRANTY, ABBOTT MEDICAL DISCLAIMS ANY REPRESENTATION OR WARRANTY OF ANY KIND, EXPRESS OR IMPLIED.

See the Terms and Conditions of Sale for further information.

[™] Indicates a trademark of the Abbott group of companies.

[‡] Indicates a third-party trademark, which is property of its respective owner.

Pat. <http://www.abbott.com/patents>

© 2023 Abbott. All Rights Reserved.

EN: English

EnSite Precision™ Surface Electrode Kit

Instructions for Use

Description

The EnSite Precision™ Surface Electrode Kit is used to transmit information from the patient to the EnSite Precision™ Cardiac Mapping System for signal conditioning and data inference. The kit contains the following single-use, non-sterile devices for use in conjunction with the EnSite Precision™ Cardiac Mapping System:

Device	Model	Quantity	Description
EnSite™ System Reference Electrode	SYS-REF-V1	One (1)	A reference electrode enclosed between adhesive foam and a layer of silver carbon film and conductive hydrogel with an integrated cable. The EnSite™ System Reference Electrode serves as a system electrical reference for impedance and ECG measurements.
3M® Red Dot® Radiolucent Monitoring Electrodes	2570-5	Two (2) Packs - Ten (10) total ECG electrodes	Self-adhesive, non-sterile, single use disposable ECG electrodes which include a silver/silver chloride sensing element and conductive sticky gel.
EnSite™ Locating Electrodes	SURF-ELEC-2PRS	One (1)	The EnSite™ Locating Electrodes consist of the following: • EnSite Precision™ Surface Electrodes (Six (6)) – An electrode enclosed between adhesive foam and a layer of silver carbon film and conductive hydrogel. The electrodes are color-coded and have a placement diagram printed on the surface. The electrodes serve as diagnostic tools that impart low-intensity alternating currents generated by the EnSite Precision™ system to the patient's thorax and sense the resultant electrical potentials which supports the acquisition of impedance signals used for locating the position of electrodes on conventional EP catheters. • Patient Reference Sensor (PRS) Patches (Two (2)) - A foam adhesive patch that is intended to hold a PRS in a stable position on the patient during an EnSite Precision™ Cardiac Mapping System procedure.

Indications for Use

Device	Indications for Use
EnSite Precision™ Surface Electrode Kit	The EnSite Precision™ Surface Electrode Kit is indicated for use with the EnSite Precision™ Cardiac Mapping System in accordance with the EnSite Precision™ Cardiac Mapping System indication for use.
EnSite™ System Reference Electrode	The EnSite™ System Reference Electrode is indicated for use with the EnSite™ Cardiac Mapping System in accordance with the EnSite™ Cardiac Mapping System indication for use.
3M® Red Dot® Radiolucent Monitoring Electrodes	The 3M® Red Dot® 2560 and 2570 Monitoring Electrodes are intended for use in ECG Monitoring. The 2560 and 2570 electrodes can be used in all ECG applications where standard ECG monitoring electrodes are used. These electrodes can be used for short term and long term (up to 5 days) monitoring.
EnSite™ Locating Electrodes	The EnSite™ Locating Electrodes are indicated for use with the EnSite Precision™ Cardiac Mapping System in accordance with the EnSite Precision™ Cardiac Mapping System indication for use.

Clinical Benefits

As supplementary devices, these products contribute to the realization of the intended clinical benefit of the primary, compatible EnSite™ Cardiac Mapping Systems for which they support.

Contraindications

There are no known contraindications.

Precautions

- Intended for single use only. Device integrity will be compromised by any reuse, which may compromise patient safety and system performance.
- Store in a cool, dry, dark place. Inspect the packaging and surface electrode(s) prior to use.
- Do not attempt to use this product prior to completely reading and understanding these Instructions for Use.
- Do not alter device(s) in any way.
- This device must be operated by, or under the supervision of, an electrophysiologist trained in the operation of this device and supported by other qualified personnel trained in the field of cardiac EP.
- Do not warm surface electrode(s) before positioning on the patient.
- Do not use surface electrodes if the package seal is broken, the conductive adhesive is dry, the "use by" date has passed or the gel column is dislodged from the electrode.
- Before applying the surface electrode(s), ensure that the application sites are hair-free, clean, non-irritated, and dry, to prevent receiving low quality signals from body surface electrodes. Dryness is particularly important if the preparation agent is flammable or a potential skin irritant.

- It is recommended to avoid removing the liner until the surface electrode is ready to be placed on the patient. Once the liner has been removed, do not attempt to stick the liner back on the adhesive side of a surface electrode. Doing so may result in the liner becoming stuck to the surface electrode.
- To minimize the risk of electric shock, make sure that electrodes and lead connectors do not contact one another or contact ground. Use care when attaching the ECG electrodes to lead connectors.
- For proper skin management and to minimize skin irritation, do not abrade skin more than one time.
- During surgical procedures, place ECG electrodes as far as possible from the electrosurgical site to minimize RF current flow through the electrode site. Otherwise burns could result at the ECG electrode site.
- Discontinue use of this product if skin sensitivity is observed.
- When removing the EnSite™ System Reference Electrode from the patient's abdomen, remove the electrode from the patient before disconnecting from the EnSite™ system hardware.
- To reduce likelihood of skin damage when removing, slowly pull back surface electrode(s) at a low angle with one hand while holding the skin taut with the other hand.
- Do not attempt to relocate the device(s) after initial application.
- Replace if self-adhesive electrodes if they no longer stick firmly to skin.
- The 3M® Red Dot® Monitoring Electrode was determined to be MR-Conditional via non-clinical testing. For information on the use of a 3M® Red Dot® monitoring electrode during an MRI procedure please see below.

Limitations

- Do not place the surface electrode(s) over devices such as pacemakers or ICDs.
- Device(s) should only be connected to the EnSite Precision™ Cardiac Mapping System.

Adverse Events

Skin and subcutaneous tissue disorders including irritation, damage, sensitivity, burns, ecchymosis, inflammation, and allergic skin reaction.

MRI Safety Information

Important Notes

- 3M® MR Conditional EKG/ECG Electrodes have not been tested for MRI-related issues with any MR Conditional monitoring system. Therefore, the acceptable use of 3M® MR Conditional EKG/ECG Electrodes when connected to an MR Conditional monitoring system is unknown.
- Disconnected 3M® MR Conditional Electrodes can be worn safely during an MRI examination according to the conditions of use presented in the MRI Safety Information labeling. Any deviation from those specific conditions may result in injury (e.g., burn) to the patient.

MR Conditional

The 3M® Red Dot® Monitoring Electrode with Foam Tape and Sticky Gel 2570 series was determined to be MR Conditional. Non-clinical testing demonstrated that the 3M® Red Dot® Monitoring Electrode with Foam Tape and Sticky Gel 2570 series is MR Conditional. A patient with these electrodes can be scanned safely immediately after placement under the following conditions:

- Static magnetic field of 3-Tesla or less
- Maximum spatial gradient magnetic field of 4,000-gauss/cm (40-T/m)
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 4-W/kg for 15 minutes of scanning (i.e., per pulse sequence) in the Normal Operating Mode

Under the scan conditions defined, the 3M EKG/ECG electrode is expected to produce a maximum temperature rise of 2.5°C after 15 minutes of continuous scanning (i.e., per pulse sequence).

In non-clinical testing, the image artifact caused by the 3M EKG/ECG electrode extends approximately 3-mm from this device when imaged using a gradient echo pulse sequence and a 3-Tesla MR system.

NOTE: The EnSite™ Locating Electrodes and the EnSite™ System Reference Electrode have not been determined to be MR-Conditional.

Directions for Use

NOTE: Please consult the EnSite Precision™ Cardiac Mapping System IFU before using this kit.

Skin Preparation

NOTE: When prepping the skin sites for PRS patch placement, do not abrade skin.

NOTE: 3M® Red Dot® electrodes have a built-in abrader on the back of the liner that can be used to abrade skin.

Good skin preparation is essential to all electrophysiology studies. To help minimize noise related to surface electrode placement, perform these steps at application site(s) before placing the device(s) on the patient's skin:

1. Trim excess hair.
2. Gently abrade the skin using a gauze pad or like material or gel prep (such as Omni Prep) that contains small particles of gritty material. Make sure all the gel material is removed from the skin.
3. Clean skin surface with soap and water. Avoid cleaning with isopropyl alcohol unless the skin is excessively oily or covered with lotion.
4. Ensure the skin is completely dry and that all alcohol has evaporated.

Placement Guidelines

When placing the EnSite Precision™ Surface Electrodes on the patient, it is important to place them in the correct position, not only in relationship to each other, but also in relationship to other surface electrodes.

It is recommended to avoid opening the package and removing the liner until the surface electrode is ready to be placed on the patient. Once the liner has been removed, do not attempt to stick the liner back on the adhesive side of the surface electrode. Doing so may result in the liner becoming stuck to the surface electrode.


Front (Anterior) and Back (Posterior) Electrodes

The precordial lead V2 ECG electrode placement at the 4th intercostal space and the left sternal border guides the placement of the anterior EnSite Precision™ Surface Electrode. If both are placed correctly, the anterior EnSite Precision™ Surface Electrode will be offset slightly to the left of the patient's sternum. It is important that the posterior EnSite Precision™ Surface Electrode be placed directly opposite of the anterior EnSite Precision™ Surface Electrode. As a result, it should be offset slightly to the left of the patient's spine. This placement will help maintain appropriate anteroposterior (AP) orientation. Deviation from the recommended location will skew the AP orientation in the map display.

Right and Left Electrodes

The right and left side EnSite Precision™ Surface Electrodes must also be directly across from each other and orthogonal to the AP surface electrode arrangement. Any other placement may skew side-to-side orientation.

Figure 1. Caudal view of thorax - Correct (left) and incorrect (right) placement

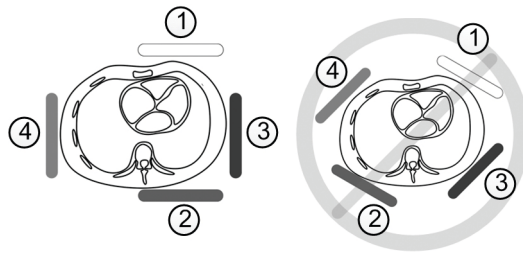EnSite Precision™ Surface Electrodes:

1. Front (white)
2. Back (brown)
3. Left (red)
4. Right (orange)

Neck and Left Leg Electrodes


The neck and left leg EnSite Precision™ Surface Electrodes should be vertically aligned to maintain an orthogonal relationship to the other four EnSite Precision™ Surface Electrodes. The neck EnSite Precision™ Surface Electrode can easily be placed at the midline. Placement of the left leg EnSite Precision™ Surface Electrode on the inner thigh is recommended to achieve as near of a midline location as possible. Any other placement may skew the vertical orientation.

NOTES:

- None of the EnSite Precision™ Surface Electrodes should overlap each other or overlap any other surface electrodes.
- Edge-to-edge contact with ECG electrodes should not create interference.
- Avoid direct contact of any EnSite Precision™ Surface Electrodes with an RF (radiofrequency) return patch.

Surface Electrode Placement

Figure 2. Device placement

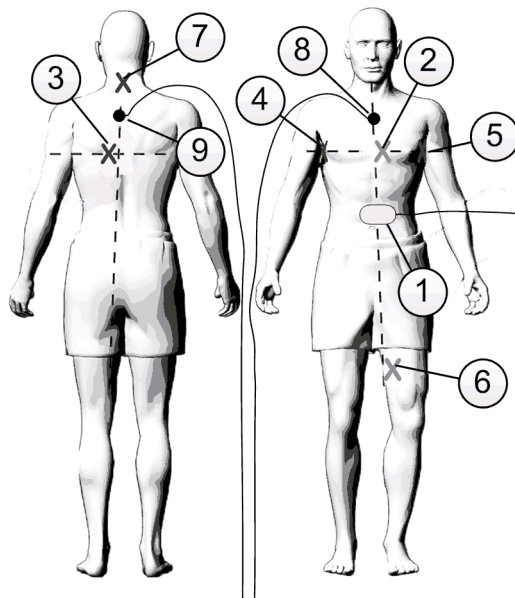1. EnSite™ System Reference Electrode

EnSite Precision™ Surface Electrodes:

2. Front (white)
3. Back (brown)
4. Right (orange)
5. Left (red)
6. Left Leg (blue)
7. Neck (green)

PRS Patches:

8. Front PRS patch (PRS-A)
9. Back PRS patch (PRS-P)

1. Place the EnSite™ System Reference Electrode on the patient's abdomen as shown in the figure above.
2. Place the V2 ECG electrode at the 4th intercostal space, left sternal edge.

NOTE: Avoid pressing on the center of the ECG electrode. Smooth the ECG electrodes down using a circular motion.

3. Place each EnSite Precision™ Surface Electrode on the patient per the illustration printed on the foam of the EnSite Precision™ Surface Electrode. See the device placement diagram above for proper placement. The orientation of the EnSite Precision™ Surface Electrode leads is not important.

NOTES:

- Grasp an area away from the lead wire when removing the liner from EnSite Precision™ Surface Electrodes.
 - It is recommended to place the front EnSite Precision™ Surface Electrode first.
 - Press firmly on the center of the EnSite Precision™ Surface Electrodes with a circular motion to ensure proper contact.
4. Place the remaining ECG electrodes on the patient.
 5. Connect each surface electrode lead to EnSite Precision™ system hardware as instructed in the EnSite Precision™ Cardiac Mapping System IFU.

NOTE: The EnSite™ System Reference Electrode must be the first patient electrode connected to the EnSite™ system hardware at the beginning of a study and the last patient electrode to be disconnected at the end of the study.

PRS Patch Placement

1. Remove the liner from the underside of the PRS patch, exposing the adhesive.
2. Place a PRS patch on the front and back of the patient as shown in the figure above.

NOTES:

- Avoid removing and repositioning the PRS patch.
 - The PRS patch must be placed with the PRS wire channel pointed toward the patient's head.
3. For each of the two PRS patches, remove the liner from the front of the PRS patch, exposing the adhesive. This is where the PRS will be placed.

PRS Patch Removal

1. Peel off each PRS patch from the body using the white foam pull tab, leaving the PRS attached.

2. Release the PRS from the PRS patch by grasping and pulling the white foam pull tab and the colored pull tab in opposite directions.
3. Remove the PRS and dispose of the PRS patch.

Surface Electrode Removal

To minimize skin trauma during ECG electrode removal, pull up slowly and gently from the edge of the ECG electrode.

When removing EnSite Precision™ Surface Electrodes and the EnSite™ System Reference Electrode, slowly pull back at a low angle with one hand while holding the skin taut with the other hand.

NOTE: Soap and water may aid in surface electrode removal.

Disposal

- The IFU is recyclable.
- After use, device(s) and packaging should be appropriately classified for disposal, e.g. electrical equipment, non-hazardous waste, etc. and carefully disposed of in compliance with facility procedures and applicable laws and regulations. Discard any unused components after the procedure.

Reporting Device Incidents

The user and/or patient should report any serious incident that has occurred in relation to this product to the manufacturer and the competent authority of the Member State in which the user and/or patient is established. Manufacturer contact information is located on the back cover of this IFU.

Limited Warranty

Abbott Medical warrants that its products shall be free from defects in materials and workmanship under normal use. This warranty does not exceed the "Expiration" date stated on any product labeling. The authorized uses and approved methods of use of each of our products is set forth in the related "Instructions for Use" that accompany each product. Abbott Medical disclaims any responsibility and liability for the use of its products in a manner that has not been authorized or approved. Abbott Medical's liability under this warranty is limited to replacing its products. The foregoing warranty excludes and is in lieu of all other warranties whether expressed or implied including, but not limited to, any implied warranties of merchantability or fitness for a particular purpose. Abbott Medical disclaims any liability for any incidental or consequential loss, damage, or expense directly or indirectly arising from the use of this product, other than as expressly provided by specific law. Abbott Medical neither assumes nor authorizes any other person to assume for it any other or additional liability for loss, damage, or expense in connection with this product. For more details please review complete Abbott Medical warranty policy available from Abbott Medical or on the back of an Abbott Medical invoice.

Symbols

Symbol	Description
	Assembler of Procedure Pack
Surface Electrode Kit	Surface Electrode Kit
Locating Electrodes	Locating Electrodes
System Reference Electrode	System Reference Electrode
Radiolucent Monitoring Electrode	Radiolucent Monitoring Electrode
MD	Medical Device
UDI	Unique device identification number
	Importer
	Not made with natural rubber latex.
	"The Green Dot" trademark on packaging means that for such packaging, a financial contribution has been paid to a national packaging recovery company that has been set up in accordance with the principles defined in European Directive No. 94/62 and the corresponding national law.
	Affixed to this device in accordance with European Council Directive 2012/19/EU. This directive calls for separate collection and disposal of electrical and electronic equipment. Sorting such waste and removing it from other forms of waste lessens the contribution of potentially toxic substances into municipal disposal systems and into the larger ecosystem.
CE	European conformity, affixed in accordance with the relevant provisions of European Council Regulation 2017/745. Hereby, Abbott Medical declares that this device is in compliance with the relevant provisions of this regulation.
UK CA	United Kingdom Conformity Assessed
UKRP	Responsible person in the United Kingdom

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3M Deutschland GmbH