# Operating Instructions

## Modular Hand Instruments for Laparoscopy

**Manufacturer:** DIMEDA Instrumente GmbH  
Gänsäcker 54 + 58  
D-78532 Tuttlingen  
T +49 (0) 7462-9461-13  
F +49 (0) 7462-9461-33  
Email: [info@dimeda.de](mailto:info@dimeda.de)  
[www.dimeda.de](http://www.dimeda.de)

**Made in Germany**

### Products:

- LAP-Dissecting scissors with HF connector
- LAP-Forceps with HF connector

By purchasing this instrument, you receive a quality product whose proper handling and use are presented below. In order to keep risks to patients and users as low as possible, we ask you to follow the instructions carefully. The application, disinfection, cleaning and sterilization of instruments must be carried out by trained specialists.

### General Information

1. Inspect the instrument for cleanliness, functionality and for damage after each cleaning and disinfection and prior to each use. These may include warped, cracked, worn or fractured components, as well as damage to the insulation.
2. The LAP instruments can bear a traction force of up to 50N (approx. 5,0 kg). If these forces are exceeded, damage to the handles, electrodes, jaw parts, tie rods and joints may occur.
3. Do not use damaged or defective instruments. Replace damaged items immediately with original spare parts.

To prevent damage to the electrodes, insert the instrument through the trocar.

Visually check the instrument before each use for warped, broken, cracked, worn or fractured components. Exclude damaged instruments immediately.

These instructions apply to all items with the number range 92.010.21A - 92.200.16. The following articles are excluded:

- All articles with the length >450mm.

### Intended Purpose

Laparoscopy is an optical examination and diagnostic procedure in which the abdominal cavity and the organs contained therein (liver, gall bladder, spleen, uterus, etc.) are mirrored with an endoscope. During a laparoscopy, in addition to the diagnosis, surgical procedures and biopsies are also made possible whereby surgical micro-instruments are utilised. These instruments are inserted through a trocar into the abdominal cavity.

In accordance to procedures, micro-surgical (MIS instruments) laparoscopic instruments with HF connectors are utilised to manipulate, grip and cut tissues and organs. DIMEDA Instrumente GmbH's separable laparoscopic instruments with HF connectors are suitable when using mono-polar HF current. If indicated, the instruments can be implemented for the specific coagulation of tissue and vessels.

### Indications

- General abdominal surgery
- Surgical procedures of the bile ducts
- Abdominal and oesophageal surgical procedures
- Small intestinal surgical procedures
- Colonic surgical procedures
- Hernia repair procedures

DIMEDA Instrumente GmbH's separable laparoscopic instruments with HF connectors are suitable for using monopolar HF current.

## Modular Hand Instruments for Laparoscopy

In accordance to procedures, micro-surgical laparoscopic forceps with HF connectors (MIC instruments) can be used to grasp and dissect tissue, organs and blood vessels. When used in conjunction with HF current, the instruments can be simultaneously used to dissect and coagulate tissue and vessels.

The different working lengths of the products are as follows:

- Ø3mm, Ø3.5mm, Ø5mm and Ø10mm instruments with a working length of 220 and 250mm working length are normally used in pediatrics.
- Ø3mm, Ø3.5mm, Ø5mm and Ø10mm instruments with a working length of 330mm working length are normally used for normal weight patients
- Ø3mm, Ø3.5mm, Ø5mm and Ø10mm instruments with a working length of 450mm are normally used in patients with obesity.

## Contra-indications

- The products are solely designed for use with mono-polar HF current. The use of products with bipolar HF current is contra-indicative and can cause serious injury to the patient, users or third parties.
- The products must not be used in the vicinity of flammable or explosive gases or other fire-sensitive media
- The products are not suitable for use in single port laparoscopy. The close proximity to other metallic instruments during these interventions can lead to the unwanted dissipation of the current (antenna coupling) and thus to serious complications such as burns away from the operating field when using monopolar HF current.
- The products should only be made with careful consideration of the suitability of laparoscopic procedures using monopolar instruments and are contraindicated for all other applications
- The products should solely be used "under direct vision". Coagulation should only take place if the instrument is in the field of view of the surgeon. This reduces the risk of unintentional contact with other metallic instruments.

## Hazards

- Injury to nerves, blood vessels and organs
- Circulatory disorders (abdominal inflation due to the build-up of Carbon Dioxide, and in particular, patients with existing pulmonary disorders)
- Instrument breakage caused by the exertion of excessive force by the operator (extension of operating times and rarely, by the prolonged entrapment of fragments) or by stress corrosion cracking caused by faulty reprocessing
- Electrical burns
- Infections caused by incorrectly reprocessed instruments
- Capacitive coupling (with 2 conductors flooded with current, consequential error: can lead to temperature development)

HF current can damage cardiac pacemakers. Before the procedure, Consult cardiologists.

## Combination with other products / instruments

Components of DIMEDA sets and DIMEDA individual products are compatible with each other. The use of DIMEDA single products / sets with Third-party products are not permitted.

## Modular Hand Instruments for Laparoscopy

The Dimeda instruments can be connected to the following HF generators with the appropriate cables:

| Manufacturer | Generator            | Electric wire |
|--------------|----------------------|---------------|
| Aesculap     | GN640                | GN202         |
| Olympus      | ESG-400              | A0392         |
| KLS Martin   | ME 102               | 80-332-03-04  |
| Berchtold    | Elektrotom 80 / 80 B | 90.280.31     |

### Handling

1. The instruments must not be placed under excessive force whilst twisting or levering, as this may result in damage or breakage to the parts of the instrument.
2. When used as intended, the HF-device can cause sparks. There is a risk of injury from the ignition or explosion of combustible gases. The safety information contained within the operating instructions of the HF-device must be observed.
3. Regulate the output power of the HF-device in accordance to the procedure. Take into consideration clinical experience and references.
4. Maintain a clean contact surface area on the HF-device during the operation. Wipe away dried tissue residue or body fluids with a damp swab.
5. The insulation of the instrument is designed for a maximum rated alternating voltage of 2 kVp. Based on the operating instructions of the HF-device, check the ratio of the maximum rated alternating voltage to the HF modus and the configured dosage.

To avoid unintentional HF burns, only use the instrument if the specified maximum rated voltage for the instrument is equal to or greater than the maximum set output voltage of the HF-device.

The electrodes of the activated instrument must be in the field of vision of the user. Before activating the HF-device, ensure that electrodes of the instrument do not come into contact with any conductive accessories or any conductive liquids. In the case of endoscopic procedures, where contact with active instruments cannot be excluded, use insulated accessories. Prior to each use, visually inspect the insulation for damage and changes to the surface.

Immediately replace damaged instruments and items with original spare parts. Switch off the "automatic power-on" modus on the HF-device.

**Always observe the Operating Instructions of the HF-device**

### Assembly

The instrument is to be assembled as follows:

1. a)  $\varnothing 5$  mm und  $\varnothing 10$  mm instrument inserts:  
Feed the instrument insert into the guide tube up to the thread and then tighten the screw. The two parts are thus joined together. b)  $\varnothing 3$  mm attachments:  
These instrument attachments are not separable, i.e. the guide tube and the instrument insert are an inseparable, connected unit.
2. In the next step, either the separable  $\varnothing 5$  mm and  $\varnothing 10$  mm instrument inserts and guide tube or the inseparable  $\varnothing 3$  mm attachments are to be connected to the modular handle. Here, the jaws must be **completely closed** and the handle **completely open**, so that the receptacles for the instrument attachments are exposed to a maximum. Thereafter, mount the orb on the instrument attachment into the bore hole of the receiver on the handle.
3. The hinged instrument attachment will now be secured by closing the handle.

# Operating Instructions

## Modular Hand Instruments for Laparoscopy

4. Tightly screw and attach the black union nut on the handle to the instrument attachment. The instrument is ready for use.

### Disassembly

The instrument is to be **disassembled** as follows:

1. With the handles closed, unscrew the black union nut on the instrument attachment.
2. Remove the orb on the instrument attachment from the bore hole of the receiver on the handle.
3. Unscrew the instrument insert ( $\varnothing 3,5/5/10$  mm) within the guide tube and remove.
4. Hereby, the instrument is completely dismantled.

### Cleaning, sterilization and maintenance

Due to the product design and the materials used, 50 reprocessing cycles can be carried out for the instrument and "inserts firmly attached to the shaft" and 100 reprocessing cycles for the handles.

Defective products must have gone through the entire reprocessing process before being returned for use.

### Preparation for processing

The instruments must be disassembled or opened for reprocessing.

### On-site preparation

Remove coarse dirt from the instruments directly after each use. Do not use a fixating agent or hot water ( $>40^{\circ}\text{C}$ ) as this can cause residues to fuse, which may influence the result of the reprocessing process.

### Transportation

Safe storage in a closed container and transport of the instruments to the reprocessing point to avoid damage to the instruments and contamination of the environment.

### Manual pre-cleaning

1. Disassemble the instruments
2. Flush the stem with a water gun for 10 seconds (static pressure 3.8 bar). The flush adapter 5100-905 is used to flush the shaft.
3. The instruments for 5 min. soak in cold water
4. Brush the instruments with a soft brush until all visible residue is removed.
5. Flush all openings and hard-to-reach areas, bores and threads with a water gun for 5 seconds (static pressure 3.8 bar). Flush the cavities with a water gun for 10 seconds. The flush adapter 5100-905 is used to flush the shaft.
6. The instruments 10 min. Place in an ultrasonic bath for a long time and sonicate in the cleaning solution (0.5% neodisher MediClean forte,  $40^{\circ}\text{C}$ ).
7. Flush the cavities with a water gun for 10 seconds (static pressure 3.8 bar). The flush adapter 5100-905 is used to flush the shaft.

### Manual cleaning

1. Disassemble the instruments
2. Flush the stem with a water gun for 15 seconds (static pressure 3.8 bar). The flush adapter 30-5100-905 is used to flush the shaft.
3. The instruments for 5 min. place in the cleaning solution (0.5% neodisher MediClean forte,  $40^{\circ}\text{C}$ ).
4. Brush the instruments with a soft brush until all visible residue is removed.
5. Flush all openings and hard-to-reach areas, bores and threads with a water gun for 5 seconds (static pressure 3.8 bar). Flush the shaft with a water gun for 15 seconds. The flush adapter 30-5100-905 is used to flush the shaft.

# Operating Instructions

## Modular Hand Instruments for Laparoscopy

6. The instruments 9 min. Place in an ultrasonic bath for a long time and sonicate in the cleaning solution (0.5% neodisher MediClean forte, 40 ° C).
7. Flush the stem with a water gun for 60 seconds (static pressure 3.8 bar). The flush adapter 30-5100-905 is used to flush the shaft. And then rinse the outer surfaces with cold water for 60 seconds while moving the moving parts.

### Mechanical cleaning

#### Cleaning / Sterilising equipment: G 7836 CD (Miele)

Cleaning program: Oxivario

| Step | Time (min) | Stage of Process | Reagents                                                                              | Temp. (°C) |
|------|------------|------------------|---------------------------------------------------------------------------------------|------------|
| 1    | 2          | Pre-cleaning     | Mains water                                                                           | Cold       |
| 2    |            | Drain water      |                                                                                       |            |
| 3    | 5          | Cleaning         | Mains water<br>Dosierung: 0.5% Neodisher<br>MediClean Forte (Dr. Weigert,<br>Hamburg) | 55         |
| 4    |            | Drain water      |                                                                                       |            |
| 5    | 3          | Neutralisation   | Deionized water                                                                       | 10-25      |
| 6    |            | Drain water      |                                                                                       |            |
| 7    | 2          | Rinse            | Deionized water                                                                       | 10-25      |
| 8    |            | Drain water      |                                                                                       |            |

### Disinfection

Carry out mechanical thermal disinfection taking into account the national regulations and the A0 value = 3000.

### Drying

Dry the outside of the instruments by utilising the drying cycle of the cleaner/steriliser. In addition, when necessary, instruments can be dried manually using a lint free cloth. Instrument cavities can be dried using sterile compressed air.

### Function tests, maintenance

Visual inspection for cleanliness; Assemble the instruments, care and function test according to the operating instructions.

If necessary, repeat the reprocessing process until the instrument is optically clean.

### Packaging

The instruments are double in sterilization bags packed (Wipak STERIKING flat rolls type R43 / Type R44) according to EN ISO 11607-1.

### Sterilisation

#### Sterilizer: Selectomat HP 666-1HR

Before the instrument components are prepared for sterilisation, the surface, and in particular, all the moving parts should be carefully lubricated. For this we recommend grease-free and temperature resistant silicone

# Operating Instructions

**dime**  
SURGICAL INSTR

## Modular Hand Instruments for Laparoscopy

products. This helps to keep the moving parts and threads accessible and also protects the entire instrument surface from mineral deposits, which can later lead to functional impairments. The lubrication product must be biocompatible and approved for medical devices. Please note that this instrument lubrication should be carried out routinely after each cleaning (ultrasound etc.) and before each sterilization.

Sterilisation of products with fractional pre-vacuum procedure (in accordance with EN ISO 17665) taking into consideration of the respective national requirements.

- 3 pre-vacuum phases with a pressure of at least 60 mBar
- Heating up to a sterilisation temperature of min. 132°C and max. 137°C
- Shortest exposure time: 5 min.
- Drying time: at least 10 min

If you suspect contamination with prions (CJD), prepare and dispose of the instrument.

### Storage

Store sterilised instruments in a dry, clean and dust-free environment at moderate temperatures, ranging from 5°C to 40°C.

### Additional instructions

If the chemicals and machines described here are not available and the preparation process cannot be carried out as described, it is the responsibility of the user to validate their process accordingly.

### Further information on the reprocessing of medical devices

- Internet: <http://www.a-k-i.org>
- Hygiene requirements involved in the reprocessing of flexible endoscopes and peripheral endoscopic instruments - recommendations by the Commission for Hospital Hygiene and Infection Prevention at the Robert Koch Institute (RKI), Internet: <http://www.rki.de>
- Hygiene requirements involved in the reprocessing of medical devices – recommendations by the Commission for Hospital Hygiene and Infection Prevention at the Robert Koch Institute (RKI) and the Federal Institute for Drugs and Medical Devices (BfArM) : "Hygiene Requirements for the Reprocessing of Medical Products"
- As the product is not re-sterilisable, please note: ISO 17664 Sterilisation of Medical Devices Information that is to be provided by the manufacturer for the reprocessing of re-sterilisable medical devices

### Warranty limitations

DIMEDA Instrumente GmbH guarantees to manufacture your products with the utmost care.

THIS IS THE ONLY VALID WARRANTY AND SUPERSEDES ALL OTHER DECLARATIONS OF WARRANTY.

It should be noted that due to the biological differences in individuals to be treated, no product is always completely effective under all conditions.

DIMEDA Instrumente GmbH has no influence on the applications of the product, the diagnosis of the patient and on the handling of the product outside of the company. DIMEDA Instrumente GmbH cannot guarantee effectiveness or even complication-free use of the product. Therefore, DIMEDA Instrumente GmbH assumes no liability for damages and costs. DIMEDA Instrumente GmbH will replace products that have a defect that has been accepted by DIMEDA Instrumente GmbH.

# Operating Instructions

## Modular Hand Instruments for Laparoscopy

Employees of DIMEDA Instrumente GmbH are not authorized to modify the above-mentioned conditions, to extend the liability or enter into additional product-related obligations.

Products are subject to change.

### Explanation of Symbols:

Reference number

Lot number


Follow the instructions

Attention – refer to accompanying documents

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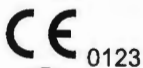Manufacturer


Date of Manufacture

Do not use if the packaging is damaged

Store dry

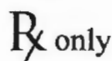To be used only by specialist personnel

Type F: galvanically insulated application part (F stands for floating) that meets the requirements for leakage currents for type B.

Type CF: galvanically insulated application part (F stands for floating), which fulfills even higher requirements for leakage currents than type B. Also suitable for direct use on the heart.

**dimed**<sup>®</sup>  
SURGICAL INSTRUMENTS  
**DIMEDA Instrumente GmbH**  
Gänsäcker 54 - 58  
D-78532 Tuttlingen, Germany  
i.A. *[Signature]*

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
To whom it may concern

STATEMENT

Herewith, «DIMEDA Instrumente GmbH», having its principal business office at Gänsäcker 54+58, 78532 Tuttlingen, Germany, represented by its General Manager, Mr. Bastian Milkau certifies that our «Operating Instructions for Modular Hand Instruments for Laparoscopy» and «Instruction Manual for Reusable Surgical Instruments» covers the following products of Class I: endoscopic forceps, grasping forceps, baskets, forceps with accessories; laryngoscopes and anascopes, proctoscopes with accessories.

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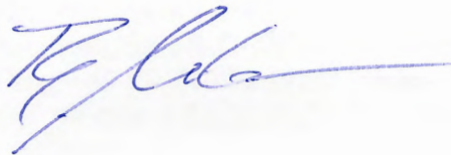[info@dimedada.de](mailto:info@dimedada.de)

## Products

This instruction manual is valid for all reusable surgical instruments, which

- are one-piece instruments
- optionally have simple joints or
- include simple moving parts
- are optionally composed of multiple replaceable parts (e.g. handle and multiple attachments)

## Important Information

*Read this instruction manual carefully before every application and keep it easily accessible for all users or the respective specialist staff.*

*Carefully read the warnings marked with this symbol. Improper use of the products may result in serious injuries to the patient, the users or third parties.*

## 1 Scope

The instruments must be used according to their intended use in the medical fields and by respectively trained and qualified staff only. The treating physician and/or user is responsible for choosing the equipment for specific applications and/or operative use, for the appropriate training and information, and for the sufficient experience regarding the handling of the equipment.

The instruments intended application duration is 60 minutes. The product must be replaced after 60 minutes.

## 2 Contraindication

The instruments must be used according to their intended use in the medical fields and by respectively trained and qualified staff only. The instruments of Dimedada GmbH may not be used in the areas of cardiothoracic surgery and neurosurgery.

## 3 Precautions and Warnings

### Attention!

The instruments are only designed for surgical use.   
 Регистрация МИ в Росздравнадзоре  
[www.nevacert.ru](http://www.nevacert.ru) | [info@nevacert.ru](mailto:info@nevacert.ru)  
 other purpose. Improper handling and

maintenance as well as misappropriated use may result in the premature wear and tear of the instruments.

### Material intolerance

Under no circumstances must the instruments be used if the user or specialist staff become aware of the patient being intolerant to the material.

### Functional Impairment

Surgical instruments corrode and become impaired in their functionality if they come into contact with aggressive substances. It is therefore necessary to observe the storage and sterilisation instructions.

### Operating Conditions

The aforementioned products require correct maintenance and care in order to guarantee that the products operate safely. In addition to this, functionality testing and a visual check should be performed prior to each application. For this reason, please pay attention to the respective chapters in this instruction manual.

### Combination with other products

Should the products be reassembled after disassembly, individual parts must not be replaced with parts from other manufacturers! If the intended purpose of the product entails certain parts being exchanged (e.g. different attachments), no parts from different manufacturers must be used! We recommend to also purchase other accessories (e.g. detergents) at Dimedada Instrumente GmbH.

### Storage

There are no specific storage requirements concerning the products. Nevertheless, we recommend storing medical products in a clean and dry environment.

### Creutzfeldt Jakob Disease

With regard to the reprocessing of medical devices that have been used on patients or suspected patients suffering from or suspected of suffering from Creutzfeldt-Jacob disease (CJD) or its variant (vCJD), the requirements specified in the corresponding appendix of the guidelines for hospital hygiene and infection prevention and the requirements specified by publications in the Federal Health Gazette must be adhered to. The medical devices that were used on this group of patients must be disposed of by incineration (European Waste Catalogue EAK 18 01 03) without risk. Dry heat, ethanol, formaldehyde and glutaraldehyde have a fixing but no inactivating effect on TSE pathogens. Of the sterilization methods available, only steam sterilization (especially 134°C, 18 minutes) has been shown to have a limited effect.

## 4 Liability and Warranty

As a manufacturer, Dimedada Instrumente GmbH is not liable for consequential damage resulting from improper use or handling. This

particularly applies to use which is not compliant with the defined intended use, or non-compliance with the instructions on preparation and sterilisation. This also applies to repairs or changes to the product which are not carried out by authorized staff of the manufacturer. These disclaimers also apply to warranty services.

## 5 Sterility

### State upon Delivery

Medical products are delivered in a non-sterile condition and need to be prepared and sterilised by the user prior to the first application and any subsequent application according to the following instructions.

## 6 Reprocessing

### Warnings

- Frequent reprocessing impairs the quality of the products.
- Any tap water to be used must comply with COUNCIL DIRECTIVE 98/83/EC of 3 November 1998 on the quality of water intended for human consumption.
- The detergents and disinfectants used in order to maintain validation are listed in this instruction for reprocessing. Should an alternative detergent or disinfectant (listed at the RKI or VAH) be used, the responsibility lies with the reprocessor.
- Reassemble disassembled products prior to sterilisation.
- The preparation may only be carried out by qualified medical personnel. Machine reprocessing must be qualified and validated by the user. The washer-disinfectors must fully comply with the requirements of DIN EN ISO 15883.

### Use Site

The first steps of a proper reprocessing take place in the operating theatre. Coarse contaminations must be removed prior to storing the instruments if possible. For this purpose, the instruments should be rinsed under cold tap water (<40°C). If this procedure is not sufficient to remove the obvious soiling, a soft plastic brush can be used to remove soiling. Whenever possible, dry removal (moistened, closed system) should be the method of choice. A drying of any residues should be avoided! Wherever possible, dry disposal is to be preferred, since with wet disposal the prolonged lying of the medical devices in solutions can lead to material damage (e.g. corrosion). Long periods of waiting until the reprocessing, for instance overnight or over the weekend, must be avoided with both types of removal (<60 minutes).

### Transport

The products must be disposed of in a dry state immediately (<60 min) after use, if possible. This means that the products have to be transported in a closed container from

# Instruction manual

For reusable surgical instruments

**dimedica**  
SURGICAL INSTRUMENTS

Valid from:

10.12.2020

Version:

5

the place of application to the purification, so that the products do not dry up.

## Preparing the Decontamination

The products must be disassembled prior to the following reprocessing steps and/or must be exposed to the following reprocessing steps in an open condition, where possible. Rinse residue must be avoided. The products must be reprocessed in appropriate screen baskets or rinsing shields (choose size according to product). The products must be positioned in the cleaning basket at a minimum clearance from one another. Avoid overlapping so that the damaging of the products during the cleaning process can be excluded.

## Pre-cleaning

Rinse the products under cold tap water suitable for use as drinking water (<40°C) until all visible contaminations are removed. Remove persistent dirt with a soft brush. Movable parts on the instrument need to be moved. Cavities, lumen, gaps and slits need to be rinsed intensively (>60 sec) with cold tap water suitable for use as drinking water (<40°C) using a water pistol (or similar). Insertion of the products into an ultrasonic bath (<40°C) with an alkaline cleaner (0.5% Neodisher® MediClean forte, Dr. Weigert), time of treatment 5 min. at a frequency of about 35 kHz. For this please observe the instructions of the cleaning agent manufacturer. Rinse instruments briefly (<15 sec) under cold water. Movable parts must be moved. Respectively rinse cavities, lumen, gaps and slits again (>30 sec.) with cold city water (<40°C) using a water pistol (or similar).

## Cleaning/disinfection

### Automated cleaning and/or disinfection process

(Miele Disinfecter G7835 CD as per ISO 15883):

- Pre-clean for 1 minute
- Drain water
- Pre-clean for 4 minutes
- Drain water
- Clean for 6 minutes at 58°C +/- 1°C using 0.5 % alkaline detergent (0,5 % Neodisher Mediclean forte)
- Drain water
- Neutralize for 3 minutes (0.1% 0,1 % NeodisherZ) with cold tap water suitable for use as drinking water <40°C
- Drain water
- Clean for 2 minutes with FD water <40°C

Please observe the special instructions of the manufacturer of the automatic cleaning machine.

### Automated Disinfection

Automated thermal disinfection in a cleaning and disinfection device taking into consideration the A<sub>0</sub> value, for instance, A<sub>0</sub> value 3000, >5 minutes at >95°C..

### Automated Drying

Automated drying in accordance with the drying operation of the cleaning and disinfection device for at least 30 minutes at 92°C +/- 2°C.

### Sterilisation

(Typ B Autoklav by Tuttnauer as per DIN EN 13060)

Sterilisation of products with a fractionated pre-vacuum method (according to DIN EN ISO 17665-1) taking into consideration the respective national requirements. The sterilisation of the products must be conducted in suitable sterilisation packaging according to DIN EN ISO 11607-1 and EN 868

The sterilisation must be completed using a fractionated pre-vacuum method with the following parameters:

134°C,

5 minutes hold time,

3 pre-vacuum cycles

Drying in vacuum for at least 20 minutes

The instruction manual of the manufacturer of the autoclave and the recommended directions for maximum loading with goods to be sterilised must be observed. The autoclave must be installed, maintained, validated and calibrated in accordance with requirements.

### Additional Information

The reprocessor is responsible for ensuring that the actual reprocessing, including the used equipment, materials and the staff involved in the reprocessing facility, achieves the desired results. This typically requires the validation and routine monitoring of the method and the equipment used.

## 7 Lifespan of the Products

Frequent reprocessing has little effect on the surgical instruments. Product life is limited by wear due to proper use and damage to the instrument. The product may no longer be used under the following conditions, among others:

Corrosion, damage, fractures, cracks, deformation, porosity, functional limitations, products with unrecognizable or missing marking. For this reason, the conditioner must comply with the relevant instructions for functional testing (section 8).

## 8 Functional testing

Check products after reprocessing and before sterilization with regard to the following aspects:

- Cleanliness
- Damage, including but not limited to signs of corrosion (rust, pitting), discoloration, deep scratches, flaking, cracking and wear.
- Proper function
- Missing or removed (ground off) part numbers.
- Do not use improperly functioning or defective and excessively worn products, products with

unrecognizable markings, m or removed (ground) part numbers.

- Check products for perfect surfaces, correct assembly a functionality. Do not use severely damaged products.
- Reassemble dismantled proc before sterilization.

## 9 Service and Repair

### Service and Repair

Do not carry out any repairs or changes the product yourself. Authorized staff of manufacturer are solely responsible for work. Should you wish to make complaints, or offer us any advice regarding products, please feel free to contact us.

### Returns

Defective or non-compliant products must be returned through the entire reprocessing process before being sent back for repairs/service.

## 10 Packaging, Storage and Disposal

Standard packaging of the products for sterilisation according to ISO 11607 and EN 868.

Store sterile products in a dry, clean, and dust-free environment, secured against damage, at moderate temperatures.

The medical products of the manufacturer should be stored and kept in single packaging, boxes or protective containers. Please handle the instruments with care during transportation, storage and reprocessing. The user and/or specialist intended for this is responsible for ensuring that the sterile state is maintained after sterilisation.

The disposal of the products, packaging well as the accessories must be performed in accordance with current rules and laws. specific instruction regarding this matter provided by the manufacturer.

## 11 Description of Symbols

|  |                                      |
|--|--------------------------------------|
|  | Attention!                           |
|  | Observe the instruction manual       |
|  | Item number                          |
|  | Lot designation                      |
|  | CE labelling                         |
|  | Indication of a non-sterile product  |
|  | Name and address of the manufacturer |

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