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PT00117507

Part No. PT00117507

COVIDIEN, COVIDIEN with logo, and Covidien logo and Positive Results for Life are U.S. and internationally registered trademarks of Covidien AG. TM* brands are trademarks of their respective owners. Other brands are trademarks of a Covidien company.

©2020 Covidien. All rights reserved.

Made in Mexico.

May be covered by U.S. Patents: www.covidien.com/patents

 Manufactured for Covidien llc, 15 Hampshire Street, Mansfield, MA 02048 USA.

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
REF SEA3700 **ULPA Filter**

Intended Use

A disposable ULPA Filter for the cleansing of surgical smoke using a RapidVac Smoke Evacuator.

Installing the Filter

1. Remove the SEA3700 ULPA Filter from the shipping box and protective wrapping.
2. Examine the filter for damage incurred during shipping or storage. If the filter is damaged, replace with a new filter before proceeding.
3. Orient the filter with the locking clip on the top edge.
4. Insert the filter into the opening on the right side of the smoke evacuator. Slide the filter into place until it is fully seated and the locking mechanism clicks into place.

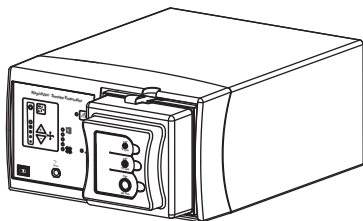Notice

If the filter is improperly installed, the smoke evacuator will not operate and the lowest filter life indicator is amber and flashing.

Caution

The filter must be kept dry and free from fluids. If needed, install a fluid-collection device to prevent the introduction of fluids into the filter.

Failure to install a fluid-collection device may cause filter blockage and electrical damage to the smoke evacuator.

Replacing the Filter

1. Turn off power to the smoke evacuator by pressing the power switch (O) on the front panel.
2. Disconnect the tubing from any reusable instrument that will not be discarded with the tubing.
3. Disconnect the tubing from the smoke evacuator filter.
4. Release the locking tab on the top of the filter by depressing it until the filter is free.
5. Slide the filter out of the smoke evacuator.

Caution

The filter and tubing capture potentially hazardous particles. When removed, handle as you would any biohazardous material. Dispose of these items with other operative waste materials according to the procedures for your institution.

6. Install the new filter as described in the previous section of these instructions.

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REF SEA3720 **Smoke Lap Tubing Kit with Valve**
1/4" x 10' (0.6 cm x 3 m)

Intended Use

The tubing with valve attaches to a trocar or cannula for the evacuation of surgical smoke in a laparoscopic procedure.

Do not use if package is opened or damaged.

Warning

This product cannot be adequately cleaned and/or sterilized by the user in order to facilitate safe reuse, and is therefore intended for single use. Attempts to clean or sterilize these devices may result in a bio-incompatibility, infection or product failure risks to the patient.

Set Up

1. Using the proper aseptic techniques, remove the sterile tubing from the packaging.
2. Connect the Luer connector on the tubing to a cannula or the stopcock on a trocar.
3. Connect the free end of the tubing to the 1/4" (6 mm) RapidVac Smoke Evacuator filter inlet port.

Notice

If fluid will be evacuated, a fluid suction container must be attached.

4. Position and attach the drape clip to hold the tubing in place.

During a Procedure

1. Select the desired level of suction on the RapidVac Smoke Evacuator by pressing the up-arrow button (+) to initiate and increase suction, or the down-arrow button (-) to reduce suction.

Caution

Use the lowest possible setting for the removal of smoke from the surgical site. The Smoke Lap Tubing Kit should be used at approximately 20% - 40% suction level (20%=1 illuminated LED, 40%=2 LEDs). Do not exceed 40% during a surgical procedure.

Extended smoke removal by the stopcock on the trocar may cause reduction of the pneumoperitoneum.

2. Depress the valve on the tubing as needed to maintain a clear field.

After a procedure

1. Turn off power to the smoke evacuator by pressing the power switch (O) on the front panel.
2. Disconnect the tubing from any reusable instrument that will not be discarded with the tubing.
3. Disconnect the tubing from the smoke evacuator filter.

Caution

The tubing assists with the capture of potentially hazardous particles. When the tubing is removed, handle as you would any biohazardous material. Dispose of tubing with other operative waste materials according to the procedures for your institution.

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REF SEA3730 Interlink Cable

Intended Use

The Interlink Cable attaches to a RapidVac Smoke Evacuator and Covidien generator to synchronize smoke evacuation with the activation of the generator. The cable is compatible with the Covidien generator models Force FX™, Force EZ™, and ForceTriad™.

Installing the Interlink Cable


1. Make sure the power to generator and smoke evacuator are off and the units are unplugged.
2. Remove the metal cover plate in the center of the generator back panel to access optional ports.
3. **Force FX only:** A replacement cover is provided for the Force FX to allow access to the RF Activation port. Cut out the ECG section on the replacement cover and secure it to the replacement port using the screws removed from the original cover plate.

Notice


Replacement covers are not provided for the Force EZ or the ForceTriad. Covers are not required for operation.

4. Plug the interlink cable into the generator RF Activation port.
5. Plug the interlink cable into the Remote Activator jack on the back panel of the smoke evacuator.


Force FX

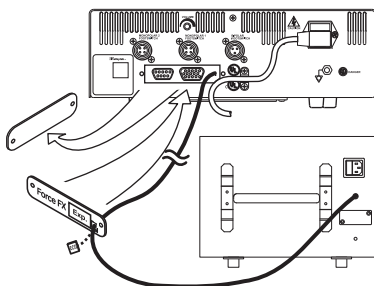Force EZ

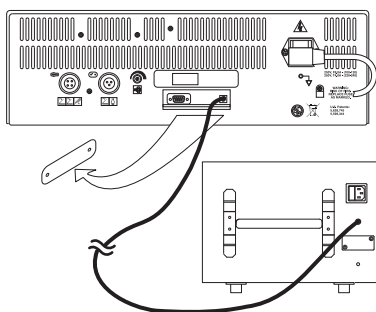ForceTriad

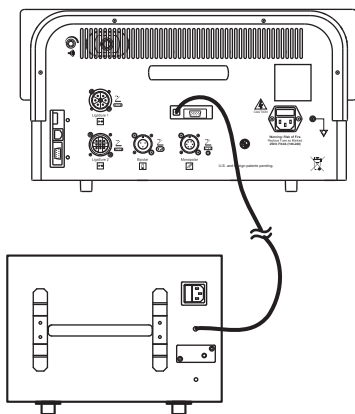Setting the Pre-Set Suction Level

The level of suction to be used when the smoke evacuator is remotely activated is selected using the RapidVac+ button.

1. Plug the generator and smoke evacuator into outlets of the proper voltage.
2. Attach hoses, surgical pencil, patient return electrode, and any other accessories to be used.
3. Turn power on to both devices by setting their respective Mains power switches to on (I).
4. When the devices have completed their self tests, set the level of suction that will be used when the smoke evacuator is activated:
 - Press the up-arrow (+) until reaching the desired suction pre-setting
 - Press the RapidVac+ button 

The pre-set suction level is now selected. The smoke evacuator continues to run and suction reduces to a low flow. Use of the surgical pencil automatically activates the RapidVac to the pre-set suction level. Deactivating the pencil returns suction to a low flow.

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ME20

Part No. 1016893

COVIDIEN, COVIDIEN with logo and the Covidien logo are U.S. and/or internationally registered trademarks of Covidien AG.

Other brands are trademarks of a Covidien company.

© 2011 Covidien.

 Manufactured for Covidien llc,
15 Hampshire Street, Mansfield, MA 02048 USA.
EC REP Covidien Ireland Limited,
IDA Business & Technology Park, Tullamore.

Made in USA. Printed in USA.

www.covidien.com

REV 02/2011

RapidVac™

Remote Switch Activator

- en** Instructions for Use
- fr** Mode d'emploi
- pt** Instruções de uso
- de** Gebrauchsanleitung
- es** Instrucciones de uso
- it** Istruzioni per l'uso
- nl** Gebruiksaanwijzing
- sv** Bruksanvisning
- ru** Инструкция по применению
- zh** 使用说明

RapidVac™

for use with the RapidVac Smoke Evacuator

REF SEA3741 **Remote Switch Activator**

Intended Use


The RapidVac Remote Switch Activator synchronizes activation of an electrosurgical generator in cut or coag mode with smoke evacuation.

Installing the Remote Switch Activator

The Remote Switch Activator (RSA) is programmed to recognize cut and coag activation currents emitted from the generator. The recognized current initiates activation of a smoke evacuator.

The following instructions explain the placement and set-up of the RSA for synchronized generator operation with smoke evacuation.

1. Select a space on the front of the generator or generator cart near the 3 prong ESU pencil cable input where the Remote Switch Activator will be mounted.
2. Clean the surface of the selected area with alcohol. Allow the alcohol to evaporate.
3. Remove the backing from the gasket on the RSA.
4. Affix the RSA to the generator or generator cart.
5. Plug the RSA connector into the Remote Activator jack (1) on the back panel of the RapidVac Smoke Evacuator.

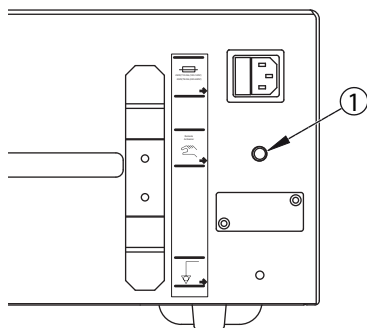Set-up


1. Plug the generator and smoke evacuator into outlets of the proper voltage.
2. Attach hoses, surgical pencil, patient return electrode, and any other accessories to be used.

Warning

When activating an electrosurgical pencil for set up, the operator must be properly grounded to avoid unintended surgical affect. A patient return electrode connected to the generator and properly adhered to the operator (i.e. on the arm opposite the hand holding the pencil) provides a safe path for electrical current from the pencil back to the generator.

3. Turn power on to both devices by turning their respective Mains power switches to on (I).
4. When the devices have completed their self tests, set the smoke evacuator suction level to 40%.
 - Press the up-arrow (+) twice. Two air flow indicator are illuminated.
 - Press the RapidVac+ button 

The suction reduces to a minimal level and the pre-set suction is set at 40%.
5. Drape the ESU wire into the RSA as shown (2). More than one ESU wire may be placed in this cradle.

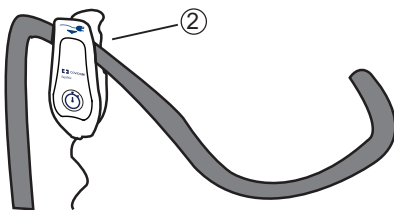Test the connectivity of the RSA by depressing the cut or coag button on the electrosurgical pencil. The smoke evacuator should start. Deactivate the pencil.

6. An optional delayed stop can be set for the smoke evacuator once deactivation of the pencil occurs. Press and hold the button located on the front of the RSA for the desired amount of delay. 

This delay can be set from 0-10 seconds. To return to delay of 0 seconds, tap the button once.

Notice

If both the cut and coag tones are not activating the smoke evacuator, disconnect the RSA cord from the smoke evacuator. Reinsert the connector and repeat step 5.

The RSA is now programmed to activate the smoke evacuator when the cut and coag tones are emitted from a generator.

Notice

If the attached smoke evacuator is to be separated from the generator, unplug the RSA from the smoke evacuator. The RSA should remain affixed to the generator to which it was programmed. The RapidVac Smoke Evacuator cannot be operated manually or with a footswitch while the RSA is plugged into the smoke evacuator.

If the RSA is moved to another generator, it must be reprogrammed using the steps in this section.

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
sv Aktiv

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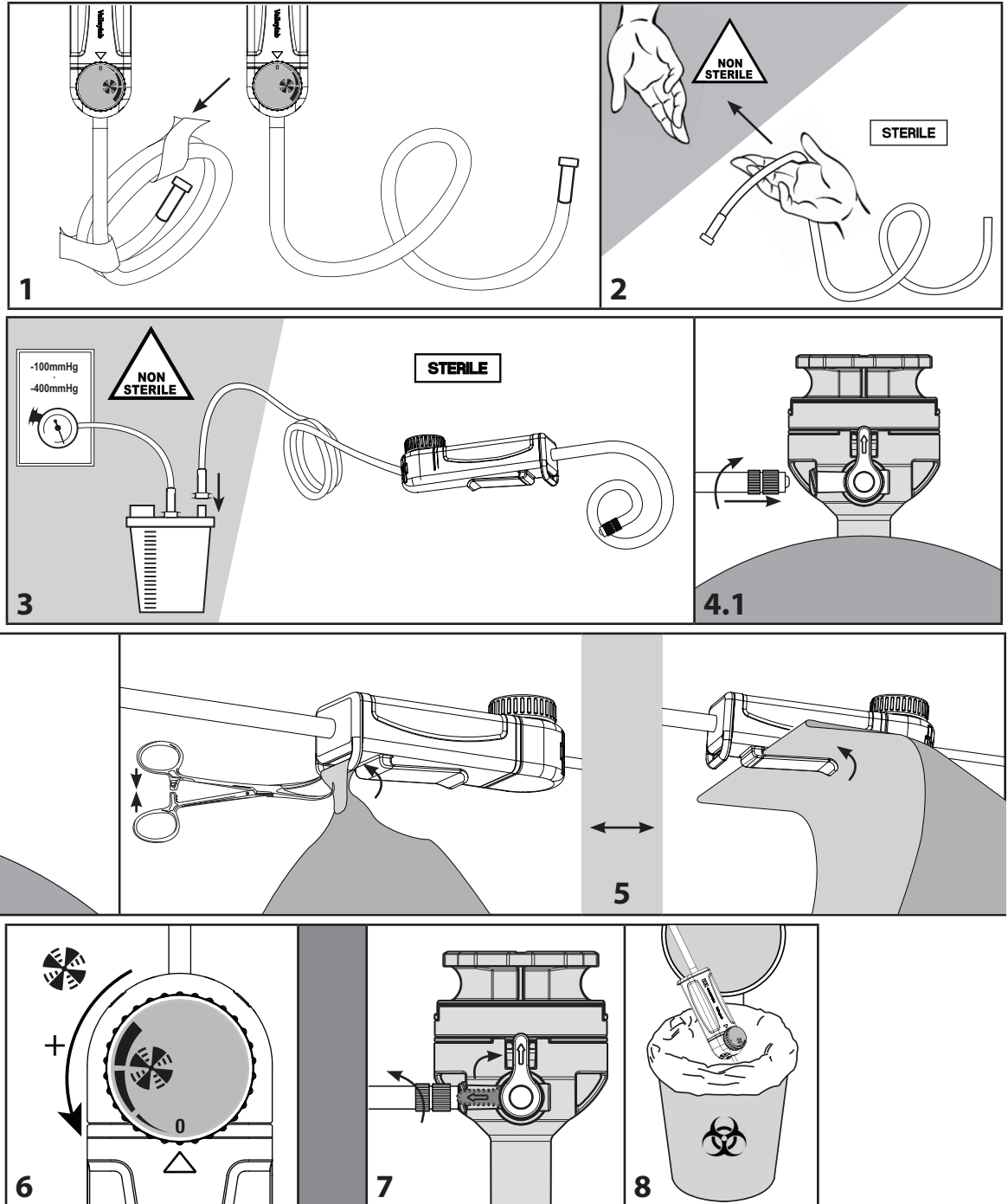Part No. PT00095524 904855 Rev B

COVIDIEN, COVIDIEN with logo, and Covidien logo and Positive Results for Life are U.S. and internationally registered trademarks of Covidien AG. Other brands are trademarks of a Covidien company, ™* brands are trademarks of their respective owner.

May be covered by U.S. Patents: www.covidien.com/patents.

©2017 Covidien

Manufactured for Covidien llc, 15 Hampshire Street, Mansfield MA 02048 USA.

EC REP Covidien Ireland Limited, IDA Business & Technology Park, Tullamore, Ireland.

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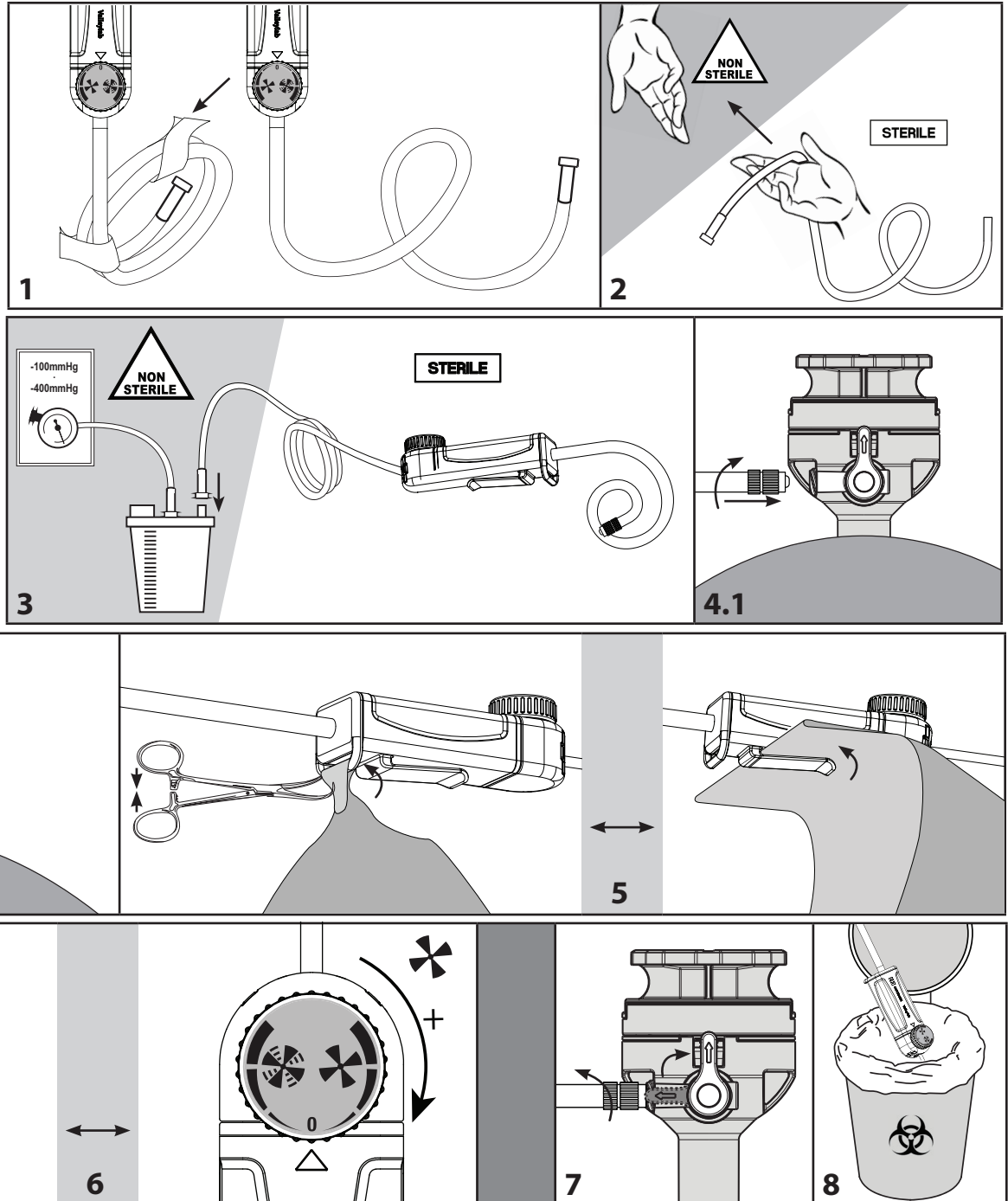
Part No. PT00100933 (OEM part no. 904856 Rev A)

COVIDIEN, COVIDIEN with logo, and Covidien logo and Positive Results for Life are U.S. and internationally registered trademarks of Covidien AG. Other brands are trademarks of a Covidien company, ™* brands are trademarks of their respective owner.

May be covered by U.S. Patents: www.covidien.com/patents.
 ©2017 Covidien

Manufactured for Covidien Inc, 15 Hampshire Street, Mansfield MA 02048 USA.
EC REP Covidien Ireland Limited, IDA Business & Technology Park, Tullamore, Ireland.
 Made in Mexico. www.covidien.com Rev 05/2019

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COVIDIEN, COVIDIEN with logo and Covidien logo are U.S. and internationally registered trademarks of Covidien AG.