Instruction for use

Medical gloves: latex surgical gloves powder free sterile First+Care (sizes: 6, 6.5, 7, 7.5, 8, 8.5, 9)

Manufacturer: Shandong Yuyuan Latex Gloves Co., Ltd. (Industrial Zone, Linqing City, Shandong Province, China)

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1. Scope.

1. 1 Purpose.

Medical gloves: latex surgical gloves powder free sterile First+Care (sizes: 6, 6.5, 7, 7.5, 8, 8.5, 9) (hereinafter Sterile Medical gloves) – sterile disposable products for general use in medical and medical and preventive medical institutions are designed to protect a health-care worker from exposure to harmful factors of the external production environment, to reduce the risk of infection of patients with microbes that are part of the resident flora of the hands of medical workers.

1.2 Use.

Sterile surgical gloves are intended for use in medical and medical-preventive medical institutions as protective equipment for surgical operations, diagnostic, medical procedures or for medical manipulations requiring compliance with sterility conditions.

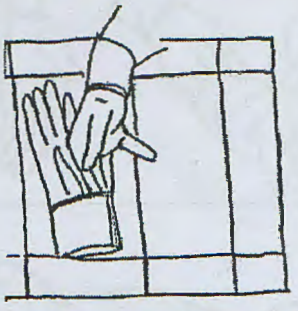
1.3 Indications for use.

- for any surgical intervention: for large and small surgical operations, for small surgical operations on the skin and other manipulations associated with penetration into the tissues under the skin, in the mucous membranes, for endoscopic manipulations, for obstetrics, and for bandaging;
- to protect a health-care worker and patients during manipulations related to access to sterile areas of the body (vessels, genitourinary tract, tracheobronchial tree), rectal, vaginal examination, contact with mucous membranes and etc.;
- for the preparation of drugs for total parenteral administration;
- gloves are used when covering a sterile table and laying on a sterile sheet of instruments, removing moisture residues, placing them in a sterile tool bag, endoscopes after sterilization;
- to protect a health-care worker and a patient when there is a risk of infection with direct contact with the patient (his/her biological matrixes) or contact with mucous membranes, damaged skin, with blood or other biological substrates that are potentially or clearly contaminated by microorganisms;

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2. Rules for use.

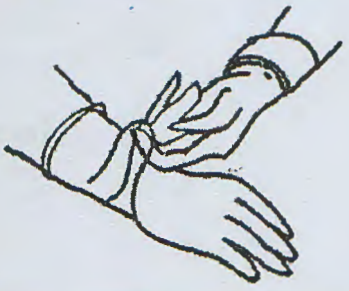
2.1 How to don the gloves

1. Take the fingers of the right hand in the glove for the lapel on the left glove, touching it from the outside, make a lapel; 2. Take the glove for the right hand for the lapel so that the fingers do not touch the inner surface of the glove;	
3. Shorten the fingers of your right hand and insert them into the glove;	
4. Put on the glove without breaking its lapel, breaking the fingers of the right hand;	
5. Take the glove for the left hand with the second, third and fourth fingers of the right hand (in the glove) for the lapel so that the fingers do not touch the inner surface of the glove;	

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6. To shorten the fingers of the left hand and insert them into the glove; Straighten the lapels on the left, then the right gloves, pulling them onto the sleeve.	

2.2 How to remove the gloves

process	justification
1. Take the fingers of the right hand in the glove for the lapel on the left glove, touching it from the outside, make a lapel;	
2. Take the fingers of the left hand in the glove for the lapel on the right glove, touching it from the outside, make a lapel;	
3. Remove the glove from the left hand, turning it inside out and holding the lapel in the right hand;	
4. With your left hand, take the glove on your right hand for the lapel from the inside and remove it by turning it inside out (the left glove was inside the right).	

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2.3 The rules for use:

- The use of medical sterile gloves does not replace hand hygiene. It is necessary to carry out hand treatment in accordance with the manipulation.
- It is necessary to select sterile medical gloves of the appropriate size.
- It is necessary to select sterile medical gloves depending on the nature of the manipulations performed - they should be easily don and not impede movement or squeeze the hands.
- When donning medical sterile gloves, no great effort should be exerted.
- When donning sterile medical gloves, the established technique should be used, that will ensure the sterility.
- Before donning the sterile medical gloves for performing medical procedures, it is necessary to prepare a workplace - ready-to-use sterile instruments and consumables should be on the manipulation table. It is necessary to ensure their availability without compromising the sterility of the hands.
- Do not use medical sterile gloves in damaged packaging. Sterility may be impaired.

2.3.1 To replace the gloves in case of:

- a visible defect;
- possible or actual puncture.
- if the glove begins to swell, stretch, stick or not seated on the hand as a result of the absorption of liquids by the material.
- in accordance with the rules adopted in the institution.
- during the transition from performing "dirty manipulations" to "clean ones", and also when changing from one patient to another one.
- every two hours of continuous operation;

3 Product description and its main specifications.

Medical gloves: latex surgical gloves powder free sterile First+Care (sizes: 6, 6.5, 7, 7.5, 8, 8.5, 9).

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Gloves' material – 100% latex
Surface: textured, ensure reliable gripping of both dry and wet tool
Form: anatomical, not ambidextrous.
Cuff options: with roller on the edge of the cuff.
The cuff fits tight.
Control Level and AQL (Acceptable Quality Limit) 1,5

Physical sizes:

Size	Length, mm	Width, mm	Cuff's thickness mm.	Fingers' thickness mm.	Palm's thickness mm
6	≥265	77±5	0.11±0.01	0.17±0.01	0.15±0.01
6,5	≥265	83±5	0.11±0.01	0.17±0.01	0.15±0.01
7	≥270	89±5	0.11±0.01	0.17±0.01	0.15±0.01
7,5	≥270	95±5	0.11±0.01	0.17±0.01	0.15±0.01
8	≥270	103±5	0.11±0.01	0.17±0.01	0.15±0.01
8,5	≥280	109±5	0.11±0.01	0.17±0.01	0.15±0.01
9	≥280	114±5	0.11±0.01	0.17±0.01	0.15±0.01

Weight of the glove

Size	Weight
6	7.0 g±3
6,5	7.5 g±3
7	8.0 g±3
7,5	8.6 g±3
8	9.5 g±3
8,5	9.9 g±3
9	10.5 g±3

Physical properties:

Force at break (N)		Elongation at break %	
Before Aging	After Aging	Before Aging	After Aging
12.5	9.5	≥700	≥550

4. Labeling.

Single package.
Sterile gloves packaging
The following data are indicated on a single package of sterile gloves:
1) name and trade mark, manufacturer

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Group package.

The group package contains a certain number of single packages of gloves of the same size, ensuring safe transportation and storage. The group packages are marked with the number of pairs of gloves and additional instructions for storage

5. Storage and transportation conditions.

5.1 Gloves are transported by all modes of transport in accordance with the transportation rules applicable for this type of transport in original factory packages - cardboard boxes in a protected from direct sunlight at a temperature of 5 ° C to 30 ° C and 65% relative humidity.

5.2 Gloves are stored in original factory group cardboard packages in a well-ventilated area away from heat sources, sunlight, direct exposure to light sources and electrical equipment, with an optimum humidity of 65%. At a storage temperature of 5 ° C to 30 ° C. - otherwise there may be difficulties due to loss of elasticity by gloves.

6. Package

The gloves are packed in packages and transport boxes.
Single package - 1 pair
Quantity in the package - 50 pairs
In the box - 10 packs * 50 pairs

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In the box there are 500 pairs

7. Shelf life

The shelf life of the product (indicated on the package and the box) is 5 years from the manufacture date indicated on the packaging.

8. Warranty

The manufacturer guarantees the compliance with all specifications and characteristics, during the entire shelf life – that is 5 years - subject to the integrity of the packaging and the compliance with the storage, transportation and use conditions in accordance with the instructions for use.

9. Special instructions

Restrictions on age and user status are determined by the nature of the work performed. The gloves should be used as indicated, in accordance with a certain size. Do not use after the expiration date. In exceptional cases in people prone to allergy to latex, an allergic reaction is possible. The employee must be familiar with this instruction and have a qualification appropriate to the nature of the work performed. The gloves are disposable, methods of care are not provided for.

10. Procedure for the recycling and disposing of gloves.

Sterile surgical gloves are disposable medical products. After use, the gloves are to be collected, temporarily stored, removed and disposed in accordance with the medical waste management rules adopted by the health care organization.

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