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1. Frantzen F, Studies on glycated haemoglobin and new analytical methods for its determination. Dr. Philos. Thesis. NTNU, Trondheim 1998. ISBN 82-7861-116-5.
2. The Diabetes Control and Complications trial Research Group, The Effect of Intensive Treatment of Diabetes on the Development and Progression of Long-Term Complications in Insulin-Dependent Diabetes Mellitus. N Engl J Med 1993; 329:977-986.
3. Sacks DB et al., Guidelines and Recommendations for Laboratory Analysis in the Diagnosis and Management of Diabetes Mellitus. Clin Chem 2002; 48:436-472.
4. Hoelzel W et al., IFCC reference system for measurement of haemoglobin A1c in human blood and the national standardization schemes in the United States, Japan and Sweden: a method-comparison study. Clin Chem 2004; 50:166-174.
5. Nathan DM et al., A1c-Derived Average Glucose Study Group. Translating the A1c Assay into estimated average glucose values. Diabetes Care 2008; 31:1473-1478.
6. http://www.ngsp.org/docs/methods.pdf
7. http://www.ngsp.org/protocol.asp
8. Standards of Medical Care in Diabetes. Recommendations of the American Diabetes Association. Diabetes Care Vol. 34. Jan 2011; 1:S62-S69.
9. Parafilm® M is a trademark of Pechiney Plastics Packaging.

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EN - Instructions For Use

Test Cartridge Storage

Refrigerated 2 to 8°C (36 to 46°F)
Quo-Lab A1C Test Cartridges are stable until the expiry date only when stored refrigerated. The expiry date is the last day of the month stated on the cartridge, cartridge tray and kit carton.

Room Temperature Storage 15 to 25°C (59 to 77°F)
Quo-Lab A1C Test Cartridges can be stored at room temperature for up to 90 days, if still within the expiry date stated on the cartridge, cartridge tray and kit carton. Write the date of removal from the refrigerator on the cartridge tray or kit carton. Avoid exposure to direct sunlight.

Test Cartridge Preparation and Integrity Check

Allow at least 50 minutes for unopened cartridges to equilibrate to the analyzer's environmental temperature before use, even if previously stored at room temperature. **Do not** use heaters to warm the cartridges.

Check that the foil seal covering the top of the cartridge is not damaged. If the foil seal is damaged then the cartridge **must not** be used.

If problems are experienced, contact your local distributor or contact customer support on +49 39203 511 414 or support@ekf-diagnostic.de

Intended Use

The Quo-Lab A1C Test is intended for the *in vitro* quantitative determination of glycated hemoglobin (HbA1c) in whole blood obtained from finger prick or from venous blood samples collected into EDTA tubes.

The Quo-Lab A1C Test is indicated in the management and treatment of diabetes mellitus and for monitoring long term glycemic control in patients diagnosed with diabetes.

The Quo-Lab A1C Test is designed for professional use only.

Summary and Explanation of the Quo-Lab A1C Test

Human red blood cells (erythrocytes) are freely permeable to glucose present in the surrounding liquid (plasma) of the blood.

During the lifetime of the erythrocytes (normally up to 120 days) exposure to blood glucose results in binding of glucose to the hemoglobin A molecule present in erythrocytes. This is referred to as glycated hemoglobin[1].

Chronic elevated blood glucose will result in tissue and organ damage[2].

Good control of blood glucose, evident as lower HbA1c values, has been proven to result in delayed onset and slower progression of these complications[2,3].

Principle of the Test

The Quo-Lab A1C Test Cartridge contains all of the reagents required for the determination of glycated hemoglobin in whole blood.

The Quo-Lab A1C Test combines the chemical binding of boronate to glycated hemoglobin with the fluorescent quenching effect that this binding exerts on a fluorescent marker bound to the boronate molecule.

After removing the foil seal, the cartridge is placed in the Quo-Lab Analyzer.

When instructed, the user pushes a mixing bead containing the fluorescent boronate conjugate into the buffer contained within the cartridge, using the blunt end of the Sample Stick, and the analyzer automatically rehydrates it.

The user then places the blood sample, contained on the Sample Stick, gently into the top of the cartridge, breaks off and removes the handle and closes the door of the analyzer to introduce the blood sample it into the buffer.

The analyzer automatically mixes the blood sample to release the hemoglobin from the erythrocytes.


The total hemoglobin concentration is determined from the initial decrease in the fluorescent signal.

The fluorescent boronate conjugate binds to the glycated hemoglobin, which is measured by monitoring a decrease in the fluorescence of the active ingredient.

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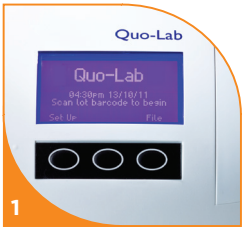The ratio of glycated hemoglobin to total hemoglobin is determined and the result is presented in up to two user selectable units: % DCCT, mmol/mol IFCC, % JDS, eAG mg/dl or eAG mmol/l.

Linear relationships have been established between these reportable units[4,5]:

mmol/mol IFCC = (% DCCT – 2.15) x 10.929
% JDS = (0.09274 x mmol/mol IFCC) + 1.724
eAG mg/dl = (28.7 x % DCCT) - 46.7
eAG mmol/l = (1.59 x % DCCT) - 2.59

eAG values are based on a correlation study linking % DCCT to the patient's average blood glucose concentration, resulting in published formula to derive the eAG[5].

Cartridge Pack Contents

REF 3101-3051 / REF 0055

25/50 Quo-Lab A1C Test Cartridges (1/2 trays of 25 Cartridges).
1 carton containing a minimum of 25/50 Sample Sticks.
1/2 Calibration Barcodes (one on each tray of 25 Cartridges, indicated by the CAL symbol).
1 Instruction For Use leaflet.

Accessories Required (not included)

Quo-Lab Analyzer
Quo-Lab A1C Control Kit
Standard Blood Collection Equipment

Active Ingredients

Fluorescent boronate conjugate
Plastic vial containing ammonium chloride buffer, lysis agent and preservative
Sample Stick treated with EDTA and surfactant

Calibration of the Quo-Lab A1C Test

The Quo-Lab A1C System is certified by the NGSP[6].

The Quo-Lab Analyzer and A1C Test Cartridges have been calibrated using samples provided by the ERL via the NGSP network.

Results obtained using the Quo-Lab A1C System are traceable to the IFCC reference method[7].

Limitations of Use

For hemoglobin concentrations outside of the range of 6.5 g/dL to 20.4 g/dL (4.0 to 12.7 mmol/l) or for HbA1c values outside the range of 4.0 to 15.0 % DCCT (20 to 140 mmol/mol IFCC), a result will not be reported and an error message will be displayed. Consult the Quo-Lab Analyzer System User Manual for further details.

eAG values may differ significantly from a patient's glucose concentration if checked at the same time.

As with all HbA1c analyses, the results obtained are dependent upon the patient's clinical history and having red blood cells with a typical lifespan.

Warnings and Precautions

For *in vitro* diagnostic use only.

Do not use Test Cartridges after the expiry date or if the Test Cartridges have not been stored in accordance with the recommendations.

Do not use the Test Cartridge if the foil seal coving the top of the cartridge has been damaged.

The Test Cartridge contains sodium azide as a preservative. The concentration is < 0.1 %, which is below the level considered to be hazardous in normal use. In case of leakage, avoid contact with eyes and skin. In case of contact, wash with plenty of water.

Do not re-use any part of the Test Cartridge.

Use Gloves – Used Test Cartridges, sampling equipment, patient samples and controls are potentially infectious. Used Test Cartridges should be disposed of immediately as clinical waste in accordance with local / national regulations.

Operate the Quo-Lab Analyzer System according to the User Manual.

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Sample Preparation and Collection

Finger Prick Blood Sample

The patient's finger must be warm, dry and clean, including being free of substances such as hand cream.

Use a single use lancet (not provided) to prick the finger.

Allow a droplet of blood to form, about the width of the Sample Stick, before sampling.

Fill the Sample Stick, avoiding air bubbles and excess sample on the outside of the tip of the Sample Stick.
Do not wipe excess blood from the Sample Stick.

Place the Sample Stick into the Test Cartridge held in the analyzer within 1 minute.

Operate the Quo-Lab Analyzer System according to the User Manual.

Venous Blood Sample

Venous samples must have been collected into EDTA tubes.

Venous samples can be stored refrigerated at 2 to 8 °C (36 to 46 °F) for up to 10 days before use and can be used without equilibration to room temperature.

Thoroughly mix the sample by inverting the tube 8 to 10 times before collecting a sample.

Place a droplet of blood, approximately 10 to 20 µl, onto a non-metallic and non-absorbent surface such as Parafilm® M[9].
Do not attempt to sample directly from the tube.

Fill the Sample Stick, avoiding air bubbles and excess sample on the outside of the tip of the Sample Stick.
Do not wipe excess blood from the Sample Stick.

Place the Sample Stick into the Test Cartridge held in the analyzer within 1 minute.

Operate the Quo-Lab Analyzer System according to the User Manual.

Expected Values

Interpret test results with careful consideration to the patient's age, ethnicity, medical history, clinical examinations and other laboratory results, see limitations of use.

Users should refer to national guidance regarding patient reference ranges, an example is as follows[8]:

The ADA recommends an HbA1c range of 4 to 5.7 % DCCT (20 to 39 mmol/mol IFCC) for non-diabetic patients.

Patients with values in the range of 5.7 to 6.0 % DCCT (39 to 42 mmol/mol IFCC) are considered to be at risk of developing diabetes and those in the range of 6.0 to 6.5 % DCCT (42 to 48 mmol/mol IFCC) are considered to be at high risk of developing diabetes.

The ADA recommends a goal of < 7 % DCCT (53 mmol/mol IFCC) for effective management of diabetes and to minimize long term complications.

Quality Control

The Quo-Lab A1C Control Kit must be used with Quo-Lab A1C Test Cartridges to ensure that the system is operating correctly.

Only Quo-Lab A1C Controls should be used with the Quo-Lab A1C Test.

Please read the IFU included with each Quo-Lab A1C Control Kit to ensure that the correct procedure is followed.

Quality control should be run:

- With each new lot of Test Cartridges.
- With each new shipment of Test Cartridges.
- At any time there is concern that the result may be wrong, the Test Cartridges may not have been stored correctly, the user is unfamiliar with the test procedure or there is a concern that a user may be performing the test incorrectly.

Use the Quo-Lab A1C Control Kit in accordance with the quality standards set by your local organization or laboratory.

If the Quo-Lab A1C controls do not give the results within the range stated in the control kit IFU, contact your local distributor or contact customer support on +49 39203 511 414 or support@ekf-diagnostic.de

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Performance Characteristics

Method Comparison

Quo-Lab A1C results were compared with a laboratory HPLC system giving a linear regression of r = 0.989 (y = 0.978x + 0.278, n = 46).

Precision

Blood samples with high (12.8 %) and low (5.3 %) A1C values were tested in duplicate every morning and afternoon for twenty days on a single Quo-Lab Analyzer. The total imprecision for each level was < 3 % CV respectively.

Linearity and Test Range

High and low pooled blood samples were used to create eleven linearity standards, which were tested on the Quo-Lab A1C Test to establish the linear range. The Quo-Lab A1C Test has a linear working range of 4 to 15 % DCCT (20 to 140 mmol/mol IFCC).

If a result is outside of this range then the analyzer will display an error message. Refer to the Quo-Lab Analyzer System User Manual for further details.

Interfering Substances

No affect from bilirubin (20 mg/dl), cholesterol (500 mg/dl), creatinine (30 mg/dl), triglycerides (1,500mg/dl), uric acid (20 mg/dl), acetaminophen (20 mg/dl), ascorbic acid (3 mg/dl), dopamine (13 mg/dl), ibuprofen (40 mg/dl), salicylate (50 mg/dl), tetracycline (4mg/dl), tolazamide (100mg/dl) and tolbutamide (100 mg/dl) was found.

Hemoglobin Variants

The Quo-Lab A1C Test was found to be unaffected by labile glycated hemoglobin and the following hemoglobin variants: Hb AS, Hb AC, Hb AD, Hb AE, Hb AJ, Hb SS, Hb CC, Hb SC, and elevated fetal hemoglobin (up to 6 %).

It is possible that other substances and/or factors not listed above may interfere and cause false results.

Glossary of Terms

EDTA	Ethylenediaminetetraacetic Acid
ADA	American Diabetes Association
DCCT	Diabetes Control & Complications Trial
ERL	European Reference Laboratory
NGSP	National Glycohemoglobin Standardization Program
IFCC	International Federation of Clinical Chemistry
IFU	Instructions For Use
eAG	Estimated Average Glucose
HPLC	High Pressure Liquid Chromatography

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How to Run a Quo-Lab A1C Test

Refer to images 1 to 12 on the inside of the back cover.

Step 1

Ensure that the analyzer is ready to run a test. Refer to the Quo-Lab Analyzer System User Manual for details.

Step 2

Remove a tray of cartridges from the carton, open the lid and place it next to the analyzer. Allow at least 50 minutes for unopened cartridges to equilibrate to the analyzer's environmental temperature before use, even if previously stored at room temperature.

Step 3

The analyzer prompts the user to scan the calibration barcode to start the test. The barcode is printed on the lid of the cartridge tray, indicated by the CAL symbol.

Place the barcode scanner over the calibration barcode and press the barcode scanner button, until you hear a beep. 'Insert New Cartridge' will be displayed on the analyzer screen along with the lot number scanned.

It is critical that the lot number displayed on the screen matches the lot number stated on the cartridges, otherwise incorrect results could be reported.

If the system returns to the home screen before a cartridge is inserted, rescan the barcode to restart the test sequence.

Step 4

Remove a cartridge from the tray. **Do not** handle the lower part of the cartridge containing the liquid.

The cartridge must be free from condensation before use. **Do not** wipe condensation from the cartridge.

Do not use the cartridge if the foil covering the top of the cartridge is damaged.

Step 5

Carefully remove the foil lid of the cartridge, ensuring that the reagent bead remains within the top of the cartridge. Once the foil lid has been removed, the cartridge must be used within 1 minute.

Step 6

Insert the cartridge into the analyzer. The analyzer will automatically detect the cartridge and will the display the message 'Cartridge Inserted. Running…'.

Do not remove the cartridge from the analyzer until the test is complete.

Do not close the door of the analyzer until instructed to do so in Step 11.

Step 7

When prompted to 'Insert Reagent', use the blunt end of the Sample Stick to push the reagent bead into the vial. Use one finger to hold the blue section of the top of the cartridge in place as the Sample Stick is removed.

There is a time limit of 60 seconds to insert the reagent bead once prompted. If the time runs out the test will be aborted. Remove the cartridge to reset the analyzer.

After the reagent bead has been added the analyzer will display the message 'Rehydrating Reagent'. This process takes approximately 50 seconds. Whilst the system is rehydrating the reagent, prepare the blood sample.

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Step 8

For finger prick sample collection, the patient's finger must be warm, dry and clean, including being free of substances such as hand cream. Use a single use lancet (not provided) on the patient's finger to get a droplet of blood about the same width as the Sample Stick.

For venous samples previously collected into EDTA tubes, ensure that the sample is thoroughly mixed and use standard sampling equipment to place a droplet of blood, approximately 10 to 20 µl, onto a non-metallic and non-absorbent surface such as Parafilm® M[9].
Do not attempt to sample directly from the tube.

Step 9

Lightly touch the pointed tip of the Sample Stick to the top surface of the droplet of blood, as per image 9. The Sample Stick will take up the blood. The blood must completely fill the slot.

Note: Avoid taking up air bubbles and excess sample on the outside of the tip of the Sample Stick. **Do not** wipe excess blood from the Sample Stick.

Step 10

When the analyzer displays the message 'Insert Sample and Close Door', gently insert the Sample Stick, containing the blood sample, so that half of the tip is in the top of the cartridge. Bend the Sample Stick handle towards the front of the analyzer to break it off.

Do not push the Sample Stick fully into the cartridge otherwise the analyzer will report an error message. Refer to the Quo-Lab Analyzer System User Manual for further details.

Step 11

Close the analyzer door fully, using a single smooth movement, to complete the introduction of the Sample Stick into the cartridge. There is a time limit of 60 seconds to insert the Sample Stick and close the door once prompted. If the time runs out the test will be aborted. Remove the cartridge to reset the analyzer.

The used cartridge must be disposed of as clinical waste, as per local regulatory requirements. The analyzer will display a progress bar whilst the test is running.

You will be prompted to scan the Operator and / or Patient ID if these features have been activated. Refer to the Quo-Lab Analyzer System User Manual for further details.

Step 12

When the test has finished the result will be displayed on the screen, until the door is opened. Remove the cartridge and dispose of it as clinical waste, as per local regulatory requirements. The analyzer will be ready to run another test.

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