KARL STORZ — ENDOSKOPE

KARL STORZ SCB control NEO System

- OR1™ control NEO 20 097120-1
- SCB control NEO Software 20 090001S431

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1 Welcome

1.1 Important note

Thank you for your expression of confidence in the KARL STORZ brand name. This product is the result of our wealth of experience and attention to detail.

This instruction manual is designed to help you operate and use the KARL **S**TORZ **C**ommunication **B**us (SCB) control NEO system.

The KARL STORZ SCB control NEO system is supplied as a unit comprising complementary hardware and software. Any changes to the hardware and software which are not authorized by KARL STORZ will result in a loss of approval for medical use and also void the warranty.

All steps necessary for operating the equipment are clearly explained. We thus ask that you read this manual carefully and keep it available for ready reference.

NOTE: *Setup and installation of the OR1™ control NEO system (hardware) is described in the separate installation manual (see index 3). Please read this carefully.*

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Content

1 Welcome
1.1 Important note ... IV

2 General information
2.1 Description of the KARL STORZ SCB control NEO system ... 4
2.2 Scope of supply ... 6

3 Safety instructions
3.1 Explanation of warnings and cautions ... 7
3.2 Safety instructions for the OR1™ control NEO ... 9
3.3 Intended use ... 11
3.3.1 KARL STORZ SCB control NEO system ... 11
3.4 User qualification ... 12
3.5 Safety precautions at the site of installation ... 12
3.6 Safety precautions when operating the KARL STORZ SCB control NEO system ... 13

4 SCB control NEO software user interface
4.1 Optimal use of the SCB control NEO software ... 14
4.2 Information area ... 16
4.2.1 Information line ... 16
4.2.1.1 Information ... 16
4.2.1.2 Warnings ... 17
4.2.2 Date/Time ... 19
4.3 Navigation area ... 20
4.3.1 Explanation of buttons in the navigation area ... 20
4.3.2 Main menu ... 21
4.3.3 Workflow bar ... 21
4.4 Application window/content area ... 22
4.4.1 Application header (menu bar) ... 23
4.4.1.1 Presets ... 23
4.4.1.2 Options ... 25
4.4.1.3 Menu “Help” ... 27

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2 General information

2.1 Description of the KARL STORZ SCB control NEO system

The KARL STORZ SCB control NEO system is based on a communication interface between the connected KARL STORZ devices which, in conjunction with the KARL STORZ OR1™ control NEO with graphical user interface (see section 4), allows the parameters of the connected SCB devices to be displayed centrally and controlled remotely, even from the sterile area.

For you as the user, the following major requirements are fulfilled:

- safety and reliability
- simple operation
- hot plug-and-play*
- user and application-specific Presets
- can also be used from the sterile area (touch screen: with sterile protective cover)
- central remote control
- modular structure
- possibility of later extensions
- simplified documentation (only in combination with AIDA compact NEO or KARL STORZ OR1™ FUSION)

***WARNING:** Connection of two active KARL STORZ SCB control NEO systems can result in temporary loss of function.*

***NOTE:** An active KARL STORZ SCB control NEO system consists of at least one of the following activated devices: TELECAM SL camera control unit, TRICAM® SL camera control unit, or TRICAM® SL II camera control unit, IMAGE 1 HUB™ or OR1™ control NEO respectively as a single device or together with other SCB devices.*

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Function of the SCB control NEO software

With the KARL STORZ OR1™ control NEO, parameters of the SCB devices connected to the KARL STORZ SCB control NEO system can be displayed and controlled remotely.

Operation is by means of a graphic user interface, and commands are entered via a maximum of two touch-sensitive LCD monitors* (touch screens). The graphic user interface is designed in such a way that the entries can be with a finger without any aids.

* The SCB control NEO software release 20090001S431 supports the following touch screen models:

- 20090337 KARL STORZ touch screen 19"
- 20090437 KARL STORZ touch screen 19"

For use in the sterile area, the LCD monitor must be fitted with a disposable sterile cover. Outside the sterile area, a keyboard or mouse can optionally be used.

The following basic functions are possible:

- Remote control of the main functions of the SCB devices
- Display of all relevant SCB device parameters
- User and application-specific “Presets” (preprograming) for certain SCB devices. In general, the desired Preset is created in advance and offline, i.e., without the configured SCB devices actually being connected or switched on (exception: ARTHROPUMP® POWER, AUTOCON® II 400 and KLS Martin maxium® ME402).
- The SCB control NEO software allows a Preset to be linked to a predefined type of surgery. The type of surgery determines which widgets are shown in the overview. This overview enables you to check the most important SCB device parameters at a glance and make any changes if necessary.

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2.2 Scope of supply

Unpacking the equipment

Check for missing items and evidence of shipping damage. Immediately file any complaints with the manufacturer or supplier.

If possible, retain the original packing materials for later use; these can come in handy if the control computer has to be transported.*

* Only if a KARL STORZ SCB control NEO system PC set is purchased.

Basic equipment

The basic equipment varies depending on the set purchased.

You will find a description of the available set(s) in the installation manual 96206262 DER on the KARL STORZ SCB control NEO system (see index 3 in this folder).

Spare parts and recommended accessories

ⓘ ***NOTE:*** *Spare parts and recommended accessories can be found in section 6 in this manual.*

For detailed information please contact:
KARL STORZ SE & Co. KG Postfach
230
78603 TUTTLINGEN/GERMANY

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3 Safety instructions

Safety instructions are measures intended to protect the user and patients from the risks which could arise in conjunction with use of the system.

3.1 Explanation of warnings and cautions

Please read this manual carefully and follow its instructions precisely. The words Warning, Caution and Note convey special meanings. Wherever they are used in this manual, the accompanying text should be carefully reviewed to ensure the safe and effective operation of the device. To make the words stand out more clearly, they are accompanied by a pictogram.

WARNING: *A Warning indicates that the personal safety of the patient or physician may be involved.*

CAUTION: *A Caution indicates that particular maintenance procedures or safety precautions must be followed to avoid possible damage to the device.*

NOTE: *A Note indicates special information about operating the device or clarifies important information.*

WARNING: *Read this instruction manual thoroughly and be familiar with its contents prior to using the KARL STORZ SCB control NEO system. Pay particular attention to the safety instructions section to avoid putting your patients, personnel and yourself at risk.*

WARNING: *The instructions and interface specifications for medical devices and/or system components used in combination must be observed precisely.*

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Safety instructions

WARNING: *Combinations of medical devices are only assured to be safe if*

- *they are identified as such in the respective instruction manuals or*
- *the intended purpose and interface specifications of the devices used in combination permit this (cf. IEC 60601-1-1). Furthermore, all configurations must comply with the system standard IEC 60601-1-1. Anyone who connects additional equipment to the signal input part or signal output part of the system components configures a medical system, and is therefore responsible for ensuring that the applicable system requirements are met. If questions arise, please consult the technical service department or your local representative.*

WARNING: *Risk of explosion if the SCB devices are used in the presence of flammable anesthetics.*

WARNING: *To avoid the risk of electrical shock, do not open line-powered SCB devices! Refer servicing to qualified personnel.*

WARNING: *The electrical installation of the operating room must comply with the applicable IEC standards.*

CAUTION: *Avoid allowing fluids to enter the housing of the SCB devices. Do not store liquids on or above the SCB devices.*

CAUTION: *Only operate the SCB devices with the voltage indicated on the rating plate. Check whether the voltage selection switch is set to your line voltage.*

CAUTION: *When replacing fuses, use only fuses of the correct rating.*

NOTE: *Follow local governing ordinances and recycling plans regarding disposal or recycling of device components.*

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3.2 Safety instructions for the OR1™ control NEO

WARNING: *This device is not intended for use in potentially explosive atmospheres.*

WARNING: *Risk of explosion in the presence of flammable anesthetics.*

WARNING: *Danger of electrical shock! Do not open the housing of the OR1™ control NEO! Refer servicing to the manufacturer KARL STORZ or personnel authorized by the manufacturer only. Removal of covers by unauthorized personnel will void the equipment's warranty.*

WARNING: *Grounding reliability for the KARL STORZ OR1™ control NEO can only be achieved when the equipment is connected to the 'Hospital Only' or 'Hospital Grade' outlet (i.e. approved for use in an operating room environment). Routinely inspect the electrical plug and cord. Do not use if inspection reveals damage.*

WARNING: *The electrical installations in the operating room in which the unit is connected and operated must comply with the applicable IEC standards.*

WARNING: *The instructions for use and the interface specifications of the medical devices and/or system components used in combination must be observed precisely. Accessory equipment connected to the analog and digital interfaces of the unit must be set up and tested according to the relevant IEC standards, i.e., IEC 60601-1 for medical electrical equipment and also, for example, IEC 60950 for data processing equipment. Only signals with the specified values may be connected to the signal inputs and outputs.*

WARNING: *All system components (OR1™ control NEO, keyboard, etc.) are not parts to be used on the patient and must therefore not come into contact with the patient or be placed in the vicinity of the patient. For safety reasons, when the system is in use the output sockets of the OR1™ control NEO and the patient must not be touched at the same time.*

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Safety instructions

WARNING: *Test this device prior to each surgical procedure to ensure that it functions correctly.*

WARNING: *If the unit does not function reliably or safely, the manufacturer recommends that the unit is not operated, and, if necessary, is sent for servicing.*

WARNING: *Always unplug the KARL STORZ OR1™ control NEO from the mains before doing any maintenance work on it.*

WARNING: *An exchangeable 3.0 V lithium battery for the "onboard" clock and CMOS RAM is provided with the unit. There is danger of explosion if the lithium battery is replaced incorrectly. Used batteries may only be replaced by the same battery type. Refer all service to authorized personnel. Dispose of old batteries in an environmentally friendly way: take used batteries to official collection points in your area or return to the supplier (also see "KARL STORZ MEDICAL PCs" instruction sheet).*

WARNING *for touch screen: Even if no picture can be seen on the SCB touch screen, the touch screen functions may still be active. To avoid unintentional control commands, either the SCB control NEO software must be exited or the OR1™ control NEO switched off. Please note that in the case of a fault, the touch screen functions may still be operational irrespective of the image displayed. In this case, operations should be performed directly on the corresponding devices, and the SCB control NEO system turned off.*
To clean the touch screen the SCB control NEO software must be exited and the OR1™ control NEO switched off.

CAUTION: *Do not store liquids on or above the unit. Make certain to prevent any penetration of liquid into the housing.*

CAUTION: *Only operate the unit in the voltage range stated on the rating plate.*

NOTE: *Any damage to the system components resulting from incorrect operation is not covered by the manufacturer's warranty.*

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Before using the system, it is imperative that you be acquainted with how the KARL STORZ SCB control NEO system operates and is controlled.

3.3 Intended use

3.3.1 KARL STORZ SCB control NEO system

The KARL STORZ SCB control NEO system is designed for networked control of the SCB devices and accessories approved by KARL STORZ for this purpose. These SCB devices are indicated by the type number extension "-1" on the KARL STORZ SCB control NEO plate and/or the "SCB logo" in front of the device name.

Depending on the configuration (see examples of configurations in the installation manual), the KARL STORZ SCB control NEO system consists of at least one SCB input/output device for controlling and/or displaying the devices connected to the SCB.

Use of the KARL STORZ SCB control NEO system in fields other than those indicated above is not allowed for safety reasons.

The KARL STORZ SCB control NEO system may only be operated with accessories which have been designated by KARL STORZ as suitable for the KARL STORZ SCB control NEO system.

Unauthorized conversions or modifications to the KARL STORZ SCB control NEO system are not allowed for safety reasons.

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3.4 User qualification

The KARL STORZ SCB control NEO system may only be used by physicians and medical assistants who have a suitable specialized qualification and who have been instructed in the use of the system.

Requirements

You should be familiar with the use of the devices connected to the KARL STORZ SCB control NEO system (e.g., HAMOU ENDOMAT® model 26 331020-1, THERMOFLATOR® model 26 432020-1, etc.).

NOTE: *This instruction manual does not substitute the instruction manuals of the devices connected.*

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
3.5 Safety precautions at the site of installation

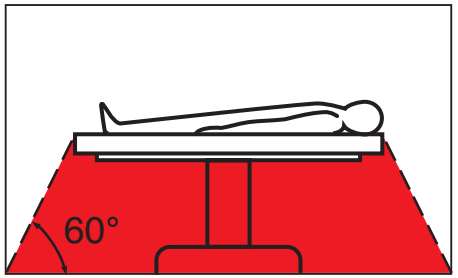The KARL STORZ OR1™ control NEO may only be used in medical rooms whose electrical systems have been installed in accordance with applicable national regulations.

The unit is equipped with a connector for attaching a ground line. It should be connected in compliance with the applicable national regulations.

It is not intended for use in hazardous zones. This means, for example:

When using easily combustible and explosive inhalation anesthetics or mixtures thereof, the unit must not be operated inside the hazard zone shown in the diagram. This also applies for easily combustible and explosive chemicals, e.g., skin disinfectants and fast-acting surface disinfectants.

Furthermore, please observe the instructions regarding place of installation in the instruction manuals of the individual devices connected to the SCB.

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3.6 Safety precautions when operating the KARL STORZ SCB control NEO system*

It is the user's responsibility to make sure the complete system and all devices connected to the SCB are safe and operate properly before using the KARL STORZ SCB control NEO system*.

During treatment in conjunction with the KARL STORZ SCB control NEO system the patient must be treated and kept under observation with the usual medical care. This includes, in particular, monitoring sterile application conditions where required by the type of intervention.

*(with or without KARL STORZ OR1™ control NEO)

The KARL STORZ OR1™ control NEO has been tested and found to comply with the EMC limits as per the Medical Device Directive 93/42/EEC (EN 55011 Class A and EN 60601-1-2:2001).

Please read the Electromagnetic Compatibility Information in the installation manual.

The ground line should be installed by a qualified electrician.

CAUTION: *If more than one device is connected to the KARL STORZ SCB control NEO system, please ensure that the max. permissible value for the housing leakage current is not exceeded. If necessary, fit an isolation transformer.*

NOTE: *Only applies if the entire system is connected to a carriage and* ***one*** *portable multiple socket outlet.*

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4 SCB control NEO software user interface

4.1 Optimal use of the SCB control NEO software

The basic functions described in the installation manual can be optimally used by you to save time and work:

1. The relevant SCB devices and the parameters associated with these devices can be grouped into a Preset according to the corresponding type of surgery and surgeon. This Preset is opened prior to the start of the operation, thereby setting the correct parameters on all the necessary SCB devices.
2. During the operation, all the important device parameters can be shown in a clear manner and, if necessary, changed directly using the so-called “widgets” (compact menus).
3. Depending on the situation, those specific SCB device operating functions which are not supported by the “widgets” (compact menus) are carried out centrally via the realistic view of the SCB device. To this end, an SCB device is chosen from the device list and can be operated via the device view displayed as well as via the front panel of the SCB device.
4. Operating functions can be controlled both from the sterile area (touch screen with sterile cover) and from the nurses’ station (touch screen, mouse or keyboard).

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User interface

The interface is divided into three main areas:

①	Information area
②	Application window
③	Navigation area

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4.2 Information area

4.2.1 Information line

SCB device information and warnings are shown in the information line in all applications of the KARL STORZ OR1 NEO™ system according to their priority. The details on the information line in the system manual must therefore also be noted.

Information messages are shown system-wide in the header. If several messages need to be displayed simultaneously, these are repeatedly shown one after the other. In the case of warning messages, the warning which was received last is always shown until it is acknowledged, the SCB device retracts this warning, or a warning with a higher priority status is received.

The display of the device messages as of SCB control NEO software version 20 0900 01-41 differs from the display in the camera image. The display of the device warnings is explained in the corresponding chapter for every device.

4.2.1.1 Information

Information messages are shown in white font in the header. The information messages are hidden automatically once the cause thereof has been remedied.

In addition to the display in the information line, the SCB devices which display messages in the information line are shown with a white frame in the device list.

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
- Medium (yellow, flashes slowly): The problem should be remedied at the next opportunity.

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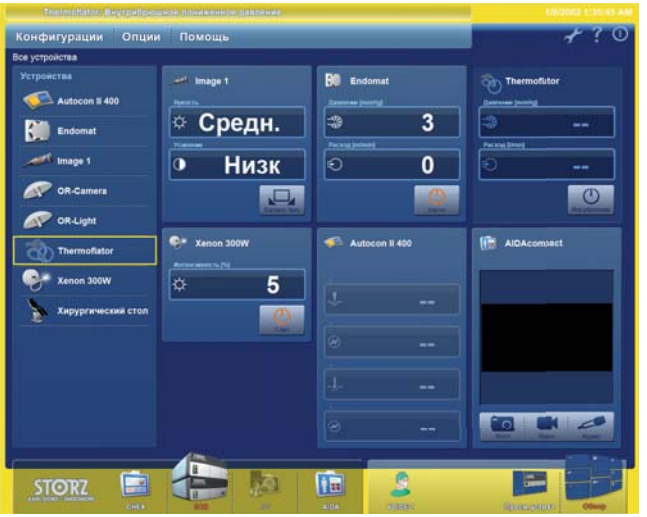
- Low (turquoise, does not flash): You should remedy this problem after the operation for example.

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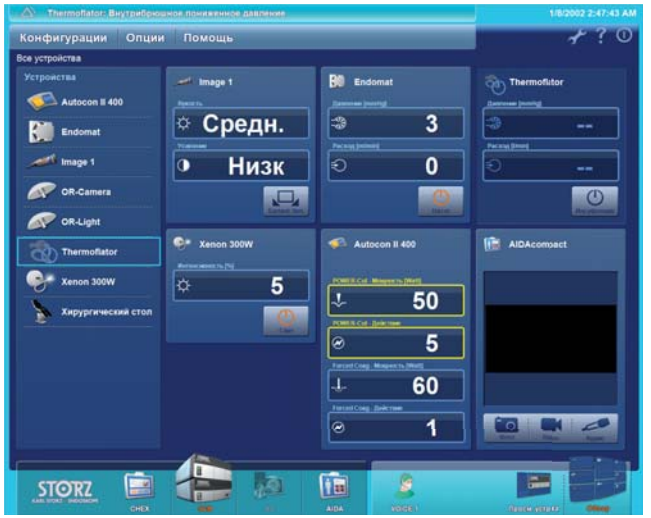The priority level is also shown in the device list by a flashing box in the corresponding color.

 Button to confirm a warning

 Button to mute a warning

The device determines whether the information shown is to be acknowledged using the “Confirmation button”. The option of silencing the acoustic signal on the device is also given by the device via the “Mute On/Off” button as well as the option of displaying further assistance via the “Info” button.

4.2.2 Date/Time

Please observe the information in the system manual.

The Date/Time can be set in the SCB control NEO software via the menu item Options/Maintenance Manager…/Start other KARL STORZ programs…/“Regional Settings” and “Date Time”.

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4.3 Navigation area


More detailed information can be found in the system manual. The navigation area is divided into the main menu and the workflow bar.

①	Main menu
②	Workflow bar

4.3.1 Explanation of buttons in the navigation area

The buttons of active applications are enlarged on the main menu. If an application is not available, the button is shown in gray or not displayed. On touching the buttons a change in the color of the font indicates to the user which buttons will become active when “released” (see 4.4.3.1).

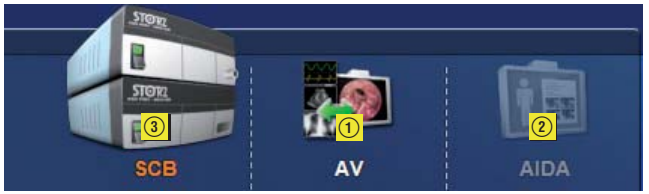①	Application inactive – button displayed normally
②	Application not ready – button shown in gray
③	Application active – button enlarged
	Application not available – no button shown

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
4.3.2 Main menu

More detailed information can be found in the system manual. In the main menu you will find the buttons to switch between the individual applications.

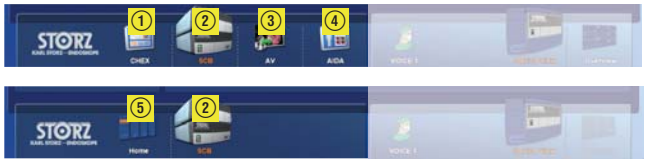① Button: OR1™ CHEX or an application specified by AIDA.
② Button: SCB control NEO
③ Button: AV system NEO
④ Button: AIDA compact NEO
⑤ Button: KARL STORZ OR1™ FUSION

NOTE: *The type and number of available applications may differ from what is shown depending on your system configuration. The displayed applications can be selected via the "Maintenance Manager.../Add/Remove device ...".*

NOTE: *A detailed description of each application can be found in the respective instructions manuals.*

4.3.3 Workflow bar

You can access the functions within the SCB control NEO software via the workflow bar:

- Realistic view of the device
- Overview
- Optional: VOICE1®

The overview provides a compilation of the most important device parameters in the so-called "widgets" (compact menus) depending on the preprogrammed type of surgery. This allows the device parameters to be monitored at a glance and, if necessary, changed.

The realistic view enables you to operate the device in exactly the same way as from the front panel of the device. As such, complete functionality of the SCB device is granted without you having to learn another operation mode.

The SCB control NEO software can be optionally enhanced with the KARL STORZ voice control VOICE1®. The workflow button VOICE1® then enables the user interface to be activated using VOICE1®. A detailed description of VOICE1® can be found in the respective instruction manual.

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

4.4 Application window/content area

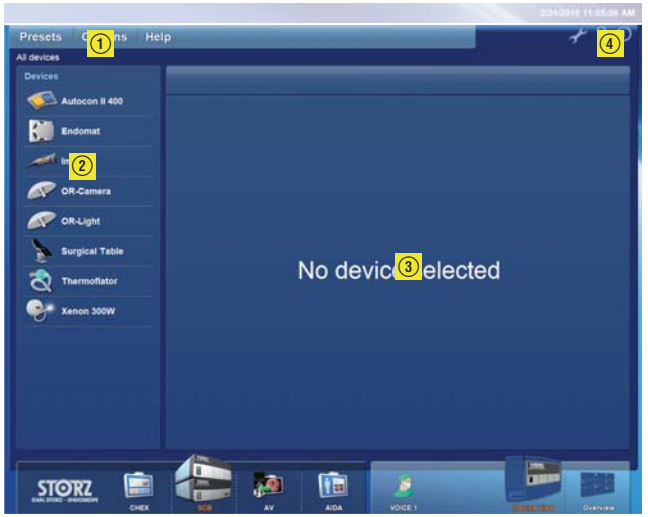The application window is divided into four “fixed” areas:

①	Application header/main menu
②	Device list/main application window
③	Device window or “widgets” (compact menu)/main application window
④	General functions/main menu

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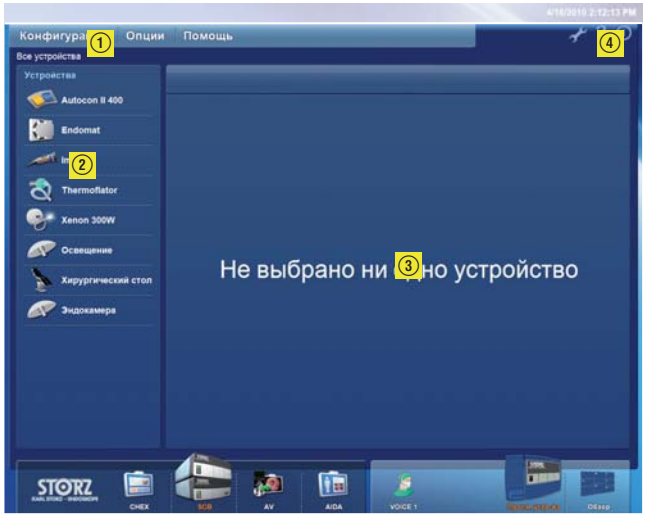4.4.1 Application header (menu bar)

The header is divided into several areas:

①	Menu Pulldown menu Presets, Pulldown menu Options and Pulldown Help
②	Name/Info on current Preset

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
4.4.1.1 Presets

The Preset submenu is used to manage the various compilations of SCB devices with their device parameters.

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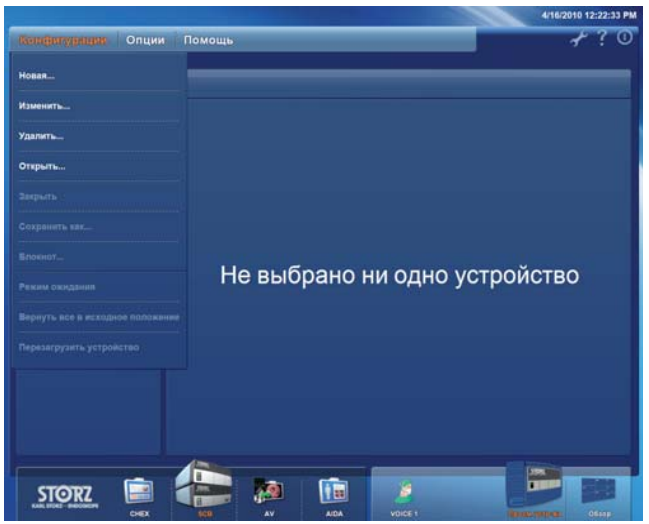The submenus mean:

Entry/ Submenu	Explanation
New	Create new Presets. (can be protected with a password)
Change...	Change and open existing Presets. (can be protected with a password)
Delete...	Delete an existing Preset. (can be protected with a password)
Open...	Select and open an existing Preset.
Close	Close the current Preset. The system returns to the default Preset “All devices”.
Save as...	The current device settings are saved under the same or a new Preset. (Can be protected with a password). In addition to the device parameters of the open Preset, the current device parameters of the SCB device shown in the realistic view are also saved. (Only available in the realistic view)
Notepad	The user can save notes on the created Presets on the electronic notepad. These are saved together with the Preset.
Standby	The connected devices are set to “Standby”.
Reset all	All devices are reset to the original settings of the currently open Preset. Only active in the overview.
Reset device	The currently called up device is reset to the original settings of the currently open Preset. Only active in the realistic view.

NOTE: *In the header directly under “Presets” menu the name of the current Preset is displayed as follows “Name/Description (optional)/ Type of surgery”.*

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4.4.1.2 Options


The “Options” menu is used to set application settings. Further information and explanations can be found under “Modal Dialogs”.

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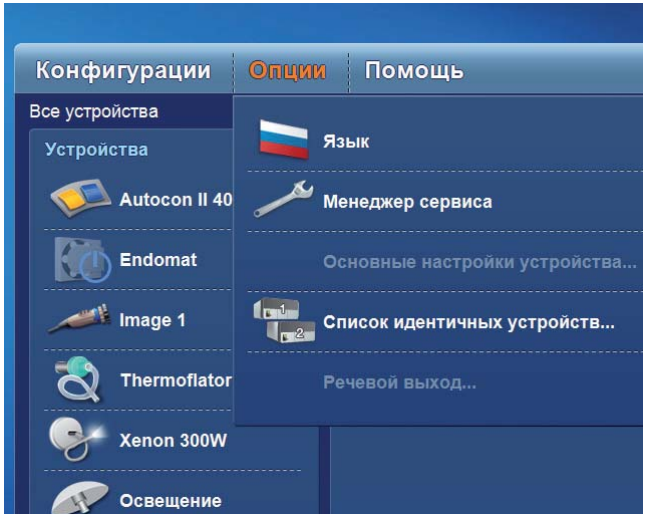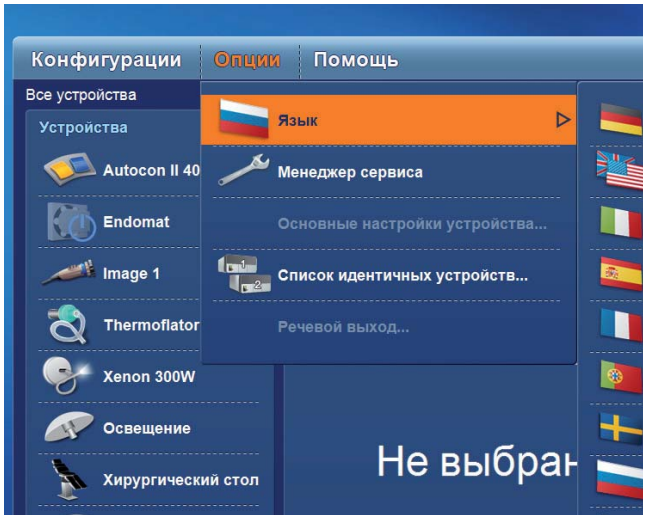
SCB control NEO software user interface

The submenus mean:

Entry/ Submenu	Explanation
Language	To set the language of the user interface. The names of the available languages are shown in the respective language.
Maintenance Manager…	To start the “Maintenance Manager” which offers the following options: • Add or delete SCB devices for possible selection in Presets • Management of aliases (replacement of device names) • Export of Presets • Import of Presets • SCB control NEO software update • Calibration of the KARL STORZ touch screen • Management of passwords • Calling up of programs for setting the system
Device setup	The menu item “Device Setup” can only be selected if the SCB device selected from the device list supports this functionality. The SCB device appears in the device list, the status of the SCB device is displayed with a superimposed icon and a blue frame. An additional dialog box is shown below the realistic view. The values shown on the realistic view are “Device Setup parameters”.
Identical device list...	The SCB control NEO system offers the option of operating several identical devices of a device type in the SCB network and controlling these centrally via the SCB control NEO software. To be able to distinguish between the devices, they are each assigned an individual distinctive number. A maximum of 3 devices each can be controlled from the following devices via the SCB control NEO software. AUTOCON® II 400, COVIDIEN ForceTriad™, Electronic ENDOFLATOR®, IMAGE1 HUB™, powerLED 175, THERMOFLATOR®, TRICAM® SL II, XENON 175 and XENON 300

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Entry/ Submenu	Explanation
Settings for speech output options, e.g., volume	This menu item opens a window in which settings for speech output can be made. Amongst other things, the speech output can be switched on and off here. (The menu item is gray if the currently selected language is not supported).
Log In / Log Out...	If password protection has been configured, this menu item is shown. A dialog box opens which enables the user to log in/log out for access to protected functions.

4.4.1.3 Menu “Help”

- Quick guide
- Info... (“About” Dialog)

The “Quick guide” menu item displays a quick guide/brief description for the KARL STORZ OR1 NEO™ system on the screen.

Information on “Operation” can be found under 4.4.5.7.1. The “Info... (“About” Dialog)” menu item provides you with information about the system (see 4.4.5.7.2).

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
4.4.2.3 Exit

- Closes the application and shuts down or restarts the OR1™ control NEO computer (depending on the setting “SCB control NEO Exit Options”).

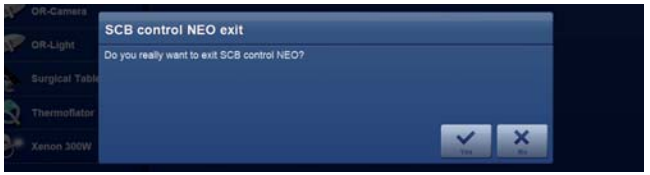• Switch off

- Select the “Exit” button in the header. Confirm the prompt.
- According to the setting in the “Maintenance Manager/Start other KARL STORZ programs/Exit Option” the OR1™ control NEO is restarted or shut down (also see 4.4.5.1 Exit SCB)
 Default setting: Shutdown

WARNING: *If the computer is turned off at the On/Off button of the OR1™ control NEO, this may impair functionality. Always switch off the computer using the button.*

• Restart

The SCB control NEO application is restarted when the OR1™ control NEO is switched on again at the ON/OFF button.

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4.4.3 Device list

The device list gives the user the option of noting the status of the devices used in the Preset at a glance.

If more than 12 SCB devices are connected, their device names (i.e., from SCB device 12 onwards) are no longer displayed. The device names which are not visible can be selected using the two arrow buttons at the bottom of the left column. Tap the button briefly to move sideward to the next device; keep the button depressed for longer to scroll through the devices automatically (Scroll/Autorepeat).

ⓘ ***NOTE:*** *In the default Preset "All Devices" which is included in all devices, devices which are not switched on are not shown in the device list. With this Preset, the devices are installed without values.*

4.4.3.1 Device status display

Font:

White: The device is switched on and is currently not shown on the right-hand side.

Gray: The device is not switched on or cannot be operated in the current status.

Orange: This device is shown as active on the right-hand side.

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Icons:

Normal: The device is switched on and ready for use.

Standby: The device is in Standby mode.
Pumps and insufflators: The pumps are switched off.
Light sources: The emission of light is interrupted.

Switched off: The device is included in the Preset but is not yet switched on.

Device Setup: The device is in Device Setup mode. The device setup parameters can be shown and changed in this mode. The device is thus not in operation mode.

Frame color:

Gray: The device is in the device setup.
White: The device contains information which is shown in the information line (see 4.2.1.1 “Information”).
Turquoise: The device has a low-priority error (see 4.2.1.2 “Warnings”).
Yellow: The device has a medium-priority error (see 4.2.1.2 “Warnings”).
Red: The device has a high-priority error (see 4.2.1.2 “Warnings”).

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4.4.4 Device window or “Widgets” (compact menu)

4.4.4.1 Device view (realistic view)

A realistic view of the device which is selected in the device list is shown.

If the device is in the device setup, in some cases an additional setup dialog box is displayed in which certain service settings are to be made.

WARNING: *The correct working of the graphic user interface must be checked on the active device. To this end, activate at least one function using SCB control NEO and check that the real SCB device reacts correctly.*

The realistic view of the devices is shown as soon as a Preset is opened, the user selects the device view in the workflow bar or the name in a widget is selected.

The device view of a specific device is shown by clicking on the device name. The “widgets” (compact menus) available for the selected Preset can be activated with the open Preset via the “Overview” button in the workflow bar.

WARNING: *Devices whose presettings are saved in the selected Preset are initialized with this.*

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
In the left-hand column all the preprogrammed devices are run for the activated Preset. The Preset description can be found above the device list (Name / Description (optional) / Type of surgery).

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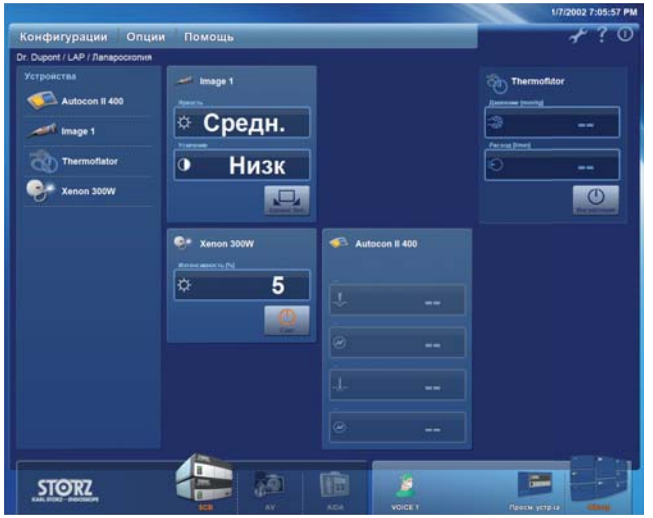4.4.4.2 “Widgets” (compact menu – Overview)

The displayed widgets give the user the option of checking the functions offered by the devices at a glance. The following functions are supported:

Change the camera setting, set the light intensity, control the irrigation and suction pumps, control the insufflation devices, set the speed of the motor-operated devices, operate documentation devices, set HF devices, etc.

A widget is available for each of these functions. The respective device is displayed in the widget header. Click on the name to switch to the realistic view of the device shown. The widget view is accessed via the Overview button in the workflow bar.

Various devices, e.g., UNIDRIVE® ENT, NEURO and OMFS support several widgets, e.g., the shaver and pump widgets.

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
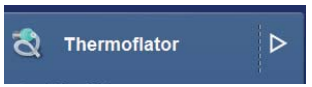
If, on the other hand, several devices support a function, an arrow button is visible on the right-hand side of the respective widget for scrolling.

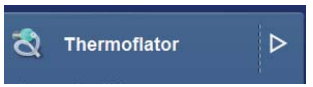The arrangement of the widgets is deliberately fixed so that the functions are always in the same place for the user. Hereby the selected type of surgery under “Create Preset” determines which widgets are shown and which not.

Arrangement of widgets

Camera systems	Pumps	Insufflators
Light sources	HF devices	Documentation
Motor-operated systems		

The following table shows which functions are available in the various surgery types:

	Camera systems	Light sources	Documentation	Pumps	Motor-operated systems	HF devices	Insufflators
ENT	X	X	X	X	X	X	
Laparoscopy	X	X	X	X	X	X	X
Gynecology	X	X	X	X	X	X	X
Urology	X	X	X	X	X	X	X
Arthroscopy, orthopedics	X	X	X	X	X	X	
Cardiovascular surgery	X	X	X	X		X	X
Neurology	X	X	X	X	X	X	
Plastic surgery	X	X	X	X		X	

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A special display and input field is used for a compact display. This shows the actual value in the operating mode.

As soon as this is clicked on and activated, the displayed set value can be set. When entering the set value, the relevant buttons for changing the value are displayed at the sides of the input field and the value is shown in orange. Continuous values additionally show the corresponding value range via a level display. If no more changes are made in the set value display for approx. two seconds, the display or input field switches back to the actual value display.

If a device is set to standby mode, this is shown by a change of color in the widget and the standby icon is added to the device icon. The display or input fields can still be used for entering the set value.

NOTE: *"Widgets" (device compact menus) cannot be explicitly selected for a device.*

Widgets are shown when:
- *the device has been added to the Preset,*
- *there is a widget for the selected device function and*
- *the selected area of specialization supports the widget.*

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
The following widgets are supported:

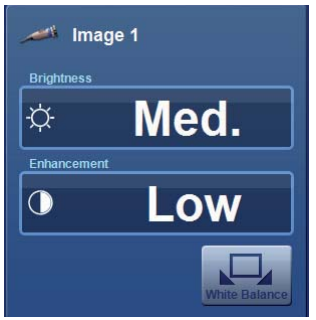- Endoscopic camera
 Brightness setting (Low – Medium – High)
 Enhancement setting (Low – Medium – High – Filter A/B)
 Carry out white balancing

ⓘ ***NOTE:*** *The setting parameters of the IMAGE 1 H3-M deviate from this (see corresponding section).*


- Light source
 Brightness setting (0…100%; increment of 20%)
 Switch between standby and “active”.

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
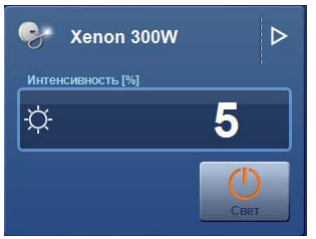
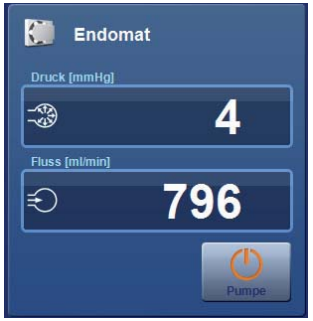
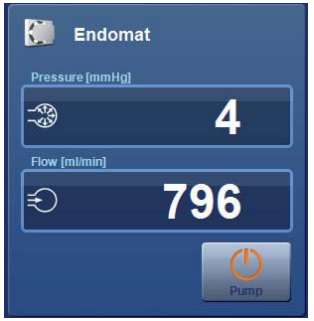
- Pump
 Pressure setting (depending on the device)
 Flow setting (depending on the device)
 Switch between standby and “active”

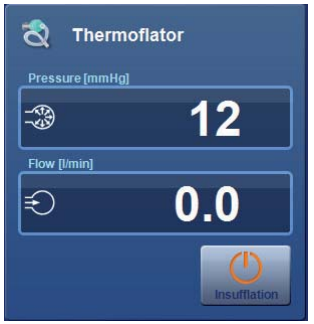- Insufflator
 Pressure setting (depending on the device)
 Flow setting (depending on the device)
 Switch between standby and “active”

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

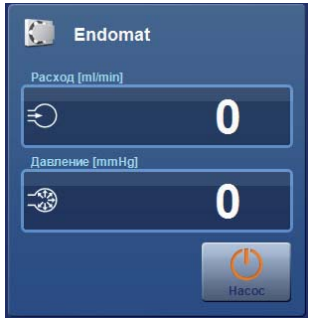
- Motor-operated systems
 Speed setting (depending on the device)
 Display of direction of rotation

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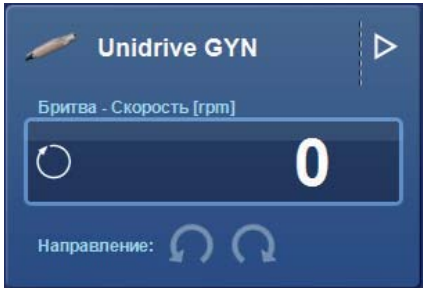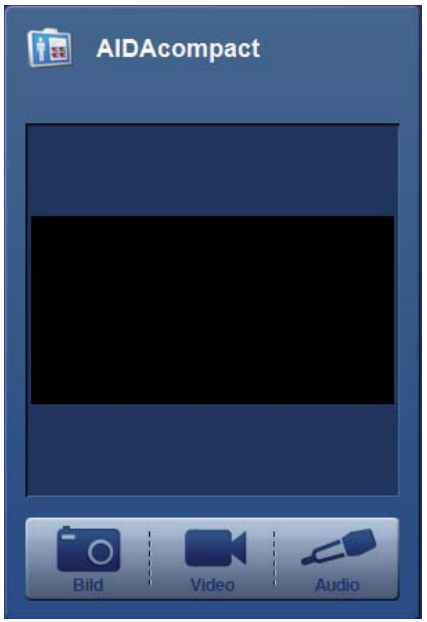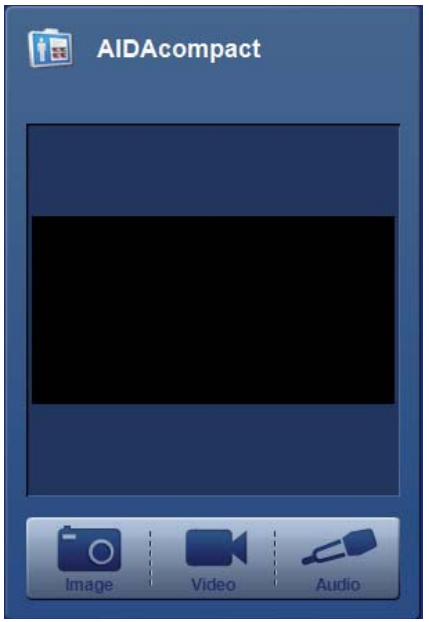- AIDA compact
 Start image recording
 Switch between start/stop video recording
 Switch between start/stop audio recording

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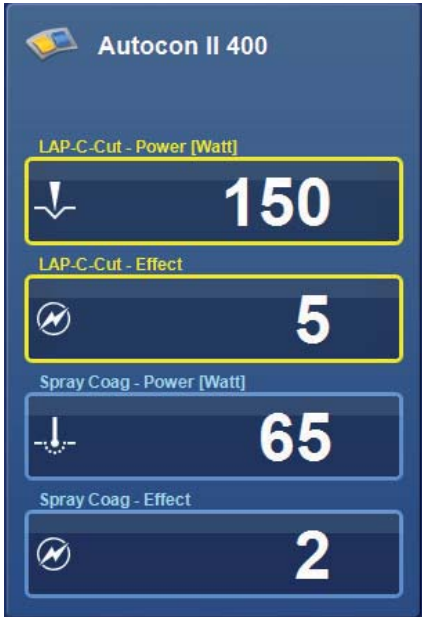
- HF device
 Display: Mode (Cut, Coag, Standby)

 “Cut Power”:
 Watt setting (depending on the device and mode)

 “Cut Effect”:
 Effect setting (depending on the device)

 “Coag Power”: Watt setting
 (depending on the device and mode)

 “Coag Effect”: Effect setting (depending on the device)
 Switch between standby and “active”

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
You can activate the (device) detail display by clicking on the device name in the left-hand column or the “Device view” button in the workflow bar.

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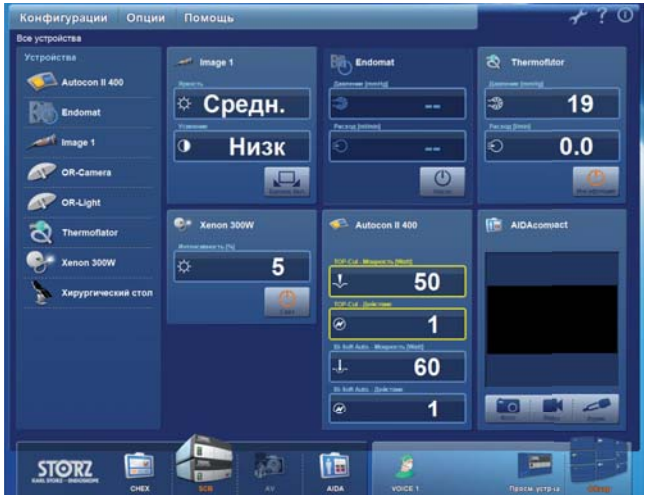
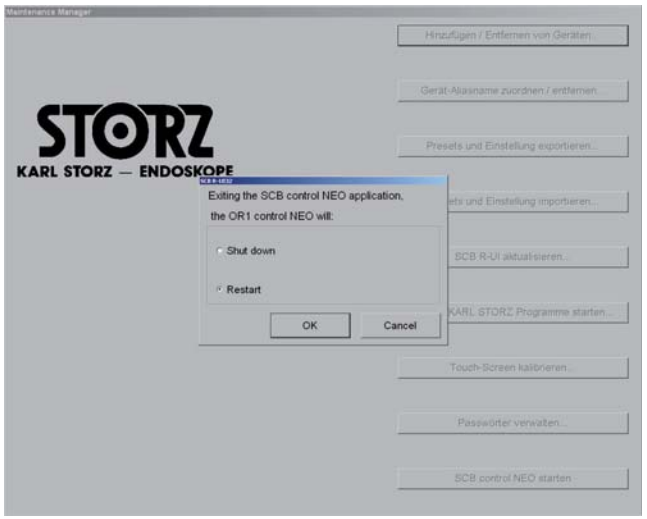
4.4.5 Modal dialogs

4.4.5.1 Exit SCB

This dialog box is opened when the “Exit” button is called up. See **4.4.2.3 Exit**

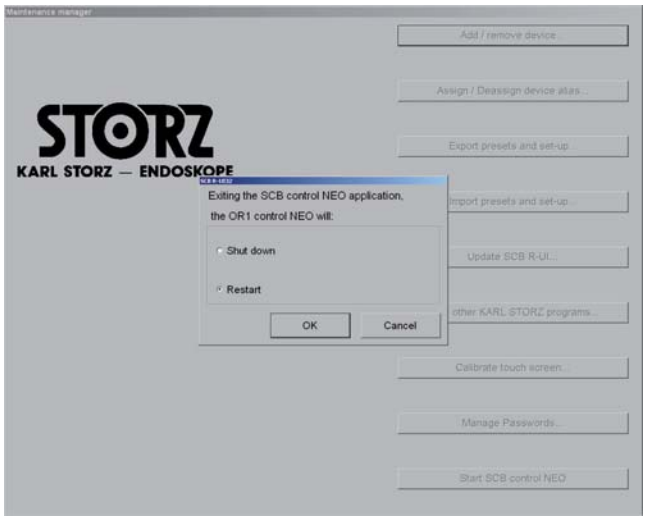
- The exit procedure is determined via the Maintenance Manager… / Start other KARL STORZ programs… / “SCB control NEO Exit Options”. One of the following options must be selected:
 – Shutdown: The OR1™ control NEO is shutdown.
 – Restart: The OR1™ control NEO is restarted.

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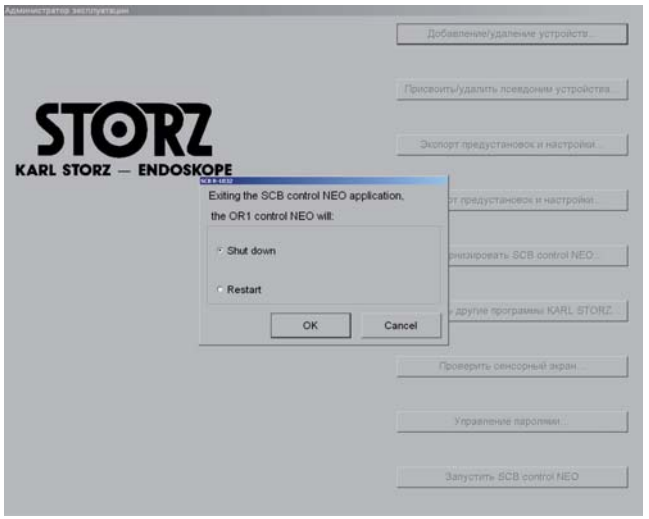On-screen keyboard

Function

With the on-screen keyboard you can enter text and all necessary commands and characters using the touch screen function. The on-screen keyboard is displayed automatically when you click on a text input field. The layout of the on-screen keyboard adapts automatically to the currently set country-specific hardware keyboard. Select “Regional Settings” to adapt the country setting of the operating system of your keyboard if necessary (Regional Settings – Input). Text is entered by actuating either the buttons on the on-screen keyboard or the hardware keyboard.

SCB control NEO software text entry using the on-screen keyboard

The on-screen keyboard can be used to enter or change text in the text input fields of the SCB control NEO software (e. g., user name of a Preset). The on-screen keyboard essentially looks the same as a normal hardware keyboard but without a number pad.

The on-screen keyboard is activated by touching the text input field. The content to be edited appears in the on-screen keyboard. In the title window the current context (e. g., “New Preset”) and below that the name of the text input field (e. g., “Name”, ...) is shown. The length of the useable characters is adapted to the text input field.

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Description of the buttons

The buttons of the on-screen keyboard function in the same way as those of a hardware keyboard with the following exceptions:

- The buttons “OK” or “Return” close the on-screen keyboard. The input characters are entered into the text input field.
- The buttons “Cancel” or “Esc” close the on-screen keyboard. An input character will not be entered. The original text in the input field remains unchanged.
- The active buttons “Caps Lock”, “Shift”, “Ctrl”, “Alt” and “Alt Gr” are shown in orange. The same button must be pressed again to deactivate.
- The “Shift” button is an exception, the character which is next entered is input as a capital letter and the orange font then disappears.
- Upon activation of the “Caps Lock” button the characters which are then entered are input as capital letters. When the “Caps Lock” button is pressed again the characters which are then entered are “small”.

NOTE: *The selection of elements on the user interface via the keyboard is not supported.*

NOTE: *The buttons F1 … F12, SysRq, Scroll Lock (scrolling), Pause, Num Lock and the special Windows® buttons are not needed to enter text and are thus not displayed.*

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

The special TELECAM SL/TRICAM® SL/TRICAM® SL II/IMAGE 1 HUB™ mode of the on-screen keyboard enables you to run commands and characters for the cameras. The command and character functions are available just as with the camera hardware keyboard and can be referred to in the instruction manuals of the respective devices.

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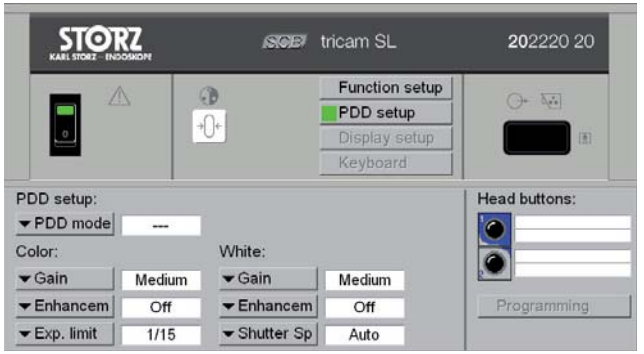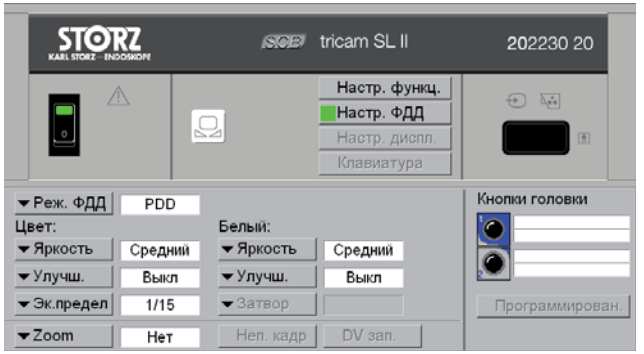
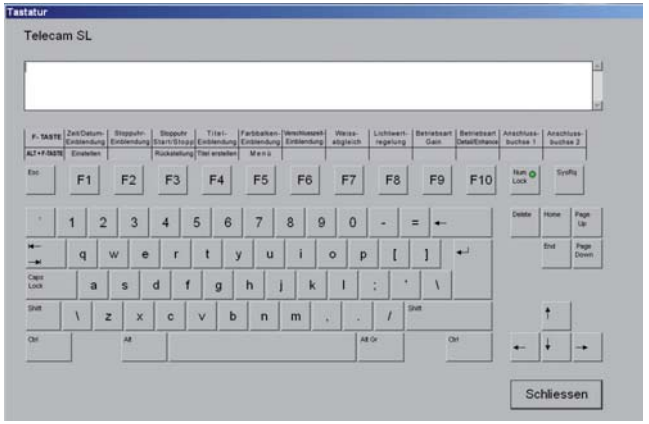
Keyboard in TELECAM SL/TRICAM® SL/TRICAM® SL II/ IMAGE 1 HUB™ and IMAGE 1 H3-M mode

The on-screen keyboard offers a special input mode for the KARL STORZ video devices TELECAM SL, TRICAM® SL/SL II, IMAGE 1 HUB™ and IMAGE 1 H3-M with which control characters and text (insertion of title in the video image) can be entered. The on-screen keyboard is activated by actuating the “Keyboard” button in the realistic device view. The commands and characters which you enter are continually shown.

Commands and characters can be repeatedly entered. As well as via the open on-screen keyboard, commands and characters can also be entered via the hardware keyboard of the OR1™ control NEO. All buttons are supported.

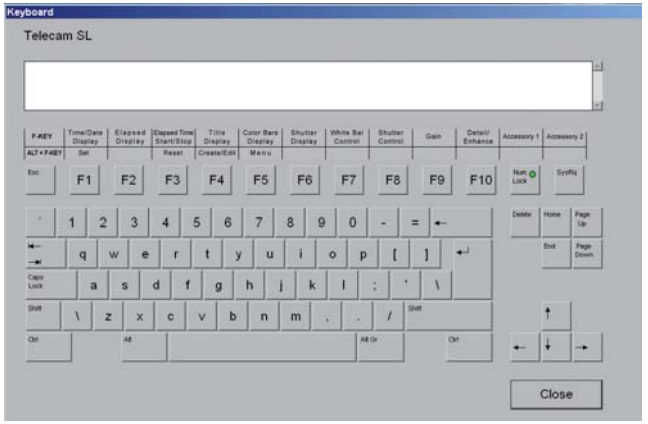Special command buttons


A description of the command buttons can be found in the instruction manuals for the TELECAM SL, TRICAM® SL, TRICAM® SL II, IMAGE 1 HUB™ and IMAGE 1 H3-M.

The following command buttons can be executed:
F1 … F10, NumLock, SysRq, Alt-F1 … Alt-F10, Alt-a, Alt-s, Alt-t, Alt-x, Alt-z, Alt-Home, Ctrl-8, Ctrl-s, Esc, Alt-l, Alt-p, Ctrl-F1 … Ctrl-F6, Ctrl-Alt-F7, Ctrl-a, Ctrl-c, Ctrl-d, Ctrl-e, Ctrl-f, Ctrl-i, Ctrl-k, Ctrl-l, Ctrl-n, Ctrl-t, Ctrl-u, Ctrl-w, Ctrl-x, Ctrl-y, Ctrl-z, Ctrl-Home.

Клавиатура в режиме TELECAM SL/TRICAM® SL/TRICAM® SL II/ IMAGE 1 HUB™ и IMAGE 1 H3-M Modus

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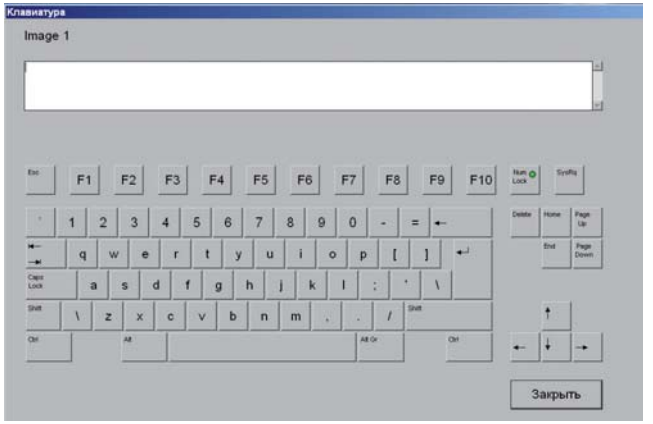Description of the buttons

- The “Close” button closes the on-screen keyboard.
- The buttons “Caps Lock” and “Shift” include an additional LED display which shows whether the respective button is currently active. The same button must be pressed again to deactivate.
- The “Shift” button is an exception, the character which is next entered is input as a capital letter and the LED display then goes out.
- Upon activation of the “Caps Lock” button the characters which are then entered are input as capital letters. When the “Caps Lock” button is pressed again the characters which are then entered are “small”.

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
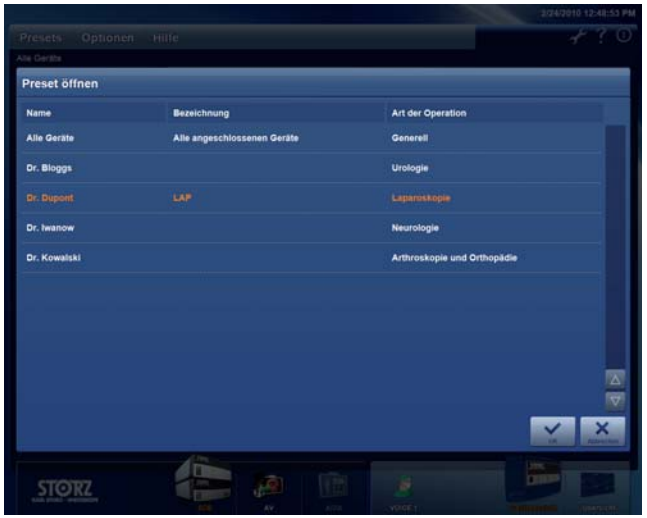
4.4.5.3 “Open…”/“Delete…” Presets

“Open…”

The menu item “Open…” is used to open/start a saved Preset.

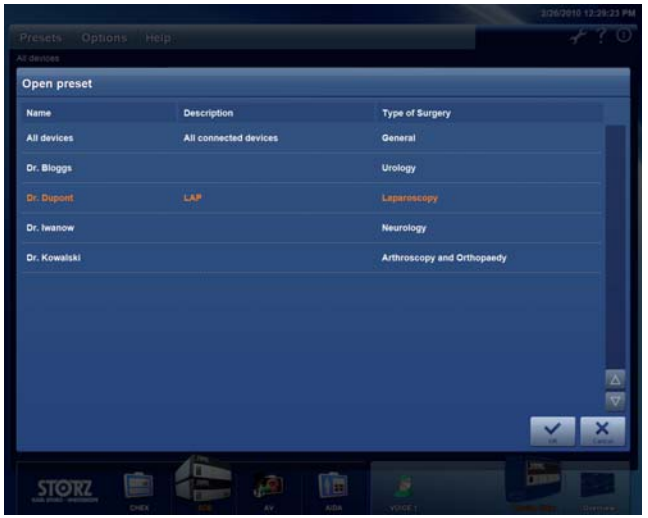You firstly see a three-column list with name, optional description and type of surgery of the existing Presets.

Select the desired Preset. Confirm your selection with the “OK” button.

The selected Preset is started straight away.

WARNING: *All SCB devices are initialized with the saved data of the selected Preset.*

WARNING: *Check the plausibility of the device set values before starting the operation.*

WARNING: *For SCB NT devices the selected procedure/program settings on the device can also be included in the selected Preset depending on the setting.*

Following successful initialization you will see the device list with the SCB devices specified in the selected Preset.
In the event of an error, an error dialog box (see page 60) is displayed.

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The Preset name is displayed below the menu. Once the devices have been initialized with the Preset data, the device list is shown on the “left” and the device view of the activated device is shown on the “right”. If possible, the first device is switched to active automatically otherwise the message “No device selected” appears. The user can activate a device manually.

If a device in a Preset is not available because it is switched off or not available, the device icon and name in the device list are shown in gray.

NOTE: *Up to 250 Presets can be saved. As of 13 Preset entries, a scroll bar (below) appears which allows you to move sideward.*

NOTE: *If, when importing configurations from older versions (see 4.4.5.8.4 Import Presets and settings), Presets are found with an ambigious name, “#1”, “#2”, “#3”, etc. are added to such names in order to make then unambigious. The names may also have to be abbreviated to comply with the maximum length for Preset names.*

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
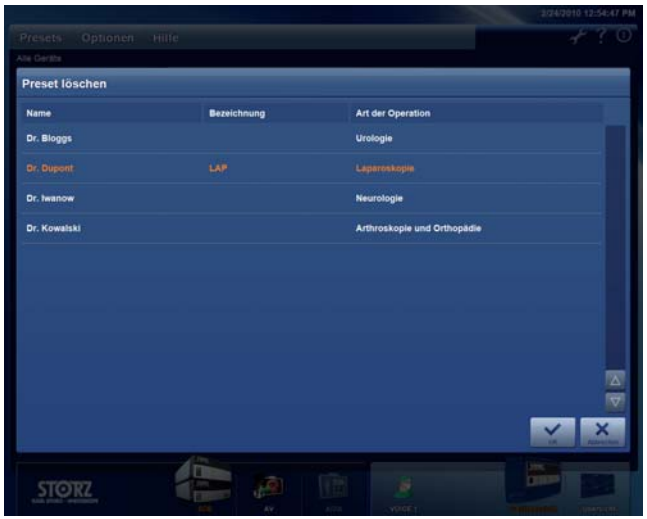
“Delete…”


The menu item “Delete...” is used to delete an existing Preset.

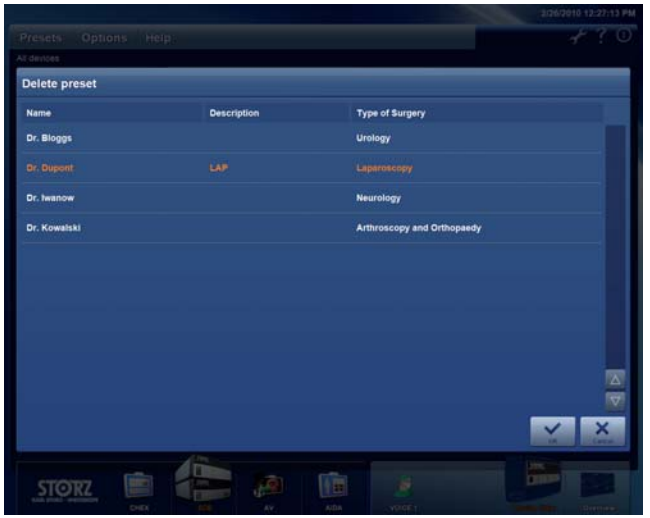
The system shows you a three-column list with name, description (optional) and type of surgery of the existing Presets (see menu item “Open…”).
Select the desired Preset by selecting the line.
Confirm deletion of your selection with the button “OK”, the Delete Preset dialog box opens and, after confirmation, the selected Preset is deleted.
The process can be stopped at any time by pressing the “Cancel” button.

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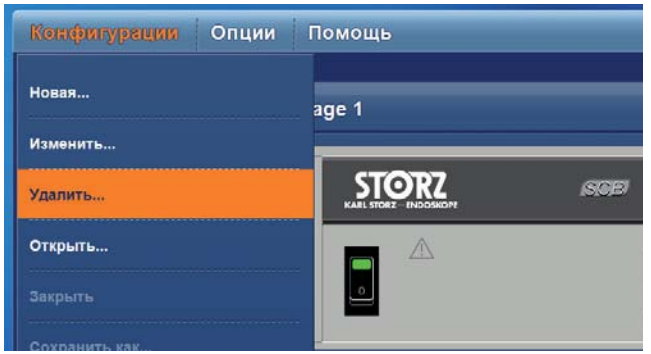
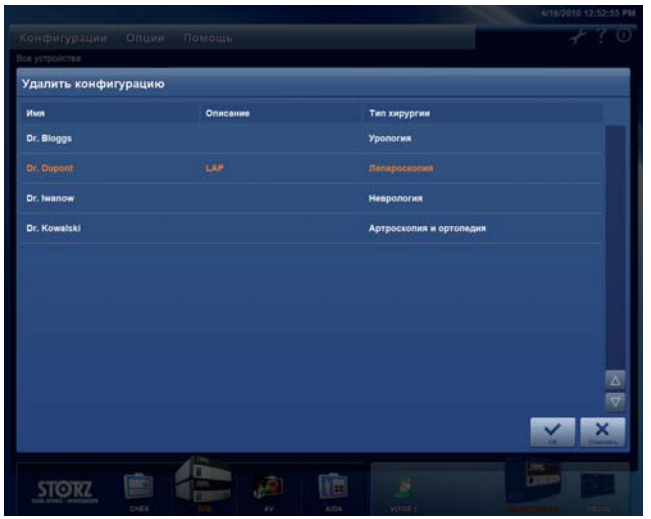
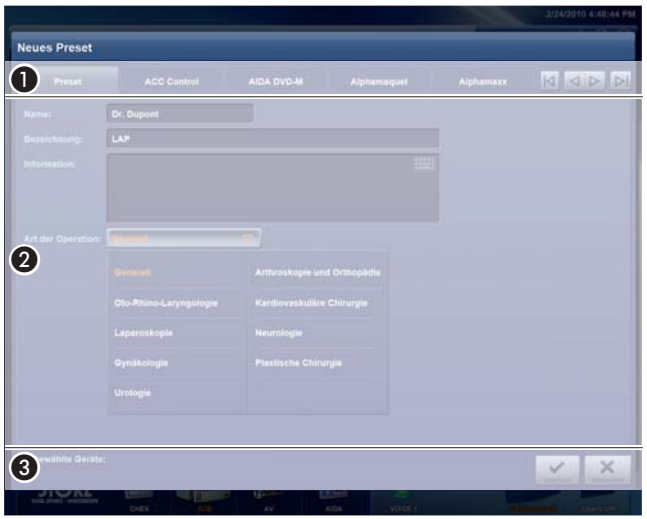
4.4.5.4 “New…”/“Change…” Presets

“New…”

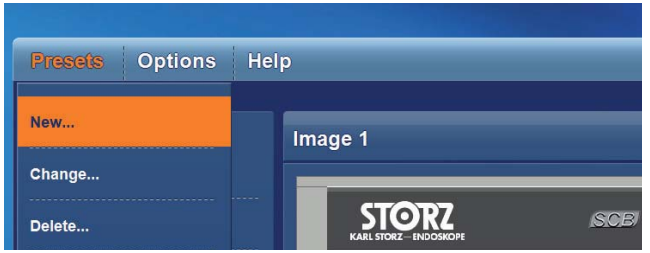The menu item “New…” is used to create new Presets

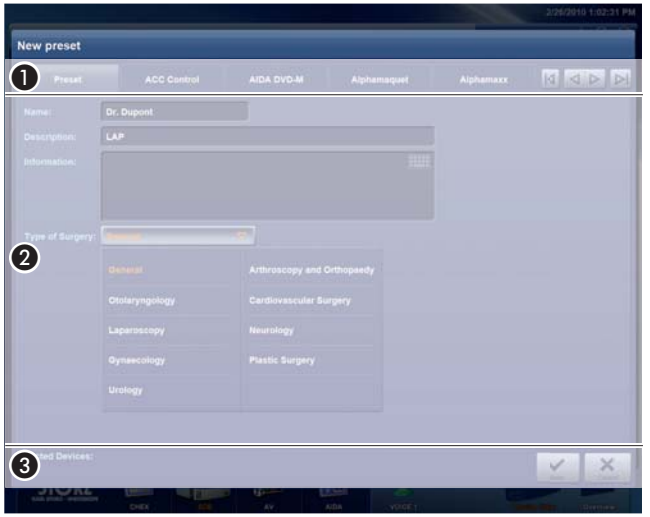An input window with the Preset dialog box appears:

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Tab line

With the help of the tabs, select the devices which you wish to include in the new Preset.

In the info and save field the devices added to the Preset are shown in the bottom left corner.

Scrolling between pages is done with the inner two arrow buttons “◁” and “▷” in the tab line. The two outer arrow buttons “|◁” and “▷|” allow you to jump to the beginning or end of the list. The arrow buttons are shown in “gray” when the beginning or end of the available tabs has been reached.

Info and save field


Pressing “Save” saves and starts the new Preset. The presettings you have set on the devices are then transferred to the devices shown under “Selected devices”. At the same time the device settings for the devices to be synchronized are added to the Preset (see input field).

The configuration process can be stopped at any time by pressing the “Cancel” button. The synchronization dialog box is called up by pressing Save.

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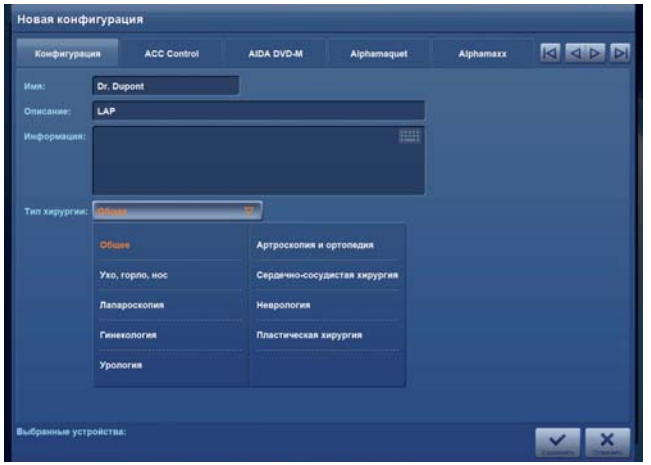Save / Cancel


Input field

The following information must be entered for the whole Preset in the first tab.

- Name
- Description (optional)
- Information (more detailed text, optional)
- Type of surgery (general, ...)

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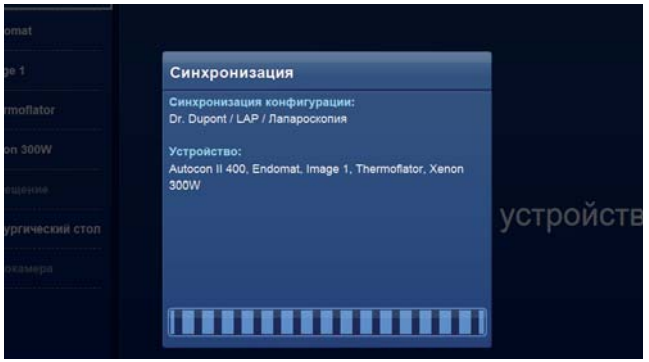The entries can be input using the on-screen keyboard (also see page 42 ff.) or, in the case of surgery types, selected from a pull-down menu. The text “Name/Description (optional)/Type of surgery” appears in the header (directly below “Presets”) after a Preset has been opened.

ⓘ ***NOTE:*** *If no “Name” has been entered for a new Preset an error message appears.*

ⓘ ***NOTE:*** *The name must be unambigious. The Preset name can have up to 23 characters. Two or more Presets with the same name cannot be saved.*

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NOTE: *The display of the selected device function groups in the "Widget display" (compact menu) depends on the type of surgery. All devices are available with "General".*

NOTE: *"Widgets" (device compact menus) cannot be explicitly selected for a device. Widgets are displayed if:*

- the device was added to the Preset.
- there is a widget for the selected device function and
- the selected type of surgery supports the widget.

As soon as the general Preset details are entered in the first tab, the corresponding devices can be selected, displayed and added to the Preset via the tab line. There are three types of devices to differentiate.

Devices which can only be added to one Preset

Devices which do not support the setting of Preset data have no further control elements.

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Devices which can be added to a Preset and preconfigured

Select the tab labeled with the desired device in order to make the necessary settings (in the example above: ENDOMAT®). In the field below the tab the display of the selected device can be seen.

Add the device to the Preset by pressing the “Add device to Preset” button.

NOTE: *Once the device has been added to the Preset, the tab lettering turns orange.*

NOTE: *Settings can only be performed once the corresponding device has been added to the Preset (“Add device to Preset” checkbox is selected).*

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

Devices which can be saved to a Preset and whose set values can be copied to the Preset upon saving (ARTHROPUMP® POWER, AUTOCON® II 400 and KLS Martin maxium® ME402)

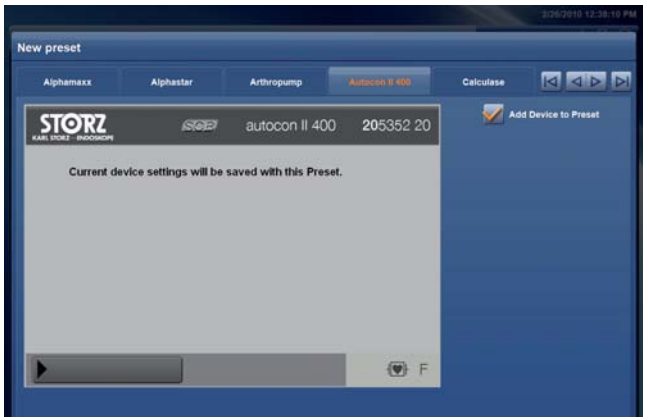- These devices include a field which describes the individual steps in order to load the desired values from the device to the Preset.

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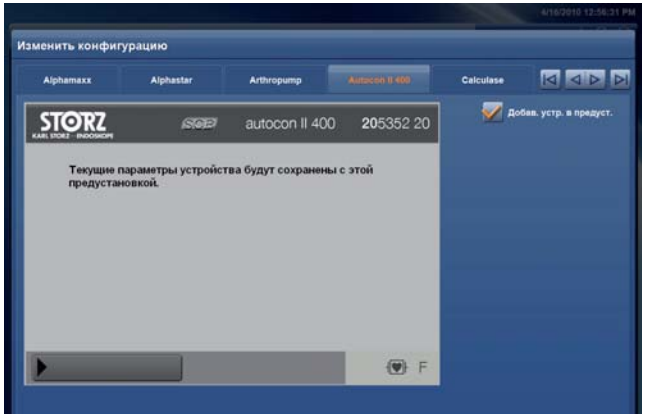
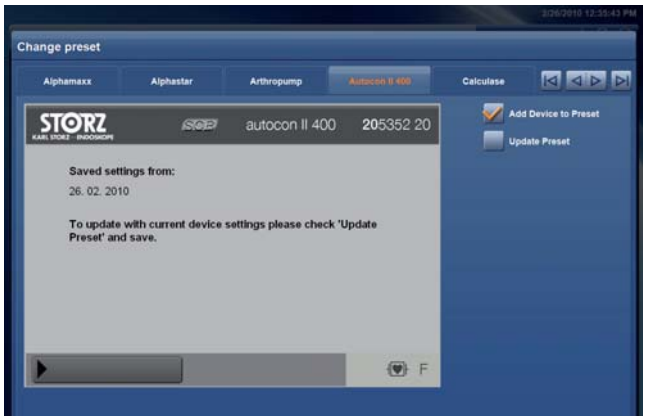
- Once the values have been copied into the Preset once, on the right-hand side the button “Update Preset” is displayed. Synchronization is started automatically when a device is added to a Preset. This must then be started manually by the user.

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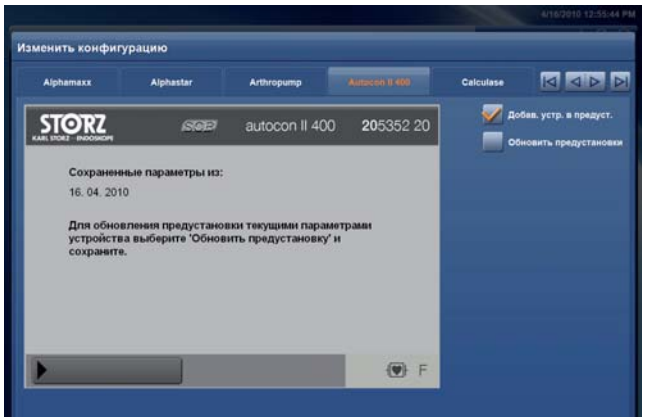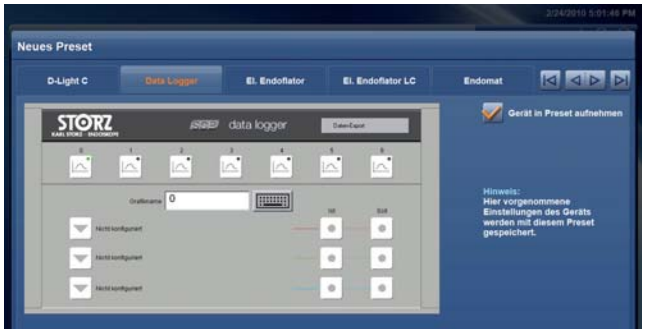
Special case: Devices which support data loggers

Once the data logger has been selected, the storable log data of the selected devices are displayed and can be selected.

If the selected device supports the data logger, the parameters to be logged are shown on the right-hand side. Each individual value can be freely selected for logging. Further information can be found in the corresponding section “Data logger”.

In the “**data logging**” field, you can choose which values are to be logged by pressing the corresponding button.

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
NOTE: *If no devices were assigned to the new Preset, an error message is displayed. In this case, the Preset cannot be saved.*

NOTE: *By pressing the button "Add" again the device is removed from the Preset. Previously saved presettings are hereby deleted.*

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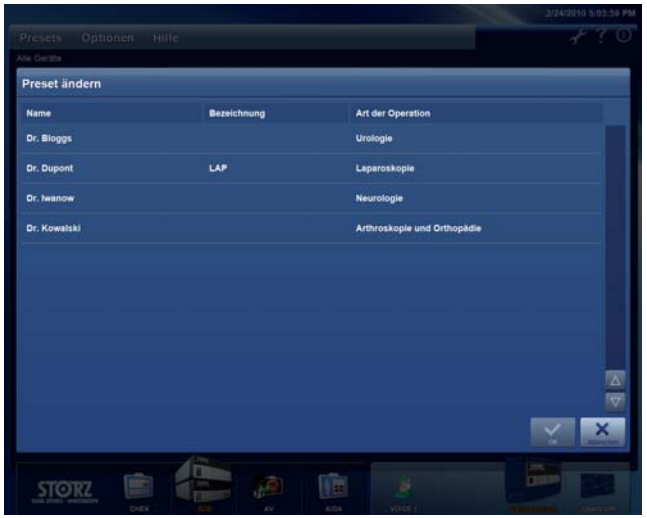
"Change…"

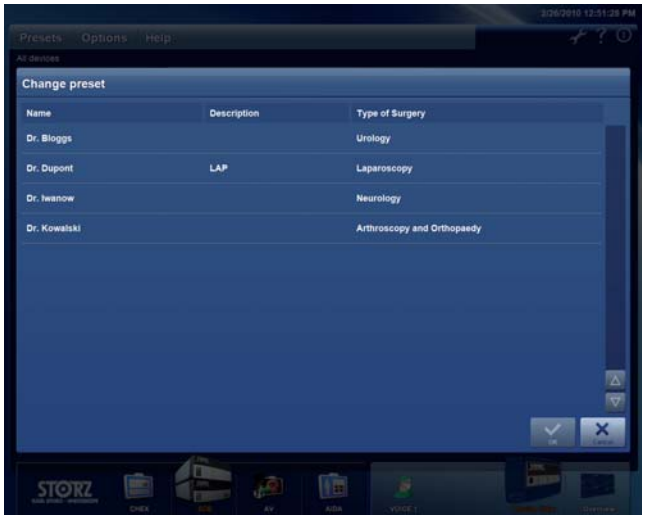The menu item "Change…" is used to change existing Presets.


The system shows you a three-column list with name, optional description and type of surgery of the existing Presets.

Select the desired Preset.

Confirm your selection with the "OK" button. The process can be stopped at any time by pressing the "Cancel" button.

The selected Preset is called up for changing as described in the following section.

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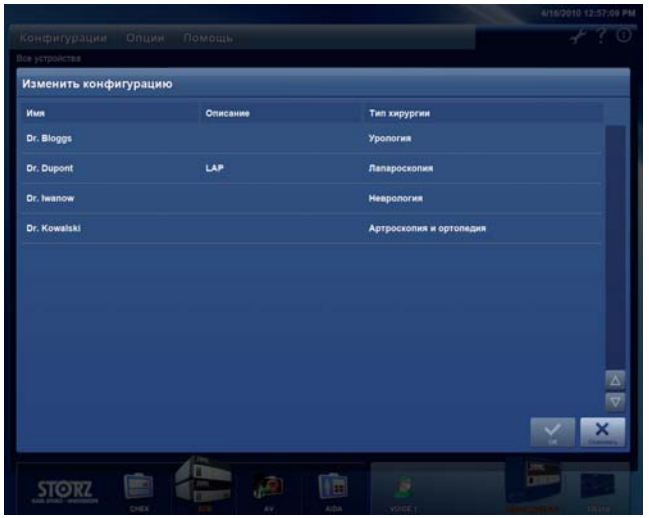Select the devices on the tab.

For devices whose options/settings can be set directly in the Preset, the image of the device is visible in the “Change Preset” entry screen. Change the values accordingly as on the device.

For devices whose options/settings are taken directly from the devices, a symbolic representation of the device is shown. When selecting the device the button “Update Preset” is displayed on the right-hand side and on the left-hand side the corresponding information for further action.

Start the synchronization of the required presettings on the device in the Preset according to the instructions displayed in the dialog box.

NOTE: *A highlighted “Update Preset” button indicates that data from this device will be saved in the Preset upon pressing “Save”.*

The process can be stopped at any time by pressing the “Cancel” button.

Pressing “Save” saves and starts the whole Preset. The presettings you have set using the devices are then transferred to the devices. At the same time the device settings for the devices to be synchronized are added to the Preset.

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
“Synchronization Error” dialog box

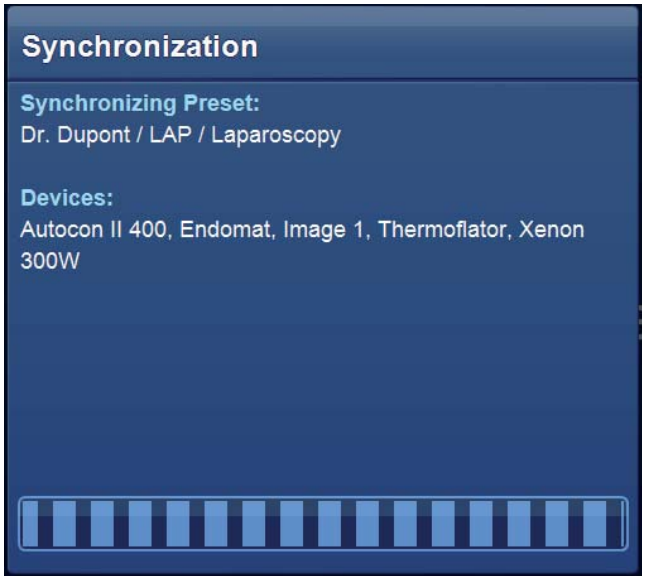
The synchronization dialog box provides information on the status of device initialization with the selected Preset values. It is displayed as soon as a Preset is opened, created or changed.

To start with, all devices in the Preset are listed. The initialized devices are then removed from the list. If all devices have been synchronized, the dialog box disappears once the last device has been removed from the list.

If all devices have not been removed from the list within 45 seconds, an error dialog box is emitted indicating the device which could not be synchronized.

The other information must be noted for troubleshooting.

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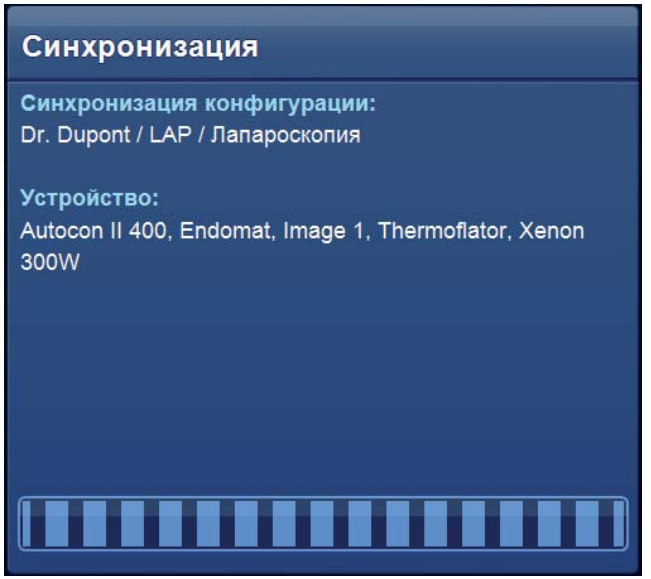
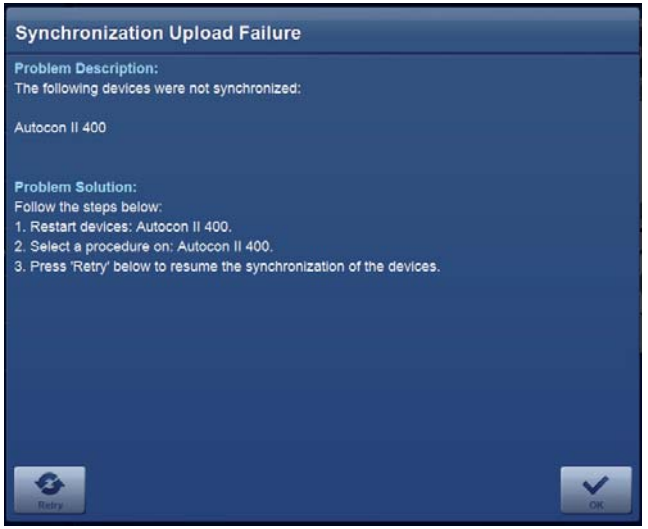
The user can restart the entire process using the “Retry” button or confirm the dialog box and manually set the necessary Preset values on the non-initialized device.

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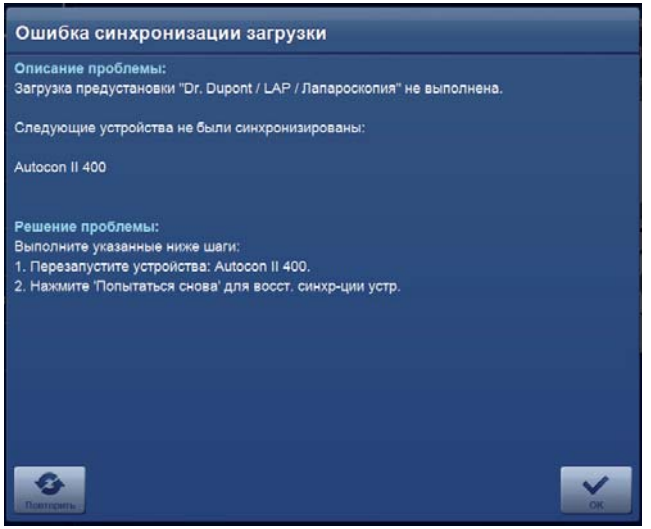4.4.5.5 “Save as…”

The menu item “Save as…” is used to save the device settings under the current or a new Preset name. The menu item “Save as” is only active when a Preset has been opened and a realistic device view is visible.

The active device in the device list (orange) assumes the visible device settings. For example, THERMOFLATOR®: If the THERMOFLATOR® is in intermittent mode, only these values are assumed even if the values in normal mode have been reset. All other existing configuration settings remain unchanged.

The synchronization dialog provides information on the progress of the saving procedure for the respective device. It is displayed as soon as a Preset is saved.

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

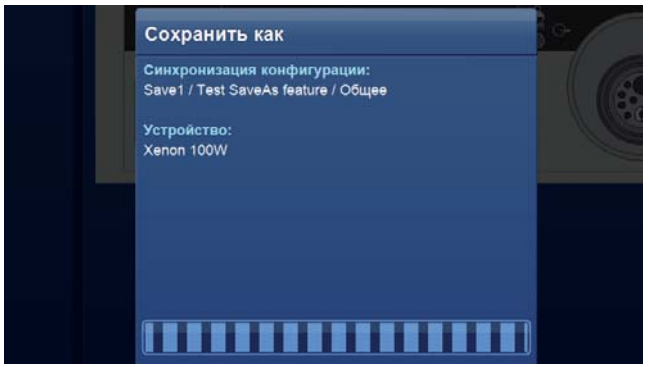
If the saving procedure is unsuccessful, a failure dialog is displayed.

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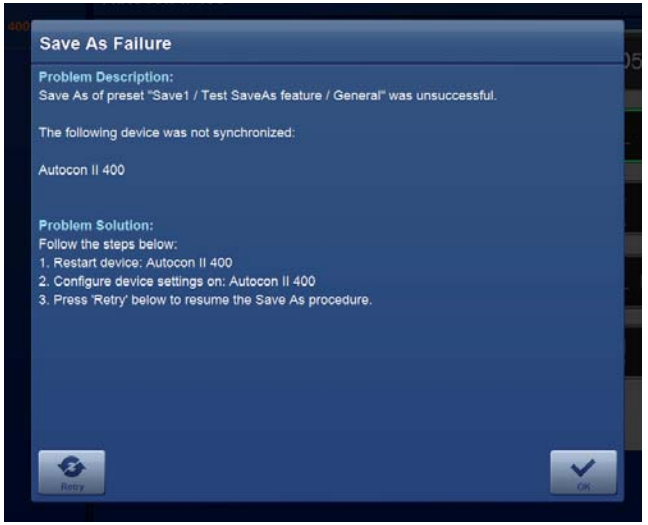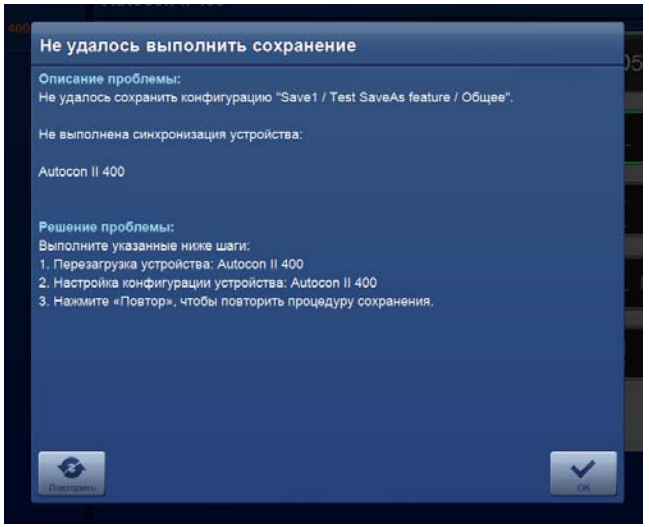
“Save as…” current Preset

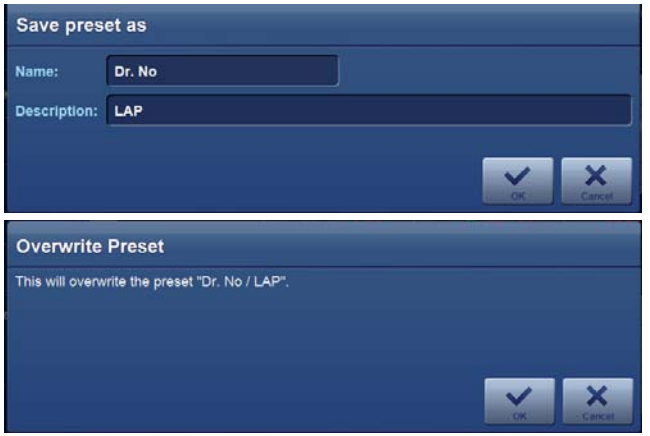If you wish to save the new device settings under the current Preset you must not change the Preset name (“Name”, “Description” and “Type of surgery”). Press the “OK” button to save the currently visible device settings under the current Preset. With the “Cancel” button, the procedure is discontinued without accepting any changes.

“Save as…” new Preset

If you wish to save the new device settings under a new Preset name you must change the Preset name (“Name”, “Description” and “Type of surgery”). With the “OK” button the currently visible device settings are saved under the new Preset name which you have entered (“Name, Description and Type of surgery”). With the “Cancel” button the procedure is discontinued without accepting any changes.

NOTE: *Clicking on the text “Name”, “Description”, or “Type of surgery” in the “Save Preset as” window will open the on-screen keyboard automatically..*

NOTE: *Since the menu item “Save as...” is used to create new Presets, it can be executed even if the device displayed does not support the function as such. Devices whose realistic views are displayed in the Preset are supported.*

NOTE: *As not all devices support this function, the menu item may also be inactive.*

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
4.4.5.6 Notepad


The user can save notes on the created Presets on an electronic notepad. These are saved together with the Preset. As such, this menu item can only be selected with an open Preset. Notes can also be entered when creating a Preset.

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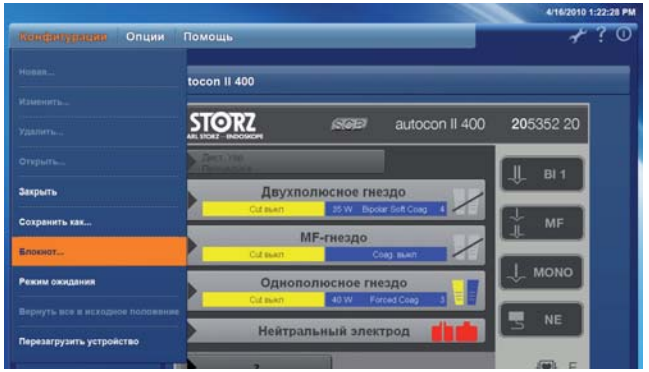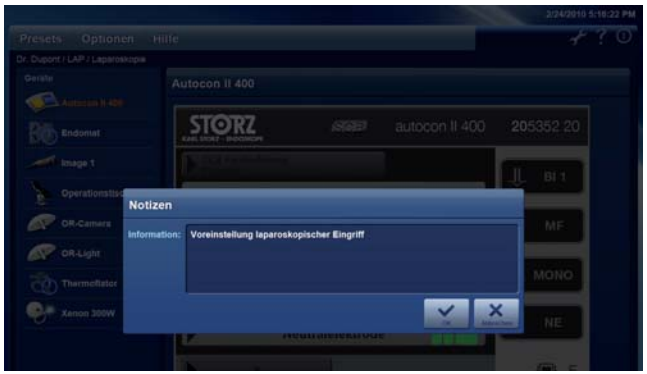
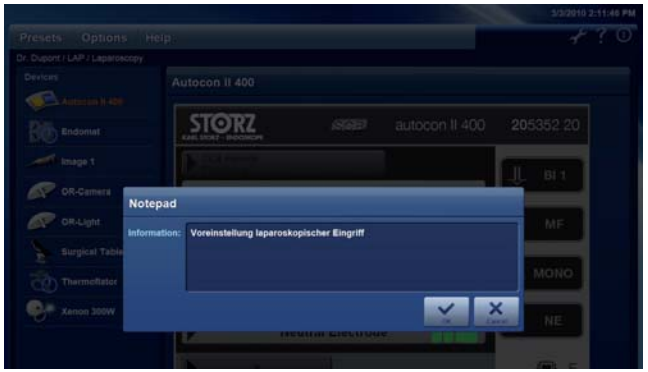
Once a Preset has been opened, the notes, assuming there are any, are shown automatically once and can also be expanded upon straight away.

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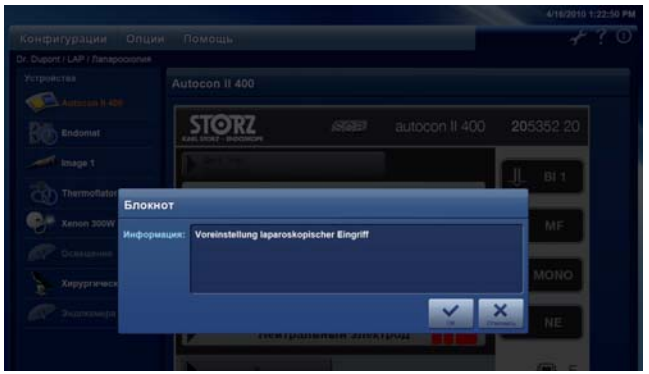
4.4.5.7 Help

4.4.5.7.1 Quick guide

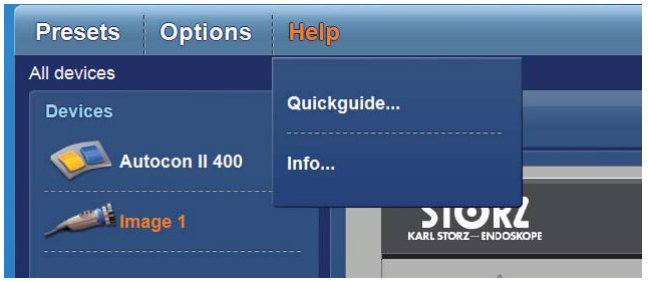The “Quick guide…” menu item provides you with a brief description of the KARL STORZ OR1 NEO™ system. This quick guide can also be called up by clicking the button in the general functions (see p. 27).


NOTE: *The quick guide appears in the selected language (if available), otherwise it is in “English”.*

You can adapt the page view to fit your screen using the “Page” and “Page width” buttons.

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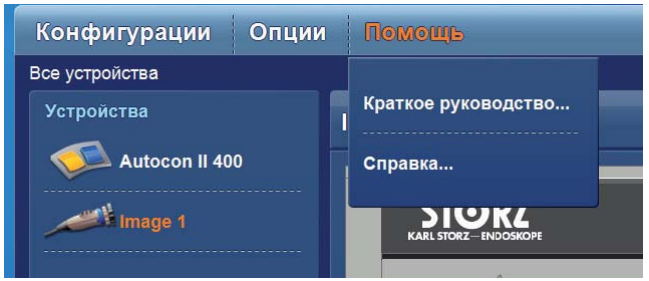Browsing in the quick guide

You can browse forward/backward in the quick guide using the “touch function”. You can also browse in the document by “dragging” the bar.

It is also possible to turn the pages with the inner arrow buttons “◁” and “▷”. The two outer arrow buttons “|◁” and “▷|” allow you to jump to the beginning or the end of the document. The arrow buttons are shown in “gray” when the beginning or end of the document has been reached.

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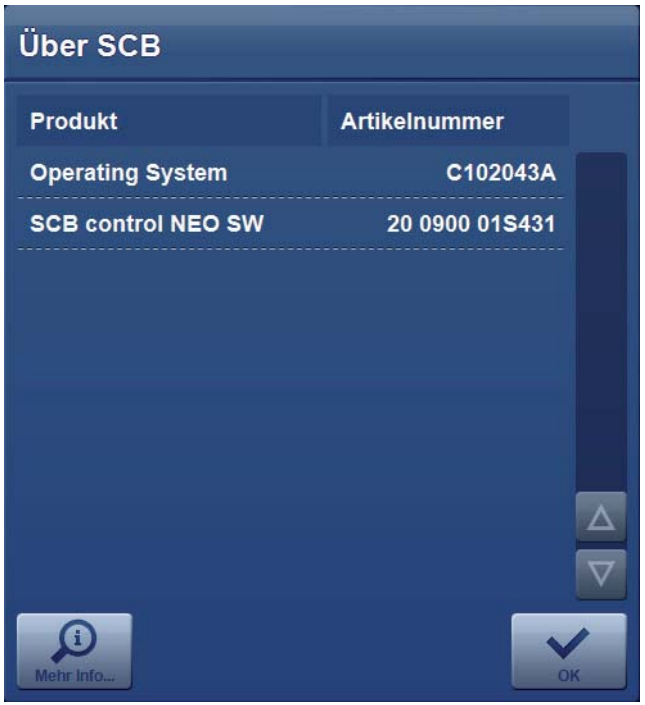
4.4.5.7.2 Info... (More Info)


Clicking on “Info...” shows information about the current software version.

“More Info” provides more detailed information on the individual system components.

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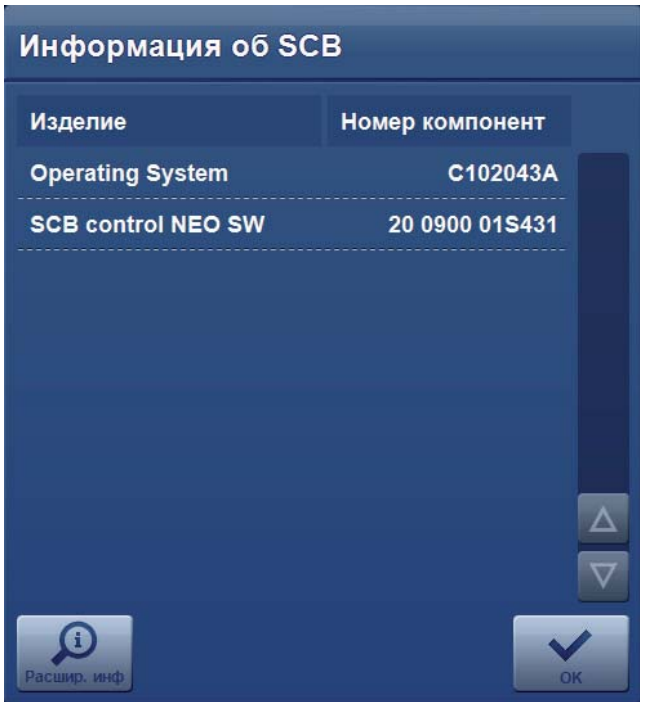
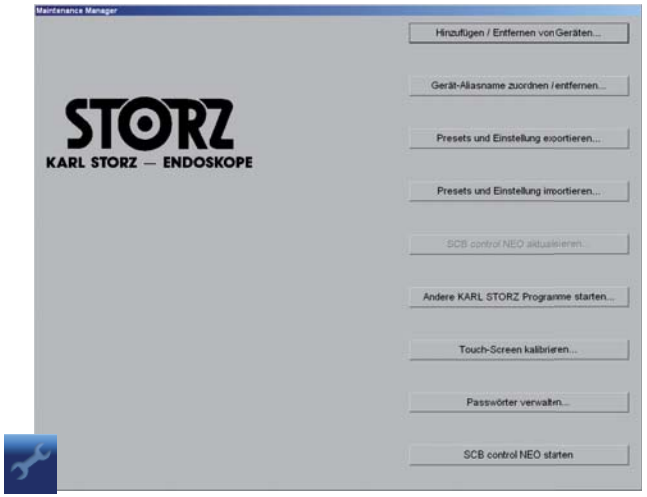
4.4.5.8 “Maintenance Manager”

Select the menu item “Options – Maintenance Manager…” or the button in the general functions to start the “Maintenance Manager”.

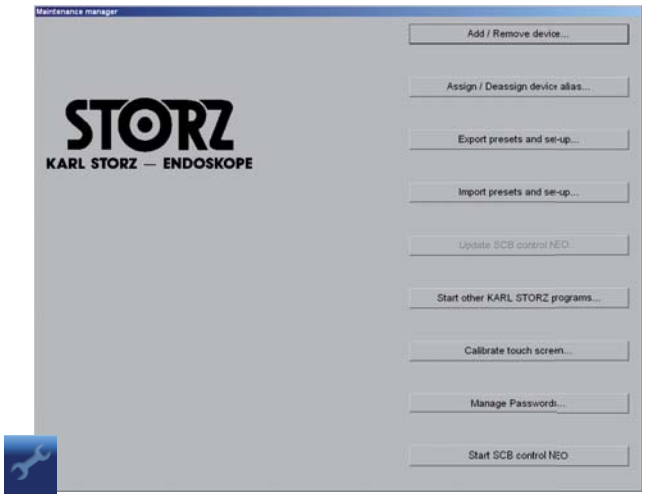Upon starting the “Maintenance Manager” the following menu items are shown:

- Add/remove device
- Assign/Deassign device alias
- Export Presets and set-up
- Import Presets and set-up
- Update SCB control NEO
- Start other KARL STORZ programs
- Calibrate touch screen
- Manage Passwords
- Start SCB control NEO

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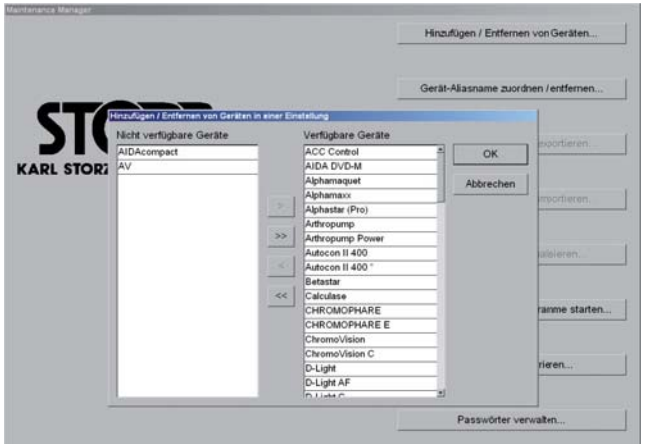
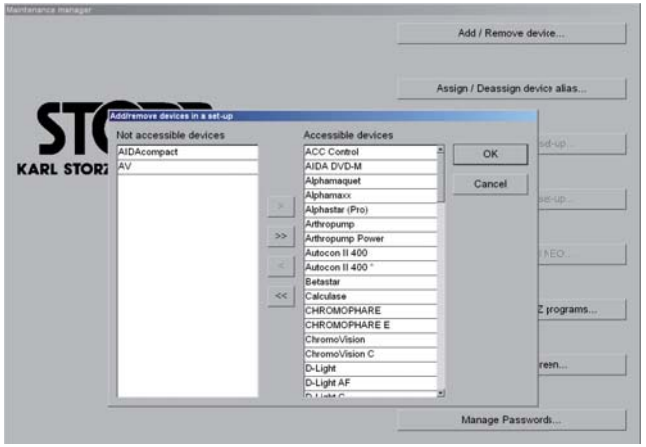
4.4.5.8.1 Add/remove an SCB device to/from the set-up

On selecting the menu item “Add/remove device” two lists appear on the screen:

- A list of “Not available devices” with SCB devices which are not accessible in Presets (left column)
- A list of “Available devices” with accessible SCB devices (right column)
- 4 Buttons: “**>>**” and “**<<**” (all devices), “**>**” and “**<**” (individual device)


To add/remove an/all SCB device(s), click on the relevant button.

Example: Add SCB device

Select SCB device: e. g., HYSTEROFLATOR®.
Clicking on the button “**>**” adds the HYSTEROFLATOR® to the “accessible” devices. Proceed in the same way to remove an SCB device.

- Click on “OK” to carry out the change.
- Click on “Cancel” to reject the change.

NOTE: *The AIDA compact or AV devices cannot be displayed together with the OR1™ FUSION device in the list of accessible devices.*

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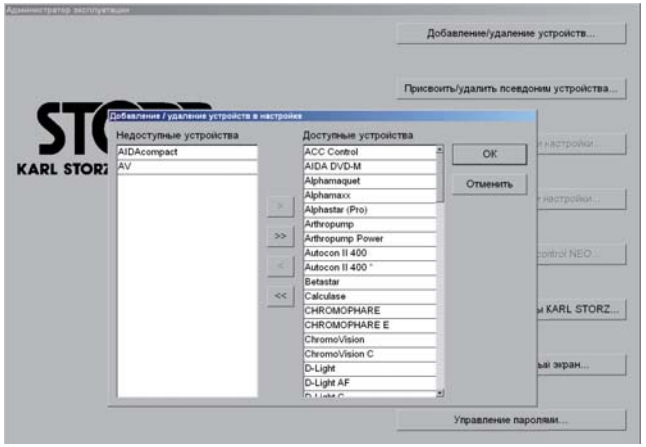
4.4.5.8.2 Assign/Deassign SCB device alias

This function can be used to assign an SCB device a general descriptive name, e.g., you could change the name THERMOFLATOR® to Insufflator. It is only possible to rename certain SCB devices which are available in the list. The names available for renaming are predefined by the system. A separate list of aliases is available for each of the supported system languages. Renaming an SCB device in one language has no influence on the other languages. Once a name has been assigned it is no longer available for other SCB devices as two or more SCB devices with the same name would lead to confusion.

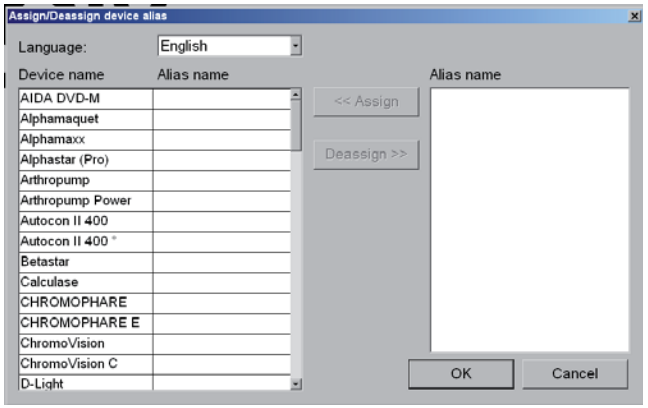On selecting the menu item “Assign/Deassign device alias” three lists appear on the screen:

- A selection list for the language
 All the aliases assigned to SCB devices only apply to the language selected here.
- A list of available SCB devices to which an alias can be assigned (left column).
- A list of aliases which can be assigned to the selected SCB device (right column).
- 2 Buttons: “**<<** Assign” and “Deassign **>>**”

To assign/remove an alias click on the relevant button.

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Example: Rename SCB device

- Select SCB device: e. g., “Printer UP51”
- Select alias: e. g., “Printer”
- Clicking on the button “Assign” moves the alias name “Printer” to the selected device in the column “Aliases”. This name is then no longer available for this SCB device or other SCB devices of the same type.

Click on “OK” to carry out the change.

Click on “Cancel” to reject the change.

The alias assigned here is then used instead of the original device name in all sections of the SCB control NEO software program.

NOTE: *The video overlay of the IMAGE 1 HUB™ also uses the alias.*

NOTE: *With IMAGE 1 HUB™ the camera control unit (CCU) must be equipped with software version 470206-0106, the S1/S3 camera head with version 470200-0104 or higher so that errors, alarms and warnings can be displayed with the relevant distinctive numbers and any device aliases.*

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4.4.5.8.3 Export Presets and settings

Presets can only be exported when a USB memory stick is plugged in. The menu item “Export Presets and settings” can then be selected.

- If an exported Preset file is already on the USB memory stick, the user must confirm that this file may be overwritten.
- The Preset file is exported to the USB memory stick.

NOTE: *The Presets and settings should be exported regularly as a precaution.*

4.4.5.8.4 Import Presets and settings

Presets can only be imported when a USB memory stick is plugged in. The menu item “Import Presets and settings” can then be selected.

- It is verified that all selected SCB devices are available.
- It is verified that the Preset name is unique. If the Preset name already exists, #1 will be added (Preset name#1).
- If any test fails, a Preset will not be added.

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4.4.5.8.5 Update SCB control NEO (software update)

The menu item “Update SCB control NEO” can only be selected if a KARL STORZ update CD is in the DVD-ROM drive.

- It is verified that the SCB control NEO installation program is available.
- The SCB control NEO installation program starts.
- The “Maintenance Manager” closes automatically.

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
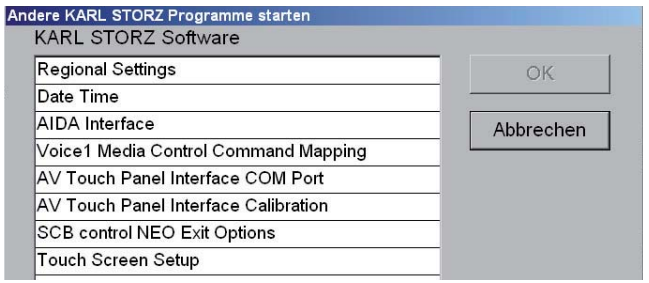
4.4.5.8.6 Start other KARL STORZ programs

After selecting the menu item “Start other KARL STORZ programs”:

- If there is an SCB program CD in the DVD-ROM drive, the programs on the program CD are listed. If necessary, select the desired program and start by pressing “OK”.

- If there is no SCB program CD in the DVD-ROM drive, the message “No CD found” appears. Press “OK” to start the local programs. When you press the “Cancel” button, the “Maintenance Manager” appears again. After confirming with “OK” the following selection list appears.

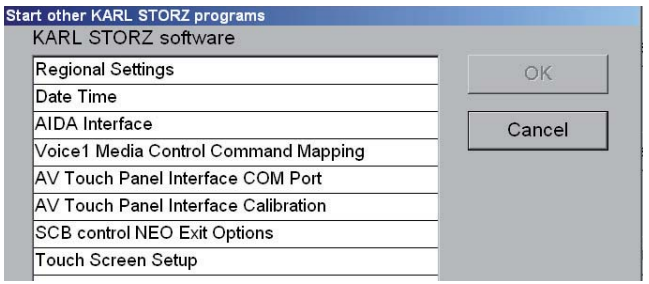Regional Settings

- Select “Regional Settings” to adapt the country setting of the operating system of your keyboard if necessary (Regional Settings – Input).

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Date Time

- Select “Date Time” to set the date and time if necessary.

AIDA Interface

- Select “AIDA Interface” to check the changeover function of the AIDA interface changeover device.

AV Touch Panel Interface ComPort

- Select “AV Touch Panel Interface ComPort” to select the ComPort for AV panel connection.

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SCB control NEO Exit Options

- Select “SCB control NEO Exit Options” to select between the following options:
 – Shutdown: The OR1™ control NEO is shutdown.
 – Restart: The OR1™ control NEO is restarted.

Touch screen setup

- Select the support of a second touch screen (“Two SCB Touch Screens”) if a second monitor is connected (see installation manual, section 1.2 “Controls, displays, connections and their uses”, No. ㉒ “Serial interface COM 2 for touch screen function of the second touch screen”).

ⓘ ***NOTE:*** *A second touch screen can only be used in connection with AIDA Interface II 20 0921 32.*

- Select the support of only one touch screen (“One SCB Touch Screen”) if no second monitor is connected.

ⓘ ***NOTE:*** *The user guidance for the programs Regional Settings, Date Time, AIDA Interface, AV Touch Panel Interface ComPort, AV Touch Panel Interface Calibration, SCB control NEO Exit Options and Touch Screen Setup is only in English.*

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
4.4.5.8.7 Calibrate touch screen

The touch screens are to be calibrated prior to use. The number of used touch screens and touch screens to be calibrated must be noted under the above point “Touch Screens Setup”.

4.4.5.8.8 Start SCB control NEO software

After this menu item has been selected, the SCB control NEO software is started.

4.4.5.8.9 Manage passwords

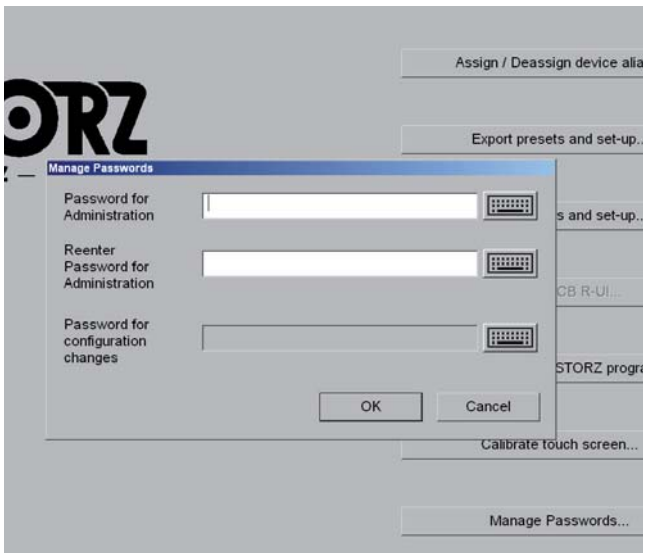This dialog box is used to enter passwords in order to protect access to certain critical functions should this be desired by the user. There are 2 different protected access levels:

- Access to password management and thus the Maintenance Manager
- Access for changing Presets and general data (new, change, delete, identical device list and notepad)

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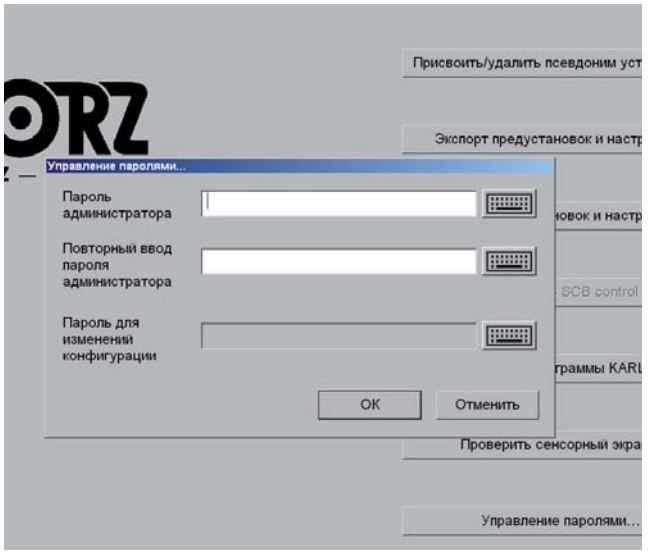
The password for the administrator must be entered twice to rule out the possibility of typing errors.

Once an administrator password has been entered, a password for changing Presets and data can be entered. If the dialog box is exited by pressing OK a warning appears reminding the user to note the passwords entered.

NOTE: *Please note the passwords you have entered and keep these in a safe place. Should you lose the passwords which you entered, you will no longer have access to the protected functions. It will then only be possible for KARL STORZ staff to access these protected functions.*

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4.4.5.9 Device setup

Device setup – general


Some SCB devices allow the device setup parameters to be set. These are mostly values which do not have to be set under normal operating conditons but rather affect the device's setup. The possible device-specific settings are explained in Section 4.5 Devices.

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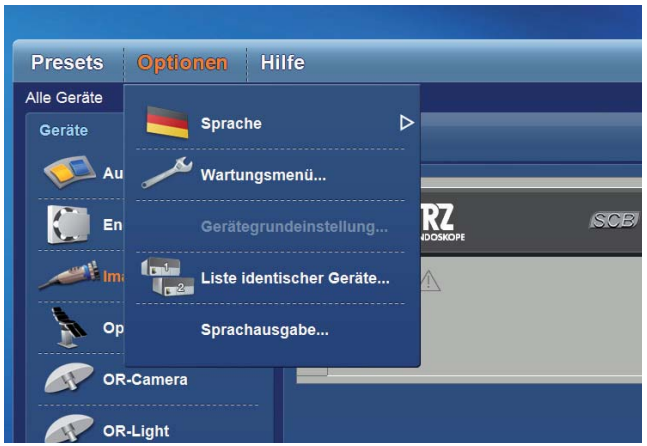
4.4.5.10 Identical devices on SCB

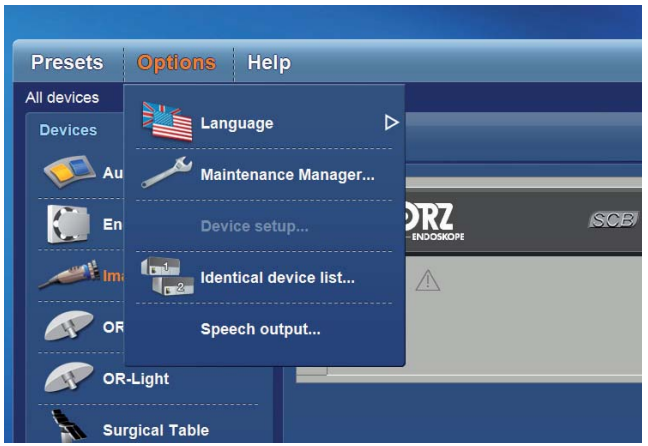
The SCB control NEO software offers the option of operating several identical devices of a device type in the SCB network and controlling these centrally via the SCB control NEO software. To be able to distinguish between the devices, they are each assigned an individual distinctive number.

A maximum of 3 devices each can be controlled from the following devices via the SCB control NEO software.

AUTOCON® II 400, COVIDIEN ForceTriad™, Electronic ENDOFLATOR®, IMAGE1 HUB™, PowerLED 175, THERMOFLATOR®, TRICAM® SL II, XENON 175 and XENON 300

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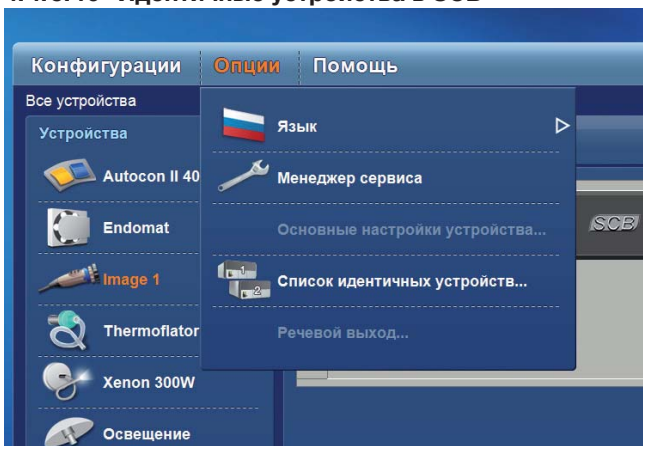AUTOCON® II 400, COVIDIEN ForceTriad™, Electronic ENDOFLATOR®, IMAGE1 HUB™, PowerLED 175, THERMOFLATOR®, TRICAM® SL II, XENON 175 и XENON 300

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Assignment of “identical” SCB devices

Requirements: In order to operate two or more identical SCB devices via the SCB control NEO software, these SCB devices must be connected with the SCB control NEO system and switched on.

This window shows all the “Identical devices” currently linked to the SCB control NEO system and switched on, together with device names and serial numbers (S/N).

It serves to assign the identical SCB devices a unique distinct number.

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Click on a device in the identical device list.
The selected line is shown in “orange”. The selected device gives the following signal:

IMAGE 1 HUB™:
SCB indicator on video overlay flashes green.

XENON 175:
Xenon “waves” by switching to and from the standby mode.*

AUTOCON® II 400:
SCB indicator on the device display flashes the “#” symbol.


PowerLED 175:
SCB indicator display shows the device’s distinctive number.

*If there is no SCB indicator on the front panel of the device, the device switches periodically to and fro between “Standby” and “Not on Standby”.

COVIDIEN ForceTriad™, Electronic ENDOFLATOR®, THERMOFLATOR® and XENON 175/300:
These SCB devices do not have an SCB indicator. Please attach the corresponding sticker (“#1”, “#2” or “#3”) next to the SCB symbol on the front plate after you have assigned the relevant number.

You will find the serial number of the SCB device on the rating plate.

Choose a distinctive number for the selected SCB device.

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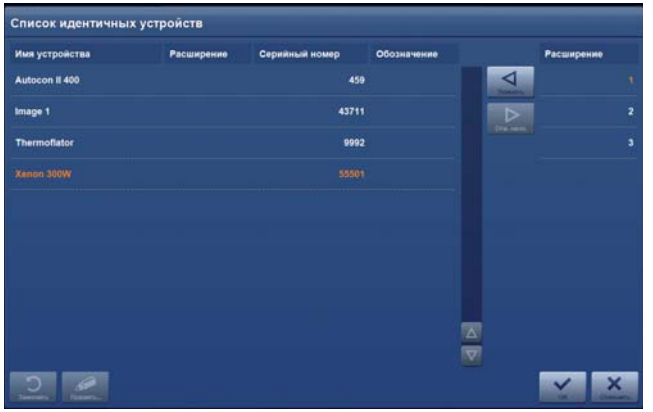NOTE: *The distinctive numbers "-1", "-2" and "-3" can only be assigned in ascending order.*

Example: The distinctive number "-3" can only be assigned to an SCB device if the numbers "-1", and "-2" have already been assigned to other SCB devices. This means if two identical SCB devices are connected to the KARL STORZ SCB control NEO system, only the number "-2" can be assigned as the highest.

Press the button "◁ Assign" to assign the desired distinctive number to the selected SCB device. A "-1", "-2" or"-3" is added to the device name.

Deassigning a distinctive number

Select the SCB device and then press the button "▷ Remove".

NOTE: *The numbers can only be desassigned in descending order.*

Example: The number "-1" assigned to a SCB device can only be deassigned if the numbers "-2", and "-3" are not assigned. If necessary, these must be deassigned first.

COVIDIEN ForceTriad™, Electronic ENDOFLATOR®, THERMOFLATOR® and XENON 175/300:
Also remove the respective sticker ("#1", "#2" or "#3") from the front panel of the device.

NOTE: *The device names of several identical SCB devices are not visible until assignment of a number in the SCB control NEO device list and in the Presets.*

In the list on the left-hand side all SCB devices with device names and S/N are displayed which can be used in parallel as identical SCB devices.

Once assigned, a distinctive number applies to the entire SCB control NEO as well as to all configurations.

NOTE: *If a distinctive number is deassigned from a SCB device, it is also removed from all Presets in which it was used. A message appears and the command is executed following confirmation by the user.*

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Replacing identical SCB devices

It is possible to replace failed SCB devices. If you wish to use the same distinctive number for the new SCB device as for the old SCB device, proceed as follows:

NOTE: *SCB devices can only be replaced if an additional SCB device of the same type has been connected to the SCB control NEO system and no distinctive number has yet been assigned to this in the SCB control NEO software.*

- Select the SCB device you wish to replace in the "Identical device list"
- Press the "Replace" button

The dialog box shown below is displayed:

A list shows the SCB devices connected to the KARL STORZ SCB control NEO system which are available to replace the device that has failed.

Click on a device line in the device list. The selected line is shown in "bold". The selected device makes itself noted as above:

Press the "OK" button to replace the SCB device with the selected SCB device.

Press the "Cancel" button to discontinue the procedure. The SCB device is not replaced.

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The types of messages can be set for which speech output is required to take place can be set under “Speech output”.

The “Speech output…” button switches the entire speech output on or off.

The “Alarms, Warnings, Errors” button allows speech output to be switched on or off for alarms, warnings and errors. If this function is deactivated, no alarms, warnings or errors are given for devices via the speech output. This does not affect the messages on the monitor.

The “Devices turned on/off” button is used to specify whether the turning on and off of devices via speech output is to be signalled e. g., “XENON 175 has been turned on”. This setting also only affects the speech output.

The volume of the speech output can be adjusted via the speech output volume field. The +/– buttons are used to increase or reduce the volume, next to the % display this is shown by a bar.

Press the “Cancel” button to reject all settings and close the window.

Press the “OK” button to accept all settings and close the window. These settings are accepted and loaded automatically the next time the program is started.

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

4.4.5.12 Log In/Log Out dialog box

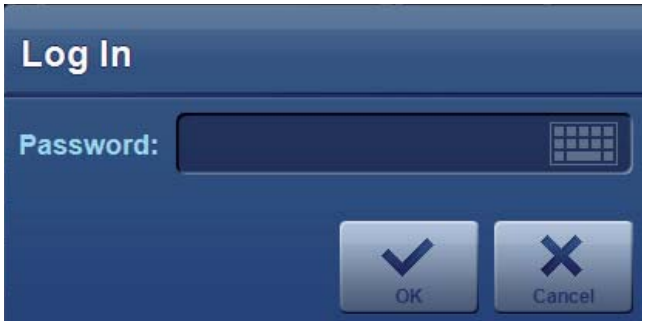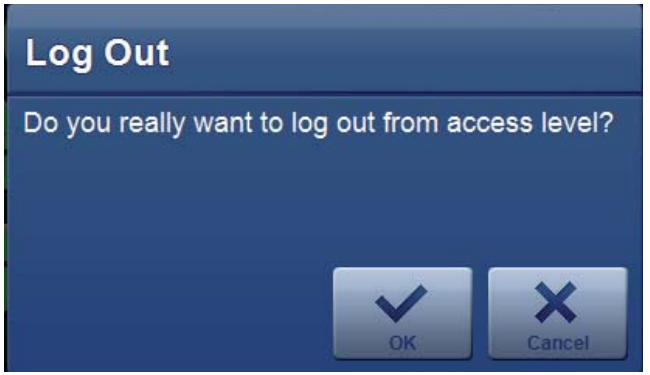The Log In/Log Out dialog box is used either to enter a password to activate the password-protected menu items (Log In) or to cancel activation (Log Out).

To log in, enter the password in the input field. When you press the OK button, registration is carried out and the access rights for this password are granted.

If the password entered is not correct, a warning is given. You can then correct the password.

If you have already logged in you can log out using the menu item Log In/Log Out. Press the OK button to log out and you will then only have access to unprotected functions.

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
4.4.5.13 IR code dialogs

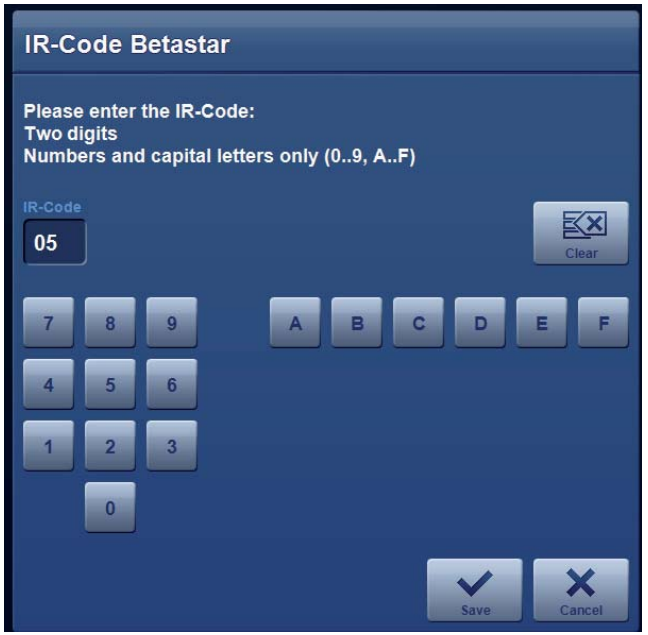The dialog shows the current context in the title: in this case “IR code” with the device name (Betastar).

In addition to a short explanation, the current IR code is shown in double digits (single-digit codes are prefixed with a zero). Depending on the manufacturer, the IR system code is either a double-digit or triple-digit hexadecimal figure.

The value is set with the numerical keypad. However, before you can enter new values, the value must be cleared with the “Clear” button. The “Save” button is then temporarily deactivated until at least one figure is entered.

Initial zeros are automatically added by the software if necessary when the dialog is exited by pressing the “Save” button.

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4.5 Devices

Complete table of all controllable devices

Art.-Nr./ Art. no./ № изделия	Gerät/Device/Устройство	SCB control NEO software Gerätename/device name/ имя устройства
20 7117 20-1	UNIDRIVE® NEURO	Unidrive NEURO
20 7118 20-1	UNIDRIVE® OMFS	Unidrive OMFS
Medien Geräte/Media devices/Медиаустройства		
20 0906 32	ACC Control®	ACC Control
20 0910 32	Media Control®	Media Control
20 0914 32	SCB OR1™ Telestration	Telestration
Dokumentation/Documentation/Документирование		
20 2045 20-1xx	AIDA DVD-M	AIDA DVD-M
Lasergeräte/Laser devices/Лазерные устройства		
27 7501 20-1	CALCULASE	Calculase
HF-Geräte/HF devices/ВЧ устройства		
20 5352 20-11	AUTOCON® II 400	Autocon II 400
20 0903 90	VALLEYLAB	Force FX
20 0912 90	KLS MARTIN maxium®	Maxium ME 402
20 0915 90	ConMed System 2450™	System 2450
20 0915 90	ConMed System 5000™	System 5000
20 0913 22	COVIDIEN ForceTriad™	COVIDIEN ForceTriad™
20 0913 32	ERBE VIO 300D	VIO
OP-Tische/OR tables/Операционные столы		
20 0908 32	MAQUET OR-Table	Alphamaquet
20 0908 32	MAQUET OR-Table	Alphamaxx
20 0908 32	MAQUET OR-Table	Alphastar (Pro)
20 0908 32	MAQUET OR-Table	Betastar
20 0928 32	MAQUET OR-Table II	Magnus
20 0907 32	Trumpf OR-Table	JUPITER
20 0907 32	Trumpf OR-Table	MARS
20 0907 32	Trumpf OR-Table	MARSEndouro
20 0907 32	Trumpf OR-Table	TruSystem
20 0919 90	BERCHTOLD OPERON®	OPERON B 710
20 0919 90	BERCHTOLD OPERON®	OPERON B 810
20 0916 90	OPT OR-Table	OPT OR-Table
20 0907 32	Trumpf OR-Table	SATURN Select
20 0905 90	STERIS® OR-Table 3085 SP	Steris OR-Table
20 0907 32	Trumpf OR-Table	TITAN
20 0925 22	Universal OR-Table	OR-Table

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4.5.1 Camera systems

TELECAM SL, TRICAM® SL and TRICAM® SL II

The main window of TELECAM SL, TRICAM® SL and TRICAM® SL II is basically the same as the view of the front panel of the device and contains a white balance button which activates the white balance function.

In addition, the following buttons are included for opening and closing associated windows:

- “Function setup” button
- “Display setup” button

NOTE: *If a new window is activated the previously activated window is closed automatically.*

Function setup

With the “Function setup” button you can open and close the window for the function settings.

NOTE: *The function setup button is only enabled if there is no PDD camera head connected to the camera.*

Exposure/shutter speed

- Range selection [s]: “Auto”, 1/60, 1/125, 1/250, 1/500, 1/1000, 1/2000, 1/4000, 1/10000
- By pressing the “Shutter Sp” button you switch to the next value. The current setting is always displayed to the right of the keypad.

Gain

- Range selection: “Low”, “Medium”, “High”, “Peak”

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- By pressing the “Gain” button you switch to the next value. The current setting is always displayed to the right of the keypad.

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- Activation of the color test image.

Camera head buttons

- This field emulates the blue and silver button on the camera head. Briefly pressing the button activates the primary camera head button function (“I:”), while pressing for longer (>2 s) activates the secondary camera head button function (“II:”).
- The current camera head button functions are displayed to the right of the keypad with “I:” and “II:”.

Programming

- Call up the video image menu to program the camera head buttons.

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
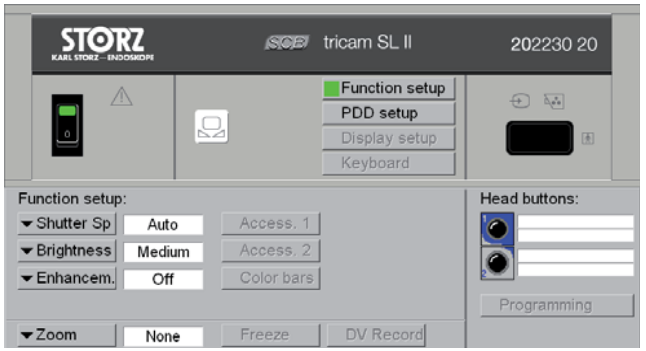
Zoom, still image and DV recording (TRICAM® SL II only)

The windows for function setup and PDD setup display the “Zoom”, “Freeze” (still image) and “DV Record” buttons.

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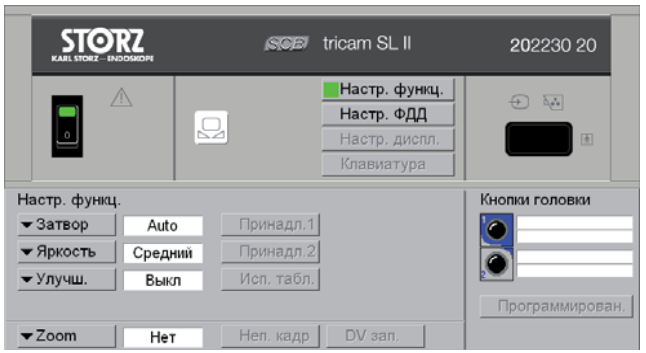Zoom

- Range selection: “none”, 1.25x, 1.5x, 1.75x, 2x
- By pressing the “Zoom” button you switch to the next value. The current setting is always displayed to the right of the keypad.

Freeze (still image)

- Activation/deactivation of the still image

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
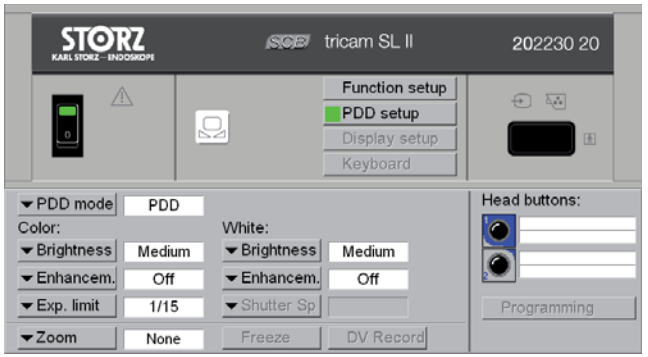
PDD setup

With the “PDD setup” button you can open and close the window for PDD settings. In the PDD setup window, additional set values are displayed for a camera with connected PDD camera head.

NOTE: *The PDD setup button is only enabled if there is a PDD camera head connected to the camera.*

PDD mode

- Range selection: “---”, “ALA”, “AF”
 Range selection with D-LIGHT or D-LIGHT AF:
 “---”, “ALA Uro”, “AF Pulmo”.
 Range selection with D-LIGHT P: “AF”, “CG”
 Range selection with D-LIGHT C and TRICAM® SL II:
 “PDD”, “AF”, “NEURO”, “HYP”
- By pressing the button you switch to the next value and the current setting is always displayed to the right of the keypad.

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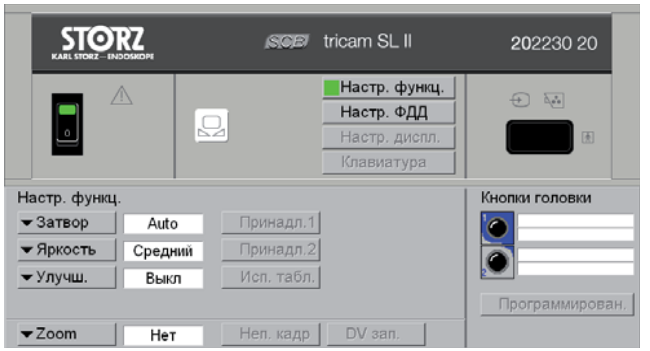The buttons and display fields under “Color:” are only activated when the camera is in the color light mode. The buttons and display fields under “White:” are only activated when the camera is in the white light mode. In addition to the “Brightness” and “Enhancem.” buttons there is the “Exp. limit” button in the color light mode.

Brightness

- Range selection: “Low”, “Medium”, “High”, “Peak” (TRICAM® SL II only)
- By pressing the “Brightness” button you switch to the next value. The current setting is always displayed to the right of the keypad.

Enhancement

- Range selection: “Off”, “Low”, “High”, “Filter A”, “Filter B”
- By pressing the “Enhancem.” button you switch to the next value. The current setting is always displayed to the right of the keypad.

Exposure limit – exposure time upper limit

Setting of the maximum exposure time in automatic mode

- Range selection [s]: 1/60, 1/30, 1/15, 1/8, 1/4, 1/2, 1, 2
- By pressing the “Exp. limit” button you switch to the next value. The current setting is always displayed to the right of the keypad.

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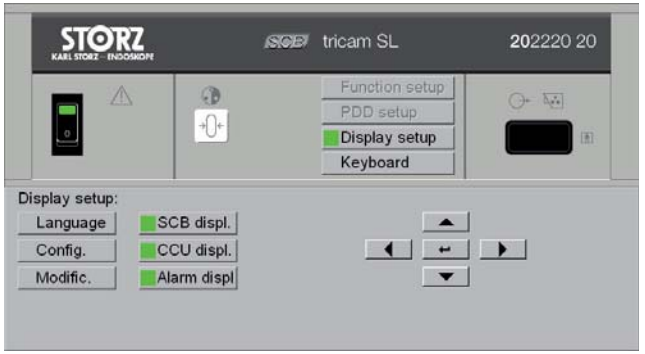
Display setup

With the “Display setup” button you can open and close the window for display settings.

Language*

- Call up the menu in the video image for selection of language.

Configuration*

- Call up the menu to configure the SCB overlay, any existing configuration will be deleted.

Modification*

- Call up the video image menu to modify the SCB overlay. The configuration already saved will be modified.

ⓘ ***NOTE:** The emulated cursor buttons are used to navigate through the menu called up in the video image. Alternatively, operation can also be performed with the video hardware keyboard. To close/exit the respective menu press the RETURN button.*

SCB display

- Activation (on/off) of the SCB overlay (button combination CTRL-ALT-F6).

CCU display


- Activation (on/off) of the CCU overlay (button combination ALT-F6).

Alarm display

- Activation/deactivation (on/off) of the SCB alarm display on the CCU monitor.

Cursor key block

- The key block emulates the cursor buttons and the Return button of the CCU’s hardware keyboard.

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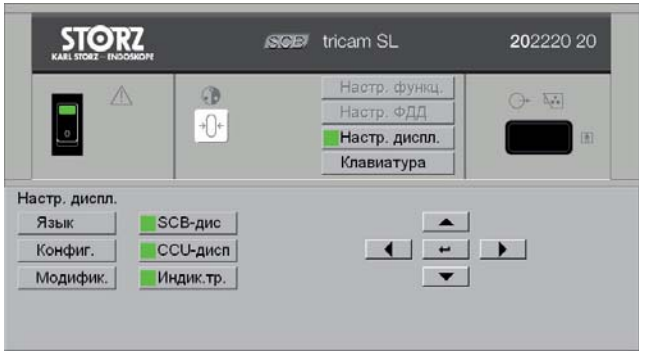Brightness

- Range selection: “Low”, “Medium”, “High”, “Peak”
- By pressing the “Brightness” button you switch to the next value. The current setting is always displayed to the right of the keypad.

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- Activation of the color test image.

Camera head buttons

- This field emulates the menu button (middle) as well as the upper (left) and lower (right) buttons on the camera head. Pressing the button activates the corresponding camera head button function.
- The current camera head button functions are displayed to the right of the keypad.

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- By pressing the “Insert” button you switch to the next value. The current setting is always displayed to the right of the keypad.

Swap picture

- Selection: “IMAGE 1”, “EXTERN”, “IMG1+ext.”, “EXT.+img1”
- Meaning:
 IMAGE 1: The camera image is displayed.
 EXTERN: The external image is displayed.
 IMG1+ext: The camera image is displayed and the external source is shown at the bottom right.
 EXT.+img1: The external image is displayed and the camera image is shown at the bottom right.
- Press the “Swap Pict.” button to switch to the next choice. The current setting is always displayed to the right of the keypad.

External source

- Selection: “SDI”, “S-Video”, “Composite”, “Auto”
- By pressing the “Source” button you switch to the next value. The current setting is always displayed to the right of the keypad.

NOTE: *The “Picture in Picture” button is only enabled if the IMAGE 1 HUB™ supports the “Picture in Picture” functionality.*

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
IMAGE 1 H3-M

The main window of IMAGE 1 H3-M is basically the same as the view of the front panel of the device.

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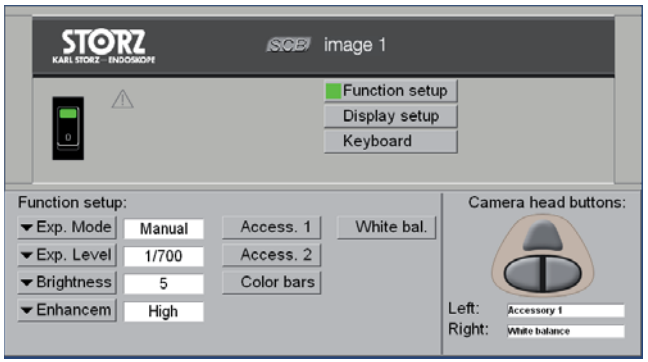Exposure mode

- Range selection: “Auto”, “Auto+Peak”, “Auto+Small”, “Auto+Med.”, “Manual”
- By pressing the “Exp. Mode” button you switch to the next value. The current setting is always displayed to the right of the key.

Exposure level

- Range selection: Automatic: “Min”, 1, 2, ..., 15, “Max” Manual: “1/60”, “1/60 +2”, “1/60 +4”, “1/60 +8”, “1/60 +Max”, “1/17000”, …, “1/85”
- By pressing the “Exp. Level” button you switch to the next value. The current setting is always displayed to the right of the key.

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NOTE: *The exposure level depends on the set frequency of the camera head. At a set frequency of 50 Hz, "1/50" is shown, at a set frequency of 60 Hz, "1/60" is shown.*

Brightness

- Range selection: "1" … "5"
- By pressing the "Brightness" button you switch to the next value. The current setting is always displayed to the right of the key.

Enhancement

- Range selection: "Off", "Low", "1" … "5", "High"
- By pressing the "Enhancem." button you switch to the next value. The current setting is always displayed to the right of the keypad.

Accessories 1, 2


- With these buttons, the accessory functions are activated according to the wiring of the accessory outputs.

Color bars (test image)

- Activation of the color test image.

Camera head buttons

- This field emulates the menu button (above) as well as the left and right buttons on the camera head. Pressing the button activates the corresponding camera head button function.
- The current functions of the left and right camera head buttons are shown below the keypad.

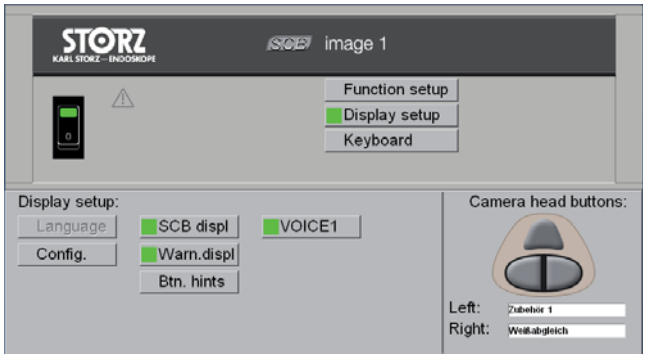Display setup

With the "Display setup" button you can open and close the window for display settings.

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Button information

- Activation/deactivation (on/off) of the button information on the CCU monitor.

VOICE1® (IMAGE 1 HUB™ only)

- Activation/deactivation of the speech tree overlay on the IMAGE 1 H3-M.

NOTE: *This button will be supported in a future SCB control NEO software version.*

Camera head buttons (see page 94)

Widget (compact menu overview)

The IMAGE 1 H3-M can be operated via both the device view and the compact view in the overview.

NOTE: *For information on the remote control of the IMAGE 1 H3-M via widgets, please see Section 4.4.4.2.*

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Device messages

Message	Priority
Temperature is too high	low
*Lamp defective. Lamp must be replaced	low

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NOTE: *For safety reasons, the light intensity can only be preset in the range 0% – 100% (indicated in green) within the Preset dialog box. The “Boost” and “General purpose” ranges are not accessible during presetting in the Preset.*

NOTE: *If the “Reset Device” or “Reset All” menu items are selected while the light source is in “Standard” mode, the pump is switched off and the light intensity and pump level are set to the preconfigured values.*

Device messages

Message	Priority
Temperature is too high	low

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4.5.3 Insufflators

Electronic ENDOFLATOR® *
HYSTEROFLATOR® **
THERMOFLATOR® ***
Read the instruction manual for the respective device carefully.
All the settings must be made as described in the instruction manual and the displays are to be interpreted accordingly.

Device messages

Message	Priority
*Intraabdominal pressure exceeds setpoint	low
**Intrauterine pressure exceeds setpoint	low
***Insufflation pressure exceeds setpoint	low

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4.5.4 Pumps

ARTHROPUMP®
ENDOMAT®
ENDOMAT® LC / ENDOMAT® LC (RD)
EQUIMAT®
Read the instruction manual for the respective device carefully.
All the settings must be made as described in the instruction manual and the displays are to be interpreted accordingly.

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HYSTEROMAT® II

Read the instruction manual for the respective device carefully.

All the settings must be made as described in the instruction manual and the displays are to be interpreted accordingly. The following are exceptions:

A button has been added to the realistic view which is not present on the actual device. Press the button for >3 s and the optional external pressure sensor is set to “zero”. With the HYSTEROMAT® II 263315 20-1 this function can be activated via two-button operation.

Device messages

Message	Priority
External pressure sensor defective	low
Volume difference exceeds setpoint	low

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ARTHROPUMP® POWER

Read the instruction manual for the respective device carefully. All the settings must be made as described in the instruction manual and the displays are to be interpreted accordingly. The following are exceptions:

NOTE: *To ensure safe operation with the ARTHROPUMP® POWER touch panel, the user must choose a so-called procedure at startup. If no procedure has been selected on the device, the device will report to the system with an error message, but cannot be operated.*

Choose a procedure on the device in order to operate the device via the SCB control NEO software.

The following dialogs are supported by the SCB control NEO software:

- main screen
- fluid supply monitor

No other dialogs or functions are supported.

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
NOTE: *An ARTHROPUMP® POWER which is connected to the OR1 control NEO computer is synchronized with the following SCB system settings.*

- Language
- Date
- Time

"Widgets" (compact menu overview)

The ARTHROPUMP® POWER can be operated via both the device view and the compact view in the overview.

NOTE: *For information on the remote control of the ARTHROPUMP® POWER via widgets, please see Section 4.4.4.2.*

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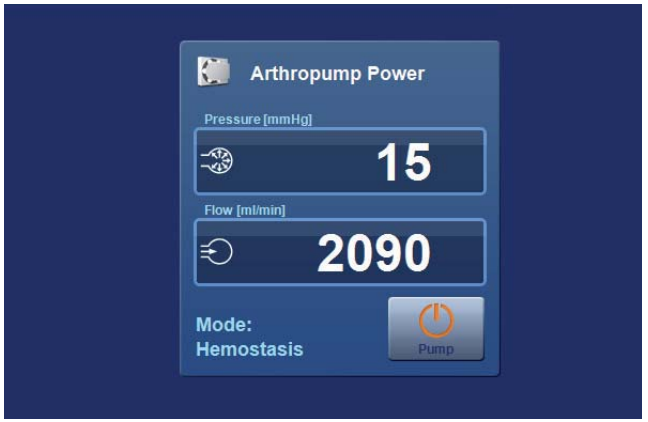- Pressure [mmHg]
- Flow [ml/min]
- Mode (hemostasis or lavage).

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The KARL STORZ POWERSHAVER system S1 (consisting of Power Unit S1, shaver handpiece and footswitch) supports the configuration of speed for all combinations of handpieces and blade classes for all four speed memories.

When creating a configuration, all the settings are assigned the default settings given by the mean value of the speed range. The set speed range is transferred automatically to a buffer when another handpiece, blade class or memory speed is selected. This is adopted on confirmation of the configuration. Otherwise the values entered are rejected.

The oscillation speeds are saved in the device and cannot be set in the Preset.

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The KARL STORZ POWERSHAVER system S2, consisting of Power Unit S2, shaver handpiece (drill, saw and footswitch) supports the presetting of speed for all combinations of handpieces and blade classes for all four speed memories.


When creating a configuration, default values are assigned to all settings.

You have the following options:

Sockets

- Upper sockets range selection:
 Drill On/Off
 Saw On/Off
- Press the “Drill” button to activate the drill setting for each socket. This is indicated with a green rectangle on the button and in the display (drill symbol with border).
- Press the “Saw” button to activate the saw setting for each socket. This is indicated with a green rectangle on the button and in the display (saw symbol with border and unit in [%]).
- Lower sockets range selection:
 Shaver On/Off
 Oscillation On/Off

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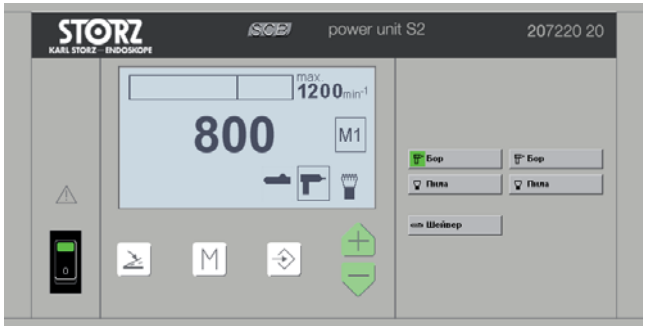- Press the “Shaver” button to activate the shaver settings. This is indicated by a green rectangle on the button, the shaver symbol with border, the enabled handpiece buttons and the enabled oscillation button.
- Press the “Oscillation” button to activate the shaver setting for oscillation. Otherwise the settings are the same as for shaver operation.
- Handpiece range selection:
 High Speed On/Off
 Standard On/Off
 Small Joint On/Off
- Press the “High Speed” button to activate the shaver setting for the high speed handpiece and enable the blade class buttons.
- Press the “Standard” button to activate the shaver setting for the standard handpiece and enable the blade class buttons.
- Press the “Small Joint” button to activate the shaver setting for the small joint handpiece. The blade class buttons remain disabled as the small joint handpiece does not support different blades.

Select the required memory location using the “M” button just as on the device. Then select the instrument as described above and set it using the “Up” and “Down” buttons. Save the changes by pressing the “Enter” button.

IMPORTANT: *When a Preset is opened, the device is loaded with these settings, which are applied according to the activated instrument and the last saved memory location.*

NOTE: *For the speed settings, please refer to the Power Unit S2 instruction manual.*

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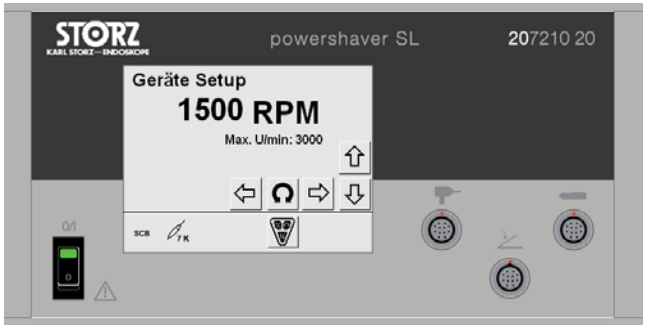
Operation of the devices' own settings is not supported. This functionality is already supported across the whole system by the Presetmanagement of the SCB control NEO which, in this context, extends the range of applications of the POWERSHAVER SL. The specific device settings, such as switching the language and the Service mode with the respective password are not supported by the SCB control NEO; these are made, instead, directly on the POWERSHAVER SL. With SCB control NEO operation, these functions cannot be called up and as such are also not shown in the respective dialog boxes.

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
Creating a configuration corresponds to handling the real device with the exception of selecting the application handpiece.

You can choose this directly by selecting the operating mode “SHAVER”, “MICRO SHAVER” or “DRILLCUT-X”. The selected dialog box for the application handpiece is then displayed. Selecting the application handpiece or selecting the “Operating Mode” button activates the selected operating mode and grants direct access to the device view.

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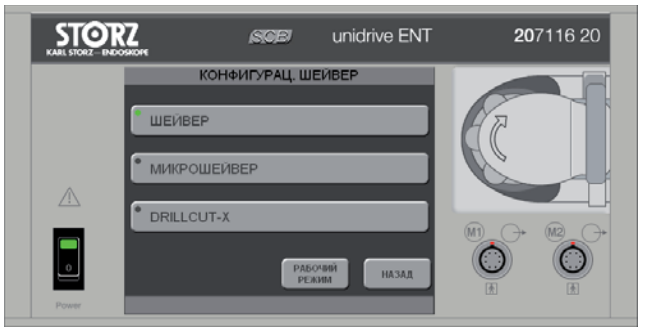Otherwise return to the operating mode selection dialog box without setting new values by pressing the “Back” button.

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The “SINUS BURR” operating mode is no longer displayed in the new device versions of the UNIDRIVE® for the “SHAVER” and “MICRO SHAVER” application handpieces. This operating mode is still displayed in the SCB control NEO software.

The following dialog boxes and functions are supported:

- Work screen:
- Presets: Operating mode selection
- Presets: Application handpiece selection

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4.5.6 Media devices

SCB ACC Control® 200906 32

Read the instruction manual for the device carefully.

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KARL STORZ SCB ACC Control® is the communication module for modern operating rooms. Data transmission, telephoning, live video surgical operations – KARL STORZ OR1™ control NEO opens the door to the information age directly from the OR using the simplest key available: A simple touch with your finger tip on the touch screen monitor is all it takes.

Requirements: An analog telephone connection and a video crossbar (see video crossbar range selection on p. 134).

The SCB control NEO user interface SCB ACC Control® enables operation of connected peripheral devices, such as telephone, accessories and video switching.

WARNING: *Devices which are in direct contact with the patient or which could endanger the health of the patient may not be operated with the SCB ACC Control®.*

SCB ACC Control®

SCB ACC Control® display

The available media devices are displayed in the main window.

As the SCB ACC Control® media system is customer-configured, only the buttons for those functions on the SCB ACC Control® user interface are enabled which can be activated via the “Options/Device Setup/General” menu item.

NOTE: *In the “All Devices” Preset the “Direct Call” button cannot be selected.*

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Telephone function window

Use the telephone function window to establish a telephone connection to another telephone participant or receive incoming calls.


Symbols used:

- Establish a connection
- Disconnect connection
- Clear incorrectly entered numbers

The telephone number entered via the keypad is shown on the left in the number field. The last number can be deleted with the ⌫ button. If the ⌫ button is held down for longer than 2 seconds, all numbers are deleted.

Establish a connection

In order to establish a connection press the ✆ button. The flashing green status display on the ✆ button indicates that the connection is being established. As soon as you are connected, the status display will remain illuminated in green.

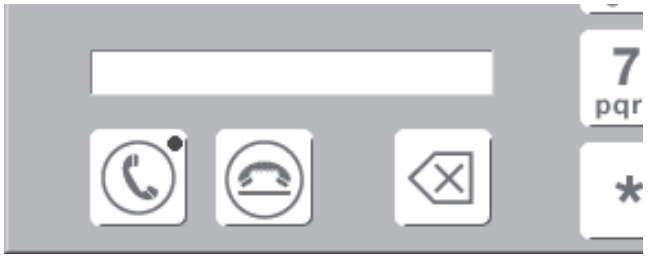If any inputs are made via the telephone buttons during a call, they are sent directly to the telephone system. This is useful if the established telephone connection expects such entries. These additional entries are each shown between two dashes.

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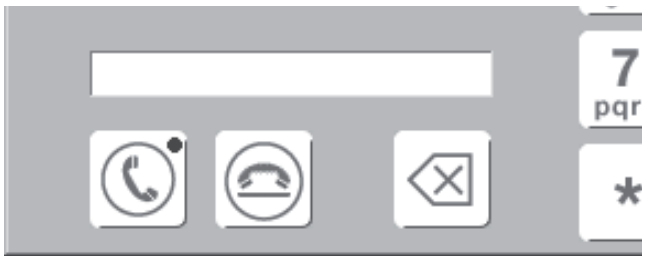
Receive call

A flashing green status display on the button indicates an incoming call. As soon as the button is pressed, the connection is established and the status display remains illuminated in green.

Disconnect connection

Press the button to disconnect the call.
The status display is no longer illuminated.

NOTE: *Sound is output via the headphones of a headset or via the external audio system.*

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
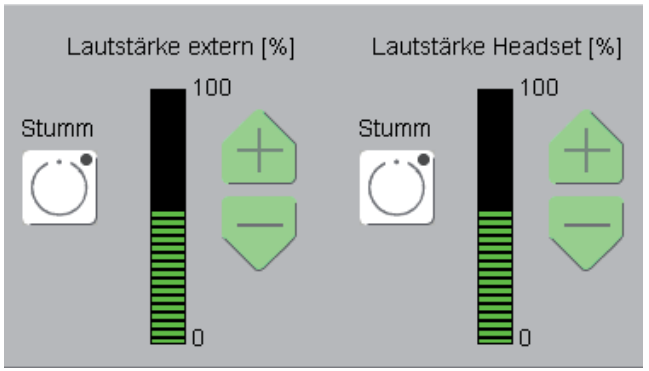
Volume control

In order to increase the quality of comprehension, it is possible to adapt the volume of the headset headphones and the external audio output to the telephone connection.

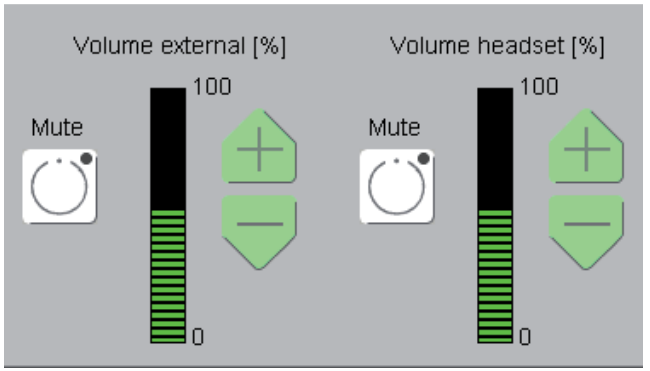
These adjustments to the volume can be made directly in the **telephone** function window. This applies both for the headphones of the headset and for the audio channel connected externally by means of the “External Volume” buttons.

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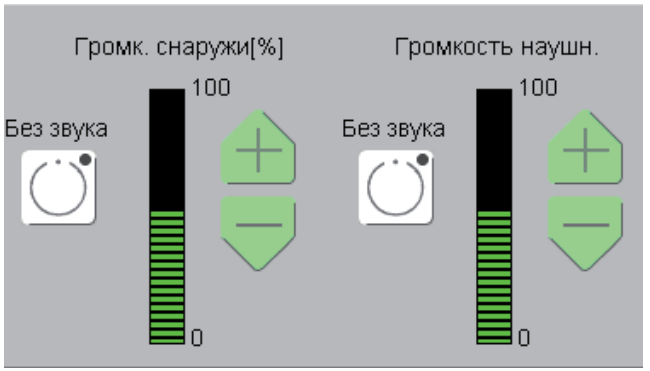The volume is thus increased or decreased using the corresponding buttons. Ths status display shows whether the mute switch (mute) is activated or deactivated.

The normal values for the volume are set via the "Options/Device Setup/Telephone" menu item.

After disconnecting the call, the different volume levels are reset to the values set in the Device Setup dialog box.

NOTE: *If the mute button is pressed in "Headset Volume" not only the headphone but also the microphone of the headset is switched to "mute".*

NOTE: *If mute is enabled the volume cannot be changed.*

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
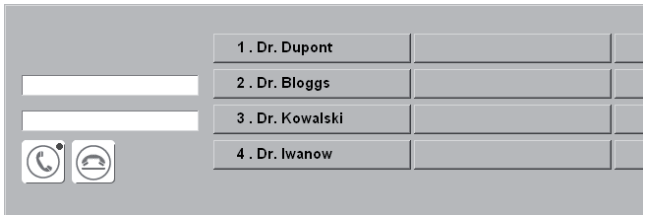
Direct call function window

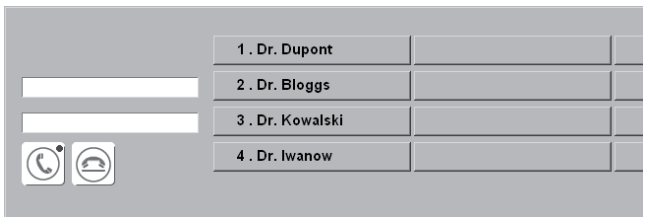A pre-assigned telephone connection can be selected and established in the direct call function window. Telephone connections must have previously been created in the SCB control NEO software via a Preset (see direct call configuration).

NOTE: *The same telephone is operated with the telephone and direct call functions. The status display and function on the ✆ button is therefore the same in both windows.*

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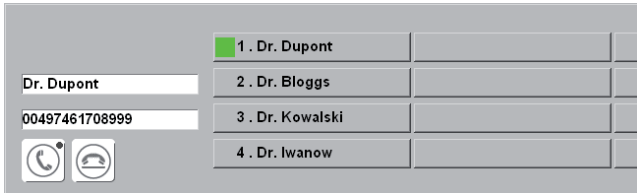
If a button for a pre-assigned telephone participant is pressed, this lights up green. The name and number of the participant are shown on the left.

NOTE: *The volume is set in the* **telephone** *function window.*

Establish a connection

In order to establish a connection, press the button. A flashing green status display on the button indicates that the connection is being established. As soon as you are connected, the status display will remain illuminated in green.

Disconnect connection

Press the button to disconnect the connection. The status display is no longer illuminated.

NOTE: *If during configuration no name or number is saved for a direct call button, then the button is disabled.*

Accessory function window

Depending on the site facilities and customers' specifications, additional external devices, such as lamps or blinds can be controlled via the **accessories** function window. The functions are activated via the displayed buttons.

In principle, a distinction is made between switch, button, pulse transmitter and pulse.

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Switch

In a switch, the relay contacts are continuously switched from one state to another. The switching status is indicated by a rectangular green status display on the left of the button.

The status display will light up green, for example, when a lamp connected to it has been turned on, and will turn off as soon as the lamp is turned off using the switch.

Button

In a button, the relay contacts remain switched as long as the button is pressed. The functionality of the button is indicated by button symbols. The following three functions are distinguished:

- Increase (▲),
- Decrease (▼) and
- Button actuation (no symbol).

Pulse transmitter

The functionality correponds to that of a button, however the hold and interval of the relay contacts must be set. If the button is kept depressed, the relay contact is closed for the hold time and then opened for the interval time. If the button is still depressed, the process keeps repeating itself until the button is released.

Pulse

The functionality of a pulse corresponds to that of a button, however, the hold time of the relay contacts can be set. Even if the button is kept depressed, the relay contact only remains activated for the set hold time.

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
Video Routing function window

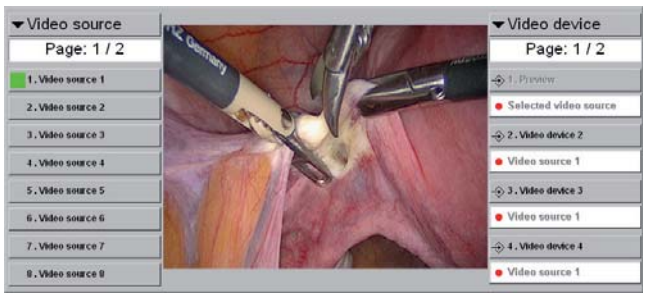
The main image of the video routing function window can be reached via the Video Routing button.

In the Video Routing function window, video sources and video (output) devices can be associated. The window is divided into three sections. On the left, the available video sources are shown, in the middle the preview window and on the right the video (output) devices.

Use the “Video source” button to scroll through the available video sources. The maximum number of video sources and the resulting number of pages to scroll through is determined by the crossbar used. If the OR1™ Frame Grabber Card RGBHV AV Control 20 040476 has been installed, a video preview of the selected video source is displayed in the middle of the dialog box.

Use the “Video device” button to scroll through the available video devices. The maximum number of video devices is also determined by the crossbar used.

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
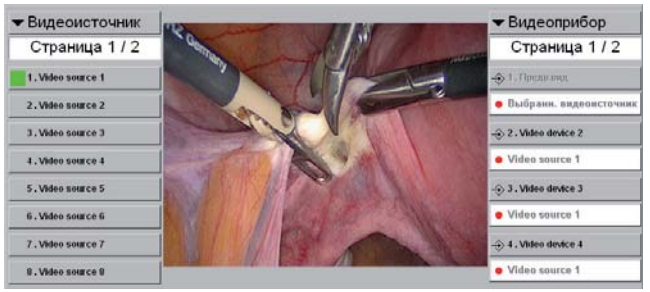
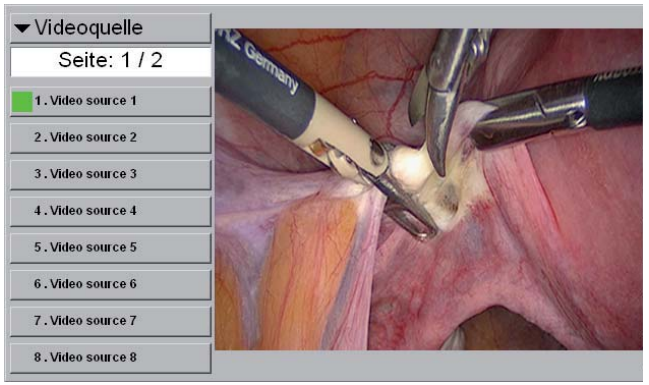
Assign a video source to a video device

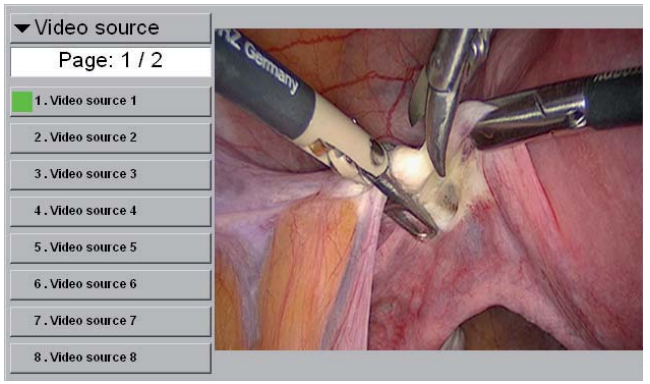To assign a video source to a video device, first select the video source. The selected video source is indicated by a green rectangle on the button. Then press the button of the video (output) device. As soon as the connection has been established, the video source label is shown under the video device.

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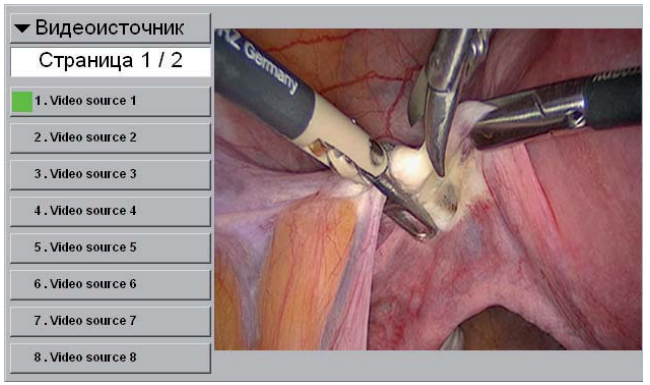
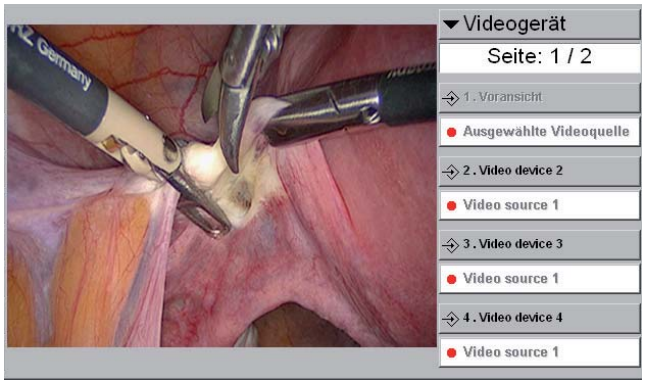
Deassign a video source from a video device

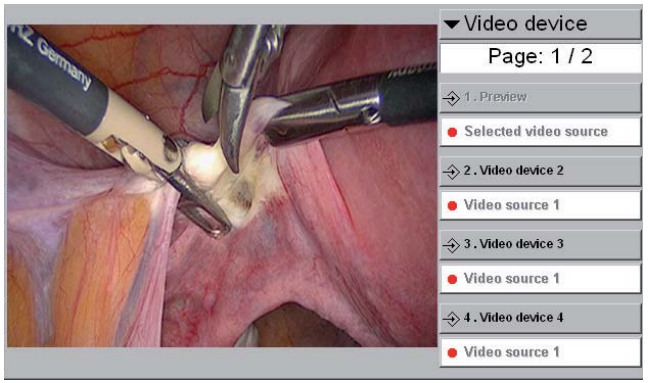To deassign a video device from a video source, press the button of the video device for more than two seconds. Once the connection has been deleted, the video source in the field under the video device disappears.

NOTE: *The designation of the available video sources and video devices as well as the video crossbar used is further described below in the section “Device Setup Video Routing”.*

NOTE: *Depending on the hardware in your OR1™ control NEO, a video preview of the currently selected video source is supported. You will find the corresponding hardware components under “Spare parts, recommended accessories”.*
If the video preview is not supported, a blue rectangle is displayed in the middle. If the preview is supported, a device is always reserved in accordance with the setting in the device setup and therefore cannot be assigned to any other video source.

NOTE: *The deletion of a connection is not supported by the following video crossbars: Gefen 8x8, 16x16*

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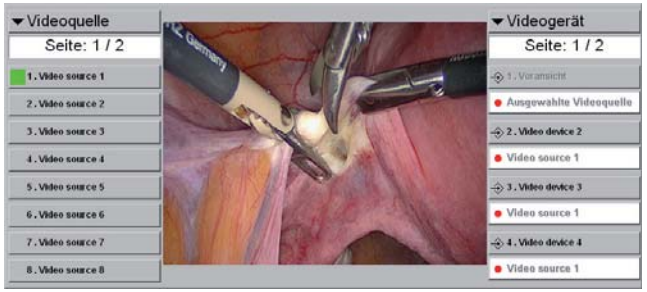
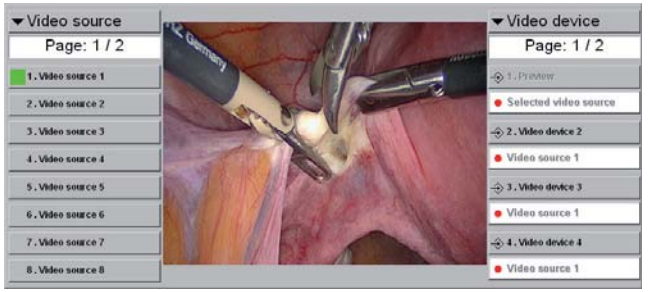
If the video preview is supported in your system*, the current video image of the selected video source will be displayed in the middle of the Video Routing function window.

* The OR1™ Frame Grabber Card RGBHV AV Control 20 0404 76 must be installed.

As a rule, the input that receives a valid video signal is selected automatically on the Frame Grabber Card in the OR1™ control NEO.

For quality reasons, we recommend using the S-Video input.

NOTE: *If the image displayed is not "clean", check whether the video format of the video sources is the same as the set video format (PAL or NTSC). The video format can be changed under the menu item "Options/Device Setup".*

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
Device Setup

The video format can be set in the Device Setup dialog box. See description on p. 134.

Configure SCB ACC Control®

Certain settings of the SCB ACC Control® device can be preset via the SCB control NEO configuration. Presets are only possible for the Direct Call and Video Routing functions.

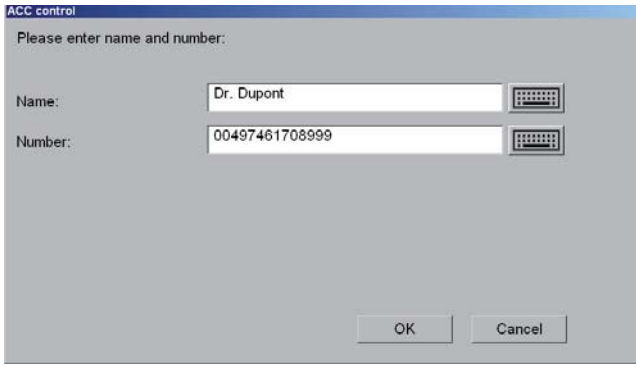Direct Call configuration

In the Direct Call configuration window, the name and telephone number of a Direct Call button are configured. Press a button to display the above input window.

Touching one of the keyboard symbols calls up the on-screen keyboard and allows entries to be made directly via the touch screen. Only the numbers 0-9 as well as the asterisk (*) and hash sign (#) are permitted when entering numbers. You can enter a maximum of 20 characters.

Confirm your entries by pressing the OK button. The window is closed and the entries are saved. The window is also closed when the Cancel button is pressed. However, the entry is not stored and the entry process is canceled.

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
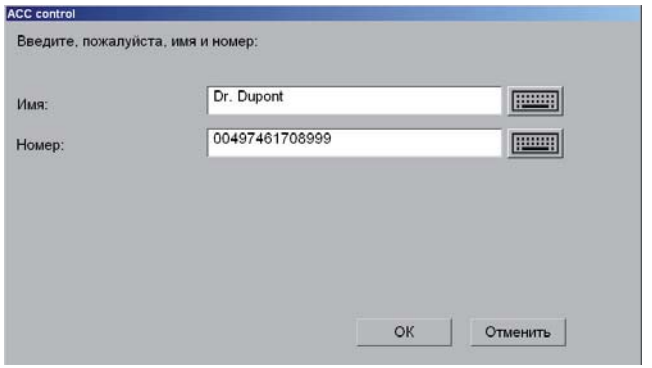
Video routing configuration

Use the **Video Routing** configuration window to preset which video source is to be routed to which video device. Operation of this window is described in the **Video Routing** function window section.

NOTE: *The designations for the video sources and video devices are entered via the* ***Options/Device Setup/Video Routing*** *menu item.*

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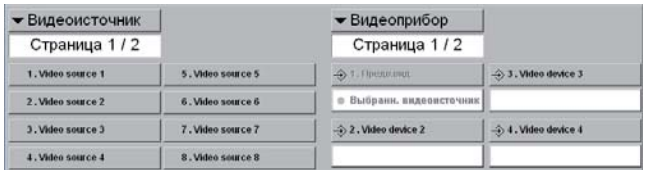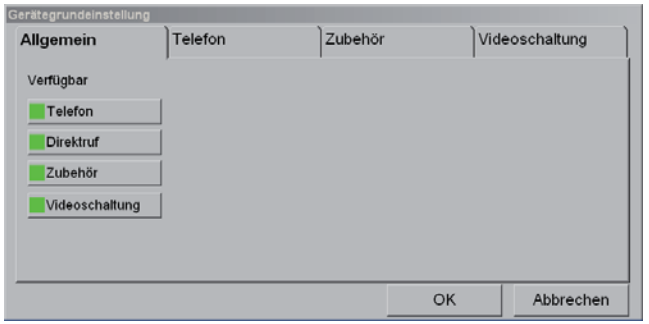
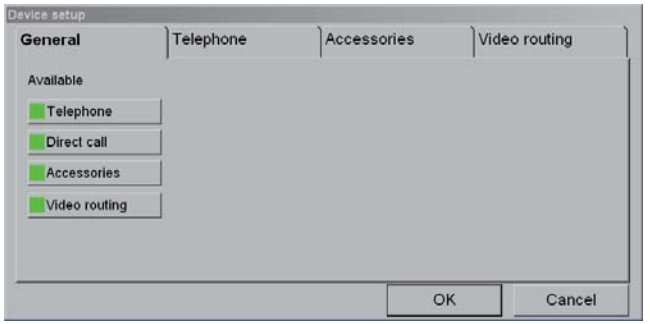
The **Device Setup** dialog box contains the specific SCB ACC Control® setup values. These are divided into four groups which can be selected via the corresponding tab. Press the OK button to accept the changes made, or press the "Cancel" button to reject them. To call up the dialog box, click on **Device Setup** in the "Options" menu. Using the **Telephone**, **Direct Call**, **Accessory** and **Video Routing** buttons, the buttons for the function windows (see "SCB ACC Control®" display) can be disabled or enabled.

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A green rectangle on the left edge of the button indicates that the function window is enabled (SCB ACC Control® display).

The respective controls below the individual tabs are disabled in the Device Setup dialog box if the allocated function window is disabled. The controls in the "Telephone" tab are disabled if the "Direct Call" and "Telephone" function windows are deactivated.

Telephone: This button can be used to disable or enable the "Telephone" function window. The green rectangle on the left side of the button indicates that the "Telephone" function window and the controls in the "Telephone" Device Setup are enabled. The controls in the "Telephone" Device Setup are only disabled when the "Direct Call" and "Telephone" buttons are deactivated.

Direct call: This button can be used to disable or enable the "Direct Call" function window. The green rectangle on the left side of the button indicates that the "Direct Call" function window and the controls in the "Telephone" Device Setup are enabled. The controls in the "Telephone" Device Setup are only disabled when the "Direct Call" and "Telephone" buttons are deactivated.

Accessories: This button can be used to disable or enable the "Accessories" function window. The green rectangle on the left side of the button indicates that the "Accessories" function window and the controls in the "Accessories" Device Setup are enabled.

Video routing: This button can be used to disable or enable the "Video Routing" function window. The green rectangle on the left side of the button indicates that the "Video Routing" function window and the controls in the "Video Routing" Device Setup are enabled.

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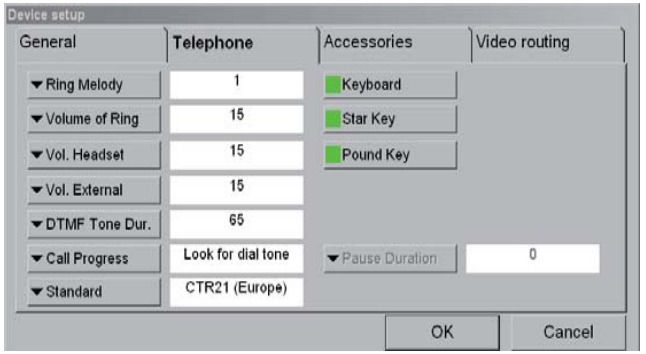
Device Setup – Telephone

With the help of the **Telephone** Device Setup it is possible to make adaptations for the telephone system used. Country-specific and telephone system settings must be taken into consideration. The corresponding details can be found in the instruction manual of the telephone system used.

Ring melody

- Range selection: 1…5
- Pressing the "**Ring Melody**" button changes to the next melody and plays this. The ring melody is output via the internal device speaker of the SCB ACC Control®, the external audio channel and the headset.

Ring volume


- Range selection: 0…31; "**0**" no ring melody is audible; "**31**" the ring melody is output at maximum volume.
- Pressing the **"Volume of Ring"** button increases the volume in increments of "1". After the maximum volume is reached (31) the setting begins again at "**0**".

Headset volume (volume of the headset)

- Range selection: 0…31; "**0**" the other participant is not audible via the headset; "**31**" the other participant is output at maximum volume via the headset.
- Pressing the **Headset Volume** button increases the volume of the headset in increments of "1". After the maximum volume is reached (31) the setting begins again at "**0**".

External volume (volume of the external audio channel)

- Range selection: 0…31; "**0**" the other participant is not audible via the external audio channel; "**31**" the other participant is output at maximum volume on the external audio channel.

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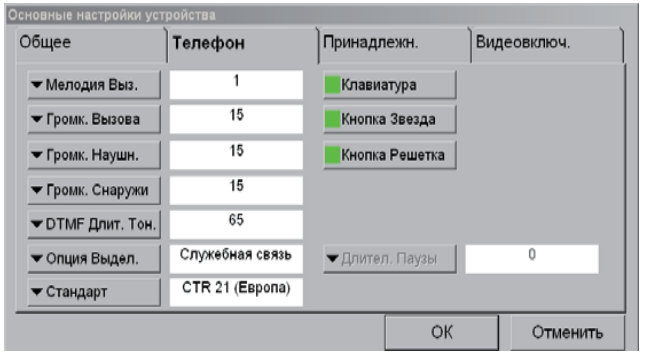- Pressing the "External Volume" button increases the volume of the external audio channel in increments of "1". After the maximum volume (31) is reached the setting begins again at "0".

DTMF Tone Duration (duration of the DTMF tone signal)

- Range selection [ms]: 65…255
- Pressing the "DTMF Tone Duration" button sets the duration of the DTMF tone signal to the value displayed. The increment of the value displayed increases if the button is held down for a longer period.

Call Progress

- Range selection: Look for dial tone, disregard dial tone
- Pressing the "Call Progress" button switches between "Look for dial tone" and "Disregard dial tone connection". If a "Disregard dial tone" connection is selected, the "Pause Duration" button is enabled.

Pause Duration (Duration between picking up the receiver and sending the telephone number)

- Range selection [s]: 1…3;
- Pressing the "Pause Duration" button switches to the next value. After the maximum duration is reached, the setting begins again at "1". This button is only enabled if a "Disregard dial tone" connection has been selected in Call Progress.

Standard (establishes the telephone standard)

- Range selection: FCC, CTR21 (Europe)
- Pressing the "Standard" button switches between the two standards FCC and CTR21 (Europe). The values can be found in the relevant documents of the telephone system.

NOTE: *Changes to the country-specific values (Standard) and to the type of telephone connection (Call Progress) cause the telephone to stop functioning. You must* ***always*** *write down the corresponding values (e. g., in the ready reference guide).*

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Keyboard

- Range selection: Keyboard with alphabet, keyboard without alphabet
- Pressing the “Keyboard” button changes the appearance of the keyboard. You can choose between the keyboard which displays the alphabet below the digit and the display with digits only.

Asterisk button (*)

- Range selection: ON = “*” is displayed,
 OFF= “*” is not displayed.
- Pressing the “*” button renders the asterisk button visible or invisible on the keyboard.

Hash button (#)

- Range selection: ON = “#” is displayed,
 OFF= “#” is not displayed.
- Pressing the “#” button renders the hash button visible or invisible on the keyboard.

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
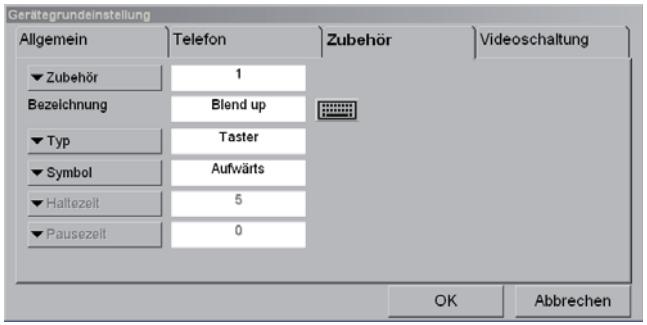
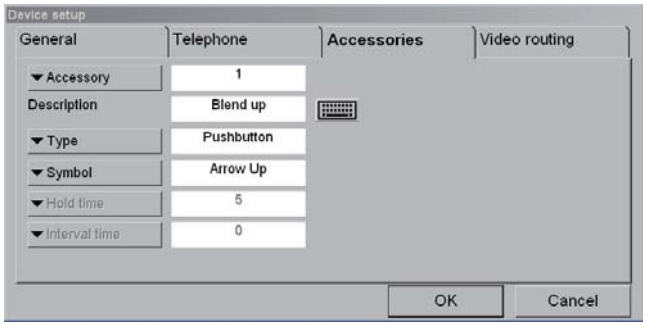
Device Setup – Accessories

The labeling and function of the individual accessory buttons can be set in the Accessory dialog box. The corresponding accessory is selected via the Accessory button. The numbering relates to the Accessory function window and can be viewed in columns from top to bottom and from left to right. Selecting the corresponding accessory button renders the set functions visible.

***WARNING:** Devices which are in direct contact with the patient or which could endanger the health of the patient may not be operated with the SCB ACC Control®.*

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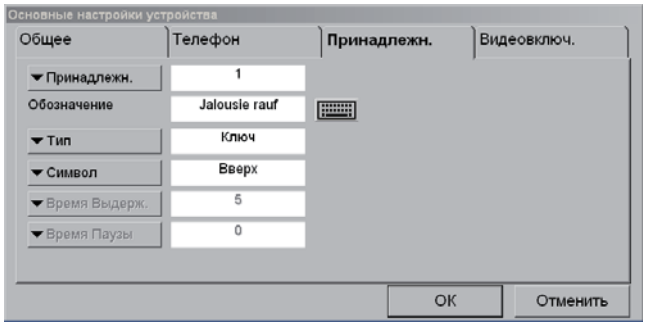Accessories

- Range selection: 1…16;
- Pressing the “Accessory” button switches to the next accessory channel. The description, type, symbol, hold time and interval time of the accessory channel are shown in the lines beneath.
- The description of the accessory channel is changed via the input keyboard (on-screen keyboard) which is displayed by clicking on the labeling field or via the keyboard button. After the last accessory channel is reached, the setting begins again at “1”. The numbering is always in columns from top to bottom and from left to right.

Type

- Range selection: Switch, button, pulse, pulse transmitter
- Pressing the “Type” button switches to the next accessory channel type. Depending on the type selected, the symbol, hold time and interval time buttons are enabled or disabled. The individual types and their time functioning have already been described in the Accessory function dialog box section.

Symbol

- Range selection: No symbol, up arrow, down arrow
- Pressing the “Symbol” button switches to the next symbol. The selected symbol is displayed on the front right of the Accessory channel button.

NOTE: *Symbol is enabled for button, pulse and pulse transmitter.*

Hold time

- Range selection: 1…255
- Pressing the “Hold Time” button sets the duration of the hold time to the value displayed. The time functioning has already been described in the Accessory function dialog box section. The increment of the value displayed increases if the button is held down for a longer period.

NOTE: *Hold time is enabled for button.*

Interval time

- Range selection: 0…255
- Pressing the “Interval Time” button sets the duration of the interval time to the value displayed. The time functioning has already been described in the Accessory function dialog box section.
- The increment of the value displayed increases if the button is held down for a longer period.

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
NOTE: *Interval time is enabled for button and pulse.*

NOTE: *When operating the Universal OR Light, a function of the "Accessory" function window can also be triggered via the light's room control. For further information, please refer to the "Universal OR Light/OR Light Cam" section from page 180 onwards.*

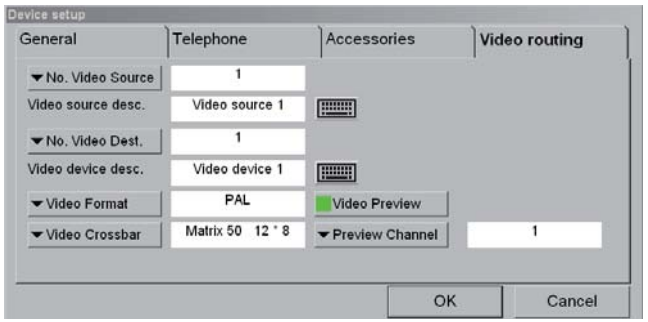Device Setup – Video Routing

Labeling of the individual video sources and video devices can be changed in the video routing dialog box.

No. Video Source

- Range selection: 1…12 or 1…32 depending on the crossbar used
- Pressing the "No. Video Source" button switches to the next video source. The respective description of the video source is displayed in the line underneath. The description of the video source is changed via the input keyboard (on-screen keyboard) which is displayed by clicking on the labeling field or via the keyboard button. After the last video source is reached, the setting begins again at "1". The numbering is always shown in columns from top to bottom and from left to right.

No. Video Destination

- Range selection: 1…8 or 1…16 depending on the crossbar used
- Pressing the "No. Video Destination" button switches to the next video device. The respective description of the video device is displayed in the line underneath. The description of the video device is changed via the input keyboard (on-screen keyboard) which is displayed by clicking on the labeling field or via the keyboard button. After the last video device is reached, the setting begins again at "1". The numbering is always shown in columns from top to bottom and from left to right.

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Video Format

- Range selection: PAL, NTSC
- Pressing the "**Video Format**" button switches between the two video formats, PAL and NTSC. The values can be found in the relevant documents of the video sources used.

NOTE: *Video Format establishes the video format used by a video card integrated in a SCB control NEO computer.*

Video Crossbar (establishes the video crossbar used)

- Range selection: Gefen 8x8, Gefen 16x16, MAV+ 8x8, MAV+ 12x8, MAV+ 16x16, MATRIX 50 12x8, MATRIX 50 12x4, MATRIX 50 8x8, MATRIX 50 8x4, DXP44 4x4, DXP88 8x8
- Pressing the "**Video Crossbar**" button switches between the different video crossbars. The values can be found in the relevant documents of the video crossbars used. In principle, all video crossbars which the Extron protocol supports can be used.

NOTE: *If the selected video crossbar has more inputs and/or outputs than the connected video crossbar or an incorrect type has been selected, the SCB control NEO software displays a corresponding error message.*

Video Preview

- Range selection: On = Preview
Off = No Preview
- Pressing the "**Video Preview**" button either enables or disables the "**Preview Channel**" button. If the "**Preview is ON**" the video signal which is in the "**Preview Channel**" is always shown in the Preview window. If the Preview function is not active, no image is shown in the Preview window. (Black or blue rectangle).

Preview Channel

- Range selection: Number corresponds with the video devices which have been selected with the video crossbar.
- Pressing the button selects the channel (video device) on the video crossbar which is used for the preview and which is to be connected to the OR1™ control NEO.

Telestration 20 091432

Read the instruction manual for the interface module and carrier device carefully.

All the settings must be made as described in the instruction manual and the displays are to be interpreted accordingly.

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The AIDA DVD-M is displayed as a pseudorealistic device. The following control functions can be executed.

Still image, video recording and speech recording:

- Press the “Still image” button to start still image capture. At the same time speech recording is started automatically depending on the device setting. Active still image capture is indicated by the green LED in the button.
- Press the “Video recording” button to start video recording. An active recording is indicated by the green LED in the button.
- Press the “Voice recording” button to start voice recording. The button is enabled as soon as a video recording starts. An active voice recording is indicated by the green LED in the button.

NOTE: *The warning “Not in capture mode” indicates that the device is not in the capture dialog box. End the action on the device and switch back to the capture dialog box.*
The warning “Video signal” indicates that no video source is connected. Connect a suitable video source to the device.

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
NOTE: *To ensure safe operation with the AUTOCON® II 400 touch panel, the user must choose a so-called procedure at startup. If no procedure has been selected on the device, the device will indicate this with an error message, but cannot be operated. Choose a procedure on the device in order to operate the device via the SCB control NEO software.*

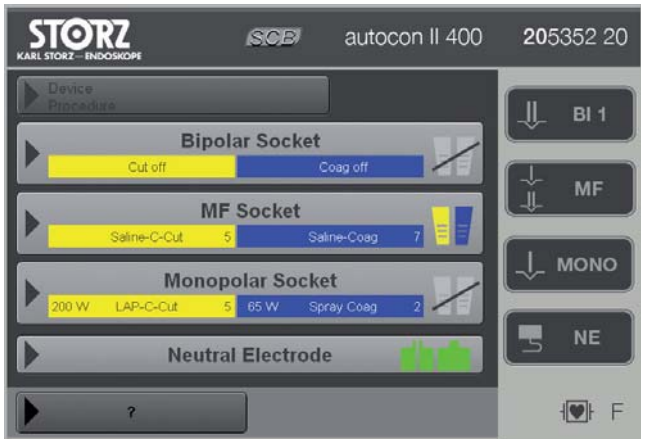For an optimum view and operation, the content of the touch panel is displayed in a full screen format. Operation is similar to operation directly on the device, and should be performed as described in the AUTOCON® II 400 instruction manual.

NOTE: *Devices from the first generation (manufactured before 11.2011, or up to serial no. xx03508) are shown in the SCB control NEO software as "Autocon II 400 °".*

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The SCB control NEO software does not support the “procedure management” of the device itself. This function is already supported by the Preset management of the SCB control NEO software throughout the system and therefore extends the AUTOCON® II 400 scope of application. The SCB control NEO software also does not support “Device Setup” but this function can still be carried out directly on the AUTOCON® II 400. In the SCB control NEO software it is not possible to call up these functions, and thus these are not shown in the corresponding dialog boxes.

If these functions are selected directly via the AUTOCON® II 400, the text depicted above is displayed in the SCB control NEO software.

In such cases, the user can immediately change to the main screen by pressing the “Back” button at the top left.

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- Help texts
- Back to the procedure list
- Procedure list

Device setup: Selection options with the respective selectable menu items:

- User settings
- Application settings
- Mode settings
- Date and time
- Device configuration
- System log
- Change password
- Service programs

Procedure management with the corresponding menu items:

- Overwrite procedure
- Delete procedure
- Save as new procedure
- Rename procedure
- Create procedure level B

NOTE: *When setting values in a Preset, please observe section 4.4.5.4 "New…/Change"… Presets.*

NOTE: *An AUTOCON® II 400 which is connected to the OR1™ control NEO computer is synchronized with the following SCB system settings.*

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Force FX

Read the instruction manual for the interface device/carrier device and interface module carefully.*

* The special SCB control NEO device functions are described in the instruction manual of the device.

Device messages

Message	Priority
Activation not possible	low
Internal error	low
RS232 failure	low

ForceTriad™

Read the instruction manuals for the carrier device and interface module carefully.

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The pseudorealistic view of the COVIDIEN ForceTriad™ energy platform contains controls and display elements (and thus functions) from the real ForceTriad™ energy platform (with the exception of the ON/OFF switch and system tray, including functions selected via the system tray). Operation via the SCB control NEO software is largely identical to operation via the device and can be looked up in the COVIDIEN ForceTriad™ instruction manual.

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The screen arrangement in this view is as follows:

- The upper line contains the left and middle device screens.
- The lower line contains the right device screen.

NOTE: *The KARL STORZ SCB control NEO software can be used to operate up to three COVIDIEN ForceTriad™ energy platforms.*

NOTE: *The COVIDIEN ForceTriad™ energy platform can be operated via the KARL STORZ SCB control NEO software and directly via the device itself simultaneously, provided that no functions are mutually exclusive.*

NOTE: *The SCB control NEO software does not support an additional acoustic signal itself when a button is pressed, this is only sounded by the COVIDIEN ForceTriad™ energy platform.*

NOTE: *A COVIDIEN ForceTriad™ energy platform connected with the KARL STORZ SCB control NEO is synchronized with the language setting of the system.*

Dialog boxes and functions which are not supported:

In principle, all settings must be performed in accordance with the COVIDIEN instruction manuals.

The functions stored in the system tray can only be performed directly via the COVIDIEN ForceTriad™ energy platform.

NOTE: *By touching one of the three device screens in the device view the user can change to the main screen.*

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
The device messages listed below are displayed by the following pop-up:

- Demo mode
- Software upgrade
- Calibration required
- Error deactivated

If a serious error occurs, the following display is shown in the pseudorealistic view.

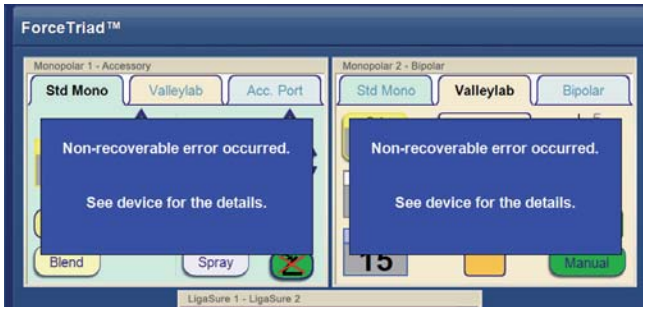Preset

A preset can be added to the COVIDIEN ForceTriad™ energy platform, but the settings of the ForceTriad™ device cannot be preset.

“Widgets” (compact menu – Overview)

The COVIDIEN ForceTriad™ energy platform can be operated via both the device view and the compact view in the overview.

An “info widget” is shown when you change to the compact view.

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NOTE: *The "info widget" is shown when:*

- *an instrument is connected to (or removed from) the COVIDIEN ForceTriad™ energy platform;*
- *a pop-up note concerning the settings shown in the widget is displayed;*
- *you change from the device view to the compact view;*
- *the menu is shown on the unit.*

When an instrument is activated, the corresponding settings are displayed in the widget.

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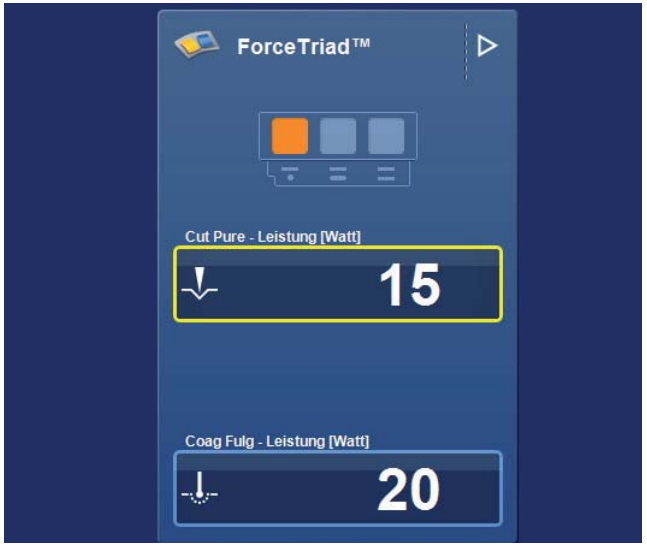
The following settings can be displayed and edited in the widget:

- Standard Monopolar Mode (Power Cut/Coag)
- Valleylab Mode (Power Cut/Valleylab/Coag)
- Accessory connector (Power Cut/Coag)
- Bipolar electrode (Power)
- LigaSure 1 and LigaSure 2 (degree of desiccation)
- Bipolar resection (Effect Cut/Coag)

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An orange screen display of the device view in the widget indicates the screen to which the current display belongs.

NOTE: *For information on the remote control of the COVIDIEN ForceTriad™ energy platform via widgets, please see Section 4.4.4.2.*

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KLS Martin maxium® ME402

Read the instruction manuals on the carrier device and interface module carefully.

The pseudorealistic view of the KLS Martin maxium® electrosurgical generator contains almost all the buttons of the real KLS Martin maxium® device. In the view the rotation regulator has been replaced with the buttons ▲ "Increment/up", ▼ "Decrement/down" and ● "Confirm". Moreover, the "MENU" button is not shown on the view as this function is not supported by KARL STORZ SCB control NEO. The button assignment is identical to the real device.

NOTE: *A maximum of one KLS Martin maxium® electrosurgical generator can be operated via the KARL STORZ SCB control NEO at any one time.*

NOTE: *The KLS Martin maxium® electrosurgical generator can be operated via the KARL STORZ SCB control NEO and directly via the device itself simultaneously provided that no functions are mutually exclusive.*

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Dialog boxes and functions which are not supported

In principle, all settings as per the KLS Martin instruction manual are to be carried out. However, the KARL STORZ control NEO software does not support operation of specific device settings, i.e., maxium® setup, display of maxium® information functions and maxium® program management.

These functions are to be executed directly on the KLS Martin maxium® electrosurgical generator.

With KARL STORZ SCB control NEO operation, these functions are not displayed and as such are also not shown in the respective dialog boxes.

If program management or setup functions are selected on the KLS Martin maxium® electrosurgical generator, the text shown below is displayed on the pseudorealistic view. In such cases, the user can switch immediately to the main screen by touching the text field. In such cases, the user can immediately change to the main screen by pressing the “Cancel” button.

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
If the service function has been selected on the KLS Martin maxium® electrosurgical generator, the text shown below is displayed on the pseudorealistic view. In such cases, operation is then only possible via the device.

All dialog boxes which can only be displayed on the KLS Martin maxium® electrosurgical generator via the “MENU” button, e.g., programs, setup, audio settings, service, etc., are not supported by the SCB control NEO software.

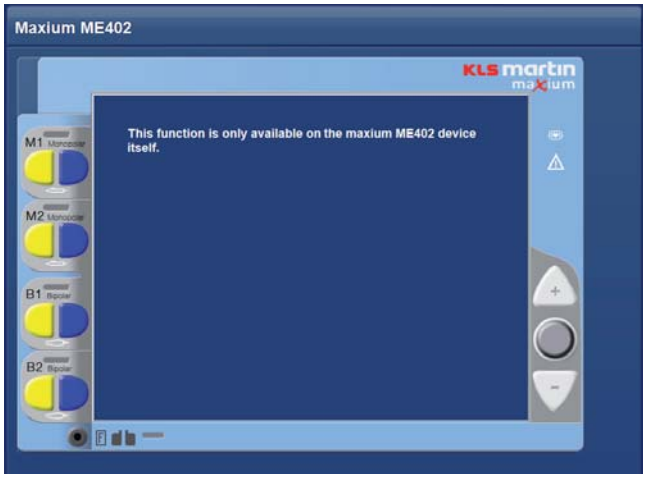Creating a preset

The KLS Martin maxium® electrosurgery generator can be added to the preset when creating a preset within KARL STORZ SCB control NEO.

The KARL STORZ SCB control NEO offers the possibility of saving the set device settings on the OR1™ control NEO in the form of a preset and then loading them onto the KLS Martin maxium® electrosurgery generator later.

NOTE: *If the downloading/saving of the device settings fails, they will be automatically downloaded and saved the next time the preset is opened.*

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
If no device settings are saved, the dialog box shown above is displayed upon opening the Preset. This is indicated to the user by the fact that all the channel elements in the preview dialog box are displayed as “deselected” (gray). On confirming this dialog box, the current device settings are saved on the OR1™ control NEO.

If a saved Preset is to be saved on the KLS Martin maxium® electrosurgical generator, all the channel settings are shown in the preview dialog box (as shown in the image below) for inspection by the user. The user can now overwrite the Preset on the OR1™ control NEO by activating the “Update Preset” checkbox.

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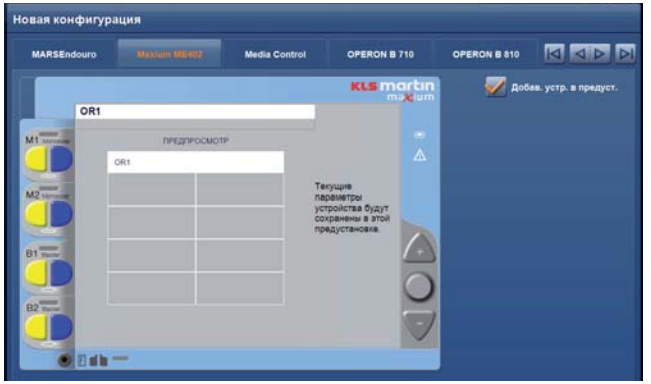System 2450™/System 5000™

Read the instruction manuals on the carrier device and interface module carefully. Using the interface module **20**091590 a ConMed System 2450™ and/or a ConMed System 5000™ can be connected to the SCB control NEO system.

In principle, all settings as per the ConMed instruction manual are to be carried out.

Using ConMed system 2450™

The pseudorealistic view of the ConMed system 2450™ contains all buttons (and therefore functions) of the ConMed system 2450™ control panel. Operation via the SCB control NEO software is identical to the original and can be looked up in the ConMed system 2450™ instruction manual.

NOTE: *Alongside operation with the SCB control NEO software, the ConMed system 2450™ can, of course, still be controlled via the device's control panel. This means the "Original" has a higher priority than the SCB control NEO software.*

NOTE: *An additional acoustic symbol is only sounded on the ConMed system 2450™ device and not within the SCB control NEO software when a button is pressed.*

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Using ConMed system 5000™

The pseudorealistic view of the ConMed system 5000™ contains all buttons (and therefore functions) of the ConMed system 5000™ control panel. Operation via the SCB control NEO software is largely identical to operation via the device's control panel and can be looked up in the ConMed system 5000™ instruction manual.

NOTE: *In accordance with the ConMed system 5000™ instruction manual, the power can be adjusted using the manually controlled accessories. This control (via the accessories) is not shown in the SCB control NEO software.*

Alongside operation with the SCB control NEO software, the ConMed system 5000™ can, of course, still be controlled via the device's control panel. This means the "Original" has a higher priority than the SCB control NEO software.

NOTE: *An additional acoustic symbol is only sounded on the ConMed system 5000™ device itself and not within the SCB control NEO software when a button is pressed.*

The ConMed ESU system 2450™ and/or ConMed system 5000™ can be added to a Preset; all settings, with minor restrictions, can be preset.

The restrictions that cannot be preset are as follows:

For the ConMed system 2450™:

- The volume of the device's loudspeaker.

For the ConMed system 5000™:

- The settings are always based on the selected program in the device. Only one program position (the one most recently selected in the device) can be preset.
- The buttons for selection and saving a program cannot be selected.

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VIO 300 D

Read the instruction manual for the interface device carefully.

The interface device 20 091332 enables the HF device VIO 300 D from ERBE to be connected to the SCB.

In principle, all settings as per the ERBE instruction manuals are to be carried out. However, the SCB control NEO software does not support operation of specific device settings, i.e., VIO setup, display of VIO information and instruction functions, the selection of the NEO Natal electrode, the ReMode function, and VIO program management.


These functions must still be carried out directly on the VIO.

Within the SCB control NEO software, it is not possible to call up these functions, and thus these are not shown in the corresponding dialog boxes. Should these functions be selected on the VIO device, the text shown above is displayed.

As administration of the device's own programs can be replaced by the SCB control NEO presets (see Menu "Presets"), the following points must be observed.

- Arrangement of the individual sockets
- Presets: APC VIO switchover
- Presets: Presetting of the instrument parameters
- Presets: Selection of whether NeoNatal or ReMode function is required (see note on next page).

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


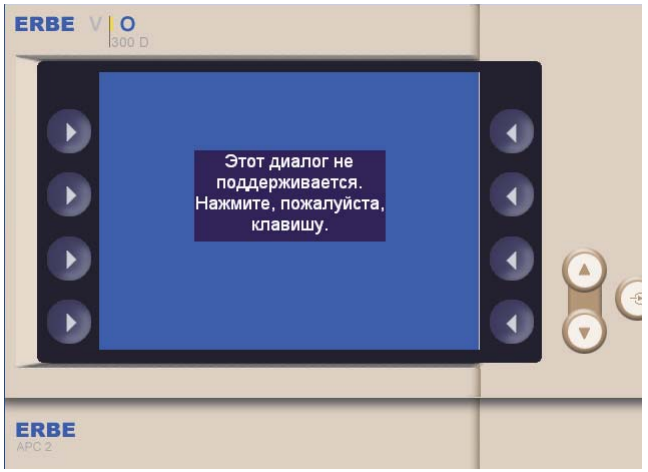
NOTE: *The NEO Natal electrode and/or the remote function are not supported via the preset. This means that if these functions are used, the device, although loaded, is not initialized.*

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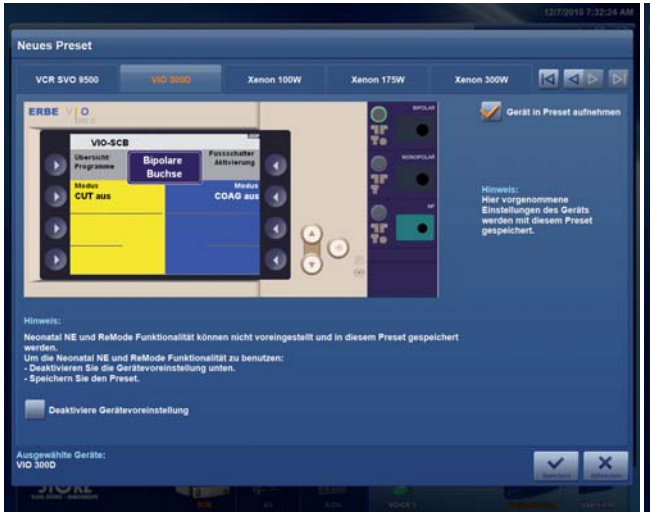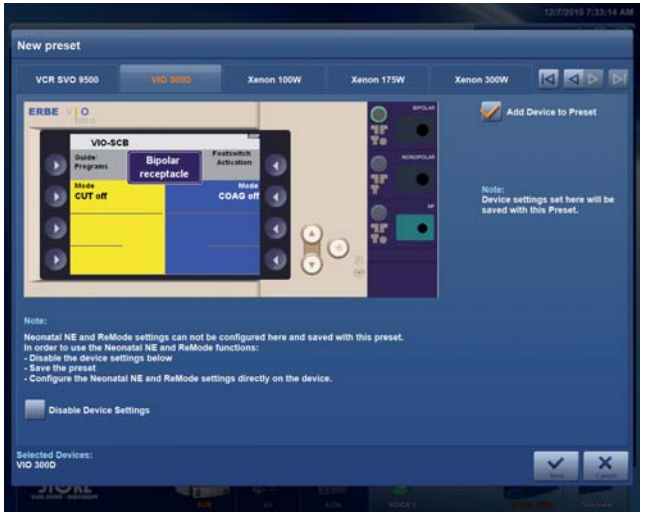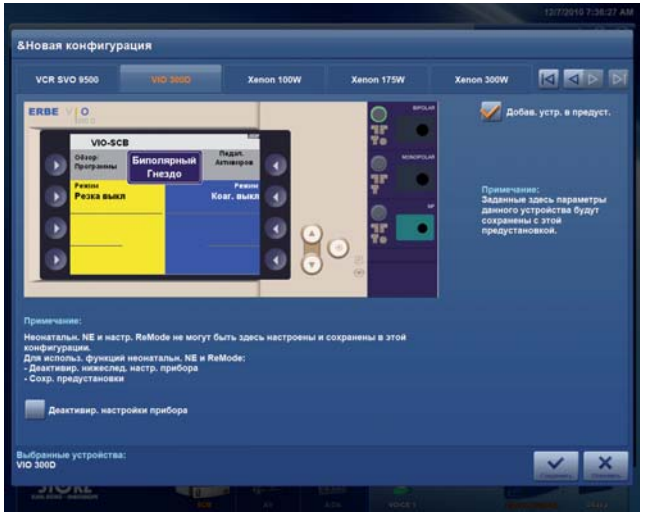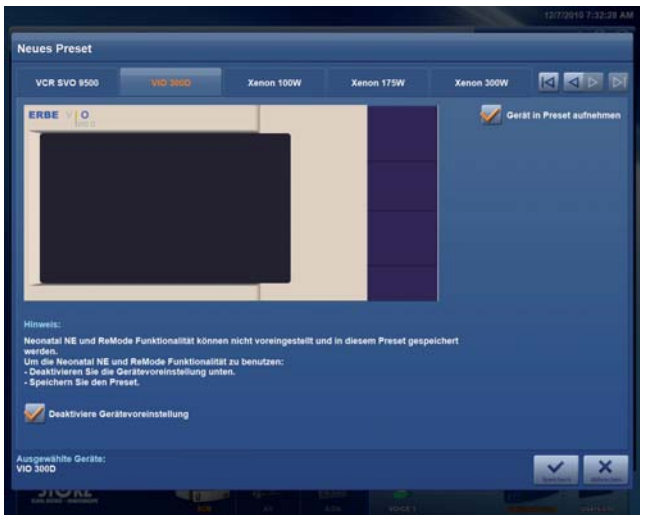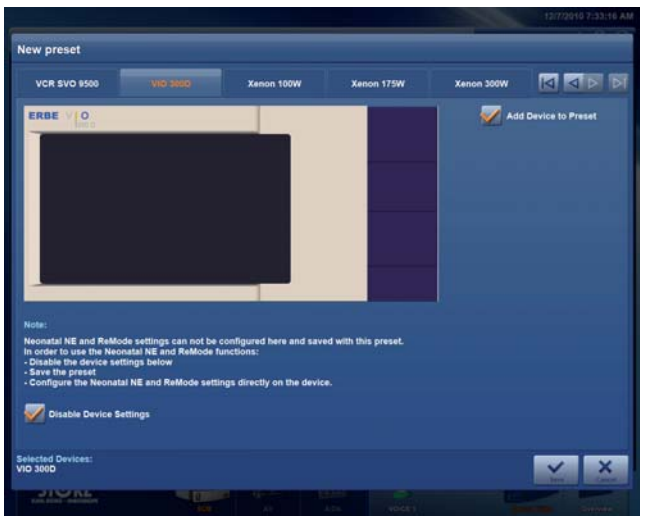
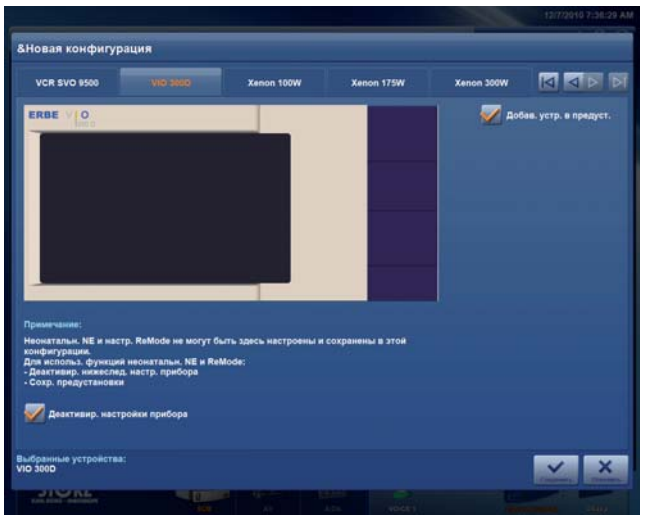
Arrangement of the individual sockets

For system-related reasons, arrangement of the individual sockets within the VIO device may differ from device to device (factory-set).

For this reason, the device must be displayed in the device list at least once before a Preset is created.

Hereby the respective arrangement of the sockets of the SCB control NEO software is adopted and saved. If the arrangement of the sockets on the device changes and therefore the stored Preset is no longer usable, the system automatically resets the default values.

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Preset: APC VIO switchover

As the view of the APC device in the Preset dialog box can no longer be attached under the VIO device, the switchover to the APC unit's sockets is carried out via an additional selection button. For this purpose, the existing VIO sockets are hidden and the APC sockets are displayed. The method of operation is the same. The VIO sockets are then displayed again via another selection button.

NOTE: *The values in the display window always refer to the active socket. Switchover from APC to VIO or vice versa does not change the status of the active socket.*

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
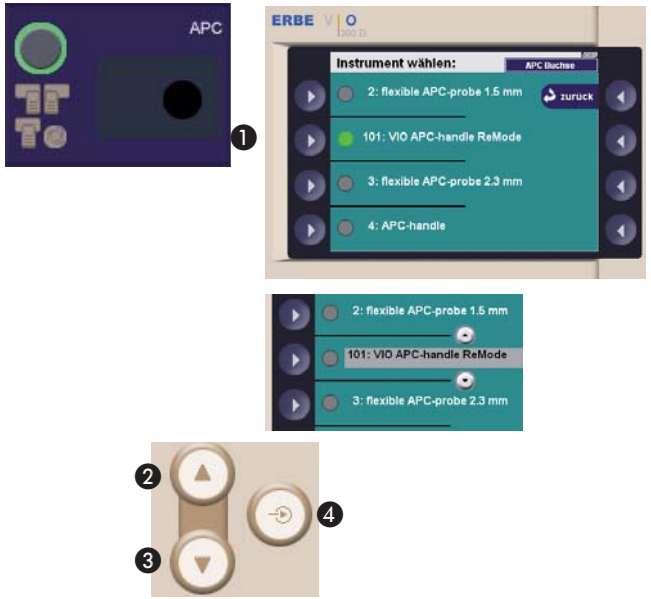
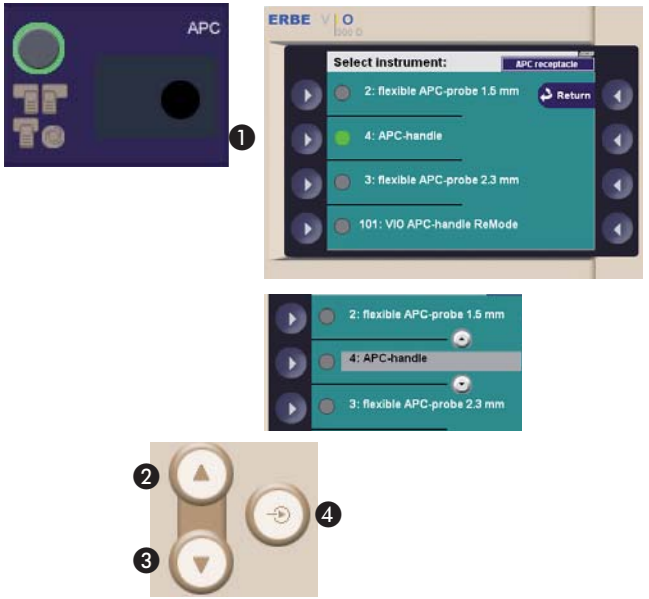
Preset: Presetting of the instrument parameters

Presettings of up to 4 different instruments can be stored in the Preset for each MF or APC socket.

If you touch on the area of the socket with a green background, the corresponding instrument dialog box is opened. If no presetting has been stored for an instrument, all the entries are displayed with “No Instrument”. If you press the selection button ❶, the instrument list is displayed. You can display the desired instrument with the Up/Down buttons ❷/❸ and activate it with the Enter button ❹. Depending on the instrument type, this also implements the corresponding presettings, such as the “Default Mode”, the maximum power and the setting for whether all the modes are supported or only the “Default mode”. To change the settings of an existing instrument, select the instrument concerned via the selection button and activate it with the Enter button.

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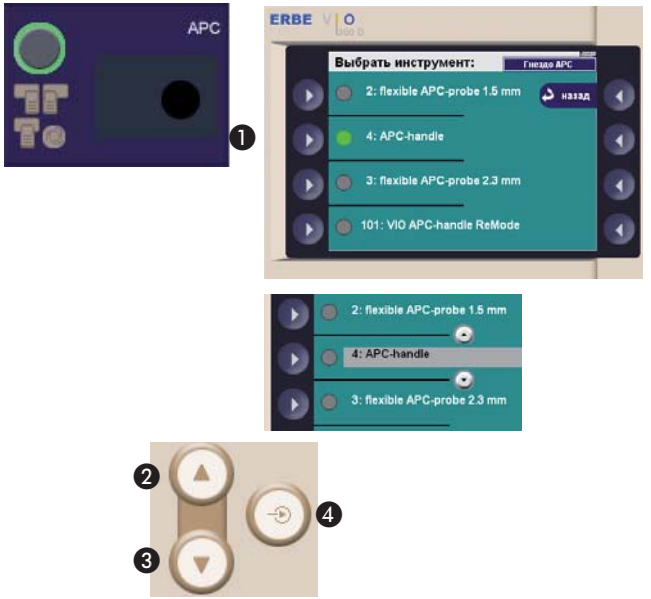NOTE: *If a new instrument from ERBE is plugged into the VIO device during operation of the SCB control NEO software with connected VIO device, the specific values are transferred from the VIO to the KARL STORZ OR1™ control NEO, with the exception of the description. For this reason, only one number may appear in some cases when an instrument is selected. The names can only be updated by performing an appropriate update.*

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- Select program
- Other functions
- Info
- Note
- Additional dialog boxes which are activated in the Expert Mode, such as: cutting duration and cutting interval with the Endo Cut I and Endo Cut Q, modulation with Coag modes, BiClamp with AUTO STOP and PRECISE APC.

Device messages

Message	Priority
VIO configuration has been changed	low
System of neutral electrode has been changed	low

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4.5.10 OR tables

Alphamaquet
Alphamaxx
JUPITER
MARS
MARS Endouro
OPERON® B 710/B 810

Read the instruction manual for the respective interface device/carrier device and interface module carefully.*

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Device messages

Message	Priority
Command sequence number erroneous.	high
Device type error	high
Command received without SCB node enable.	high
Command received without SCB control NEO software enable.	high
Device type error	high
IR drive error	low
EEPROM error	low
Communication error	low
Error	low
Please check connection to all touch screens.	low
SESAM error	low

OPERON device messages

Message	Priority
Device type error	high
Command sequence number erroneous.	high
Command without SCB enable	high
Command received without SCB node enable.	high

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When the IR field is touched, the IR code dialog opens. The IR code can be changed here. The operation of the IR code dialog is described in Chapter 4.4.5.13.

Adjusting the operating table

The pseudorealistic view of the OR table contains all the buttons (and therefore the functions) of the actual MAQUET remote control which are supported by the SCB control NEO system. The button assignment is identical to the actual remote control (with the exception of Magnus, see page 160).

The “Menu navigation button” can be used to select the current display in the display.

The “Changeover button [leg plate]” determines whether the left, right, or both leg plates are moved. Both plates are selected in the default setting.

NOTE: *Please refer to the MAQUET instruction manual for the precise description of operation.*

An operating table is moved for as long as a button is pressed via the touch screen. The movement ends automatically as soon as the touch screen button is released.

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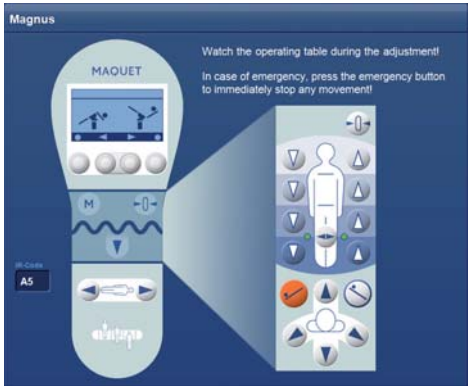
The following functions are not supported in the pseudorealistic view of the Magnus remote control:

- Only the following screens can be shown in the display: “FLEX”/“REFLEX” and “BEACH CHAIR”/“FLAT”.
- Save button
- Status display for CRANIAL/CENTRAL/CAUDAL position
- Changeover button [patient position]

BetaStar/Alphastar (Pro)/Magnus device messages

Message	Priority
Command sequence number erroneous.	high
Command received without SCB node enable.	high
Command received without SCB control NEO software enable.	high
IR drive error	low
RS232 error	low
Please check connection to all touch screens.	low

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OPT OR table OPT/90, OPT/100, Vanto

Read the instruction manuals of the carrier device and interface module carefully. Using the interface module 20091690 an OPT/90, OPT/100 or Vanto can be connected to the SCB.

In principle, all settings as per the OPT instruction manuals are to be carried out.

Using OPT Access Point

The pseudorealistic view of the OPT Access Point (OAP) is displayed if the OR table has not been switched on or registration has not been performed or no connection has been created as per the OPT instruction manual.

NOTE: *Please refer to the OPT instruction manual for the description of operation.*

NOTE: *The OPT OR table can be added to a Preset, settings cannot be preset.*

OPT/90, OPT/100, Vanto device messages

Message	Priority
Sequence command error	high
Please check connection to all touch screens.	high
General	low
Battery low	low
Incompatible OR table	low

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SATURN Select
STERIS® OR table
TITAN

Read the instruction manual for the respective interface device/ carrier device and interface module carefully.*

* The special SCB control NEO device functions are described in the instruction manual of the device/module.

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SATURN Select/TITAN device messages

Message	Priority
Command sequence number erroneous.	high
Command device type error	high
Command received without SCB node enable.	high
Command received without SCB control NEO software enable.	high
Device type error	high
IR drive error	low
EEPROM error	low
Communication error	low
Error	low
Please check connection to all touch screens.	low
SESAM error	low

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STERIS® OR Table device messages

Message	Priority
Event notification enabled	high
Device type error	high
Command without SCB enable	high
Device type error	high
Battery depleted. Please connect immediately to AC power.	middle
RS232 error	low
OR-Table communication error	low
SESAM error	low
Error	low
Please check connection to all touch screens.	low
Battery low. Please charge the battery.	low

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TruSystem 5500/7500

Read the instruction manual for the respective device carefully. All the settings must be made as described in the instruction manual and the displays are to be interpreted accordingly.

Differences in device operation:

- Touch screen on TruSystem 5500/7500 displayed on the handheld remote control via flex/reflex buttons.

By clicking on the L or R buttons, you can choose to adjust the left, right, or both leg plates. The selected leg plate is identified in the realistic view by a white marking on the end of the leg plate. Both plates are selected for adjustment in the default setting. The leg plate selection is registered for ten seconds or as long as the leg plate is moved, then the selection returns to the default setting.

Touch the IR field to open the IR code dialog. The IR code can be changed here. The operation of the IR code dialog is described in Chapter 4.4.5.13.

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Device messages

Message	Priority
Command sequence number erroneous.	high
Device type error	high
Command received without SCB node enable.	high
Command received without SCB control NEO software enable.	high
IR drive error	low
EEPROM error	low
RS232 error	low
Please note alarm message at device.	low
Please check connection to all touch screens.	low

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"Universal" OR Table

Via "OR Table", the SCB control NEO software can control different OR tables foreseen for this purpose by the respective OR table manufacturer.

The respective OR table specifies one of the three following control methods and thus defines whether the corresponding element is displayed or not.

- Table ID (unidirectional, e.g., via infrared)
- "Binding" (bidirectional, e.g., via radio, Bluetooth, etc.)
- Cable connection (bidirectional)

For control via "Table ID" (ID = identification number) the ID specified by the OR table manufacturer (range 1–32) must be entered in the Table ID window of the SCB control NEO software (see Table ID button below right).

NOTE: *For OR tables controlled unidirectionally (e.g., via infrared), the OR table manufacturer must assign the OR table an index ("Table ID"). The specification of ID nos. allows different OR tables to be remote controlled via the SCB control NEO software. Ensure that the assigned index (ID sticker) is clearly visible on the OR table.*

NOTE: *For detailed information on its accurate affixing please refer to the information manual of the OR table manufacturer. The sticker is supplied by the OR table manufacturer.*

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NOTE: *"unidirectional" OR table control:*
If the OR table is displayed in the device view, but cannot be moved, the Table ID entered might not be the same as that of the OR table. Please enter the Table ID in the KARL STORZ SCB control NEO software accordingly.

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
The "**Binding**" control takes place over three stages/statuses.
"Disconnected" status: white symbol, white text "Connect"

NOTE: *The "Connect" button is not pressed.*

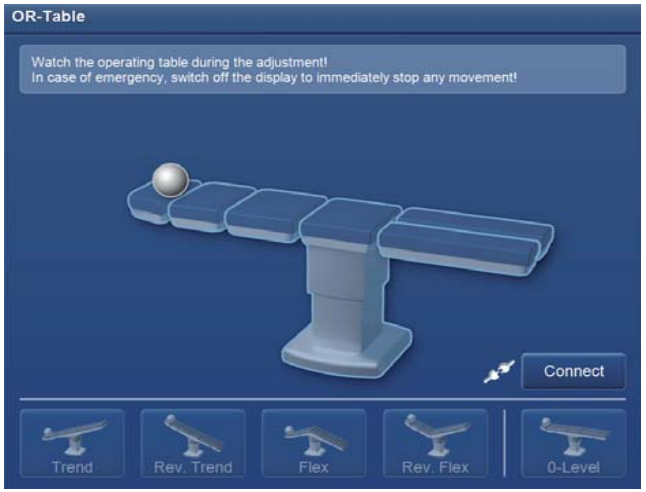"Connecting" status: orange symbol and text

NOTE: *The "Connect" button is/stays pressed.*

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The device view of the OR table supports the following functions:


- Height adjustment of the table pedestal
- Longitudinal adjustment of the table top
- Tilting of the table top
- Setting of individual panels
- Trendelenburg (Trend*)/ Anti-Trendelenburg (Rev. Trend*)
- Flex / Reflex (Rev. Flex*)
- Zero position (0-Level*)
- Patient orientation (Normal/Reverse)

* Descriptions used in the device view

The range of functions offered depends on the OR table connected.

The OR table is controlled either via elements specified by the manufacturer (OR table) (head, ..., foot/feet; see orange arrows) or via five pictograms** with predefined positions (Trend, Rev. Trend, Flex, Rev. Flex, 0-Level). **Depending on device

NOTE: *For safety reasons (danger of being distracted) the position of the individual elements is not shown graphically. As soon as the table is moved, it must be observed to avoid endangering the patient or causing collisions. The display of the OR table varies depending on the current table configuration.*

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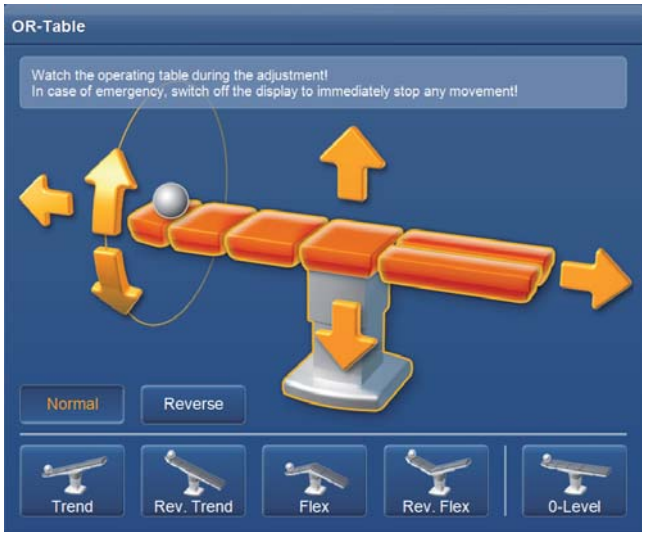NOTE: *Please refer to the instruction manual of the OR table manufacturer for the precise description of operation.*

NOTE: *An OR table is moved for as long as a button is pressed via the touch screen monitor. The movement ends automatically as soon as the touch screen button is released.*

NOTE: *Parallel to operation via the KARL STORZ SCB control NEO software the connected OR table can naturally also still be controlled with the OR table remote control or, if present, the control panel on the table pedestal. The OR table remote control or the control panel on the table pedestal have higher priority than the SCB system.*

NOTE: *For safety reasons, operating table functions cannot be activated using the mouse or keyboard.*

NOTE: *If the table moves unintentionally, activate the emergency stop (safety stop) button on the touch screen monitor.*

Activating functions

NOTE: *The figure on page 169 shows which positions/functions on the OR table can be set overall.*

General procedure

- Observe patient orientation* (Normal/Reverse)
- Select panel (element)
- Select “Movement” (not possible when choosing a predefined OR table setting via pictogram)
- Observe the reaction of the OR table – if necessary, correct setting or select new panel (element)

* If the patient is placed on the OR table “**reversed**” (head/foot area the other way around), the SCB control NEO software can take this into account via the “**Reverse**” button. The head position in the realistic view is always on the left half of the table and is thus independent of the currently set patient orientation.

NOTE: *If a table supports this patient orientation, then it will also reverse the commands accordingly if “Reverse” is set.*

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The selected element is displayed in “orange”.

ⓘ ***NOTE:*** *When choosing one of the five predefined OR table positions, the complete OR table is displayed in orange for as long as the button remains pressed- The setting is performed directly. The “ball” represents the position of the patient’s head.*

Arrows appear over the selected element showing the adjustment possibilities (up/down, left/right, tilt left/tilt right)

The display of the arrows is depending on the selection:

(triangle) for an OR table element:

(arrow) if the movement is enacted on the whole OR table.

ⓘ ***NOTE:*** *If the shoulder element is selected, the head element is marked at the same time, i.e., the two can only be moved together.*

Adjusting the OR table (element)

The adjustments are performed by touching the arrows on the touch screen monitor. The arrow then changes its appearance. The position of the OR table continues to change as long as the arrow remains pressed.

-> *and* ->

If the movement is no longer possible, the appearance also changes.

-> *and* ->

ⓘ ***NOTE:*** *The “Emergency stop (safety stop)*” function is described in the SCB/Universal OR table instruction manual.*

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
SCB control NEO displays the pseudorealistic view of the wall control panel for the CHROMOPHARE® X65 and D-series models.

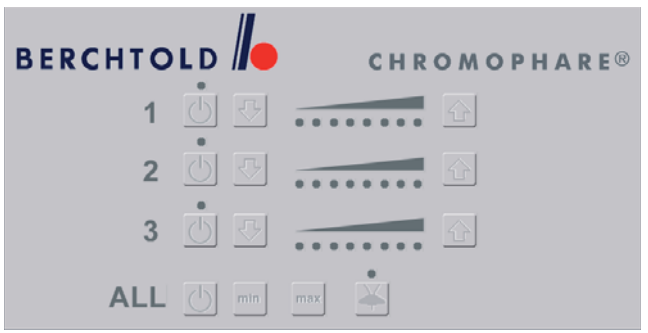The detailed view displays the light intensity and on position of the lights. You can also switch all lights on/off or switch to Endolite. Furthermore, all lights can be adjusted from minimum to maximum intensity. The detailed view varies according to the number of lights to be controlled and the availability of an Endolite.

NOTE: *With the D-series, the SCB control NEO software (in contrast to the wall control panel), displays the first light as switched on when Endolite is active. This corresponds to the light status. The intensity of the individual lights is also displayed.*

CHROMOPHARE® can be added to a Preset, settings, however, cannot be preset.

In the Preset the maximum configuration with three lights is displayed. During operation, the number of lights specified by the CHROMOPHARE® is displayed accordingly.

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
CHROMOPHARE® E-series

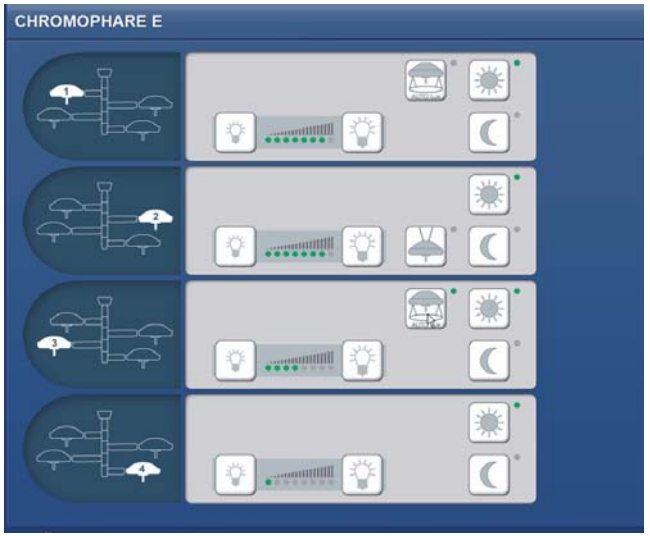Read the instruction manuals on the carrier device and interface module carefully. Using the interface module **20** 0922 90 up to four lights from the E-series (HID, BRite or LED) from Berchtold can be connected to the SCB control NEO system. In principle, all settings as per the Berchtold instruction manuals are to be carried out.


Operating OR lights

The realistic view of the Berchtold OR lights contains all buttons (and thus functions) of the individual control panels of the individual lights. Operation via the SCB control NEO software is identical to the original and can be looked up in the Berchtold instruction manual. The detailed view varies according to the number of lights to be controlled and the technology used.

NOTE: *Alongside operation with the SCB control NEO software, the Berchtold OR lights can, of course, still be controlled via the control panels of the individual illumination units. The control panels on the illumination units have a higher priority status than the SCB control NEO.*

NOTE: *To prevent the lights from being switched off unintentionally due to operation of the OFF or Sleep buttons, the button must be pressed for a longer period of time.*

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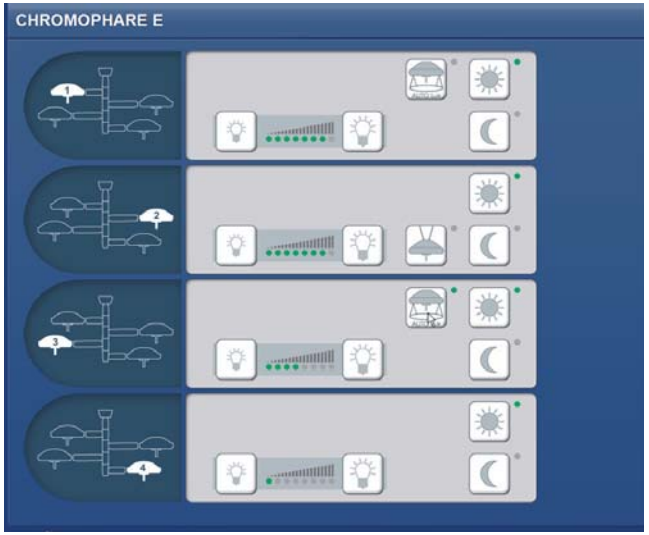HID

NOTE: *The operating panels on the lights and touch screen show the intensity in different ways during a change in the Sleep On/Off and Endo On/Off statuses. After the change, the statuses are again displayed in the same manner.*

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Device setup

The device setup dialog box allows the position and special functions of the individual lights to be set.

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If the positioning specified by the lighting system does not correspond with the individual lights on the support arm, these can be individually assigned here. The control panel of the corresponding light is moved to the desired position. A corresponding image is displayed next to the control panel.

Assignment of AutoLux

Expand the control panel to include the AutoLux option should this be supported by the light.

Assignment of OR camera ECO

Expand the control panel to include zoom functionality for the internal (ECO) camera should this be supported by the light.

Trumpf iLED

Read the instruction manuals on the carrier device and interface module carefully. Using the interface module **20**0924 90 up to 3 iLED (maximum 2 x iLED5 and 3x iLED3) lights and an internal camera (iLED5 functionality) from Trumpf can be connected to the SCB.

In principle, all settings as per the Trumpf instruction manuals are to be carried out.

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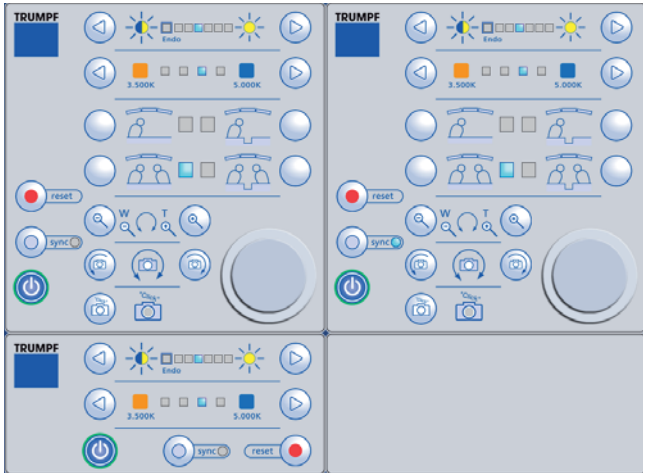
Operating OR lights

The pseudorealistic view of the iLED OR lighting system contains all the functions of both control panels (iLED3 and iLED5). Operation has been adapted to the touch screen and thus differs from the original. The individual functions are identical to the original and can be looked up in the Trumpf iLED instruction manual. The detailed view varies according to the number and type of the lights to be controlled.

OR light position display


For unique assignment of the lights and their positioning on the arm bracket a corresponding image is shown in the control panel. The arrangement is specified by the manufacturer.

Operating the internal OR camera

The internal camera is controlled via the iLED5 control panel available in the system. A distinction is not made as to the light in which the camera is located.

NOTE: *Alongside operation with the SCB control NEO software, the Trumpf iLED OR lighting system can, of course, still be controlled via the control panels of the individual illumination units. The control panels on the illumination units have a higher priority status than the SCB software.*

NOTE: *The remote interface prevents repeated switching off of the lights by disabling the ON/OFF button. A disabled ON/OFF button is shown in orange and is reactivated after 3 seconds.*

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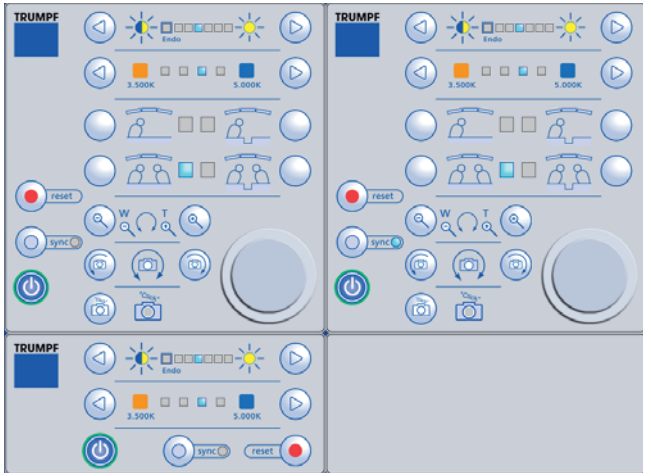KLS MARTIN marLED®

Read the instruction manuals on the carrier device and interface module carefully. Using the interface module **20**0923 90 up to four marLED® (V 10, V 16) lights from KLS MARTIN can be connected to the SCB.

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Operating OR lights

The realistic view of the marLED® OR lighting system contains all the buttons (and thus functions) of the control panel. Operation via SCB control NEO is identical to the original and can be looked up in the KLS MARTIN marLED® instruction manual.

The detailed view varies according to the number of lights to be controlled.

NOTE: *Alongside operation with the SCB control NEO software, the KLS MARTIN marLED® OR lighting system can, of course, still be controlled via the control panels of the individual illumination units. The control panels on the illumination units have a higher priority status than the SCB control NEO software.*

NOTE: *To prevent the lights from being switched off unintentionally due to operation of the OFF button, this button must be pressed for a longer period of time.*

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
Assignment of the OR light position

If the positioning specified by the lighting system does not correspond with the individual lights on the arm bracket, these can be individually assigned here. The control panel of the corresponding light is moved to the desired position.

Internal OR camera

Assignment of the internal camera and light is transferred when starting the system. Camera functionality can be selected in the control panel of this light.

KLS MARTIN marLux®

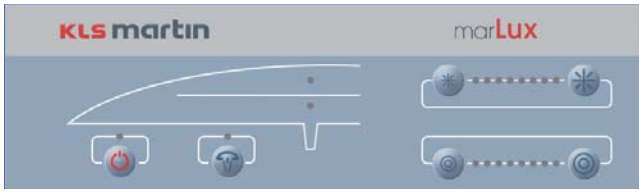Read the instruction manuals on the carrier device and interface module carefully. Using the interface module 20091490 up to four marLux® lights or a surgiCam*digital* from KLS MARTIN can be connected to the SCB.

In principle, all settings as per the KLS MARTIN instruction manual are to be carried out.

Operating OR lights

The realistic view of the marLux® OR lighting system contains all the buttons (and thus functions) of the control panel. Operation via SCB control NEO is identical to the original and can be looked up in the KLS MARTIN marLux® instruction manual. The detailed view varies according to the number of lights to be controlled.

NOTE: *Alongside operation with the SCB control NEO software, the KLS MARTIN marLux® OR lighting system can, of course, still be controlled via the control panels of the individual illumination units. The control panels on the illumination units have a higher priority status than the SCB control NEO software.*

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Universal OR Light/OR Light Cam

Using the interface module 20091790 up to four OR lights and/or an OR camera from different manufacturers can be connected to the KARL STORZ SCB control NEO system.

All settings must be made in accordance with the instruction manual supplied by the manufacturer(s) of the light and camera.

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Operating OR lights

The control panel of the universal OR light system contains functions for controlling the light intensity, switching the light on and off and (depending on which light is connected) controlling up to two special functions.

NOTE: *For safety reasons, when switching the light on and off, the ON/OFF button must be held down for at least 300 msec.*

NOTE: *For safety reasons, when switching all lights on and off, the ON/OFF button must be held down for at least 500 msec.*

NOTE: *The two special functions are preset by the light manufacturer. The detailed view varies according to the number of lights to be controlled.*

At the same time, functions for controlling all the lights, such as setting the light intensity to “low”, “medium” or “high”, switching the light on/off and the special functions are supported.

NOTE: *The intensity can only be controlled for the lights that are switched on.*

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The lights are switched off simultaneously if they were all on before. Otherwise all the lights are switched on. The two special functions operate in the same way. The function buttons for all the lights are displayed as soon as more than one light is controlled.

NOTE: *Alongside operation with the SCB control NEO software, the universal OR lighting system can, of course, still be controlled via the control panel of the individual illumination units. The control panels on the illumination units have a higher priority status than the SCB control NEO software.*

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Operating the OR camera

The control panel of the universal OR camera system contains functions for controlling enlargement, focusing, rotation and aperture, for freezing the camera image and for white balancing and, depending on the camera connected, control of up to two special functions.

NOTE: *Alongside operation with the SCB control NEO software, the connected OR camera system can, of course, still be controlled with the system on the manufacturer's control panel. The control panel on the device has a higher priority status than the SCB control NEO software.*

NOTE: *The two special functions are preset by the camera manufacturer.*

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Operating the OR light camera

The control panel combines the universal OR light system and the universal OR camera system. These systems are described in the sections “Operating OR lights” and “Operating the OR camera”.

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Operating the room control

The control panel for room control can be used in conjunction with the “Universal OR light system”. It allows the user to choose between two illumination scenarios if a KARL STORZ media unit (SCB ACC Control®) is connected which supports them.

The function and labeling of these buttons in the room control are specified in the Device Setup dialog box.

NOTE: *The buttons of the room control are activated as soon as the relevant media unit is also listed in the Preset. Otherwise the buttons are inactive, as indicated by a grey background.*

The selected media unit determines the status of the buttons.

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NOTE: *The buttons of the "All" line are only active if at least one connected light is switched on.*

The "universal OR light", the "universal OR camera" and the "universal OR light camera" can be added to a Preset, settings cannot be preset.

Device setup

WARNING: *Settings must only be carried out by skilled technical staff or KARL STORZ personnel.*

To call up the dialog box, click on "Device Setup" in the "Options" menu.


The settings made here influence the display of the corresponding control panel.

Press the "OK" button to accept the changes made, or press the "Cancel" button to reject them.

When you change between the individual tabs, a plausibility check is performed on the settings.

The green rectangle at the left end of the button indicates that the corresponding function is displayed on the control panel.

Device Setup – Camera

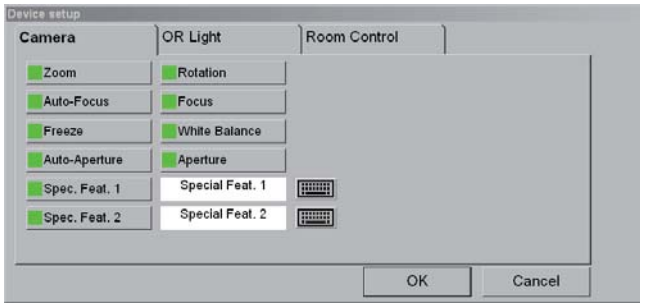In this dialog box the camera-specific settings can be made (if these are supported by the connected camera).

NOTE: *This tab is only displayed with the OR camera and OR light camera.*

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Zoom:

- Range selection:
 ON = “Zoom” buttons are displayed,
 OFF= “Zoom” buttons are not displayed.
- When the button is pressed, the “Zoom” buttons are displayed or hidden on the control panel.

Rotation:

- Range selection:
 ON = “Rotation” buttons are displayed,
 OFF= “Rotation” buttons are not displayed.
- When the button is pressed, the “Rotation” buttons are displayed or hidden on the control panel.

Autofocus:

- Range selection:
 ON = “Autofocus” button is displayed,
 OFF= “Autofocus” button is not displayed.
- When the button is pressed, the “Autofocus” button is displayed or hidden on the control panel.

Focus:

- Range selection:
 ON = “Focus” buttons are displayed,
 OFF= “Focus” buttons are not displayed.
- When the button is pressed, the manual “Focus” buttons are displayed or hidden on the control panel.

Freeze:

- Range selection:
 ON = “Freeze” button is displayed,
 OFF= “Freeze” button is not displayed.
- When the button is pressed, the “Freeze” button is displayed or hidden on the control panel.

White balance:

- Range selection:
 ON = “White Balance” button is displayed,
 OFF= “White Balance” button is not displayed.
- When the button is pressed, the “White Balance” button is displayed or hidden on the control panel.

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Auto Aperture:
- Range selection:
 ON = “Auto Aperture” button is displayed,
 OFF= “Auto Aperture” button is not displayed.
- When the button is pressed, the “Auto Aperture” button is displayed or hidden on the control panel.

Aperture:
- Range selection:
 ON = “Aperture” buttons are displayed,
 OFF= “Aperture” buttons are not displayed.
- When the button is pressed, the manual “Aperture” buttons are displayed or hidden on the control panel.

Special feat. 1:
- Range selection:
 ON = “Special feat. 1” button is displayed,
 OFF= “Special feat. 1” button is not displayed.
- When the button is pressed, the first “Special feat.” button is displayed or hidden on the control panel.

Special feat. 2:
- Range selection:
 ON = “Special feat. 2” button is displayed,
 OFF= “Special feat. 2” button is not displayed.
- When the button is pressed, the second “Special feat.” button is displayed or hidden on the control panel.

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
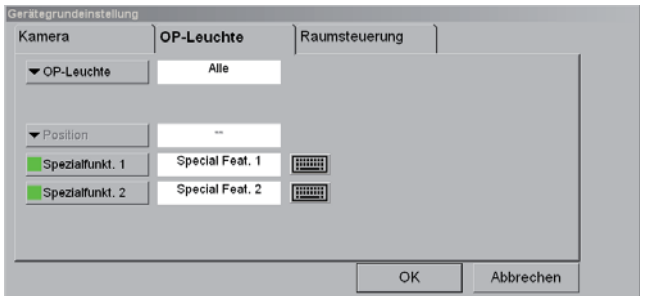
Device Setup – Light

In this dialog box the light-specific settings, e.g., the positioning of the individual lights or the special features (if these are supported by the connected light) can be made.

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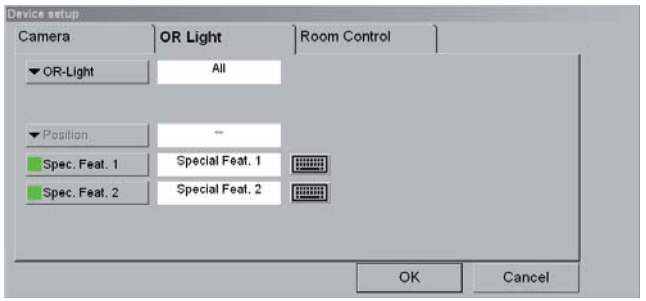OR light:

- Range selection: All, 1, 2, 3, 4
- When the button is pressed, the device changes to the light for which the position and text input are enabled.
 Once the last connected light is reached, the setting begins again at "All". If only one light is connected, this button remains deactivated.

Position:

- Range selection: 1, 2, 3, 4
- When the button is pressed, the position of the selected light on the holding arm is set. Once the last position is reached, it returns to "1".

NOTE: *The positions on the holding arm are numbered in ascending order from bottom to top. The positions must be uniquely assigned.*

Special feat. 1:

- Range selection:
 ON = "Special feat. 1" button is displayed,
 OFF= "Special feat. 1" button is not displayed.
- When the button is pressed, the first "Special feat." button is displayed or hidden in the control panel.

Special feat. 2:

- Range selection:
 ON = "Special feat. 2" button is displayed,
 OFF= "Special feat. 2" button is not displayed.
- When the button is pressed, the second "Special feat." button is displayed or hidden in the control panel.

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
Device Setup – Room control

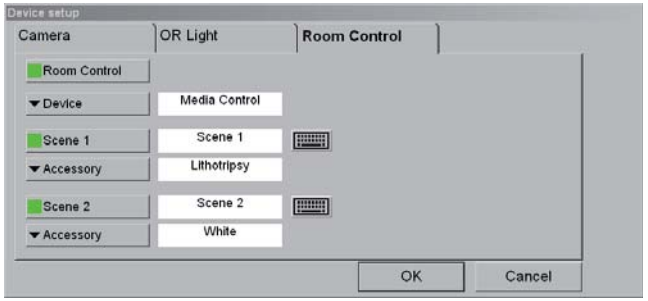NOTE: *This tab is only displayed with the "OR Light" and "OR Light Cam" devices.*

Room control:
- Range selection:
 ON = room control is displayed,
 OFF = room control is not displayed
- When the button is pressed, the "Room control" dialog box is displayed or hidden in the control panel.

Device:
- Range selection: SCB ACC Control®
- When the button is pressed, the media device that is to be used for room control is selected.

Scene 1:
- Range selection:
 ON = "Scene 1" button is displayed,
 OFF= "Scene 1" button is not displayed.
- When the button is pressed, the first "Scene" button is displayed or hidden in the control panel.

Scene 2:
- Range selection:
 ON = "Scene 2" button is displayed,
 OFF= "Scene 2" button is not displayed.
- When the button is pressed, the "Scene 2" button is displayed or hidden in the control panel.

Accessories:
- Range selection:
 Dependent on the labeling of the accessory buttons of the selected media device.
- When the "Accessory" button is pressed, the system changes to the next selectable accessory channel of the selected media device.

NOTE: *The two scenes must not both indicate the same accessories.*

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NOTE: *Depending on the manufacturer, the priorities of the device messages can deviate from the priorities in the list below.*

OR light/OR light cam (additional*)/OR camera (only) device messages**

Message	Priority
Please note technical alarm at device (light 1…4)	low
*/**Please note technical alarm at device (camera 1…4)	low
Please note technical alarm at device (light 1…4)	middle
*/**Please note technical alarm at device (camera 1…4)	middle
Please note technical alarm at device (light 1…4)	high
*/**Please note technical alarm at device (camera 1…4)	high
Please note physiological alarm at device (light 1…4)	low
*/**Please note physiological alarm at device (camera)	low
Please note physiological alarm at device (light 1…4)	middle
*/**Please note physiological alarm at device (camera)	middle
Please note physiological alarm at device (light 1…4)	high
*/**Please note physiological alarm at device (camera)	high
Control temporarily locked for safety reasons (light 1…4)	low
*/**Control temporarily locked for safety reasons (camera)	low
Control temporarily locked for safety reasons (light 1…4)	middle
*/**Control temporarily locked for safety reasons (camera)	middle
Control temporarily locked for safety reasons (light 1…4)	high
*/**Control temporarily locked for safety reasons (camera)	high

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Harmony LA OR light

Read the instruction manuals for the interface device/carrier device and interface module carefully.

Using the interface device/interface module 20 091632/20 090690 the Harmony LA OR lighting system from STERIS® can be connected to the SCB control NEO software.

In principle, all settings as per the STERIS® instruction manuals are to be carried out. However, the SCB control NEO software does not support operation of specific device settings, e.g., the camera menus 2, 3 and 4 as well as the system recording error display.

These functions must still be carried out directly at the device. With SCB control NEO software operation, these functions cannot be called up and as such are also not shown in the respective dialog boxes.

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Operating OR lights

The pseudorealistic view of the Harmony LA OR lighting system contains all buttons (and thus functions) of the real STERIS® Harmony LA Control Center which can be operated via the SCB control NEO software.

Button assignment is identical to the real Control Center of the Harmony LA OR lighting system.

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Operating the OR camera

Button assignment is identical to the real Control Center of the Harmony LA OR lighting system. A few of the lesser used camera functions of the STERIS® Harmony LA Control Center are not supported.

NOTE: *Please refer to the STERIS® instruction manual for the respective detailed description.*

NOTE: *Alongside operation with the SCB control NEO software the STERIS® Harmony LA OR lighting system can, of course, still be operated with the STERIS® Harmony LA Control Center or auxiliary control panel of the individual illumination units. The STERIS® Harmony LA Control Center and auxiliary control panels on the illumination units have a higher priority status than the SCB control NEO system.*

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The pseudorealistic view of the OR camera contains all the buttons (and therefore the functions) of the real BERCHTOLD control device which can be carried out via the SCB control NEO system. The button assignment is identical to the real control device.

NOTE: *If the BERCHTOLD ChromoVision® device is operated directly on the front panel, it cannot be operated via the pseudorealistic view at the same time. This is indicated by flashing LEDs in the pseudorealistic view.*

NOTE: *If the BERCHTOLD ChromoVision® device is operated via the pseudorealistic view, it cannot be operated on the front panel at the same time.*

NOTE: *If the OR camera control device is in the standby mode, all buttons except the OFF button are disabled and cannot be selected. It is not fully operable until the control device is switched to the operating mode by clicking on the OFF button.*

NOTE: *There is no correlation between the standby mode of the BERCHTOLD ChromoVision® device and the standby function of the SCB control NEO software. Also see 4.4.3.1 Device status display.*

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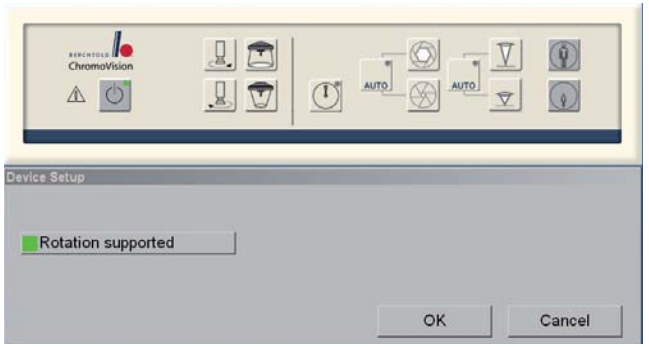
Camera rotation support by the OR camera control device or the OR camera is only optional.

The pseudorealistic view, however, supports both those cameras with motorized rotation function and those without.

If the control device or camera is to support the rotation function, this must be configured in the pseudorealistic view.


To do so, proceed as follows:

- Select the menu item “Options/Device Setup” for the ChromoVision® device activated in the device list.
- The dialog box shown below is displayed.
- Select “Rotation supported” if the OR camera control device supports motorized camera rotation.
- Exit the dialog box by clicking on “OK” if you wish to apply the changes you made.

The image shows the pseudorealistic view without the function buttons for motorized camera rotation.

ChromoVision® 1C/3C
OR camera
Read the instruction manuals on the carrier device and interface module carefully.

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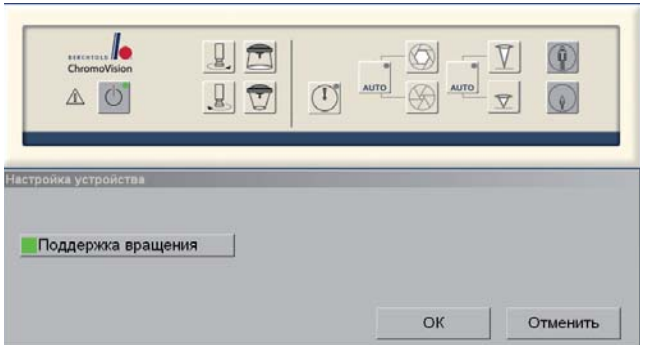
Sola camera

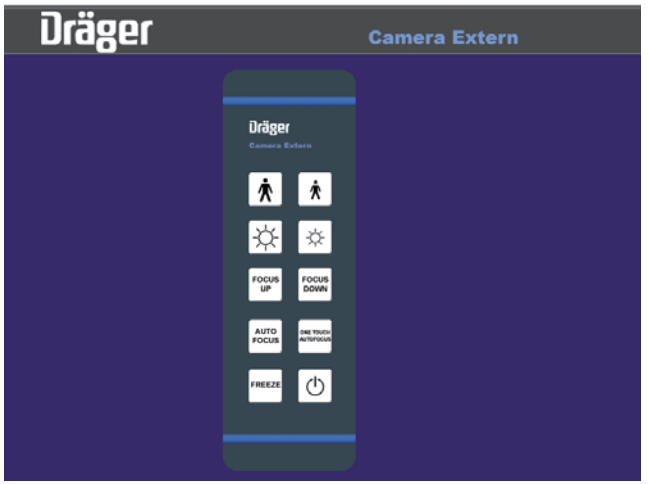Read the instruction manuals on the carrier device and interface module carefully. The camera is operated with the accompanying remote control.

NOTE: *The Sola camera is shown in the device list as soon as the respective SCB/Dräger Sola camera module is supplied with power. This also applies if there is no serial connection between the Sola camera and SCB/Dräger Sola camera module which means that control is not possible via the SCB control NEO software.*

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surgiCam*digital*

Read the instruction manuals on the carrier device and interface module carefully. Using the interface module **20**0914 90 up to four marLux® lights or a surgiCam*digital* from KLS MARTIN can be connected to the SCB.

In principle, all settings as per the KLS MARTIN instruction manuals are to be carried out.

Operating the OR camera

The realistic view of the front panel of the surgiCam*digital* OR camera system contains all buttons (and thus functions) of the KLS MARTIN surgiCam*digital* front panel which can be implemented via the SCB control NEO system.

Button assignment is identical to the real front panel. The auto-exposure, freeze and on-screen display functions are not supported. These functions must still be carried out directly on the device.

NOTE: *In accordance with the surgiCam*digital* instruction manual, the camera functions (zoom and rotate) may be executed via the front panel. This mode is not displayed in the SCB control NEO software.*

NOTE: *Alongside operation with the SCB control NEO software, the KLS MARTIN surgiCam*digital* OR camera system can, of course, still be controlled via the front panel of the KLS MARTIN surgiCam*digital*. This means the "Original" has a higher priority than the SCB control NEO software.*

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4.5.13 Videoprinters and videorecorders

Printer UP 21
Printer UP 51
VCR DSR 20
VCR SVO 9500

Read the instruction manuals for the interface device/carrier device and interface module carefully*.

*The special SCB control NEO device functions are described in the instruction manual of the device/module.

SONY video devices – General

The realistic view of the SONY devices supports the most important functions. Buttons whose functions are not supported are displayed as inactive (light gray). Functions for which there are no buttons on the unit and which are only accessible via the device's own menu are operated in the SCB control NEO software via "Options/Device Setup". This menu item is only available if menu functions are supported.

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
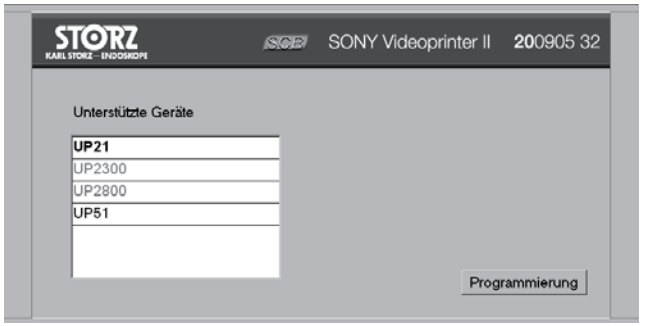
SONY Videoprinter II

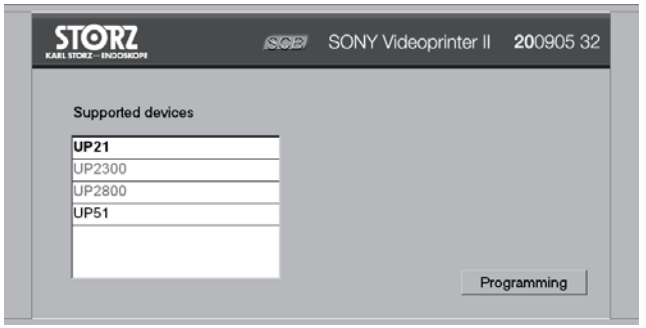The SONY videoprinter II is a virtual device which is used for programming the interface device of the same name. After the interface device is switched on, it sends the device models which it can control to the OR1™ control NEO. This transfer may take some time and is shown in a window with a progress bar. The SONY videoprinter II displays this in a list window.


NOTE: *If your device is not shown in the list, you need a newer version of the Interface Box. Contact KARL STORZ Service.*

To program the interface device, select the connected device and actuate the button “Programming”. After switching the interface device off and on (again), the device registers on the control computer.

NOTE: *If you cannot select your device and it is displayed in light gray, start the Maintenance Manager (see p. 68) and add the device to the device list.*

WARNING: *Take care to program the model you have connected exactly otherwise an incorrect view will be displayed and the device may change to an error status.*

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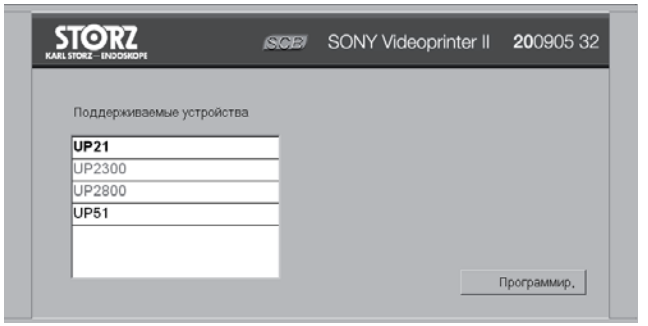If the SONY videoprinter II interface device has been programed incorrectly or the device to be controlled has been replaced by another device you must reprogram the interface device. This involves deleting the interface device programing. For this, each device that can communicate with the SONY videoprinter II interface device has a dialog box which can be opened under “Options/Device Setup”. The “Delete Device” button is used to delete the programing. After switching the interface device on and off, the device registers as SONY videoprinter II and can be programed as described above.

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4.5.14 Data logger

Function

The data logger is a “virtual device” and a component of the OR1™ control NEO software. It permits the recording and visualization of data from selected SCB devices. This data is stored, can be presented graphically on the screen and exported to a data medium or the KARL STORZ AIDA® system.

Data recording

The data logger makes it possible to record all the data which you have selected for data recording during creation of the Preset. Data is therefore only acquired when a Preset to which the data logger was added is open.

From each device which is contained in the Preset, error messages, alarms, warnings and statuses are recorded as standard. In addition, you can select function blocks for individual devices (comprising actual and/or set values) for data recording.

Over an acquisition period of 5 seconds the minimum and maximum values are determined for the individual actual and set values, the average is calculated and stored together with the errors, alarms and warnings. The total recording time is limited to 24 hours.

Each recording is stored in a file on the hard disk of the OR1™ control NEO. The data logger can manage up to 10 separate recordings.

ⓘ ***NOTE:*** *No data acquisition takes place in the Preset “All devices”.*

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
Configuration of data recording

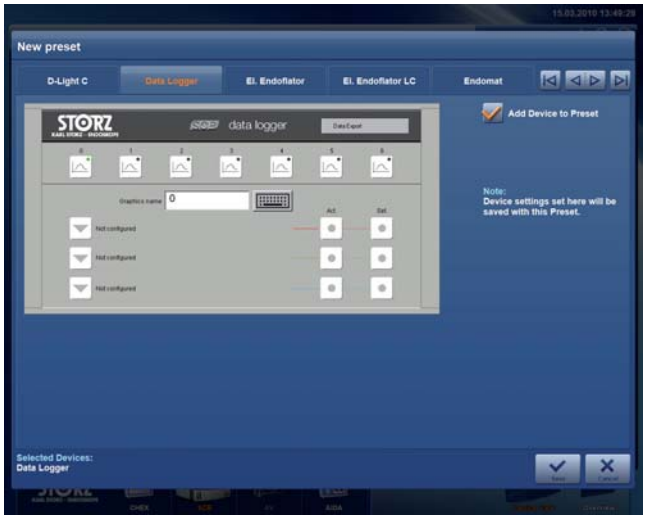You can add the data logger like an SCB device to each Preset (see “Preset - New”). You can then select function blocks for individual devices (comprising actual and/or set values) for data recording.

For each device representation there is a “Data logging” window which permits selection of the supported function blocks via checkboxes. By clicking on the respective check boxes the function blocks for data recording are added.

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NOTE: *You can select a total of up to 12 different function blocks for data recording.*

After selecting the function blocks for data recording, these trend graphs can be assigned. The trend graphs can be displayed on the screen online during a recording.

Up to seven different trend graphs are supported. Each of these trend graphs can contain up to three function blocks (comprising actual and/or set values). However, these must have identical units.

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
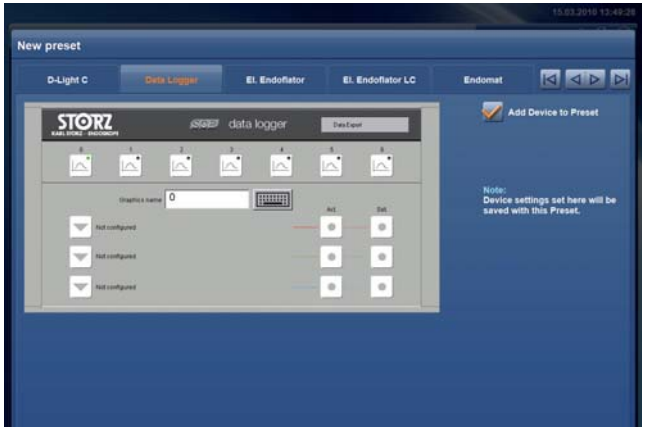
To configure a trend graph, click on the corresponding graphic button. In the input field you can assign a brief description to the trend graph.

By clicking on the first selection button several times you can then select one of the configured function blocks. This block is displayed with the name of the SCB device and the name and unit of the function block. With the aid of the two other selection buttons you can now add two further function blocks. However, you will only be shown those which have the same unit.

Generally speaking, a function block comprises an actual and a set value. By clicking on the assigned actual and set value buttons, you can set whether they are to be shown on the trend graph or not. By default, the actual and set value are always shown if they exist.

Assignment of the colors and type of line for the various function blocks is preset.

You can remove function blocks assigned to a trend graphic by clicking on the corresponding selection button until the text “Not configured” indicates that the entry has been deleted.

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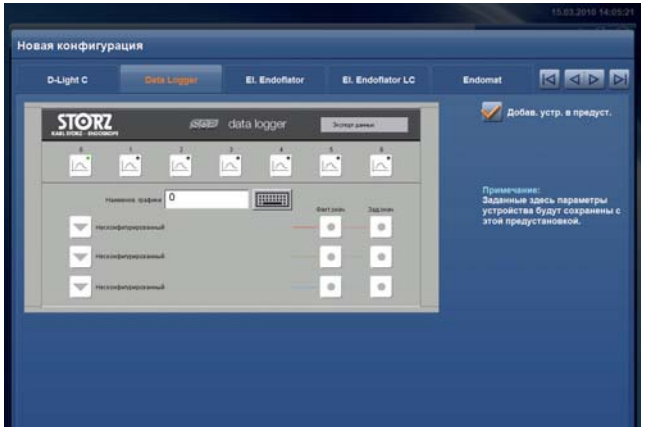Recording data

When the Preset is open you can start a new recording by clicking on the button “New”.


This automatically selects the oldest recording (out of 10). In a dialog box you will be shown the recording information (date and time, surname, first name, date of birth and patient ID) which is assigned to that recording. By clicking on the “OK” button you confirm that this recording will be overwritten.

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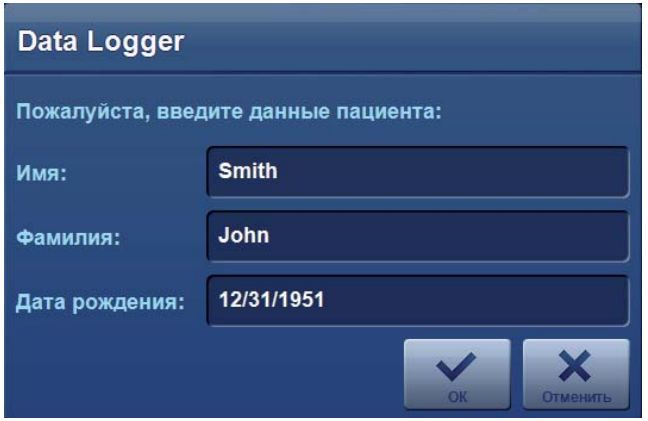If, on the other hand, you click on the “Cancel” button, no new recording will be started and the data of the old recording will remain intact.

NOTE: *If you confirm that the recording should be overwritten, the data of the old recording will be lost irretrievably. For this reason you should save the various recordings after 10 recordings at the latest by exporting them to a USB memory stick (see p. 205 “Export recording”).*

NOTE: *The total recording time is limited to 24 hours.*

By clicking on the “OK” button you confirm the entry and the new recording will be started. The new patient information will now be shown in the data logger window.

If, on the other hand, you click on the “Cancel” button, no new recording will be started and the data of the old recording will remain intact.

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Display trend graph

After starting a recording the first trend graph is displayed automatically (if it was configured). If you wish to display a different trend graph, click on the appropriate graph button.

Whilst a recording is running, the last 45 minutes will always be displayed in the trend graph. Every 5 seconds the graphic is updated with new averages. When the right-hand edge of the graphic has been reached, the graphic will shift 15 minutes.

In the legend the minimum and maximum actual values (over the entire recording) are shown. If you have only added the set value to the trend graphic, the minimum and maximum set values will be given instead.

After selection, a trend graph is always presented in the maximum range of the relevant function blocks (actual and set values). By clicking on one of the legend fields, you can scale the trend graph according to the displayed minimum and maximum values.

Click on the trend graph to return to default scaling.

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Stop recording

By clicking on the “Close” button you stop the current recording.

Continue recording

You can continue the last terminated recording if it was terminated less than an hour ago. By contrast with “New” the new data is added to the recording and existing data remains intact.

Select recording

If no recording is running, you can display the patient information of all 10 recordings one after the other by pressing the “File” button. You can then display the recording just selected again by pressing the “Open” button or overwrite it with a new recording with the “New” button.

Open recording

You can display the recording just selected by pressing the “Open” button. The graphics now show the entire course of the operation. You now export the recording either to a USB memory stick or transfer it to the KARL STORZ AIDA® system.

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
Export recording

Export to the USB stick

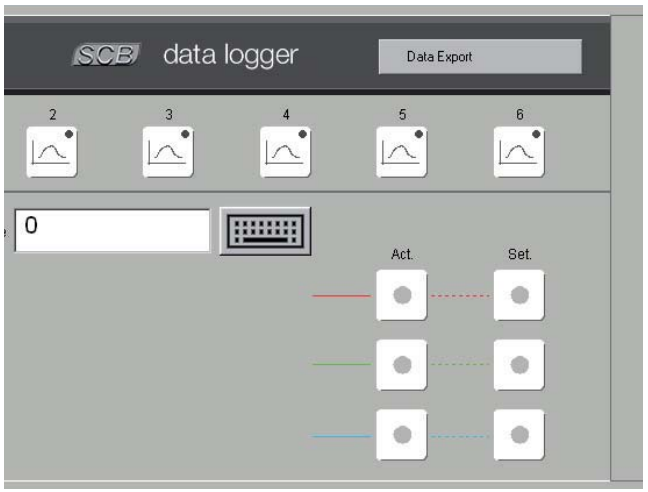When you have opened a recording, you can save it on a USB memory stick with the “Data Export” button. The data is stored there in a readable ASCII file.


Export via the network

The data logger is capable of transferring the recorded data in realtime via TCP/IP to a “third-party system”. Processing this data is not usually a basic function of this “third-party system” thus making adaptation and/or development necessary.

A separate installation manual contains a precise description of this functionality and shows how to carry out this kind of adaptation. Please contact KARL STORZ about this matter.

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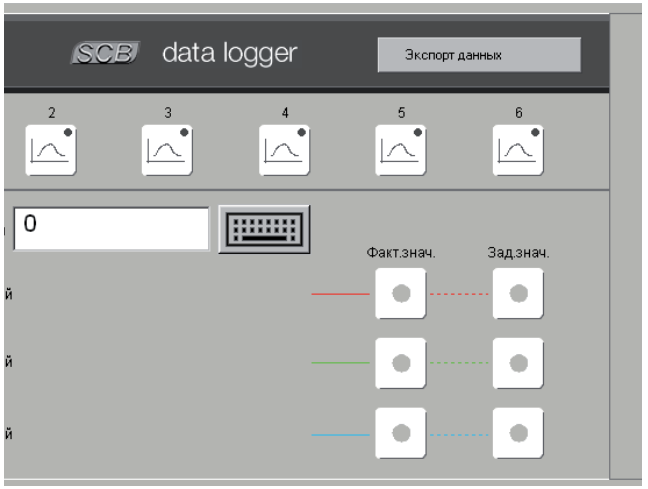4.6 Information on operation

4.6.1 KARL STORZ SCB control NEO

- If devices are removed via “Add/Remove Devices” these are also deleted from the corresponding Preset. These are not automatically added to the Preset via “Add” but have to be added manually to all Presets.
- If a password has been assigned to change Presets and the keyboard layout is subsequently changed, the password cannot be changed and, accordingly, the system cannot be unlocked.

4.6.2 OR1 NEO™

- If AIDA or Home is switched to via the main menu, but the AIDA Interface II switches back automatically to SCB control NEO (due to a lack of a keyboard connection between AIDA control NEO or KARL STORZ OR1™ FUSION <--> AIDA Interface II) no further inputs to the SCB control NEO application are accepted for safety reasons. The missing keyboard connection is to be replaced or the OR1™ control NEO is to be restarted.

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5 Maintenance

5.1 Cleaning, disinfection and care

CAUTION: *Avoid allowing fluids to enter the housing of the KARL STORZ SCB control NEO system components.*

Wipe-clean the exterior surfaces of the KARL STORZ SCB control NEO system components with a low-lint disposable cloth moistened with disinfectant. For wipe-down disinfection, we recommend using only surface disinfectants based on quaternary ammonium compounds (no solutions predominantly containing alcohol such as rapid-action disinfectants) and “disposable disinfectant methods” based on ammonium chloride. Observe all instruction manuals included with the system, e.g., for the monitors.

5.2 Servicing and repair

Faulty components of the KARL STORZ SCB control NEO system must be serviced and repaired exclusively by persons authorized by KARL STORZ; only original KARL STORZ parts may be used for repair work.

5.3 Disposal

The OR1™ control NEO has been marked in accordance with the European Directive 2002/96/EC on waste electrical and electronic equipment (WEEE). At the end of its useful operating life, dispose of the unit as electronic scrap.

Please ask either KARL STORZ SE & Co. KG, a KARL STORZ subsidiary or your specialist dealer for information on your local collection point. Within the scope of application of this Directive, KARL STORZ SE & Co. KG is responsible for the proper disposal of this unit.

According to the “Chinese Administration on the Control of Pollution caused by Electronic Information Products” (ACPEIP) this KARL STORZ product contains none of the six substances above the MCV (maximum concentration value).

That means that this device is compliant according to the Chinese regulations.

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5.4 Repair program

Usually to bridge the repair period, you will receive a unit on loan which you then return to KARL STORZ as soon as you receive the repaired unit.
In Germany you can refer repairs direct to:

KARL STORZ SE & Co. KG
Abt. Reparaturservice
Dr. Karl-Storz-Str. 34
78532 Tuttlingen, GERMANY

In other countries please contact your local KARL STORZ branch or authorized dealer.

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Important information

In order to prevent infection, it is strictly forbidden to ship contaminated medical devices. All medical devices must be decontaminated on site, to avoid contact and aerogenous infections (among staff). We reserve the right to return contaminated instruments/equipment to the sender.

Repairs, changes or expansions which are not carried out by KARL STORZ or by experts authorized by KARL STORZ will invalidate all guarantee rights. KARL STORZ gives no guarantee on the correct functioning of equipment or instruments which have been repaired by unauthorized third parties.

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5.5 Limitation of liability

KARL STORZ shall be liable for the safe operation, operational reliability, and performance of the KARL STORZ SCB control NEO system only subject to the following conditions:

- all assembly operations, system expansions, readjustments, modifications, or repairs have been performed by a person or persons duly authorized by KARL STORZ;
- all electrical installations in the room where the system is connected and operated complies with the applicable laws and standards, and
- the units have been used in accordance with their instruction manuals at all time.

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KARL STORZ SE & Co. KG
Postfach 230
78503 TUTTLINGEN

5.6 Manufacturer's warranty

For a period of two years after delivery to the end-user, we shall replace free of charge equipment that can be proven to have defects in material or workmanship. In so doing we cannot bear the cost of transportation or the risks associated with shipment. The warranty referred to in our Standard Conditions of Business shall apply. Please fill out the attached warranty card on the last page and return it as soon as possible to:

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If the equipment is opened or any repairs or modifications to the components of the OR1™ control NEO are performed by unauthorized persons, we accept no liability for its operational safety. Any such opening, repair, or modification performed during the warranty period shall void all warranty.

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Property rights

Protected in the USA by (at least one of the following) US Patent No(s) 5,788,688; 6,397,286; 6,484,221; 6,824,539.

5.7 Software User Rights

1. With respect to KARL STORZ products containing software components, the Customer has a non-exclusive License, which is limited by the meaning of the following terms and conditions, to use the programmed logic, computer programs and/or software supplied by KARL STORZ and incorporated in KARL STORZ products (“Software”).

2. This License applies for such products sold to the Customer only to the form delivered to Customer and for the sole purpose of operating such products in accordance with written instructions by KARL STORZ. Use of the Software in connection with other products bought by the Customer from KARL STORZ or from third parties is prohibited.

3. A transfer of this License by the Customer is permissible only if the Customer provides the buyer’s name to KARL STORZ, and said buyer agrees to comply with the terms and conditions regarding the usage of the software as set forth in this License. In the case of a transfer, the Customer’s License shall terminate.

4. All copyrights to the Software, as well as all modifications, enhancements, and upgrades thereto, will remain the property of KARL STORZ. In particular, the Customer may not duplicate, copy, reverse-engineer, decompile, or disassemble the Software or in any way modify the same. The Customer has no right to create derivatives of the software, and the Customer may not copy, decipher or re-create the source code of the Software.
§ 69g paragraph 2 of the Copyright remains unaffected by this. [This refers to German Copyright Law]

5. Any modification or enhancement of the Software by the Customer in violation of this Agreement is considered a violation of this License; KARL STORZ immediately shall acquire all property rights and copyrights to such modifications or enhancements. The Customer hereby in advance transfers an exclusive, free, and unlimited license to use such modifications and enhancements to KARL STORZ.

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6. KARL STORZ does not recommend use of the Software in connection with other software not supplied by KARL STORZ for the respective products. KARL STORZ is not liable for software malfunctions resulting from the use of software not supplied or recommended by KARL STORZ in connection with products supplied by KARL STORZ.

7. In the event of a failure of the Customer, or its agents, employees or representatives, to comply with any terms and conditions of the License herein granted, the License will, without termination or any further action by KARL STORZ or any other party, immediately terminate.

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5.8 Directive compliance (for **20**0900 01S431)

According to Medical Device Directive (MDD):

This medical product belongs to Class IIb

CE 0123 This medical product bears the CE mark in accordance with MDD 93/42/EEC.

NOTE: *The code number after the CE mark indicates the responsible notified body.*

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5.9 Directive compliance (for **20**0971 20-1)

According to Medical Device Directive (MDD):

This medical product belongs to Class I

CE This medical product bears the CE mark in accordance with MDD 93/42/EEC.

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WARRANTY

For a period of two years after delivery to the end-user, we shall replace free of charge equipment that can be proven to have defects in material or workmanship.
In doing so we cannot bear the cost of transportation or the risks associated with shipment. The warranty referred to in our Standard Conditions of Business shall apply.

Complete warranty card upon purchase/ delivery and return immediately to:

KARL STORZ SE & Co. KG
Postfach 230
78503 TUTTLINGENGERMANY

To be completed by the supplier/importer:

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To be filled out by instrument owner:

Return address/company stamp:

Field of application:

Type of instrument: Serial no.:

Purchase date:

Signature/Date:

Please attach sufficient postage

REPLY CARD

KARL STORZ SE & Co. KG
Postfach 230
78503 TUTTLINGEN
GERMANY

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Absender/Firmenstempel:

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Geräte-Typ: Serien-Nr.:

Kauf-Dat.:

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ANTWORTKARTE

KARL STORZ SE & Co. KG
Postfach 230
78503 TUTTLINGEN

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KARL STORZ SE & Co. KG
Dr.-Karl-Storz-Strasse 34, 78532 Tuttlingen,
Germany Postfach 230, 78503 Tuttlingen,
Germany Phone: +49 7461 708-0, Fax: +49 7461 708-105 E-Mail: info@karlstorz.de

KARL STORZ Endoskope Berlin GmbH
Ohlauer Straße 43
10999 Berlin, Germany
Phone: +49 30 30 69090, Fax: +49 30 30 19452

KARL STORZ Endoscopy Canada Ltd.
7171 Millcreek Drive, Mississauga
Ontario, L5N 3R3, Canada
Phone: +1 905 816-4500, Fax: +1 905 858-4599
E-Mail: info@karlstorz.ca

KARL STORZ Endoscopy-America, Inc.
2151 East Grand Avenue
El Segundo, CA 90245-5017, USA
Phone: +1 424 218-8100,
Toll free: 800 421-0837 (USA only)
Fax: +1 424 218-8526
Toll free: 800 321-1304 (USA only)
E-Mail: info@ksea.com

KARL STORZ Veterinary Endoscopy America, Inc.
175 Cremona Drive
Goleta, CA 93117, USA
Phone: +1 805 968-7776, Fax: +1 805 685-2588
E-Mail: info@karlstorzvet.com

KARL STORZ Endoscopia Latino-America, Inc.
815 N. W. 57th Avenue, Suite 480
Miami, FL 33126-2042, USA
Phone: +1 305 262-8980, Fax: +1 305 262-8986
E-Mail: info@ksela.com

KARL STORZ Endoscopia
Miramar Trade Center
Edificio Jerusalem, Oficina 108,
La Habana, Cuba
Phone: +53 72041097, Fax: +53 72041098

KARL STORZ Endoscopia México S.A. de C.V
Lago Constanza No 326
Col. Chapultepec Morales
D.F.C.P. 11520, México, México
Phone: +52 55 5250-5607, Fax: +52 55 55450174

KARL STORZ Marketing América Do Sul Ltda.
Rua Joaquim Floriano, nº. 413, 20º andar – Itaim Bibi,
CEP-04534-011 São Paulo, Brasil
Phone: +55 11 3526-4600, Fax: +55 11 3526-4680
E-Mail: info@karlstorz.com.br

KARL STORZ Endoscopia Argentina S.A.
Zufriategui 627 6° Piso
B1638 CAA - Vicente Lopez
Provincia de Buenos Aires, Argentina
Phone: +54 11 4718 0919, Fax: +54 11 4718 2773
E-Mail: info@karlstorz.com.ar

KARL STORZ Endoskopi Norge AS
P.O. Box 153, Rolf Olsenvei 28
N-2007 Kjeller, Norway
Phone: +47 6380 5600, Fax: +47 6380 5601
post@karlstorz.no

Subsidiaries

KARL STORZ Endoskop Sverige AB
Storsätragränd 14,
12724 Skärholmen, Sweden
Postal address: Po Box 8013,
14108 Kungens Kurva, Sweden
Phone: +46 8 50 56 4800, Fax: +46 8 50 56 4848
E-Mail: kundservice@karlstorz.se

KARL STORZ Endoscopy Suomi OY
Valimotie 1 B, 4. kerros, 00380 Helsinki, Finland
Phone: +358 96824774, Fax: +358 968247755
E-Mail: asiakaspalvelu@karlstorz.fi

KARL STORZ Endoskopi Danmark A/S
Skovlytoften 33
2840 Holte, Denmark
Phone: +45 45162600, Fax: +45 45162609
E-Mail: marketing@karlstorz.dk

KARL STORZ Endoscopy (UK) Ltd.
392 Edinburgh Avenue, Slough
Berkshire, SL1 4UF, Great Britain
Phone: +44 1753 503500, Fax: +44 1753 578124
E-Mail: customerservice@karlstorz-uk.com

KARL STORZ Endoscopie Nederland B. V.
Phone: +31 651 938 738, Fax: +31 135 302 231

KARL STORZ Endoscopy Belgium N. V.
Phone: +32 473 810 451, E-Mail: info@stopler.be

KARL STORZ Endoscopie France S. A. S.
12, rue Georges Guynemer, Quartier de l'Europe
78280 Guyancourt, France
Phone: +33 1 30484200, Fax: +33 1 30484201
E-Mail: marketing@karlstorz.fr

KARL STORZ Endoskop Austria GmbH
Landstraßer Hauptstr. 148/1/G1
1030 Wien, Austria
Phone: +43 1 71 56 0470, Fax: +43 1 71 56 0479
E-Mail: storz-austria@karlstorz.at

KARL STORZ Endoscopia Ibérica S. A.
Parque Empresarial San Fernando
Edificio Munich – Planta Baja
28830 Madrid, Spain
Phone: +34 91 6771051, Fax: +34 91 6772981
E-Mail: marketing@karlstorz.es

KARL STORZ Endoscopia Italia S. r. l.
Via dell'Artigianato, 3
37135 Verona, Italy
Phone: +39 045 8222000, Fax: +39 045 8222001
E-Mail: info@karlstorz.it

KARL STORZ Adria Eos d.o.o.
Zadarska 80
10000 Zagreb, Croatia
Phone: +385 1 6406 070, Fax: +385 1 6406 077
E-Mail: info@karlstorz.hr

KARL STORZ Endoskopija d.o.o.
Verovškova c. 60A
1000 Ljubljana, Slovenia
Phone: +386 1 620 5880, Fax: + 386 1 620 5882
E-Mail: pisarna@karlstorz.si

KARL STORZ Endoscopia Romania srl
Str. Prof. Dr. Anton Colorian, nr. 74, Sector 4
041393 Bukarest, Romania
Phone: +40 31 4250800, Fax: +40 31 4250801
E-Mail: office@karlstorz.ro

KARL STORZ Endoskope Greece M.E.P.E.
Gounari 40,
15343 Agia Paraskevi/Athens, Greece
Phone: +30 210 61 31 386, Fax: +30 210 61 31 392
E-Mail: sales@karlstorz.gr

KARL STORZ Endoskope Greece Ltd.*
Ipsilantou Str. 32
54248 Thessaloniki, Greece
Phone: +30 2310 304868 Fax: +30 2310 304862
*Repair & Service Subsidiary

KARL STORZ Industrial**
Gedik Is Merkezi B Blok
Kat 5, D 38-39, Bagdat Cad. No: 162
Maltepe Istanbul, Turkey
Phone: +90 216 442 9500, Fax: +90 216 442 9030
**Sales for Industrial Endoscopy

ООО KARL STORZ Endoskopy – WOSTOK
Derbenyevskaya nab. 7, building 4
115114 Moscow, Russia
Phone: +7 495 983 02 40, Fax: +7 495 983 02 41
E-Mail: kste-wostok@karlstorz.ru

TOV KARL STORZ Ukraine
Obolonska naberezhna, 15
building 3, office 3
04210 Kiev, Ukraine
Phone: +380 44 42668-14, -15, -19
Fax: +380 44 42668-20

KARL STORZ Endoskope
East Mediterranean and Gulf S.A.L.
Solidere – Beirut Souks
Block M, 3rd Floor
2012 3301 Beirut, Lebanon
Phone: +961 1 1999390, Fax +961 1 1999391

KARL STORZ Endoscopy (South Africa) (Pty) Ltd.
P.O. 6061
Roggebaai 8012, South Africa
Phone: +27 21 417 2600, Fax: +27 21 421 5103
E-Mail: info@karlstorz.co.za

ТОО KARL STORZ Endoskopy Kasachstan
Bokeykhan Ulica 8/2, VP – 1
010000 Astana, Respublika Kasachstan
Phone: +7 7172 57-52 16 / 57 28 49 / 57 09 34
Fax +7 7172 43 96 96
E-Mail: info@karlstorz.kz

KARL STORZ ENDOSKOPE East Mediterranean & Gulf
Building 25 (Al-Hayyan) – Units 205 – 206
Dubai Health Care City, P.O.Box 118069
Dubai - United Arab Emirates
Phone: +971 4 4471230
Fax: +971 4 4471231
E-Mail: info-gne@karlstorz-emg.com

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KARL STORZ Endoscopy India Private Ltd.
D-181, Okhla Industrial Area, Phase-1
New Delhi 110020, India
Phone: +91 11 26815445-51,
Fax: +91 11 268129 86
E-Mail: karlstorz@vsnl.com

KARL STORZ GmbH & Co. KG
Resident Representative Office
80/33 (44/19) Dang Van Ngu,
F.10–Q. Phu Nhuan
Ho Chi Minh City, Vietnam
Phone: +848 991 8442, Fax: +848 844 0320

KARL STORZ Endoscopy Australia Pty Ltd
15 Orion Road Lane Cove NSW 2066
P O Box 50 Lane Cove NSW 1595, Australia
Phone: +61 (0)2 9490 6700
Toll free: 1800 996 562 (Australia only)
Fax: +61 (0)2 9420 0695
E-Mail: info@karlstorz.au

KARL STORZ Endoscopy China Ltd.
Hong Kong Representative Office
Unit 1601, Chinachem Exchange Square
1 Hoi Wan Street, Quarry Bay,
Hong Kong, People's Republic of China
Phone: +8 52 28 65 2411, Fax: +8 52 28 65 4114
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy (Shanghai) Ltd.
Beijing Branch Company,
Room 610, China Life Tower No. 6, Chaowai Street,
Beijing, 100020, People's Republic of China
Phone: +86 10 8525 3725, Fax: +86 10 8525 3728
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Unit 3901-3904, Tower 1 Grand Gateway
No. 1 Hong Qiao Road
Shanghai, 200030, People's Republic of China
Phone: +86 21 6113-1188, Fax: +86 21 6113-1199
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Chengdu Branch Company
F-5, 24/F., Chuanxing Mansion,
No. 18 Renming Road South
Chengdu, Sichuan, 610016, People's Rep. of China
Phone: +86 28 8620-0175, Fax: +86 28 8620-0177
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd,
Shenyang Branch Company,
Rm 2225, Tower B, City Plaza,
No. 83, Zhongshan Road, Heping District,
Shenyang, Liaoning, 110001,
People's Republic of China
Phone: +86 24 6258 9911
Fax: +86 24 6258 9922
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Guangzhou Branch Company
Room 1119-20, Dongshan Plaza
69 Xianlie Road Middle, Dongshan District,
Guangzhou, Guangdong, 510095,
People's Rep. of China
Phone: +86 20 8732-1281, Fax: +86 20 8732-1286
E-Mail: info@karlstorz.com.cnk

KARL STORZ Endoscopy Asia Marketing Pte Ltd
No. 8 Commonwealth Lane #03-02
Singapore 149555, Singapore
Phone: +65 69229150, Fax: +65 69229155
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Singapore Sales Pte Ltd
No. 8 Commonwealth Lane #03-02
Singapore 149555, Singapore
Phone: +65 69229150, Fax: +65 69229155
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Korea Co. Ltd.
Jeljon Tower 1, Suite 811
17-1 Jeongja-dong
Seongnam, 463-847, Korea
Phone: +82-70-3450-7474, Fax: +82-70-8277-3299

KARL STORZ Endoscopy Taiwan Ltd.
6F-1, No. 10, Sec. 1, Beisin Rd.
Sindian, Taipei County 231, Taiwan (R.O.C.)
Phone: +886 933 014 160
Fax: +886 2 8672 6399
E-Mail: mingwang@karlstorz.tw

KARL STORZ Endoscopy Japan K. K.
Stage Bldg. 8F, 2-7-2 Fujimi
Chiyoda-ku, Tokyo 102-0071, Japan
Phone: +81 3 6380-8622, Fax: +81 (0)3 6380-8633
E-Mail: info@karlstorz.co.jp

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KARL STORZ Endoscopy Canada Ltd.
7171 Millcreek Drive, Mississauga
Ontario, L5N 3R3, Canada
Phone: +1 905 816-4500, Fax: +1 905 858-4599
E-Mail: info.canada@karlstorz.com

KARL STORZ Endoscopy-America, Inc.
2151 East Grand Avenue
El Segundo, CA 90245-5017, USA
Phone: +1 424 218-8100,
Toll free: 800 421-0837 (USA only)
Fax: +1 424 218-8526
Toll free: 800 321-1304 (USA only)
E-Mail: info@ksea.com

KARL STORZ Veterinary Endoscopy America, Inc.
175 Cremona Drive
Goleta, CA 93117, USA
Phone: +1 805 968-7776, Fax: +1 805 685-2588
E-Mail: info@karlstorzvet.com

KARL STORZ Endoscopia Latino-America, Inc.
815 N. W. 57th Avenue, Suite 480
Miami, FL 33126-2042, USA
Phone: +1 305 262-8980, Fax: +1 305 262-8986
E-Mail: info@ksela.com

KARL STORZ Endoscopia
Miramar Trade Center
Edificio Jerusalem, Oficina 108,
La Habana, Cuba
Phone: +53 72041097, Fax: +53 72041098

KARL STORZ Endoscopia México S.A. de C.V
Lago Constanza No 326 , Col. Chapultepec Morales
D.F.C.P. 11520, México, México
Phone: +52 55 1101-1520, Fax: +52 55 55450174

KARL STORZ Marketing América Do Sul Ltda.
Rua Joaquim Floriano, nº. 413, 20º andar – Itaim Bibi,
CEP-04534-011 São Paulo, Brasil
Phone: +55 11 3526-4600, Fax: +55 11 3526-4680
E-Mail: info@karlstorz.com.br

KARL STORZ Endoscopia Argentina S.A.
Zufriategui 627 6° Piso
B1638 CAA - Vicente Lopez
Provincia de Buenos Aires, Argentina
Phone: +54 11 4718 0919, Fax: +54 11 4718 2773
E-Mail: info@karlstorz.com.ar

KARL STORZ Endoskopi Norge AS
P.O. Box 153, Rolf Olsenvei 28
N-2007 Kjeller, Norway
Phone: +47 6380 5600, Fax: +47 6380 5601
post@karlstorz.no

KARL STORZ Endoskop Sverige AB
Storsätragränd 14,
12724 Skärholmen, Sweden
Postal address: Po Box 8013,
14108 Kungens Kurva, Sweden
Phone: +46 8 50 56 4800, Fax: +46 8 50 56 4848
E-Mail: kundservice@karlstorz.se

KARL STORZ Endoscopy Suomi OY
Valimotie 1 B, 4. kerros, 00380 Helsinki, Finland
Phone: +358 (0)96824774, Fax: +358 (0)968247755
E-Mail: asiakaspalvelu@karlstorz.fi

KARL STORZ SE & Co. KG
Representation Office
Kestucio g. 59
08124 Vilnius, Lithuania

KARL STORZ Endoskopi Danmark A/S
Skovlytoften 33
2840 Holte, Denmark
Phone: +45 45162600, Fax: +45 45162609
E-Mail: marketing@karlstorz.dk

KARL STORZ Endoscopy (UK) Ltd.
392 Edinburgh Avenue, Slough
Berkshire, SL1 4UF, Great Britain
Phone: +44 1753 503500, Fax: +44 1753 578124
E-Mail: customerservice@karlstorz-uk.com

KARL STORZ Endoscopie Nederland B. V.
Phone: +31 651 938 738, Fax: +31 135 302 231

KARL STORZ Endoscopy Belgium N. V.
Phone: +32 473 810 451, E-Mail: info@stopler.be

KARL STORZ Endoscopie France S. A. S.
12, rue Georges Guynemer, Quartier de l'Europe
78280 Guyancourt, France
Phone: +33 1 30484200, Fax: +33 1 30484201
E-Mail: marketing@karlstorz.fr

KARL STORZ Endoskop Austria GmbH
Landstraßer Hauptstr. 148/1/G1
1030 Wien, Austria
Phone: +43 1 71 56 0470, Fax: +43 1 71 56 0479
E-Mail: storz-austria@karlstorz.at

KARL STORZ Endoscopia Ibérica S. A.
Parque Empresarial San Fernando
Edificio Munich – Planta Baja
28830 Madrid, Spain
Phone: +34 91 6771051, Fax: +34 91 6772981
E-Mail: marketing@karlstorz.es

KARL STORZ Endoscopia Italia S. r. l.
Via dell'Artigianato, 3
37135 Verona, Italy
Phone: +39 045 8222000, Fax: +39 045 8222001
E-Mail: info@karlstorz.it

KARL STORZ Adria Eos d.o.o.
Zadarska 80
10000 Zagreb, Croatia
Phone: +385 1 6406 070, Fax: +385 1 6406 077
E-Mail: info@karlstorz.hr

KARL STORZ Endoskopija d.o.o.
Verovškova c. 60A
1000 Ljubljana, Slovenia
Phone: +386 1 620 5880, Fax: + 386 1 620 5882
E-Mail: pisarna@karlstorz.si

KARL STORZ Endoscopia Romania srl
Str. Prof. Dr. Anton Colorian, nr. 74, Sector 4
041393 Bukarest, Romania
Phone: +40 (0)31 4250800, Fax: +40 (0)31 4250801
E-Mail: office@karlstorz.ro

KARL STORZ Endoskope Greece M.E.P.E.
Gounari 40,
15343 Agia Paraskevi/Athens, Greece
Phone: +30 210 61 31 386, Fax: +30 210 61 31 392
E-Mail: sales@karlstorz.gr

KARL STORZ Endoskope Greece Ltd.*
Ipsilantou Str. 32
54248 Thessaloniki, Greece
Phone: +30 2310 304868 Fax: +30 2310 304862

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KARL STORZ Endoscopy (Shanghai) Ltd.
Shenyang Branch Company,
Rm 2225, Tower B, City Plaza,
No. 83, Zhongshan Road, Heping District,
Shenyang, Liaoning, 110001,
People's Republic of China
Phone: +86 24 6258 9911, Fax: +86 24 6258 9922
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Guangzhou Branch Company
Room 1119-20, Dongahan Plaza
69 Xianlie Road Middle, Dongahan District,
Guangzhou, Guangdong, 510095,
People's Rep. of China
Phone: +86 20 8732-1281, Fax: +86 20 8732-1286
E-Mail: info@karlstorz.com.cnk

KARL STORZ Endoscopy Asia Marketing Pte Ltd
No. 8 Commonwealth Lane #03-02
Singapore 149555, Singapore
Phone: +65 69229150, Fax: +65 69229155
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Singapore Sales Pte Ltd
No. 8 Commonwealth Lane #03-02
Singapore 149555, Singapore
Phone: +65 69229150, Fax: +65 69229155
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Korea Co. Ltd.
Jeijon Tower 1, Suite 905
16 Neuti-Ro, Bundang-Gu
Seongnam, 463-847, Korea
Phone: +82-70-4350-7474, Fax: +82-70-8277-3299
E-Mail: stephen.Lee@karlstorz.com

KARL STORZ Endoscopy Taiwan Ltd.
6F-1, No. 10, Sec. 1, Beisin Rd.
Sindian, Taipei County 231, Taiwan (R.O.C.)
Phone: +886 933 014 160
Fax: +886 2 8672 6399
E-Mail: mingwang@karlstorz.tw

KARL STORZ Endoscopy Japan K. K.
Stage Bldg. 8F, 2-7-2 Fujimi
Chiyoda-ku, Tokyo 102-0071, Japan
Phone: +81 3 6380-8622, Fax: +81 (0)3 6380-8633
E-Mail: info@karlstorz.co.jp

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KARL STORZ Endoscopy Canada Ltd.
7171 Millcreek Drive, Mississauga
Ontario, L5N 3R3, Canada
Phone: +1 905 816-4500, Fax: +1 905 858-4599
Toll free phone: 1-800-268-4880 (Canada only)
Toll free fax: 1-800-482-4196 (Canada only)
E-Mail: info.canada@karlstorz.com

KARL STORZ Endoscopy-America, Inc.
2151 East Grand Avenue
El Segundo, CA 90245-5017, USA
Phone: +1 424 218-8100, Fax: +1 424 218-8526
Toll free phone: 800 421-0837 (USA only)
Toll free fax: 800 321-1304 (USA only)
E-Mail: info@ksea.com

KARL STORZ Veterinary Endoscopy-America, Inc.
1 South Los Carneros Road
Goleta, CA 93117, USA
Phone: +1 805 968-7776, Fax: +1 805 685-2588
E-Mail: info@karlstorzvet.com

KARL STORZ Endoscopia Latino-America, Inc.
815 N. W. 57th Avenue, Suite 480
Miami, FL 33126-2042, USA
Phone: +1 305 262-8960, Fax: +1 305 262-8986
E-Mail: info@ksela.com

KARL STORZ Endoscopia México S.A. de C.V.
Edificio Atlantic, Oficina 3G
Calle D e/ 1ra y 3ra
10400 Vedado, Havanna, Cuba
Phone: +537 836 95 06, Fax: +537 836 97 76
E-Mail: kstorzcuba@gmail.com

KARL STORZ Endoscopia México S.A. de C.V.
Lago Constanza No 326,
Col. Granada Del. Miguel Hidalgo
C.P. 11520 México D.F.
Phone: +52 (55) 1101 1520
E-Mail: mx-info@karlstorz.com

KARL STORZ Marketing América Do Sul Ltda.
Rua Joaquim Floriano, nº. 413, 20º andar – Itaim Bibi,
CEP-04534-011 São Paulo, Brasil
Phone: +55 11 3526-4600, Fax: +55 11 3526-4680
E-Mail: info@karlstorz.com.br

KARL STORZ Endoscopia Argentina S.A.
Zufriategui 627 6° Piso, B1638 CAA - Vicente Lopez
Provincia de Buenos Aires, Argentina
Phone: +54 11 4718 0919, Fax: +54 11 4718 2773
E-Mail: info@karlstorz.com.ar

KARL STORZ Endoskopi Norge AS
Østensjøveien 15B
0661 Oslo, Norway
Phone: +47 6380 5600, Fax: +47 6380 5601

KARL STORZ Endoskop Sverige AB
Stonsåtragränd 14
12739 Skärholmen, Sweden
14108 Kungens Kurva, Sweden
Phone: +46 8 505 648 00
E-Mail: kundservice@karlstorz.se

KARL STORZ Endoscopy Suomi OY
Taivaltie 5
01610 Vantaa, Finland
Phone: +358 (0)96824774, Fax: +358 (0)968247755
E-Mail: asiakaspalvelu@karlstorz.fi

KARL STORZ SE & Co. KG
Representation Office
Kestucio g. 59
08124 Vilnius, Lithuania
Phone: +370 68567000
E-Mail: LT-LV-Info@karlstorz.com

KARL STORZ Endoskopi Danmark A/S
Skovlytoften 33
2840 Holte, Denmark
Phone: +45 45162600, Fax: +45 45162609
E-Mail: marketing@karlstorz.dk

KARL STORZ Endoscopy (UK) Ltd.
415 Perth Avenue, Slough
Berkshire, SL1 4TQ, United Kingdom
Phone: +44 1753 503500, Fax: +44 1753 578124
E-Mail: info-uk@karlstorz.com

KARL STORZ Endoscopie Nederland B. V.
Displayweg 2
3821 BT Amersfoort, Netherlands
Phone: +31 (0)33 4545890
E-Mail: info-nl@karlstorz.com

KARL STORZ Endoscopy Belgium N. V.
Phone: +31 (0)33 4545890
E-Mail: info-be@karlstorz.com

KARL STORZ Endoscopie France S. A. S.
12, rue Georges Guynemer, Quartier de l'Europe
78280 Guyancourt, France
Phone: +33 1 30484200, Fax: +33 1 30484201
E-Mail: marketing@karlstorz.fr

KARL STORZ Endoskop Austria GmbH
Landstraßer Hauptstr. 148/1/G1
1030 Wien, Austria
Phone: +43 1 71 56 0470, Fax: +43 1 71 56 0479
E-Mail: storz-austria@karlstorz.at

KARL STORZ Endoscopia Ibérica S. A.
Parque Empresarial San Fernando
Edificio Munich – Planta Baja
28830 Madrid, Spain
Phone: +34 91 6771051, Fax: +34 91 6772981
E-Mail: info-es@karlstorz.com

KARL STORZ Endoscopia Italia S. r. l.
Via dell'Artigianato, 3
37135 Verona, Italy
Phone: +39 045 8222000, Fax: +39 045 8222001
E-Mail: info-ita@karlstorz.com

KARL STORZ Adria Eos d.o.o.
Capraška 6
10000 Zagreb, Croatia
Phone: +385 1 6406 070, Fax: +385 1 6406 077
E-Mail: info@karlstorz.hr

KARL STORZ Endoskopija d.o.o.
Verovškova c. 60A
1000 Ljubljana, Slovenia
Phone: +386 1 620 5880, Fax: + 386 1 620 5882
E-Mail: pisarna@karlstorz.si

KARL STORZ Endoskop Polska Marketing Sp. z o.o.
ul. Bojkowska 47
44-100 Gliwice, Poland
Phone: +48 32 706 13 00, Fax: +48 32 706 13 07
E-Mail: info-pl@karlstorz.com

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KARL STORZ ENDOSKOPE
East Mediterranean and Gulf S.A.L.
Block M, 3rd Floor
Beirut Souks, Weygand Street
2012 3301 Beirut, Lebanon
Phone: +961 1 999390, Fax +961 1 999391
E-Mail: info@karlstorz-emg.com

KARL STORZ Endoscopy (South Africa) (Pty) Ltd.
P.O. 6061, Roggebaai 8012
Cape Town, South Africa
Phone: +27 21 417 2600, Fax: +27 21 421 5103
E-Mail: info@karlstorz.co.za

TOO KARL STORZ Endoskopy Kasachstan
6, Saryarka str. BC „Arman", office314
010000 Astana, Republic of Kazakhstan
Phone: +7 7172 552-549, 552-788, Fax: -444
E-Mail: info@karlstorz.kz

KARL STORZ ENDOSKOPE
East Mediterranean & Gulf (branch)
Building West Side 7A – Unit 7WA – 3008
Dubai Airport Free Zone, P.O. Box 54963
Dubai - United Arab Emirates
Phone: +971 (0)4 2958887, Fax: +971 (0)4 3205282
Service Hotline: +971 (0)4 3415882
E-Mail: info-gne@karlstorz-emg.com

KARL STORZ Endoscopy India Private Limited
11th Floor, Dr. Gopal Das Bhawan
28, Barakhamba Road
New Delhi 110001, India
Phone: +91 11 4374 3000, Fax: +91 11 4374 3010
E-Mail: corporate@karlstorz.in

KARL STORZ SE & Co. KG
Resident Representative Office
14th Floor, MPlaza Saigon
39 Le Duan, District 1
Ho Chi Minh City, Vietnam
Phone: +84 8 3823 8000, Fax: +84 8 3823 8039
E-Mail: infovietnam@karlstorz.com

KARL STORZ Endoscopy Australia Pty. Ltd .
15 Orion Road Lane Cove NSW 2066
P O Box 50 Lane Cove NSW 1595, Australia
Phone: +61 (0)2 9490 6700
Toll free: 1800 996 562 (Australia only)
Fax: +61 (0)2 9420 0695
E-Mail: info@karlstorz.au

KARL STORZ Endoscopy China Ltd., Hong Kong
Room 2503-05, 25F AXA Tower, Landmark East,
No. 100 How Ming Street, Kwun Tong, Kowloon,
Hong Kong, People's Republic of China
Phone: +852 28 65 2411, Fax: +852 28 65 4114
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy (Shanghai) Ltd.
Beijing Branch,
Room 1805-1807, Building B, 18F Beijing IFC
No. 8, Jianguomenwai Street, Chaoyang District,
Beijing, 100022, People's Republic of China
Phone: +86 10 5638188, Fax: +86 10 5638199
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Shanghai Branch
36F, Tower 1 Grand Gateway
No. 1 Hongqiao Road, Xuhui District,
Shanghai, 200030, People's Republic of China
Phone: +86 21 6113-1188, Fax: +86 21 6113-1199
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Chengdu Branch
Room 805, 8F Jin Jiang International Building
No. 1 West Linjiang Road, Wuhou District,
Chengdu, 610041, People's Republic of China
Phone: +86 28 86587977, Fax: +86 28 86587975
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Shenyang Branch,
Room 2001-2005, 20F N-MEDIA International Center,
No. 167 Youth Avenue, Shenhe District,
Shenyang, 110014, People's Republic of China
Phone: +86 24 23181118, Fax: +86 24 23181119
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Guangzhou Branch
Room 02B, 35F Teem Tower,
No. 208 Tianhe Road, Tianhe District,
Guangzhou, 510620, People's Republic of China
Phone: +86 20 8732-1281, Fax: +86 20 8732-1286
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy Asia Marketing Pte Ltd.
No. 8 Commonwealth Lane #03-02
Singapore 149555, Singapore
Phone: +65 69229150, Fax: +65 69229155
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Singapore Sales Pte Ltd
No. 8 Commonwealth Lane #03-02
Singapore 149555, Singapore
Phone: +65 69229150, Fax: +65 69229155
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Korea Co. Ltd.
I-Park Bundang1, 102dong, 8F
239 Jeongjail-Ro, Bundang-Gu
Seongnam, 13556, Korea
Phone: +82-70-4350-7474, Fax: +82-70-8277-3299
E-Mail: infokorea@karlstorz.com

KARL STORZ Endoscopy Taiwan Ltd.
12F, No. 192, Sec. 2, Chung Hsin Rd.,
Sindian District, New Taipei City, Taiwan
Phone: +886 933 014 160, Fax: +886 2 8672 6399
E-Mail: mingwang@karlstorz.tw

KARL STORZ SE & Co. KG
Representative Office Philippines
1901 Picadilly Star Bldg., 4th Avenue, BGC
Taguig City 1636, Philippines
Phone: , Fax: +63 2 31745-00
E-Mail: philippines@karlstorz.com

KARL STORZ Endoscopy Japan K. K.
Stage Bldg. 8F, 2-7-2 Fujimi
Chiyoda-ku, Tokyo 102-0071, Japan
Phone: +81 3 6380-8622, Fax: +81 3 6380-8633
E-Mail: info@karlstorz.co.jp

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KARL STORZ Endoscopy Canada Ltd.
7171 Millcreek Drive, Mississauga
Ontario, L5N 3R3, Canada
Phone: +1 905 816-4500, Fax: +1 905 858-4599
Toll free phone: 1-800-268-4880 (Canada only)
Toll free fax: 1-800-482-4198 (Canada only)
E-Mail: info.canada@karlstorz.com

KARL STORZ Endoscopy-America, Inc.
2151 East Grand Avenue
El Segundo, CA 90245-5017, USA
Phone: +1 424 218-8100, Fax: +1 424 218-8526
Toll free phone: 800 421-0837 (USA only)
Toll free fax: 800 321-1304 (USA only)
E-Mail: info@ksea.com

KARL STORZ Veterinary Endoscopy-America, Inc.
1 South Los Carneros Road
Goleta, CA 93117, USA
Phone: +1 805 968-7776, Fax: +1 805 685-2588
E-Mail: info@karlstorzvet.com

KARL STORZ Endoscopia Latino-America, Inc.
815 N. W. 57th Avenue, Suite 480
Miami, FL 33126-2042, USA
Phone: +1 305 262-8960, Fax: +1 305 262-8986
E-Mail: info@ksela.com

KARL STORZ Endoscopia México S.A. de C.V.
Edificio Atlantic, Oficina 3G
Calle D e/ 1ra y 3ra
10400 Vedado, Havanna, Cuba
Phone: +537 836 95 06, Fax: +537 836 97 76
E-Mail: kstorzcuba@gmail.com

KARL STORZ Endoscopia México S.A. de C.V.
Lago Constanza No 326,
Col. Granada Del. Miguel Hidalgo
C.P. 11520 México D.F.
Phone: +52 (55) 1101 1520
E-Mail: mx-info@karlstorz.com

KARL STORZ Marketing América Do Sul Ltda.
Rua Joaquim Floriano, nº. 413, 20º andar – Itaim Bibi,
CEP-04534-011 São Paulo, Brasil
Phone: +55 11 3526-4600, Fax: +55 11 3526-4680
E-Mail: info@karlstorz.com.br

KARL STORZ Endoscopia Argentina S.A.
Zufriategui 627 6° Piso, B1638 CAA - Vicente Lopez
Provincia de Buenos Aires, Argentina
Phone: +54 11 4718 0919, Fax: +54 11 4718 2773
E-Mail: info@karlstorz.com.ar

KARL STORZ Endoskopi Norge AS
Østensjøveien 15B
0661 Oslo, Norway
Phone: +47 6380 5600, Fax: +47 6380 5601

KARL STORZ Endoskop Sverige AB
Storsätragränd 14
12739 Skärholmen, Sweden
14108 Kungens Kurva, Sweden
Phone: +46 8 505 648 00
E-Mail: kundservice@karlstorz.se

KARL STORZ Endoscopy Suomi OY
Taivaltie 5
01610 Vantaa, Finland
Phone: +358 (0)96824774, Fax: +358 (0)968247755
E-Mail: asiakaspalvelu@karlstorz.fi

KARL STORZ SE & Co. KG
Representation Office
Kestucio g. 59
08124 Vilnius, Lithuania
Phone: +370 68567000
E-Mail: LT-LV-Info@karlstorz.com

KARL STORZ Endoskopi Danmark A/S
Skovlytoften 33
2840 Holte, Denmark
Phone: +45 45162600, Fax: +45 45162609
E-Mail: marketing@karlstorz.dk

KARL STORZ Endoscopy (UK) Ltd.
415 Perth Avenue, Slough
Berkshire, SL1 4TQ, United Kingdom
Phone: +44 1753 503500, Fax: +44 1753 578124
E-Mail: info-uk@karlstorz.com

KARL STORZ Endoscopie Nederland B. V.
Displayweg 2
3821 BT Amersfoort, Netherlands
Phone: +31 (0)33 4545890
E-Mail: info-nl@karlstorz.com

KARL STORZ Endoscopy Belgium N. V.
Phone: +31 (0)33 4545890
E-Mail: info-be@karlstorz.com

KARL STORZ Endoscopie France S. A. S.
12, rue Georges Guynemer, Quartier de l'Europe
78280 Guyancourt, France
Phone: +33 1 30484200, Fax: +33 1 30484201
E-Mail: marketing@karlstorz.fr

KARL STORZ Endoskop Austria GmbH
Landstraßer Hauptstr. 148/1/G1
1030 Wien, Austria
Phone: +43 1 71 56 0470, Fax: +43 1 71 56 0479
E-Mail: storz-austria@karlstorz.at

KARL STORZ Endoscopia Ibérica S. A.
Parque Empresarial San Fernando
Edificio Munich – Planta Baja
28830 Madrid, Spain
Phone: +34 91 6771051, Fax: +34 91 6772981
E-Mail: info-es@karlstorz.com

KARL STORZ Endoscopia Italia S. r. l.
Via dell'Artigianato, 3
37135 Verona, Italy
Phone: +39 045 8222000, Fax: +39 045 8222001
E-Mail: info-ita@karlstorz.com

KARL STORZ Adria Eos d.o.o.
Capraška 6
10000 Zagreb, Croatia
Phone: +385 1 6406 070, Fax: +385 1 6406 077
E-Mail: info@karlstorz.hr

KARL STORZ Endoskopija d.o.o.
Verovškova c. 60A
1000 Ljubljana, Slovenia
Phone: +386 1 620 5880, Fax: + 386 1 620 5882
E-Mail: pisarna@karlstorz.si

KARL STORZ Endoskop Polska Marketing Sp. z o.o.
ul. Bojkowska 47
44-100 Gliwice, Poland
Phone: +48 32 706 13 00, Fax: +48 32 706 13 07
E-Mail: info-pl@karlstorz.com

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KARL STORZ Endoscopy (Shanghai) Ltd.
Chengdu Branch
Room 805, 8F Jin Jiang International Building
No. 1 West Linjiang Road, Wuhou District,
Chengdu, 610014, People's Republic of China
Phone: +86 28 86587977, Fax: +86 28 86587975
E-Mail: info@karlstorz.com.cn

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Important information for users of KARL STORZ devices

Thank you for your expression of confidence in the KARL STORZ brand name. Like all of our other products, this product is the result of years of experience and great care in manufacture. You and your organization have decided in favor of a modern, high-quality item of equipment from KARL STORZ.

This instruction manual is intended to serve as an aid in the proper installation of the interface module* SCB/KLS MARTIN, model 20 0923 22, and it describes the differences in operation via the OR1™ control NEO. All essential details of the equipment and all actions required on your part are clearly presented and explained. We thus ask that you read this manual carefully before proceeding to work with the equipment. Keep it available for ready reference in a convenient and conspicuous location near the equipment.

***Note:** In order to use the interface module, it must be integrated in the SCB Interface control or SCB ACC Control®. This is why in the instruction manual the term »device« is frequently used.*

***Note:** The instruction manual for the SCB Interface control exists as a ring binder version. If you are using this version, please insert these module instructions in the corresponding section of the binder.*

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General warning

Users of KARL STORZ devices should be thoroughly trained in their use and the applicable medical procedures. These devices and accessories should not be used until the user has fully familiarized himself/herself with the instructions for use, assembly and care. Instruction manuals should be carefully studied and be available to the surgical team during the procedure; it is essential to follow the instructions contained in the instruction manual, paying particular attention to the cautions and warnings.

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6 Niederlassungen 18

Contents

1 General information 2
1.1 Description of the unit 2
1.2 Property rights 2

2 Safety instructions 3
2.1 Explanation of warnings and cautions 3
2.2 Warnings and cautions 4
2.3 Intended use 5
2.4 User qualification 5
2.5 Safety precautions when operating the unit 5

3 Installation and operating instructions 6
3.1 Unpacking the equipment 6
3.2 Basic equipment 6
3.3 Connection diagrams 7
3.4 Installing 8
3.5 Operating the unit 9
3.6 Operating KLS MARTIN marLED® operating lights 9
3.7 Device setup 11

4 Technical description 12
4.1 Troubleshooting 12
4.2 Technical data 13
4.3 Standard compliance 13
4.4 Directive compliance 13
4.5 Technical documentation 14

5 Spare parts, recommended accessories 15
5.1 Spare parts 15
5.2 Accessories 15

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1 General information

1.1 Description of the unit

When integrated in SCB Interface control **20**0901**2**0, or SCB ACC Control® **20**0906**3**2, the KARL STORZ Interface Module **20**0923**2**2 enables you to remotely control KLS MARTIN marLED® operating lights using the KARL STORZ SCB control NEO system. Using the SCB control NEO software, a maximum of 4 lights of V series or E series can be controlled. The lights' camera functionality is also supported. The functions that are accessible on the wall control board or the front panel can be controlled through the KARL STORZ SCB control NEO system.

Note: *Please also follow the instruction manuals for the KARL STORZ SCB control NEO system and for the KLS MARTIN surgical lights marLED® V series and E series. The KARL STORZ interface module **20**0923**2**2, together with other possible KARL STORZ interface modules, is integrated into KARL STORZ SCB Interface control **20**0901**2**0, or SCB ACC Control® **20**0906**3**2. Also follow the instruction manual for the device in question.*

1.2 Property rights

This product is protected in the USA by (at least one of the following) US Patent no(s) 5,788,688; 6,397,286; 6,484,221; 6,824,539.

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2 Safety instructions

2.1 Explanation of warnings and cautions

Please read this manual and follow its instructions carefully. The words Warning, Caution, and Note convey special meanings. Wherever they are used in this manual, they should be carefully reviewed to ensure the safe and effective operation of this product. To make the signal words stand out more clearly, they are accompanied by a pictogram.

Warning: *A Warning indicates that the personal safety of the patient or physician may be involved. Failure to observe a Warning could result in injury to the patient or physician.*

Caution: *A Caution indicates that particular service procedures or precautions must be followed to avoid possible damage to the device.*

Note: *A Note indicates special information about operating the device or clarifies important information.*

Warning: *Before putting the KLS MARTIN® Interface Module **20**0923 22 into operation, read the instruction manual thoroughly and follow the instructions carefully.*

Warning: *Combinations of medical devices are only assured to be safe if*

- *they are identified as such in the respective instruction manuals or*
- *the intended purpose and interface specifications of the devices used in combination permit this (cf. IEC 60601-1-1).*

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2.2 Warnings and cautions

Warning: *The instructions and interface specifications for medical devices and/or system components used in combination must be observed precisely.*

Warning: *Test this equipment prior to each surgical procedure to ensure that it functions correctly.*

Warning: *Before using the unit it is imperative that you be acquainted with how the unit operates and is controlled.*

Note: *This unit has been marked in accordance with the European Directive on waste electrical and electronic equipment (WEEE).*
At the end of its useful operating life, dispose of the unit as electronic scrap.
Please ask either KARL STORZ SE & Co. KG, a KARL STORZ subsidiary or your specialist dealer for information on your local collection point.
Within the scope of application of this Directive, KARL STORZ SE & Co. KG is responsible for the proper disposal of this unit.

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2.3 Intended use

The KARL STORZ interface module **20**0923 22 is intended for integration in SCB Interface control **20**0901 20, or SCB ACC Control® **20**0906 32 and enables you to remotely control KLS MARTIN marLED® light types V series and E series with the KARL STORZ SCB control NEO system.

(* is a requirement for use of the KARL STORZ interface module **20**0923 22.)

The use of this interface module in any way other than indicated above is not permitted for safety reasons.

Unauthorized conversions or modifications to the unit are not allowed for safety reasons.

The KARL STORZ interface module **20**0923 22 must only be operated with accessories designated as suitable by KARL STORZ.

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2.4 User qualification

The KARL STORZ interface module **20**0923 22 may only be used (in accordance with the specification given under »Intended use«) by physicians and medical assistants with the appropriate specialized qualifications and who have been instructed in use of the devices.

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2.5 Safety precautions when operating the unit

It is the user's responsibility to make sure that the SCB Interface control or SCB ACC Control® with integrated KARL STORZ interface module **20**0923 22 operates properly before using it.

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3 Installation and operating instructions

3.1 Unpacking the equipment

Check for missing items and evidence of shipping damage. File any complaints with the manufacturer or supplier immediately. If possible, retain the original packing materials for later use.

The interface module **20**0923**22** may have to be installed into KARL STORZ SCB Interface control **20**0901**20**, or SCB ACC Control® **20**0906**32** and tested for proper operation by KARL STORZ service personnel or by personnel authorized by KARL STORZ (see the »KARL STORZ SCB Modules Installation Instructions« in relation to this).

3.2 Basic equipment

KARL STORZ Set No.: 20092390

1 Interface module **20**0923**22**

1 RS232 connecting cable **20**0918**70**, length 200 cm

1 Instruction manual 96206226DER

1 Manufacturer's identification plate

1 »Compatible Devices« label

2 Label strips for SCB Interface control **20**0901**20**

2 Label strips for SCB ACC Control® **20**0906**32**

ⓘ ***Note:*** *Operation requires a KARL STORZ OR1™ control NEO **20**0971**20**-1 with installed SCB control NEO software **20**0900**01**-41 or higher as well as an interface board from KLS MARTIN. The interface board can be ordered directly from KLS MARTIN.*

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The KARL STORZ interface module **20**0923**22** may only be installed in the SCB Interface control **20**0901**20** or SCB ACC Control® **20**0906**32** by KARL STORZ service personnel or personnel authorized by KARL STORZ (see the “KARL STORZ SCB Modules Installation Instructions” in relation to this).

ⓘ ***Note:*** *In particular, follow all instructions in the manual for the KLS MARTIN marLED®.*

Insert the RS232 connecting cable into the appropriate connection socket of SCB Interface control **20**0901**20** or SCB ACC Control® **20**0906**32**. Connect the other end of the RS232 connecting cable to the connection of the wall control unit of KLS MARTIN marLED® - (see KLS MARTIN marLED® instruction manual).

ⓘ ***Note:*** *The KLS MARTIN marLED® light system with an optional surgiCam internal camera must meet the following system requirements in order to be remote controlled via the KARL STORZ SCB control NEO system:*

- *Lights: up to four OR lights*
 1 – 4 x V series (without internal camera)
 1 – 4 x V series (only 1 x V16/V10 can be equipped with an internal camera)
 1 – 4 x E series (without internal camera)
 1 – 4 x E series (only 1 x E15/E9i can be equipped with an internal camera)
- *Internal camera: surgiCam*
 Up to two functions of the surgiCam internal camera are supported: Zoom in/out and rotation
- *KLS MARTIN interface board:*
 Software (SW) version: FW 01-00-00

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3.5 Operating the unit

Switch on the unit (SCB Interface control 20 0901 20 or SCB ACC Control® 20 0906 32), the marLED® light system and the KARL STORZ OR1™ control NEO. Once the validated hardware and software have been checked, the system boots and automatically starts up the SCB control NEO software.

KLS MARTIN marLED® reports to the SCB control NEO software.

ⓘ ***Note:*** *Please also observe the instruction manual for the KARL STORZ SCB control NEO system.*

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3.6 Operating KLS MARTIN marLED® OR lights

The realistic depiction of the marLED® operating light shows all keys (and thereby all functions) of the real control unit that can be accessed through the SCB control NEO system. The key assignment corresponds to the real control unit.

Each light's proper position on the holding arm is indicated by a symbol in the control panel. As is the case when using the wall control unit, the "ALL" control button synchronizes all lights with one another.

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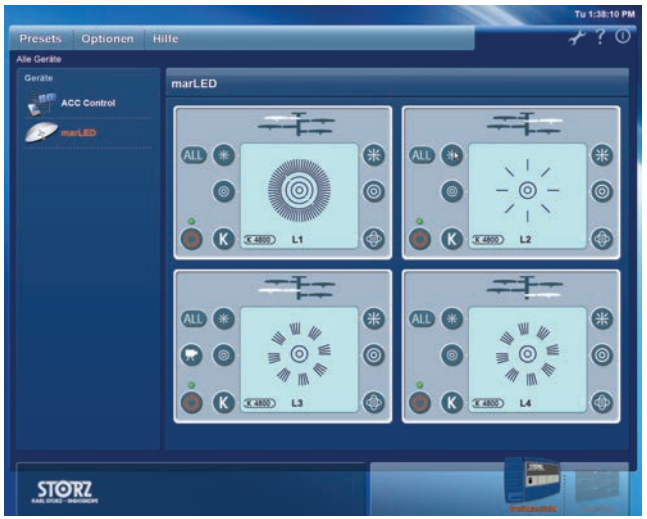
At most, 4 connected lights are displayed.

ⓘ ***Note:*** *If the light system shows the error "Cycling off the lights within 0.5 seconds", the marLED® operating lights cannot be controlled for 10 seconds. An error message appears at the same time.*

When creating a "Preset" within the SCB control NEO software, the marLED® light is grayed out since no presets can be made.

When creating a "Preset" within the SCB control NEO software, the number of lights that were most recently detected in the system is displayed. Presets cannot be made.

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Note: *Please refer to the KLS MARTIN instruction manual for the description of the operation.*

Note: *Alongside operation with the KARL STORZ SCB control NEO system, the KLS MARTIN marLED® lights can, of course, still be controlled with the control unit on the wall or the light cardanics.*

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3.7 Device setup

Warning: *Settings must only be carried out by skilled technical staff or KARL STORZ personnel.*

Perform the settings as described in the 'KARL STORZ SCB control NEO system'« instruction manual.

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4 Technical description

4.1 Troubleshooting

Warning: *Always unplug the unit before all maintenance work.*

Symptom:

- The KLS MARTIN marLED® light does not report to the KARL STORZ SCB control NEO system.

Possible causes:

- SCB connecting cable faulty or not properly inserted.
- RS232 connecting cable between MarLEDCom (board in the wall control unit or on the top hat rail) and the interface is faulty or not properly inserted.
- Control unit (SCB Interface control 20 0901 20 or ACC Control® 20 0906 32) is switched off.
- Software version of the “MarLEDCom” component is incompatible.

Remedy:

- Check SCB connecting cable and replace, if necessary.
- Check RS232 connecting cable and replace, if necessary.
- Switch on control unit.
- Request “MarLEDCom” with compatible software from KLS MARTIN.

Symptom:

- The order of the lights differs from the control fields in the SCB control NEO software.

Possible causes:

- Lights addressed in wrong sequence.

Remedy:

- KLS MARTIN will readdress the lights.

Note: *Should any problems arise which directly affect the controlled unit, please contact the unit manufacturer.*

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4.5 Technical documentation

On request, the manufacturer will provide those circuit diagrams, itemized parts listings, descriptions, sets of adjustment instructions, and other items of available documentation to suitably qualified user personnel duly authorized by the manufacturer for their use in repairing those components of the instrument that have been designated as reparable by their respective manufacturers.

Supply of such technical documentation relating to the instrument shall not be construed as constituting manufacturer's authorization of user's personnel, regardless of their level of technical training, to open or repair the instrument.

Explicitly exempted herefrom are those maintenance and repair operations described in this manual.

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5 Spare parts, recommended accessories

5.1 Spare parts

Item	Cat. no.
RS232 connecting cable, length 200 cm	20091870
Instruction manual	96206226DER

5.2 Accessories

Item	Cat. no.
RS232 connecting cable, length 500 cm	20091970
RS232 connecting cable, length 10 m	20091971
RS232 connecting cable, length 15 m	20091972

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KARL STORZ Endoscopy (Shanghai) Ltd.
Shanghai Branch
36F, Tower 1 Grand Gateway
No. 1 Hongquiao Road, Xuhui District,
Shanghai, 200030, People's Republic of China
Phone: +86 21 6113-1188, Fax: +86 21 6113-1199
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.
Chengdu Branch
Room 805, 8F Jin Jiang Intenational Building
No. 1 West Linjiang Road, Wuhou District,
Chengdu, 6100414, People's Republic of China
Phone: +86 28 86587977, Fax: +86 28 86587975
E-Mail: info@karlstorz.com.cn

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