DA

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Storage
- Store in a dry place and protect from high humidity.
- Protect from direct sun light.
- Regard temperature limits.

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EN

Description
curea P1 duo active is a sterile, hydro-active wound dressing. The surface is made of polypropylene fleece. The absorbent bed consists of cellulose composite with fluid-retaining sodium polyacrylate (SuperCore®) and a layer of activated carbon.
curea P1 duo active absorbs on both sides. curea P1 duo active absorbs exudate, blood and germs under compression.
The fluids are bound in the interior of the absorbent core.
The absorption- and retention property of curea P1 duo active allows the wound a longer resting phase between the dressing changes therapy supporting the body's own wound healing.
The product does not contain any pharmacological ingredients.

Indication
Chronic and low to heavy exuding wounds.
Acute and iatrogenic wounds.
Infected wounds or wounds vulnerable to infection.
Malodorous wounds; Ulcerated malignant tumours.

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Other notices
- The product has been sterilised with ethylene oxide.
- Latex-free

Storage
- Store in a dry place and protect from high humidity.
- Protect from direct sun light.
- Regard temperature limits.

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Storage
• Store in a dry place and protect from humidity.
• Protect from direct sun light.
• Regard temperature limits.

Other notices
• The product has been sterilised with ethylene oxide.
• Latex-free

Description
curea P1 is a sterile superabsorbent wound dressing. The surface facing the wound is made of polypropylene fleece; the absorbent pad consists of cellulose composite with fluid-retaining superabsorbent materials. The backfoil which forms as additional moisture impermeable barrier.
This allows curea P1 to absorb secretion and exudate, while maintaining no optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.
curea P2 posee un borde blando que evita la formacion de nuevos decubitos por presion. No contiene principios activos farmacologicos.

Indication
Chronic to low to heavy exuding wounds.
Acute wounds.

Contraindication
Dry wounds (signation, risk of adhesion)
Do not use for mucous membranes, eyes, sinews or bones, because deep wounds parts might dry out and take permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Application
Make sure that curea P1 is used under medical supervision and based on medical indication. The physician treating the patient will define the intervals after which the dressing has to be changed based on the actual wound situation regarding exudation and infection. We recommend to change the dressing at least after 7 days and a frequent (best daily) check of the wound area for signs of infection and of the dressing for depletion of capacity.
Ensure that a sufficient amount to the wound dressing on the skin.

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
Application
Make sure that curea P1 border is used under medical supervision and based on medical indication. The Physician treating the patient will define the interval after which the dressing has to be changed based on the actual wound situation regarding exudation and infection. We recommend changing the dressing at least after 7 days or a frequent (preferably daily) check of the wound area for signs of infection, and of the dressing for depletion of capacity.

Disposal
The wound pad can be disposed of with household waste or in accordance with the local provisions. Special guidelines in medical institutions governing the disposal (e.g. for multi-resistant MRSA or VRE) must be complied with.

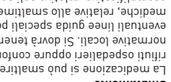Storage
• Store in a dry place and protect from humidity.
• Protect from direct sun light.
• Observe storage temperature limits.

Other notices
LATEX STERILE EO
• The product has been sterilised with ethylene oxide.
• Latex-free

Contra-indication
Dry wounds (stagnation, risk of adhesion)
Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and suffer permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Indication
Chronic and low to heavy exuding wounds.
Acute wounds

Contra-indications
Dry wounds (stagnation, risk of adhesion)
Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and suffer permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Application
curea P1 border is a sterile superabsorbent wound dressing. The surface in contact with the wound is made of polypropylene fleece; the absorbent core consists of cellulose complex with fluid-retaining superabsorbent materials. The peel layer contains a breathable but moisture impermeable polyurethane film. The adhesive edges are acrylate based.
This allows curea P1 border to absorb secre-tion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.
The wound gains a longer resting phase to support the body's own wound healing power.
curea P1 border also absorbs under pressure.
curea P1 border features a soft border to prevent the formation of renewed pressure sores. It does not contain any pharmacological ingredients.

Indications
Chronic wounds
Acute wounds

Contra-indication
Dry wounds (stagnation, risk of adhesion)
Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and suffer permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Application
curea P1 border is a sterile superabsorbent wound dressing. The surface in contact with the wound is made of polypropylene fleece; the absorbent core consists of cellulose complex with fluid-retaining superabsorbent materials. The peel layer contains a breathable but moisture impermeable polyurethane film. The adhesive edges are acrylate based.
This allows curea P1 border to absorb secre-tion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.
The wound gains a longer resting phase to support the body's own wound healing power.
curea P1 border also absorbs under pressure.
curea P1 border features a soft border to prevent the formation of renewed pressure sores. It does not contain any pharmacological ingredients.

Indications
Chronic and low to heavy exuding wounds.
Acute wounds

Contra-indications
Dry wounds (stagnation, risk of adhesion)
Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and suffer permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Application
curea P1 border is a sterile superabsorbent wound dressing. The surface in contact with the wound is made of polypropylene fleece; the absorbent core consists of cellulose complex with fluid-retaining superabsorbent materials. The peel layer contains a breathable but moisture impermeable polyurethane film. The adhesive edges are acrylate based.
This allows curea P1 border to absorb secre-tion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.
The wound gains a longer resting phase to support the body's own wound healing power.
curea P1 border also absorbs under pressure.
curea P1 border features a soft border to prevent the formation of renewed pressure sores. It does not contain any pharmacological ingredients.

Indication
Chronic and low to heavy exuding wounds.
Acute wounds

Contra-indication
Dry wounds (stagnation, risk of adhesion)
Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and suffer permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Application
curea P1 border is a sterile superabsorbent wound dressing. The surface in contact with the wound is made of polypropylene fleece; the absorbent core consists of cellulose complex with fluid-retaining superabsorbent materials. The peel layer contains a breathable but moisture impermeable polyurethane film. The adhesive edges are acrylate based.
This allows curea P1 border to absorb secre-tion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.
The wound gains a longer resting phase to support the body's own wound healing power.
curea P1 border also absorbs under pressure.
curea P1 border features a soft border to prevent the formation of renewed pressure sores. It does not contain any pharmacological ingredients.

Indications
Chronic and low to heavy exuding wounds.
Acute wounds

Contra-indications
Dry wounds (stagnation, risk of adhesion)
Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and suffer permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Application
curea P1 border is a sterile superabsorbent wound dressing. The surface in contact with the wound is made of polypropylene fleece; the absorbent core consists of cellulose complex with fluid-retaining superabsorbent materials. The peel layer contains a breathable but moisture impermeable polyurethane film. The adhesive edges are acrylate based.
This allows curea P1 border to absorb secre-tion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.
The wound gains a longer resting phase to support the body's own wound healing power.
curea P1 border also absorbs under pressure.
curea P1 border features a soft border to prevent the formation of renewed pressure sores. It does not contain any pharmacological ingredients.

Indication
Chronic and low to heavy exuding wounds.
Acute wounds

**Contra-indication

Modern wound care plus lint compress functionality now has a name: curea P1 drain

curea medical
Innovation with Care

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EN

Description
curea P1 drain is a sterile superabsorbent wound dressing. The surface is made of polypropylene fleece; the absorbent pad consists of cellulose composite with fluid-retaining superabsorbent materials. curea P1 drain absorbs on both sides.
This allows curea P1 drain to absorb secretion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.
The wound gains a longer resting phase to support the body's own wound healing power. curea P1 drain also absorbs under pressure. curea P1 drain features a soft border to prevent the formation of renewed pressure sores. It does not contain any pharmacological ingredients.

Indication
Chronic and low to heavy exuding wounds, catheters and drainage or infusion tubes.
Acute wounds

Contraindication
Dry wounds (stagnation, risk of adhesion).
Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and take permanent damage.
Caution with cavities, because the wound pad expands when absorbing fluid.

Application
Make sure that curea P1 drain is used under medical supervision and based on medical indication. The physician treating the patient will define the intervals after which the dressing has to be changed based on the actual wound situation regarding exudation and infection.
We recommend to change the dressing at least after 7 days and a frequent (best daily) check of the wound area for signs of infection and of the dressing for depletion of capacity. Make sure that the side without imprint is applied to the wound. Ensure that a sufficient amount of area is available to attach the wound dressing on the skin.

Adverse reactions
No adverse reactions are known.

Warnings
- Do not cut or tear the wound pad. Polyacrylate might be spilled, contaminating the patient. In case of polyacrylate spilled, wash the patient with clean water and the wound with a wound irrigation solution
- Dispose of damaged or expired products to avoid the risk of infection
- Designed for single use only. Do not re-use to prevent the possible risk of infection
- Do not re-sterilise

Other notices
- The product has been sterilised with ethylene oxide
- Latex-free

Storage
- Store in a dry place and protect from high humidity
- Protect from direct sun light
- Regard temperature limits

Disposal
The wound pad can be disposed of with the household waste or in accordance with the local provisions. Special guidelines in medical institutions governing the disposal (e.g. for multi-resistant MRSA or VRE) must be complied with.

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EN

Description

curea P1duo is a sterile superabsorbent wound dressing. The surface is made of polypropylene fleece; the absorbent pad consists of cellulose composite with fluid-retaining superabsorbent materials. curea P1duo absorbs on both sides.

This allows curea P1 to absorb secretion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.

Indication

Chronic and low to heavy exuding wounds.

Acute wounds

Contra-indication

Dry wounds (stagnation, risk of adhesion)

Do not use for mucous membranes, eyes, sinews or bones, because these body parts might dry out and suffer permanent damage.

Caution with cavities, because the wound pad expands when absorbing fluid.

Application

Make sure that curea P1 is used under medical supervision and based on medical indication. The physician treating the patient will define the intervals after which the dressing has to be changed based on the actual wound situation regarding exudation and infection. We recommend to change the dressing at least after 7 days and a frequent (best daily) check of the wound area for signs of infection and of the dressing for depletion of capacity.

Ensure that a sufficient amount of area is available to attach the wound dressing on the skin.

Adverse reactions

No adverse reactions are known.

Warnings

- Do not cut or tear the wound pad. Polyacrylate might be spilled, contaminating the patient. In case of polyacrylate spilled, wash the patient with clean water and the wound with a wound irrigation solution.

- Dispose of damaged or expired products to avoid the risk of infection.

- Designed for single use only. Do not re-use to prevent the possible risk of infection.

- Do not re-sterilise.

Other notices

- The product has been sterilised with ethylene oxide.

- Latex-free

Storage

- Store in a dry place and protect from high humidity.
- Protect from direct sun light.
- Regard temperature limits.

Disposal

The wound pad can be disposed of with the household waste or in accordance with the local provisions. Special guidelines in medical institutions governing the disposal (e.g. for multi-resistant MRSA or VRE) must be complied with.

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EN

Description

curea P2 is a sterile, superabsorbing wound dressing. The surface facing the wound consists of an integrated wound contact lattice made from PE, the absorbent pad consits of cellulose composite with fluid-retaining superabsorbent materials. The back contains a breathable and moisture impermeable BTBS.

This allows curea P2 to absorb secretion and exudate, while maintaining an optimal, moist wound milieu. The fluids are bound in the interior of the absorbent core without dripping or leaking out.

The wound contact lattice avoids the wound dressing to stick onto the regenerating tissue.

The wound gains a longer resting phase to support the body's own wound healing power. curea P2 also absorbs under pressure.

curea P2 features a soft border to prevent the formation of renewed pressure sores. It does not contain any pharmacological ingredients.

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Application

Make sure that curea P2 is used under medical supervision and based on medical indication. As well, the physician treating the patient will define the intervalls after which curea P2 needs to be exchanged. The maximum wear time is 7 days.

Make sure that the side without imprint is applied to the wound.

Ensure that a sufficient amount of space is available to attach curea P2 on the skin.

Adverse reactions

No adverse reactions are known.

Warnings

- Do not cut or tear the wound pad.
- Do not use if the sterile packaging is damaged or opened. Dispose of damaged or expired products.
- Designed for single use only. Do not re-use to prevent the possible risk of infection.
- Do not re-sterilise.

Other notices

- The product has been sterilised with ethylene oxide. STERILE EO
- Latex-free

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