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Rev. Nr. 04/2013-06

1. Description

curea 150 is an OR-table sheet for single use and consists of a laminate of a polypropylene nonwoven with a breathable polyethylene film, an absorbent core made of cellulose and polyacrylate covered with another layer of polypropylene nonwoven.
curea 150 can bear loads up to 150 kg, if the patient lies outstretched and the patient's weight is distributed evenly over the length of the product. curea 150 is pervious to x-rays.

2. Indication

- general use superabsorbent OR table cover
- curea 150 is intended for the transport of patients between OR-table and transport stretcher.
- highly absorbent cover to protect transport stretcher and OR-table

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curea medical
Innovation with Care

curea 150 STERILE – Instructions for use

Rev. Nr. 01/2017-08

1. Description

curea 150 is a sterile OR-table sheet for single use and consists of a laminate of a polypropylene nonwoven with a breathable polyethylene film, an absorbent core made of cellulose and polyacrylate covered with another layer of polypropylene nonwoven.
curea 150 can bear loads up to 150 kg, if the patient lies outstretched and the patient's weight is distributed evenly over the length of the product. curea 150 is pervious to x-rays.

2. Indication

- general use superabsorbent OR table cover
- curea 150 is intended for the transport of patients between OR-table and transport stretcher
- highly absorbent cover to protect transport stretcher and OR-table

3. Contraindication

- curea 150 is not a carrying sheet acc. EN 1865-1
- extremely wet conditions (e.g. outdoor use, > 1 L fluid)
- loads exceeding 150 kg and punctiform loads

4. Application

Lay the patient in the middle of the roughly structured, non-printed surface of the product. To lift the sheet, firmly grab the left and right long sides along the absorbent core. Lift smoothly and carefully, DO NOT lift abruptly. Shift the patient with appropriate promptness.

5. Adverse reactions


No adverse reactions are known. Very heavy patients (> 150 kg) cannot be lifted with curea 150, because the product might tear.

6. Warnings

curea 150 is designed for single use only. Do not re-use to prevent the possible risk of infection.

Dispose of damaged or expired products to avoid the risk of infection.

Do not re-sterilise.

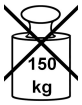Do not use for loads exceeding 150 kg, do not use damaged products

7. Other notices

curea 150 is Latex-free and phatlate-free

The product has been sterilised with ethylene oxide.

8. Storage

Store in a dry place and protect from high humidity.

Protect from direct sun light.

9. Disposal

curea 150 can be disposed of with the household waste or in accordance with the local provisions. Special guidelines in medical institutions governing the disposal (e.g. for multi-resistant MRSA or VRE) must be complied with.

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Instructions for use of Liquimat
GB
Handling
1. Place the mat on the floor, where the liquid is expected to drop. Please note that Liquimat can be placed on either side, as both sides are absorbent.
2. Let Liquimat absorb during the whole time of the surgical operation. If necessary, use additional mats.
3. After completing the surgical operation: fold the mat and dispose according to local regulations.

Safety Warning
Do not cut or tear the mat apart. Damaged mats must be disposed of. Do not use mats that are expired. The mats might get slippery when wet - do not step on them.

Made in Germany by:
curea medical GmbH
Münsterstraße 61-65
DE-48565 Steinfurt

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Professional solutions for a better working environment in the OR

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curea medical GmbH
Münsterstraße 61-65
48565 Steinfurt, Germany

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EN

Description

Liquisorb is a sterile superabsorbing skirt, designed to take up liquids, lymph or blood. The surface is made from a nonwoven and the core consists of cellulose pulp containing superabsorber.

The adhesive version can be attached to the patient or the treatment unit directly.

Indication
Absorption of large quantities of fluid during operations, invasive procedures and other treatments.

Contraindication
None

Application
Liquisorb can be combined with itself or other coverings. The layer should not prevent/obscure iatrogenic access.

Adverse reactions
There are no known side effects. An onset of allergic reactions in direct contact with patients cannot be entirely ruled out in rare cases.

Warnings

- Do not cut or tear the layer.
- Dispose of damaged or expired products to avoid the risk of infection
- Used products can be slippery – do not step on.
- The layer can be discarded with standard clinical waste or in accordance with local regulations.

Other notices

- Do not use if the sterile packaging is damaged or opened.
- Do not re-use, single use only.
- Do not re-sterilise.
- The product has been sterilised with ethylene oxide.

Storage

- Store in a dry place and protect from high humidity.
- Protect from direct sun light.
- Regard temperature limits.

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Rev. Nr. 01/2012-08

1. Description

The SafeCare patient underlay consists of a laminate of a polypropylene nonwoven with a breathable polyethylene film, an absorbent core made of cellulose and polyacrylate covered with another layer of polypropylene nonwoven on top.
Polar fluids are chemically bound by the sodiumpoylacrylate and so pollution of the bedding is prevented. The risk of bed sores is reduced.

2. Indication

- protection of bedding against pollution by urine and other fluids
- the risk of bed sores is reduced

3. Contraindication

- extremely wet conditions (e.g. outdoor use, > 1 L fluid)
- the underpad provides no lifting capacity

4. Application

Take only the number of care mats needed out of the plastic bag. Place the care mat, with the textured surface uppermost, under the patient posterior.

5. Adverse reactions

No adverse reactions are known.

6. Warning

SafeCare patient underlays are designed for single use only. Do not re-use, to prevent the possible risk of infection.

7. Other advice

SafeCare patient underlays are Latex-free

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Protect from direct sun light.

9. Disposal

SafeCare patient underlays can be disposed of with the household waste or in accordance with the local provisions. Special guidelines in medical institutions governing the disposal (e.g. for multi-resistant MRSA or VRE) must be complied with.

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