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RECOMMENDED INDICATIONS
1. Cementation of metal-based inlays, onlays, crowns and bridges
2. Cementation of resin inlays, onlays, crowns and bridges
3. Cementation of all ceramic inlays
4. Cementation of high strength (zirconia based) all ceramic crowns and bridges
5. Cementation of metal, ceramic and fiber posts

CONTRAINDICATIONS
1. Pulp capping.
2. Avoid use of this product in patients with known allergies to glass ionomer cement, methacrylate monomer or methacrylate polymer.

DIRECTIONS FOR USE

Powder / Liquid ratio (g/g)	2.0/1.0
Mixing time (sec.)	20''
Working Time (23°C / 73°F) (min., sec.)	2'30''
Finishing Time (37°C/99°F, after seating the restoration) (min., sec.)	4'30''

1. TOOTH PREPARATION
 a) Prepare tooth in usual manner. For pulp capping, use calcium hydroxide.
 b) Clean the prepared tooth with pumice and water.
 Note:
 Use GC Fuji PLUS CONDITIONER on the prepared tooth surface for 20 seconds for optimal adhesion.
 c) Rinse thoroughly with water. Remove excess moisture by blotting with a cotton pellet or gently blowing with an oil-free air. DO NOT DESICCATE. Prepared surfaces should appear moist (glistening).
2. RESTORATION PREPARATION
 Make sure that the restoration is pretreated and handled according to the manufacturer's instructions.
3. POWDER AND LIQUID DISPENSING
 a) The standard powder to liquid ratio is 2.0 g / 1.0 g. (1 level large scoop of powder to 3 drops of liquid, or 1 level small scoop of powder to 1 drop of liquid).
 b) For accurate dispensing of powder, lightly tap the bottle against the hand. Do not shake or invert.
 c) Hold the liquid bottle vertically and squeeze gently.
 d) Close bottles immediately after use.
4. MIXING
 Dispense powder and liquid onto the pad. Using the plastic spatula, add all the powder to the liquid. Mix rapidly for 20 seconds.
 Note:
 When mixing larger amounts, divide the powder into 2 equal parts. Mix the first portion with all of the liquid for 15 seconds. Incorporate the remaining portion and mix the whole thoroughly for 15-25 seconds (total 30-40 seconds).
5. CEMENTATION
 a) Coat the internal surface of the restoration with sufficient cement and seat immediately. The working time is 2 minutes 30 seconds from start of mixing at 23°C (73°F).
 Higher temperature will shorten working time.
 Note:
 For more complex procedures requiring longer working time, it is recommended to mix powder and liquid on a cold mixing slab and or use GC Fuji PLUS EWT.
 b) Maintain moderate pressure.
 c) Start removing excess cement when the excess cement feels rubbery.
 d) Finishing can be started 4 minutes 30 seconds after seating the restoration.

STORAGE
Recommended for optimal performance, store in a cool and dark place (4-25°C) (39.2-77.0°F).

SHADES
Yellow, A3 (Vita® shade), Translucent
Vita® is a registered trademark of Vita Zahnfabrik, Bad Säckingen, Germany.

PACKAGES
1. 1-1 package: 15 g powder, 8 g (7 mL) liquid, GC Fuji PLUS CONDITIONER 7 g (6.5 mL), a double-ended powder scoop, a plastic spatula, mixing pad (#23)
2. Powder 15 g with scoop
3. Liquid 8 g (7 mL)
4. GC Fuji PLUS CONDITIONER 7 g (6.5 mL)

CAUTION
1. In case of contact with oral tissue or skin, remove immediately with a sponge or cotton soaked in alcohol. Flush with water. To avoid contact, a rubber dam and/or COCOA BUTTER can be used to isolate the operation field from the skin or oral tissue.
2. In case of contact with eyes, flush immediately with water and seek medical attention.
3. DO NOT mix powder or liquid with components of other glass ionomer cements.
4. GC Fuji PLUS CONDITIONER contains 2% Ferric Chloride. It is not recommended to use GC Fuji PLUS CONDITIONER in combination with products containing Epinephrine as this may stain oral tissue.
5. Do not use GC Fuji PLUS in combination with eugenol containing materials as eugenol may hinder GC Fuji PLUS with setting or bonding properly.
6. This product is not indicated for filling or core build-up.
7. Personal protective equipment (PPE) such as gloves, face masks and safety eyewear should always be worn.
8. In rare cases the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.

Some products referenced in the present IFU may be classified as hazardous according to GHS. Always familiarize yourself with the safety data sheets available at:
http://www.gceurope.com
They can also be obtained from your supplier.

CLEANING AND DISINFECTING:
MULTI-USE DELIVERY SYSTEMS: to avoid cross-contamination between patients this device requires mid-level disinfection. Immediately after use inspect device and label for deterioration. Discard device if damaged.
DO NOT IMMERSE. Thoroughly clean device to prevent drying and accumulation of contaminants. Disinfect with a mid-level registered healthcare-grade infection control product according to regional/national guidelines.

Undesired effects- Reporting:
If you become aware of any kind of undesired effect, reaction or similar events experienced by use of this product, including those not listed in this instruction for use, please report them directly through the relevant vigilance system, by selecting the proper authority of your country accessible through the following link:
https://ec.europa.eu/growth/sectors/medical-devices/contacts_en as well as to our internal vigilance system: vigilance@gc.dental
In this way you will contribute to improve the safety of this product.

Last revised: 06/2019

CE 2797

MANUFACTURED BY
GC CORPORATION
76-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8585, Japan

DISTRIBUTED BY
GC CORPORATION
76-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8585, Japan

EU: GC EUROPE N.V.
Researchpark Haasrode-Leuven 1240, Interleuvenlaan 33, B-3001 Leuven, Belgium TEL: +32 16 74 10 00

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RECOMMENDED INDICATIONS
1. Cementation of metal-based inlays, onlays, crowns, bridges and posts.
2. Cementation of high strength (zirconia based) all ceramic crowns and bridges.

CONTRAINDICATIONS
1. Pulp capping.
2. Avoid use of this product in patients with known allergies to glass ionomer cement.

DIRECTIONS FOR USE

Powder/liquid ratio (g/g)	1.8/1.0
Mixing time (min., sec.)	20"
Working time (min., sec.) (at 23°C) (from start of mixing)	2'00"
Final finishing commencing time (after seating the restoration)	4'30"

1. TOOTH PREPARATION
 a) Prepare tooth in usual manner. For pulp capping, use calcium hydroxide.
 b) Clean the prepared teeth with pumice and water.
 c) Rinse thoroughly with water. Remove excess moisture by blotting with a cotton pellet or gently blowing with an oil-free air. DO NOT DESICCATE. Prepared surfaces should appear moist (glistening).
2. RESTORATION PREPARATION
 Make sure that the restoration is pretreated and handled according to the manufacturer's instructions.
3. POWDER AND LIQUID DISPENSING
 a) The standard powder to liquid ratio is 1.8g/1.0g (1 level scoop of powder to 2 drops of liquid).
 b) For accurate dispensing of powder, lightly tap the bottle against the hand. Do not shake or invert.
 c) Hold the liquid bottle vertically and squeeze gently.
 d) Close bottles immediately after use.

4.MIXING
 Dispense powder and liquid onto the pad. Using the plastic spatula, add all the powder to the liquid. Mix rapidly for 20 seconds.
 Note :
 When mixing larger amounts, divide the powder into two equal parts. Mix the first portion with all of liquid for 5 seconds. Incorporate the remaining portion and mix the whole thoroughly for 15 seconds (total 20 seconds).

5.CEMENTATION
 a) Coat the internal surface of the restoration with sufficient cement and seat immediately. The working time is 2 minutes from start of mixing at 23°C (73.4°F). Higher temperatures shorten working time.
 b) Maintain moderate pressure.
 c) Start removing excess cement when the excess cement feels rubbery.
 d) Finishing can be started 4 minutes 30 seconds after seating the restoration.

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SHADE
Light yellow

PACKAGES
1. 1-1 package : 35g powder, 25g (20mL) liquid, powder scoop, mixing pads (No. 20).
2. Bottle of 35g powder with scoop.
3. Bottle of 25g (20mL) liquid.

CAUTION
1. In case of contact with oral tissue or skin, remove immediately with a sponge or cotton soaked in alcohol. Flush with water. To avoid contact, a rubber dam and/or cocoa butter can be used to isolate the operation field from the skin or oral tissue.
2. In case of contact with eyes, flush immediately with water and seek medical attention.
3. DO NOT mix powder or liquid with components of other glass ionomer cements.
4. This product is not indicated for filling or core build-up.
5. Personal protective equipment (PPE) such as gloves, face masks and safety eyewear should always be worn.
6. In rare cases the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.

Some products referenced in the present IFU may be classified as hazardous according to GHS. Always familiarize yourself with the safety data sheets available at: http://www.gceurope.com
They can also be obtained from your supplier.

CLEANING AND DISINFECTING:
MULTI-USE DELIVERY SYSTEMS: to avoid cross-contamination between patients this device requires mid-level disinfection. Immediately after use inspect device and label for deterioration. Discard device if damaged.

DO NOT IMMERSE. Thoroughly clean device to prevent drying and accumulation of contaminants. Disinfect with a mid-level registered healthcare-grade infection control product according to regional/national guidelines.

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Prior to use, carefully read the instructions for use. EN

GC FujiCEM® 2

RADIOPAQUE REINFORCED GLASS IONOMER LUTING CEMENT

For use only by a dental professional in the recommended indications.

RECOMMENDED INDICATIONS

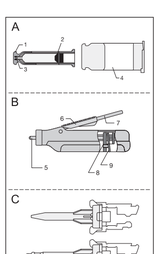
1. Cementation of metal-based inlays, onlays, crowns and bridges
2. Cementation of resin inlays, onlays, crowns and bridges
3. Cementation of all ceramic inlays
4. Cementation of high strength (e.g. zirconia based, lithium disilicate etc.) all ceramic inlays, onlays, crowns and bridges
5. Cementation of metal, ceramic and fiber posts

CONTRAINDICATIONS

1. Pulp capping.
2. In rare cases, the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.

IDENTIFICATION OF PARTS (Fig. 1)

A. Paste Pak Cartridge SL
1. Cover 2.Inner plunger 3. Cartridge tip 4. Sleeve
B. Paste Pak Dispenser
5. Piston 6. Adjusting guide 7. Lever 8. Sliding block
9. Piston release
C. GC FujiCEM 2 Mixing Tip SL (Slide & Lock system)
GC FujiCEM 2 Mixing Tip SL for Endo (with intraoral tip)

Fig. 1

DIRECTIONS FOR USE

Paste / Paste ratio (g/g)	2.0 / 1.0
Mixing time (sec.)	10"-15"
Working time (min., sec.) (at 23°C) (from start of mixing)	2'15"
Final finishing commencing time (min., sec.) (after seating the restoration)	4'30"

PASTE PAK CARTRIDGE SL LOADING

1. Make sure that the piston is completely retracted into the dispenser (Fig. 2).
Note :
To retract the piston, push the piston release and at the same time pull the sliding block (on the underside of the dispenser) with fingers.
2. Load the cartridge into the Paste Pak Dispenser (Fig .3).
Ensure the ▲ mark on the cartridge lines up with the ▲ mark on the Paste Pak Dispenser as illustrated. Then turn the cartridge fully to the right position.
3. After cartridge loading, push the piston release forward until it stops (Fig .4).
a) Push the piston release forward.
b) Do not push the upper side.
Note :
To remove the cartridge, retract the piston fully and turn the cartridge back to the one-half position.

A. CEMENTATION TECHNIQUE – AUTOMIXING WITH GC FujiCEM 2 Mixing Tip SL

1. TOOTH PREPARATION
a) Prepare tooth in usual manner. For pulp capping, use calcium hydroxide.
b) Clean the prepared tooth with pumice and water. Note:
1) **NaClO, EDTA or Clorhexidine can be used to clean the prepared tooth. Do not use substances like H_2O_2 (hydrogen peroxide) and other disinfectants as these may lower the bond strength.**
2) Use GC Fuji Plus Conditioner on the prepared tooth surface for 20 seconds for optimal adhesion.
c) Rinse thoroughly with water. Remove excess moisture by blotting with a cotton pellet or gently blowing with an oil-free air. DO NOT DESICCATE. Prepared surfaces should appear moist (glistening).
2. RESTORATION PREPARATION
Make sure that the restoration is pretreated and handled according to the manufacturers instructions.
3. DISPENSING AND AUTOMIXING
a) Remove the cartridge cover.
b) **Before using, check the two openings of the cartridge tip and confirm the pastes are at the same level to ensure even flow from the cartridge.**
Note:
Before dispensing the first time from a new cartridge or after a long interval in between use, bleed the pastes from the cartridge (Fig. 5-2).
c) Attach the GC FujiCEM 2 Mixing Tip SL on the cartridge. Push the mixing tip completely until the end. (Fig. 6-1)
Hold the sides of sleeve and slide it forward to lock the mixing tip (Fig. 6-2).
Note :
1) Make sure the GC FujiCEM 2 Mixing Tip SL is securely attached to the cartridge. If not, remove the GC FujiCEM 2 Mixing Tip SL and reattach.
2) If the material has been refrigerated, leave the cartridge in room temperature for 10 minutes.
d) Set the adjusting guide to the front of the lever (Fig. 7).
Note:
Before dispensing the material into the appliances, bleed a small amount of the mixed paste on a pad.
e) Slowly depress the lever to dispense material into the appliance.
f) After dispensing, release the lever. Slide the sleeve backward and remove theGC FujiCEM 2 Mixing Tip SL (Fig. 8). Clean the cartridge tip and replace the cartridge cover (Fig. 11).
Note :
1) Do not leave the used GC FujiCEM 2 Mixing Tip SL on the cartridge for a long period of time.
2) The amount of residual material in the cartridge can be confirmed by the position of the slide block on the underside of the dispenser. As the cartridge cylinders are translucent, the residual amount can also be confirmed by the position of the toner seen through the cylinders. When only a small amount of paste remains in the cartridge, the paste may not be dispensed from the GC FujiCEM 2 Mixing Tip SL.
3) To avoid damage, do not drop the dispenser with the cartridge and mixing tip.
4. CEMENTATION
a) Coat the internal surface of the restoration with sufficient cement and seat immediately. The working time is 2 minutes 15 seconds from start of mixing at 23°C (74°F). Higher temperature will shorten the working time and lower temperature will prolong the working time. Maintain moderate pressure.
b) Maintain moderate pressure and make sure the restoration remains in place. Start removing excess cement with a scaler when the excess cement is rubbery.
c) Finishing can be started 4 minutes 30 seconds after seating the restoration.

B. CEMENTATION TECHNIQUE – HANDMIXING

1. TOOTH PREPARATION
Refer to A. 1.
2. RESTORATION PREPARATION
Refer to A. 2.
3. DISPENSING
a) Remove the cartridge cover.
b) Depress the lever to dispense the required amounts of paste onto the mixing pad.
Note :
1) Before dispensing the first time from a new cartridge or after a long interval in between use, bleed the pastes from the cartridge (Fig. 5).
2) The adjusting guide can be moved forwards or backwards to change the quantity of material to be dispensed (Fig. 9).
1. Position of adjusting guide for one screw.
2. Position of adjusting guide for two screws.
c) While depressing the lever (Fig. 10-1), level and cut off the extruded material from the cartridge tips by moving the cartridge and dispenser to a position over the mixing pad (Fig. 10-2).
d) When the lever is released after dispensing, residual pastes at the cartridge tips are retracted into the cartridge. Clean the cartridge tip and replace the cartridge cover (Fig. 11).
Note :
The amount of residual material in the cartridge can be confirmed by the position of the sliding block on the underside of the dispenser. As the cartridge cylinders are translucent, the residual amount can also be confirmed by the position of the inner plunger seen through the cylinders.
4. MIXING
After dispensing, spread the two pastes thinly on the mixing pad using a plastic spatula. Mix thoroughly, with lapping motions, for 10 seconds. Be careful not to incorporate air bubbles.
Note :
1) Larger amount of paste requires a longer mixing time (15 seconds).
2) Do not mix GC FujiCEM 2 pastes with other materials.
5. CEMENTATION
Refer to A. 4. Refrigeration of the cartridge is recommended to extend the working time.

CLEANING / STERILIZATION OF PASTE PAK DISPENSER

For cleaning and sterilization of PASTE PAK DISPENSER, refer to the respective instructions for use.

SHADE

Light Yellow

STORAGE

Recommended for optimal performance, store at room temperature away from heat, moisture and sunlight (4-25°C).Shelf life : 2 years from date of manufacture)

PACKAGES

1. GC FujiCEM 2 Automix
Paste Pak Cartridge SL (13,3g / 7.2 mL) (2), GC FujiCEM 2 Mixing Tip SL (44)
2. GC FujiCEM 2 Refill
Paste Pack Cartridge SL (13,3g / 7.2 mL) (2) with mixing pad (No. 22).
3. GC FujiCEM 2 Automix Starter kit with GC FujiCEM 2 Dispenser
Paste Pak Cartridge SL (13,3g / 7.2 mL) (1), GC FujiCEM 2 Mixing Tip SL (15), GC FujiCEM 2 Mixing Tip SL for Endo (5), GC FujiCEM 2 Dispenser (1)
4. GC FujiCEM 2 Automix Starter kit with Paste Pak Dispenser
Paste Pak Cartridge SL (13,3g / 7.2 mL) (1), GC FujiCEM 2 Mixing Tip SL (15), GC FujiCEM 2 Mixing Tip SL for Endo (5), Paste Pak Dispenser (1)
5. GC FujiCEM 2 Mixing Tip SL, Refill pack
6. GC FujiCEM 2 Mixing Tip SL for Endo
7. GC FujiCEM 2 Mixing Tip SL for Endo, Refill pack
8. GC FujiCEM 2 Mixing Tip SL for Endo (10)
7. Paste Pak Dispenser (1)
8. GC FujiCEM 2 Dispenser (1)

CAUTION

1. In case of contact with oral tissue or skin, remove immediately with a sponge or cotton soaked in alcohol. Flush with water. To avoid contact, a rubber dam and/or cocoa butter can be used to isolate the operation field from the skin or oral tissue.
2. In case of contact with eyes, flush immediately with water and seek medical attention.
3. In rare cases metal may cause sensitivity to some people. If any such reactions are experienced, discontinue the use of Paste Pak Dispenser(nickel-plated aluminum) and refer to a physician.
4. Store Paste Pak Dispenser away from chemicals such as iodine tincture containing halogen element. Otherwise the dispenser may discolor.
5. GC Fuji PLUS CONDITIONER contains 2%Ferric Chloride. It is not recommended to use GC Fuji PLUS CONDITIONER in combination with products containing Epinephrine as this may stain oral tissue.
6. Do not use GC FujiCEM 2 in combination with eugenol containing materials as eugenol may hinder GC FujiCEM 2 from setting or bonding properly.
7. This product is not indicated for filling or core build-up.
8. The GC Mixing Tip SL (including intraoral tip) is for single use only. To prevent cross contamination between patients, do not reuse the GC FujiCEM 2 Mixing Tip SL . The GC FujiCEM 2 Mixing Tip SL cannot be sterilized in an autoclave or

9. Personal protective equipment (PPE) such as gloves, face masks and safety eyewear should always be worn.

Some products referenced in the present IFU may be classified as hazardous according to GHS.
Always familiarize yourself with the safety data sheets available at: http://www.gceurope.com
They can also be obtained from your supplier.

Undesired effects - Reporting:
If you become aware of any kind of undesired effect, reaction or similar events experienced by use of this product, including those not listed in this instruction for use, please report them directly through the relevant vigilance system, by selecting the proper authority of your country accesible through the following link:
https://ec.europa.eu/growth/sectors/medical-devices/contacts_en
as well as to our internal vigilance system: vigilance@gc.dental
In this way you will contribute to improve the safety of this product.

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GC Fuji PLUS™ CAPSULE

RADIOPAQUE REINFORCED GLASS IONOMER LUTING CEMENT IN CAPSULES

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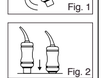
CONTRAINDICATIONS
1. Pulp capping.
2. In rare cases the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.

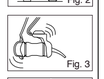
DIRECTIONS FOR USE

Powder / Liquid ratio (g/g)	0.36 / 0.18
Mixing time (sec.)	10"
Working time (min., sec.) (at 23 °C) (from start of mixing)	2'00"
Final finishing commencing time (after seating the restoration)	4'15"

Fig. 1

1. TOOTH PREPARATION
 a) Prepare tooth in usual manner. For pulp capping, use calcium hydroxide.
 b) Clean the prepared teeth with pumice and water. Note :
 Apply GC Fuji PLUS CONDITIONER for 20 seconds using a cotton pellet or sponge.
 c) Rinse thoroughly with water. Remove excess moisture by blotting with a cotton pellet or gently bowing with an oil-free air. DO NOT DESICCATE. Best results are obtained when prepared surfaces appear moist (glistening).
2. RESTORATION PREPARATION
 Make sure that the restoration is pretreated and handled according to the manufacturer's instructions.
3. CAPSULE ACTIVATION AND MIXING
 a) Before activation, shake the capsule or tap its side on a hard surface to loosen the powder (Fig. 1).
 b) To activate the capsule, push the plunger until it is fully depressed with the main body and hold it down for 2 seconds (Fig. 2).
 Note :
 1) Ensure the plunger is fully pressed to avoid the incorrect mixing ratio of powder and liquid.
 2) The capsule should be activated just before mixing and used immediately.
 c) Immediately set it into a mixer (or an amalgamator) and mix for 10 seconds (+/- 4,000 RPM) (Fig. 3).
 Note: When using a rotational mixing device, it is important not to use the centrifuge option.

Fig. 2

Fig. 3

4. DISPENSING
 a) Immediately remove the mixed capsule from the mixer and load it into the GC CAPSULE APPLIER.
 b) Make two clicks to prime the capsule (Fig. 4).
 To adjust the direction of the nozzle, hold the applier with the capsule towards you and turn the capsule body.
 c) Extrude directly into the restoration (Fig. 5).
 d) To remove the used capsule, push the applier release button. Twist the capsule and pull upwards.
5. CEMENTATION
 a) Coat the internal surface of the restoration with sufficient cement and seat immediately. The working time is 2 minutes from start of mixing at 23°C (73.4°F). Higher temperatures will shorten working time.
 Note:
 For more complex procedures requiring longer working time, it is recommended to refrigerate the capsules in the original aluminium foil.
 b) Maintain moderate pressure.
 c) Start removing excess cement when the excess cement feels rubbery (Fig. 6).
 d) Finishing can be started 4 minutes 15 seconds after seating the restoration.

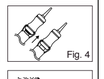
Fig. 4

Fig. 5

Fig. 6

STORAGE
Recommended for optimal performance, store the capsules in the original aluminium foil in a cool and dark place 4-25°C (39.2-77.0°F).

SHADES
Yellow
A3 (Vita® shade)
Vita® is a registered trademark of Vita Zahnfabrik, Bad Säckingen, Germany.

PACKAGE
50 capsules.
Average contents per capsule : 0.36g powder and 0.18g (0.16mL) liquid.
Minimum net volume of mixed cement per capsule : 0.19mL.

Option : GC CAPSULE APPLIER (1 piece).

CAUTION
1. In case of contact with oral tissue or skin, remove immediately with a sponge or cotton soaked in alcohol. Flush with water. To avoid contact, a rubber dam and/or cocoa butter can be used to isolate the operation field from the skin or oral tissue.
2. In case of contact with eyes, flush immediately with water and seek medical attention.
3. GC Fuji PLUS CONDITIONER contains 2% Ferric Chloride. It is not recommended to use GC Fuji PLUS CONDITIONER in combination with products containing Epinephrine as this may stain oral tissue.
4. Do not use GC Fuji PLUS CAPSULE in combination with eugenol containing materials as eugenol may hinder GC Fuji PLUS CAPSULE from setting or bonding properly.
5. This product is not indicated for filling or core build-up.
6. Personal protective equipment (PPE) such as gloves, face masks and safety eyewear should always be worn.

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Undesired effects - Reporting:

If you become aware of any kind of undesired effect, reaction or similar events experienced by use of this product, including those not listed in this instruction for use, please report them directly through the relevant vigilance system, by selecting the proper authority of your country accessible through the following link:
https://ec.europa.eu/growth/sectors/medical-devices/contacts_en

as well as to our internal vigilance system: vigilance@gc.dental

In this way you will contribute to improve the safety of this product.

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