
Brightway®

Nitrile Examination Gloves – Powder Free

Examination Gloves - instructions for use.

How to put on your gloves

1. Clean your hands with soap and water or hand sanitizer before touching the gloves.
2. Remove a glove from the box and make sure there are no obvious tears or holes.
3. Determine the thumb up side of the glove and hold the inner cuff opposing the hand. Slide in your hand into the glove working it into the glove.
4. Repeat the process with the second glove.

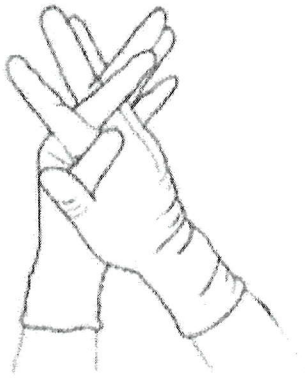How to remove gloves

1. Thoroughly clean your gloved hands with soap and water or hand sanitizer before touching the mask.
2. Avoid touching your skin when removing the gloves.
3. Pinch the glove at the cuff on the hand with you are removing the gloves from. Gently remove the glove keeping hold of the glove in the remaining gloved hand.
4. With your un-gloved hand pinch the INSIDE of the cuff at the wrist and draw the glove inside out to remove it so it enclosed the previously removed glove.
5. Place the removed gloves into the trash/waste ensuring that, where possible, the gloves are disposed where it cannot be removed from the trash/waste by a 3rd party.
6. Thoroughly clean your hands with soap and water or hand sanitizer.

If gloves have been used in an environment where there is the potential for viral contamination then gloves should be disposed of in clinical waste. If this is not possible the gloves should be disposed of in sealed double bags to avoid contamination downstream.

Brightway Holdings Sdn. Bhd. (180521-T)
Lot 1559, Jalan Ikhtiar, Batu Salah, 42100 Klang, Selangor Darul Ehsan, Malaysia. Tel:(603) 3343 1007, 3343 1094

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