Ormco™

Instructions for Use

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blūgloo™

TWO-WAY COLOR-CHANGE ADHESIVE FOR AESTHETIC BRACKETS

Instructions for Use **English**

blūgloo™ is a two-way color-change adhesive for aesthetic brackets. **blūgloo** contains a temperature-sensitive blue dye that becomes clear when it approaches body temperature and reverts to a blue color when it falls below body temperature. The blue tint makes any excess or remaining adhesive easier to see when bonding and debonding.

Tooth Preparation

- Place lip/cheek retractor and saliva ejector in the mouth.
- Apply etching solution to the bonding surfaces with applicator sponges. Dab gently and continuously for 30 seconds per tooth (one arch at a time) – do not rub. Rinse thoroughly for a minimum of 5 seconds per tooth with a forceful airway spray and suction using a high-speed evacuator. Maintain a dry field. Do not allow moisture, saliva or mucosal tissue to contact the etched enamel.
- Dry the etched enamel with clean, dry air. Inspect to ensure a uniform, dull frosty appearance.

Bonding Procedure

- Dispense small amount of Ortho Solo™ sealant on the mixing pad. Replace the cap immediately after dispensing. Apply very thin, uniform coat of Ortho Solo to prepared enamel surface using the brushes provided.
- Extrude small amount of **blūgloo** adhesive paste onto the base of the bracket pad. After application, clean the tip and replace the cap.
- Place bracket onto the tooth surface immediately after applying the adhesive.
- Position bracket, and press down firmly. A slight excess of **blūgloo** (blue tint) should be visible around the periphery of the bracket.
- Remove excess adhesive from the bracket base periphery. (The blue color will disappear when the adhesive reaches body temperature.)
- With conventional light-curing units (i.e. Optilux 501), light-cure the adhesive for 10-15 seconds for clear and 20 seconds for metal brackets. Light-cure clear brackets by shining the light directly through the bracket. Light-cure metal brackets by shining the light from as many directions as possible. For molar tubes a cure time of 30 seconds is recommended due to larger pad sizes.
- Archwires may be inserted after all brackets are cured.

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Debonding Procedure

- Follow manufacturer's recommendations to debond brackets.
- To easily see and remove any remaining **blūgloo** adhesive after debonding, cool the adhesive with short blasts of air or water to revert it back to its original blue color. Remove remaining **blūgloo** adhesive with high or low-speed carbide burs followed by sandpaper discs and a polishing cup or paste.

Note: The color change of **blūgloo** is caused by a change in temperature. It is not an indication of the degree of polymerization. **blūgloo** may revert temporarily to its blue color when the patient consumes cold liquids or foods. However, the color change reverts quickly, and should not be noticeable or objectionable under clear brackets when the adhesive is used in a thin film.

Shelf Life and Storage

- Do not expose materials to intense light or elevated temperatures.
- **blūgloo** should be used at room temperature.
- Syringe and bottle caps must be replaced immediately after use.

Note:

1. Unidose™ Delivery – for single patient use only, to prevent cross contamination between patients.
2. Unidose Composite Gun Sterilization Guidelines – Steam Autoclave: 30 Minutes @ 132°C (270°F).

⚠ Caution

1. **Etching Solution**
 Contains 37% phosphoric acid solution. Avoid contact with skin, eyes, and soft tissue. In case of contact with skin or eyes, flush immediately with water, get medical attention for eyes. Do not take internally.
2. **Sealant and Paste**
 Avoid prolonged or repeated contact with oral tissues and skin. Wash skin with soap and water after contact. Uncured materials may cause contact dermatitis. If skin irritation occurs, discontinue use. Use in a well-ventilated area. Keep containers closed when not in use. Wear gloves when handling this material.

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(Side 1)

Step 1

Apply **very thin** coat of Ortho Solo™ to the prepared tooth.


Tooth Preparation

Etch approximately 30 seconds per tooth (one arch at a time).

Rinse thoroughly for a minimum of 5 seconds per tooth with a forceful air/water spray and suction using a high-speed evacuator.

Dry etched enamel with clean, dry air. Inspect to ensure a uniform, dull, frosty appearance.

Step 2

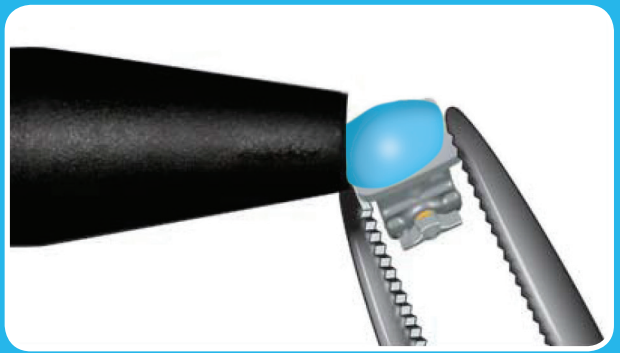Extrude small amount of blūgloo adhesive paste onto the bracket pad. Wear gloves when handling this product.

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(Side 2)

Step 3

(A) Immediately place bracket,
(B) position,
(C) press firmly,
(D) remove excess adhesive.

Note: Finalize position and begin light curing within 1 minute. Blue color will disappear when adhesive approaches body temperature.

Step 4


Light Cure:
10-15 seconds for clear brackets
20 seconds for metal brackets
30 seconds for molars due to increased pad size

Cleanup Feature:

(A) To easily see and remove any remaining blūgloo adhesive after debonding, cool the adhesive with short blasts of air or water to revert it back to its original blue color.


(B) Remove remaining blūgloo adhesive with high or low-speed carbide burrs followed by sand-paper discs and a polishing cup or paste.

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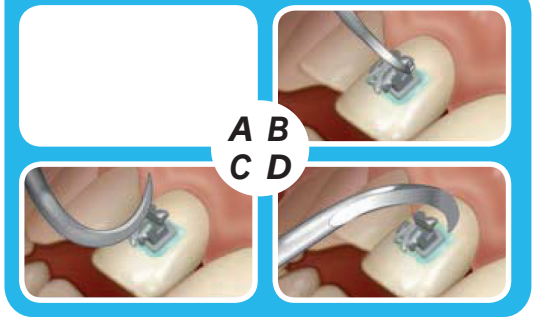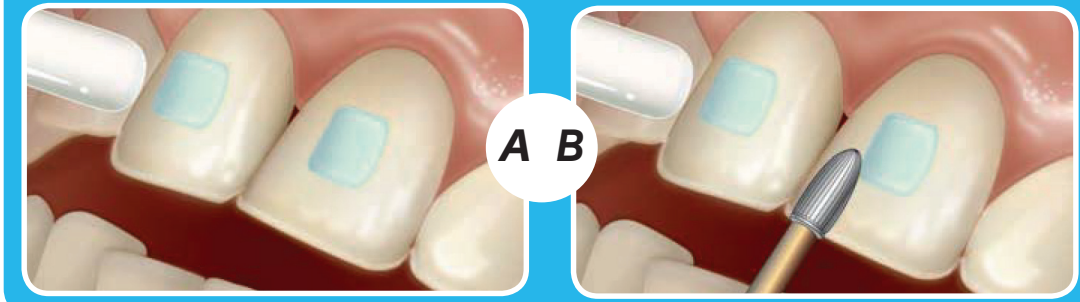ENGLISH

grēngloo™
Two –Way Color-Change Adhesive
For Metal Brackets

Instructions for Use

grēngloo ™ is a two-way color change adhesive for metal brackets. **grēngloo** contains a temperature-sensitive dye that lightens when it approaches body temperature and reverts to a green color when it falls below body temperature. The green color makes any excess or remaining adhesive easier to see when bonding and debonding.

Contra indications

Due to its dark (green) color when below body temperature this adhesive is not intended for use with clear brackets.

Tooth preparation

- Place lip/cheek retractor and salvia ejector in the mouth.
- Applying etching solution to the bonding surfaces with applicator sponges. Dab gently and continuously for 30 seconds per tooth (one arch at a time)- do not rub. Rinse thoroughly for a minimum of 5 seconds per tooth (10 seconds per tooth for gel etch) with a forceful airway spray and suction using a high-speed evacuator.
- Dry the etched enamel with clean, dry air. Maintain a dry field. Do not allow moisture, saliva or mucosal tissue to contact the etched enamel.
- Inspect to ensure a uniform, dull frosty appearance.

Bonding Procedure

- Dispense small amount of **Ortho Solo**™ sealant on the mixing pad. Replace the cap immediately after dispensing. Apply very thin, uniform coat of **Ortho Solo** to prepared enamel surface using the brushes provided.
- Extrude a small amount of **grēngloo** adhesive paste onto the base of the bracket pad. After application, clean the tip and replace the cap.
- Place bracket on the tooth surface immediately after applying the adhesive.
 - Note: All light-cured adhesives have a finite working time. The amount of working time depends on the lighting conditions (ambient and operatory) within each office.

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- Position bracket, and press down firmly. A slight excess of **grēngloo** (green tint) should be visible around the periphery of the bracket.
- Remove excess adhesive from the bracket base periphery.
- Cure bracket on tooth making sure to work the tip of light around periphery of bracket while curing. Cure times are as follows:
 - Conventional LED curing light (e.g. L.E.Demetron) - 20 seconds for metal brackets/ 30 seconds for molar tubes
 - Conventional Halogen curing light (e.g. Opitlux 501) -5 seconds for metal brackets/ 10 seconds for molar tubes
 - Note: Longer cure times may be required with curing lights from other manufacturers.
- Archwires maybe inserted after all brackets are cured.

Note: The color change of **grēngloo** is caused by a change in temperature. It is not an indication of the degree of polymerization.

Storage
- Do not expose materials to intense light or elevated temperatures.
- **grēngloo** should be stored at ambient temperature.

Note:
1. Unidose™ Delivery – for single patient use only, to prevent cross contamination between patients.
2. Unidose Composite Gun Sterilization Guidelines – Steam Autoclave: 30 Minutes @ 132°C (270°F).

Caution
1. Etching Solution – contains 37% phosphoric acid solution. Avoid contact with skin, eyes, and soft tissue. In case of contact with skin or eyes, flush immediately with water, get medical attention for eyes. Do not take internally.
2. Sealant and Paste – avoid prolonged contact with oral tissues and skin. Wash skin with soap and water after contact. Uncured materials may cause contact dermatitis. If skin irritation occurs, discontinue use. Use in a well-ventilated area. Keep containers closed when not in use. Wear gloves when handling this product.

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Step 2

Extrude small amount of grēngloo adhesive paste onto the bracket pad. Wear gloves when handling this product.

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Note: Finalize position and begin light curing within 2 minutes. Green color will disappear when adhesive approaches body temperature.

Step 4 – Light Cure

Conventional LED light:
10 seconds for clear brackets,
20 seconds for metal brackets and
30 seconds for molar tubes.


Halogen light:
5 seconds for all bracket types,
10 seconds for molar tubes.

Cleanup Feature

(A) To see remaining grengloo adhesive after debonding, cool the adhesive with short blasts of air or water to revert it back to its original green color.


(B) Remove remaining grengloo adhesive with high- or low-speed carbide burrs followed by sand-paper discs and a polishing cup or paste.

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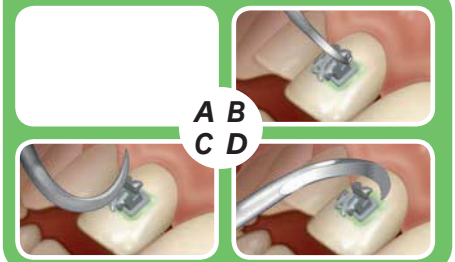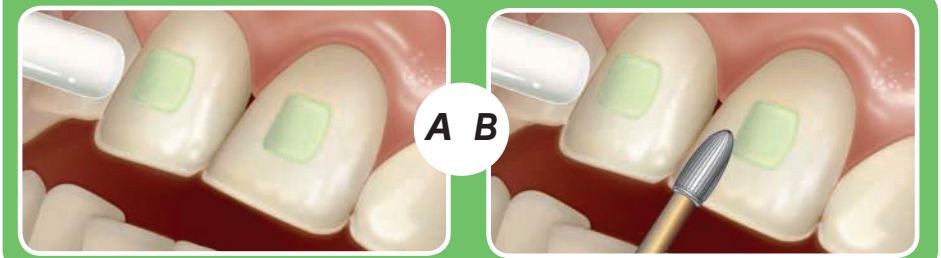ENGLISH

Tooth Preparation

- Place the lip/cheek retractor and saliva ejector in the mouth. Prophy the teeth thoroughly with pumice and water. Rinse with water while using a high-speed evacuation system. Dry the teeth with an air syringe.
- Apply the etching solution to the bonding surfaces with applicator sponges. Gently dab the etching solution: do not rub. A continuous dabbing of fresh solution for 30 seconds is most effective. Rinse thoroughly with a forceful air-water spray and suction using a high-speed evacuator. Maintain a dry field. Do not allow the patient to rinse or permit contact of the etched enamel with saliva or mucosal tissue.
- Dry the etched enamel with clean, dry air. The etched enamel must now have a uniform, dull, frosty appearance.

Bonding Procedure

1. Dispense a small amount of sealant on the mixing pad. Apply Ortho Solo™ sealant to teeth using the brushes provided. A thin, uniform coating should extend from the gingival crest and cover all of the etched enamel. Replace the cap immediately after dispensing sealant.
2. Extrude a small amount of Enlight adhesive paste onto the base of the bracket pad. After application, clean the tip and replace cap.
3. Place the bracket onto the tooth surface immediately after applying adhesive.
4. Adjust the bracket to its final position and press firmly.
5. Remove the excess adhesive from the periphery of the bracket base.

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6. With Demetron's Optilux 501, expose the adhesive for 10 seconds for metal brackets and 5 seconds for ceramic brackets. Shining the light directly through the bracket bonds ceramic brackets. Divide the curing time and shine the light on the mesial and distal edges of the bracket to bond metal brackets.
7. Active archwires may be placed immediately after curing the last bracket.

Shelf Life and Storage

- Do not expose materials to intense light or elevated temperatures.
- Enlight should be used at ambient temperature.

⚠ Caution

1. Etching Solution:
 Contains 37% phosphoric acid solution. Avoid contact with skin, eyes, and soft tissue. In case of contact with skin or eyes, flush immediately with water, get medical attention for eyes. Do not take internally.
2. Sealant and Paste:
 Avoid prolonged or repeated contact with oral tissues and skin. Wash skin with soap and water after contact. Uncured materials may cause contact dermatitis. If skin irritation occurs, discontinue use. Use in a well ventilated area. Keep containers closed when not in use.

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


ENLIGHT™

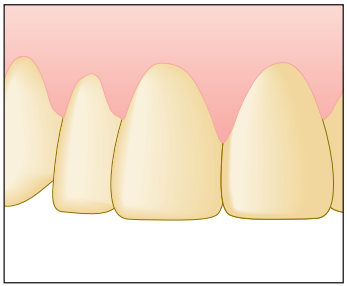1. Arch before enamel preparation.

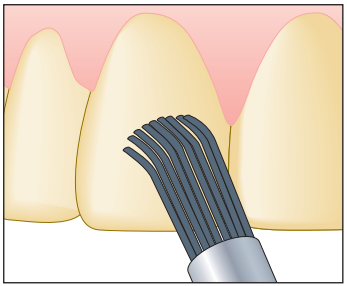2. Dab conditioner* on teeth for 30 seconds

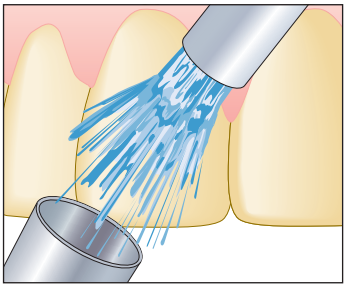
3. Rinse conditioner* off teeth (5 seconds per tooth).


*Caution: Conditioner contains phosphoric acid, avoid contact with skin, eyes and soft tissue. In case of contact with skin or eyes, flush immediately with water, get medical attention for eyes. Do not take internally.

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
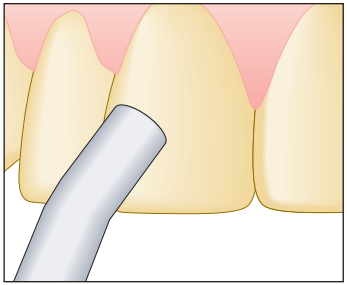
4. Air dry enamel until dull and frosty

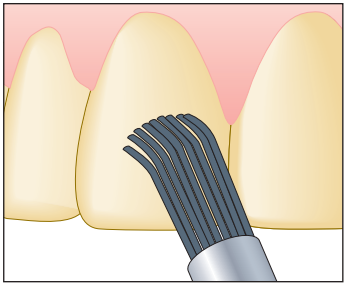
5. Apply Ortho Solo Bond Enhancing Primer over conditioned enamel.

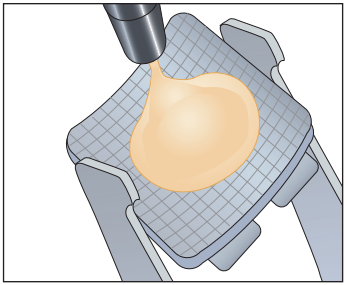
6. Syringe Enlight paste on bracket base.

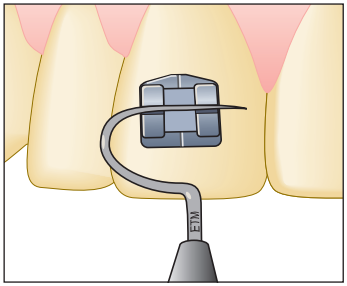7. Position bracket on tooth, clean excess and press firmly.

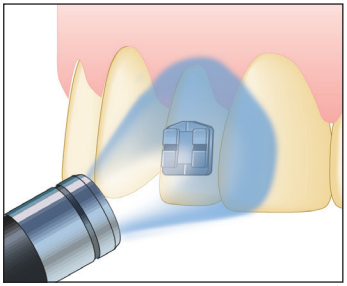8. Light cure the adhesive –
Demetron Optilux 501:
8-10 sec. metal, 5-sec. ceramic.
Conventional light:
30 sec. metal, 10-sec. ceramic.

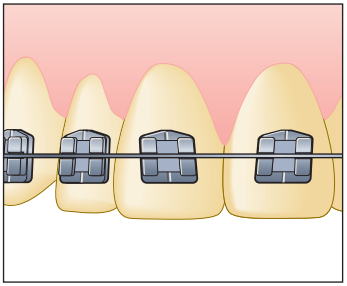9. Place archwire immediately after curing.

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English

Universal Bonding Primer for use with Chemical and Light Cure Adhesive

Instructions

Tooth Preparation

• Place the lip/cheek retractor and saliva ejector in the mouth. Prophy the teeth thoroughly with pumice and water. Rinse with water while using a high-speed evacuation system. Dry the teeth with an air syringe.

• Apply the etching solution to the bonding surfaces with applicator sponges. Gently dab the etching solution: do not rub. A continuous dabbing of fresh solution for 30 seconds is most effective. Rinse thoroughly with a forceful air-water spray and suction using a high-speed evacuator. Maintain a dry field. Do not allow the patient to rinse or permit contact of the etched enamel with saliva or mucosal tissue.

• Dry the etched enamel with clean, dry air. The etched enamel must now have a uniform, dull, frosty appearance.

Bonding Procedure

• Dispense a small amount of primer on the mixing pad. Apply Ortho Solo to teeth using the brushes provided. A thin, uniform coating should extend from the gingival crest and cover all of the etched enamel. **No air dry or cure step is necessary.** Replace the cap immediately after dispensing primer.

Note: Ortho Solo can be used with chemical and light cured adhesive systems. With light cure and 2-paste chemically cured systems, Ortho Solo takes the place of the primer. With chemically cured products like Ormco's System 1+, Ortho Solo should be used in addition to the activator.

Adhesive Bonding: Follow the manufacturer's instructions for proper bonding procedures.

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prolonged or repeated contact with oral tissues and skin. Wash skin with soap and water after contact. Uncured materials may cause contact dermatitis. If skin irritation occurs, discontinue use. Use in a well ventilated area. Keep containers closed when not in use.

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


OrthoSOLO™

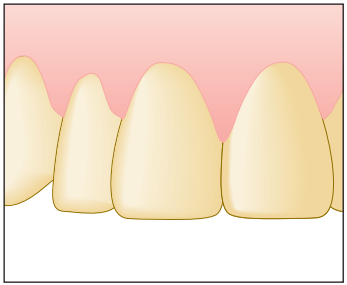1. Arch before enamel preparation.

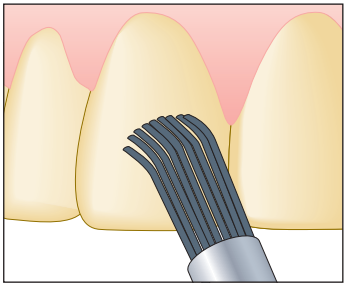
2. Dab conditioner on teeth for 30 seconds

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
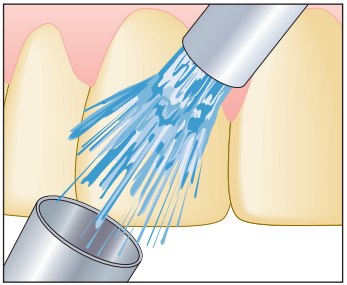
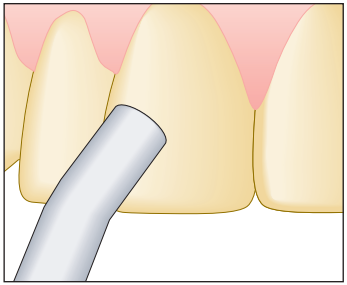
5. Dispense 1 drop of OrthoSolo into mixing well (no air dry or cure step necessary).

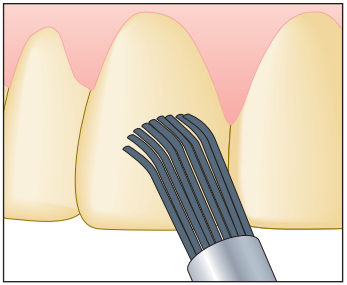6. Apply OrthoSolo over conditioned enamel (no air dry or cure step necessary).

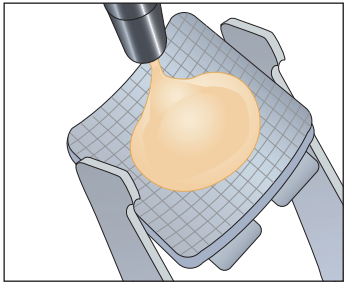7. Follow manufacturer's instructions for proper bonding (OrthoSolo takes the place of the regular sealant with all systems except one-step chemical cure).

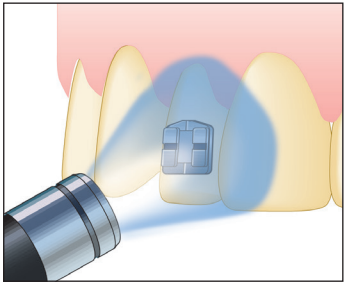8. If light curing with a Demetron 501, cure the adhesive 5 sec./clear brackets or 10 sec./metal brackets. Follow manufacturer's instructions for other light cure units.

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