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ver.01 from 20.02.2020

Instructions for the use of nonsterile wound plasters

Direction for use:
1. Make sure that the corn and the area around it are clean and dry.
2. Remove protective strips from the patch without touching the adhesive surface.
3. Place the patch firmly on the corn.
4. Leave the patch in place until it starts to peel off. The patch lasts up to 48 hours.
5. For best results, use until the skin is completely healed.
6. Replace the patch once it turns white, as this indicates it is completely filled.

Models:
0308-010 FootExpert Hydrocolloid plaster, 6 pcs. Size: 2.8x4.6 cm (oval).
0308-011 FootExpert Hydrocolloid plaster, 6 pcs. Size: 3.7x5.5 cm (oval).
0308-012 FootExpert Hydrocolloid plaster, 8 pcs. Size: 2.2x4.1 cm (oval).
0308-013 FootExpert Hydrocolloid plaster, 6 pcs. Size: 2.2x4.1 cm. (oval) – 3 pcs, 2.9x6.0 cm (rectan[illegible] pcs..

Storage conditions: Keep in a dry place at 6 to 25°.

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ver.01 from 20.02.2020

Instructions for the use of nonsterile corn plasters

Direction for use: Steam and thoroughly dry the foot. Remove the protective film from the patch. Place the pad with the active substance on a dry corn and press the patch for fixation. Repeat this procedure daily until the corn is completely removed. Replace daily. Do not reuse.

Models:
0108-011 Master Uni Corn plaster for removal of dry corns, 6 pcs. Size: 6,5x2,0 cm.

Storage conditions: Keep in a dry place at 5° to 30°C and relative air humidity of 65% away from direct sunlight, away from sources of heat and away from children's reach.

Shelf live: 5 years

Manufacturer: PharmLine s.r.o., Dobronicka 1257, Kunratice, 148 00 Praha 4, Czech Republic. Production place: No. 668 ChanHua Road, Fotang Town Industrial Functional Area, 322002 Yiwu City, Zhejiang Province, PEOPLE'S REPUBLIC OF CHINA

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ver.01 from 20.02.2020

Instructions for the use of Non-Sterile Capsicum Plasters

Directions: open the package, remove the protective paper strip, fix the patch. Effective for 24-48 hou
in the area of application must be clean and dry.

Models:
0107-001 Master Uni Capsicum plaster, 6x10 cm.
0107-002 Master Uni Capsicum plaster, 10x18 cm.

Storage conditions: Keep in a dry place at 5° to 30°C and relative air humidity of 65% away from direct sunlight, away from sources of heat and away from children's reach.

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ver.01 from 20.02.2020

Instructions for the use of sterile wound plasters

Direction for use: wash hands, clean and dry wound and surrounding area thoroughly, without touching a sorptior pad with fingers, position an adhesive plaster over wound and press down the edges firmly without creasing. Change daily. Do not reuse.


Models:
0108-002 Sterile wound plasters Master Uni Travel Kit. Polymeric base, 72x19 mm, 20 pcs.
0108-003 Sterile wound plasters Master Uni Sensitive. Non-woven base, 72x19 mm, 20 pcs.
0108-005 Sterile wound plasters Master Uni Elastic. Textile base, 72x19 mm, 56x19 mm, Ø22 mm, 20 pcs.
0108-007 Sterile wound plasters Master Uni Invisible. Polymeric base, 72x19 mm, 20 pcs.
0108-008 Sterile wound plasters Master Uni Kids. Polymeric base, 56x19 mm, 20 pcs.
0208-001 Sterile wound plasters Teneris Fun Kids. Polymeric base, 56x19 mm, 20 pcs.
0208-002 Sterile wound plasters Teneris Comfort Flex. Polymeric base, 76x19 mm, 20 pcs.
0208-003 Sterile wound plasters Teneris Sensitive. Non-woven base, 76x19 mm, 20 pcs.
0208-005 Sterile wound plasters Teneris Elastic. Textile base, 76x19 mm, 20 pcs.
0208-006 Sterile wound plasters Teneris Universal. Polymeric base, 76x19 mm, 20 pcs.
0208-022 Sterile wound plasters Teneris Aqua Silver. Polymeric base, 76x19 mm, 20 pcs.
0211-001 Sterile wound plasters with pad Teneris. Non-woven base. Size: 6x7 cm.
0211-002 Sterile wound plasters with pad Teneris. Non-woven base. Size: 8x10 cm.
0211-003 Sterile wound plasters with pad Teneris. Non-woven base. Size: 8x15 cm.
0211-004 Sterile wound plasters with pad Teneris. Non-woven base. Size: 10x10 cm.
0211-005 Sterile wound plasters with pad Teneris. Non-woven base. Size: 10x15 cm.
0211-006 Sterile wound plasters with pad Teneris. Non-woven base. Size: 10x20cm.
0211-007 Sterile wound plasters with pad Teneris. Non-woven base. Size: 10x25cm.
0211-008 Sterile wound plasters with pad Teneris. Non-woven base. Size: 10x30 cm.
0211-029 Sterile wound plasters with pad Teneris. Non-woven base, 10 pcs. Size: 6x7cm.
0211-030 Sterile wound plasters with pad Teneris. Non-woven base, 10 pcs. Size: 8x10cm.
0211-032 Sterile wound plasters with pad Teneris. Non-woven base, 10pcs. Size: 7x5cm.
0211-033 Sterile wound plasters with pad Teneris. Non-woven base, 10pcs. Size: 8x15cm.

Storage conditions: Keep in a dry place at 4 to 25° and relative air humidity of 65% away from direct sunlight, away from sources of heat and away from children's reach.

Shelf live: 5 years

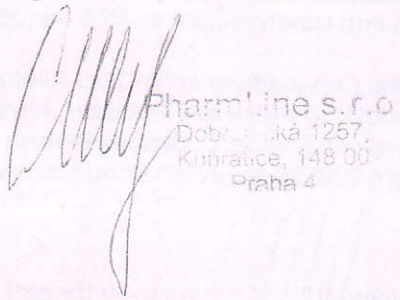Manufacturer: PharmLine s.r.o., Dobronicka 1257, Kunratice, 148 00 Praha 4, Czech Republic. Production place: Zhejiang Hongyu Medical Commodity Co., Ltd., No. 668 ChanHua Road, Fotang Town Industrial Functional Area, 322002 Yiwu City, Zhejiang Province, PEOPLE'S REPUBLIC OF CHINA

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