COVIDIEN™


SurgiSleeve™

Wound Protector

SurgiSleeve™

Wound Protector with Retraction Ring

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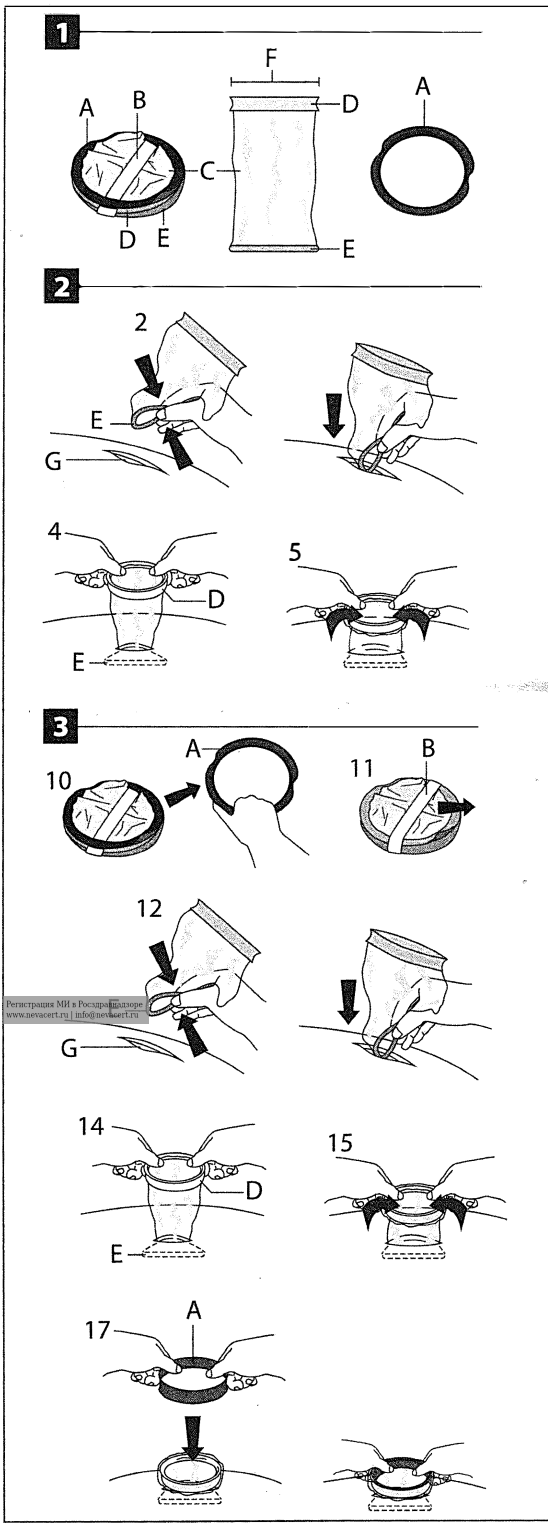EN

BEFORE USING PRODUCT, READ THE FOLLOWING INFORMATION THOROUGHLY.

IMPORTANT!

This booklet is designed to assist in using this product. It is not a reference to surgical techniques.

This device was designed, tested and manufactured for single patient use only. Reuse or reprocessing of this device may lead to its failure and subsequent patient injury. Reprocessing and/or resterilization of this device may create the risk of contamination, patient infection and/or device malfunction. Do not reuse, reprocess or resterilize this device.

DESCRIPTION

The SurgiSleeve™ wound protector is a one piece design consisting of a sheath film with a ring on both ends. The retraction ring is a separate piece, which is available on the large and extra large sizes only.

INDICATIONS

The SurgiSleeve™ wound protector is indicated for use to access the abdominal cavity during surgery through an atraumatically retracted incision, deliver maximum exposure of the abdominal cavity with minimum incision size, and protect against wound contamination during laparoscopic and open surgery. Additionally, the small size wound protector is indicated for use to access the thoracic cavity during cardiac and general surgical procedures through an atraumatically retracted incision. The extra small wound protector is also indicated for use to access the thoracic cavity and other soft tissue retraction during cardiac and general surgical procedures through an atraumatically retracted incision.

CONTRAINDICATIONS

Use of the SurgiSleeve™ wound protector is contraindicated when, in the judgment of the physician, use of this device would be contrary to the best interest of the patient.

WARNINGS AND PRECAUTIONS

Use care when introducing or removing sharp-edged or sharp-angled endoscopic instruments to minimize the potential of inadvertent damage to the SurgiSleeve™ wound protector.

This device is provided STERILE and is intended for use in a SINGLE procedure only. DISCARD AFTER USE. DO NOT RESTERILIZE.

Anesthesiologists should carefully monitor for hypotension when a retractor with an inner-ring is used in lower abdominal surgery.

1 SCHEMATIC VIEW

A) RETRACTION RING (available for large and extra large sizes only)

B) SHIPPING WEDGE

C) SLEEVE

D) BLUE PROXIMAL RING

E) GRAY DISTAL RING

F) SurgiSleeve™ WOUND PROTECTOR

INSTRUCTIONS FOR USE

2 SurgiSleeve™ wound protector without retraction ring

1. Prepare the abdominal or thoracic cavity for the SurgiSleeve™ wound protector by making an incision adequate to accommodate the wound protector.

2. Insert the gray distal ring of the SurgiSleeve™ wound protector into the incision.

G) INCISION

3. Examine and sweep the incision area to ensure that tissue is not caught between the ring and the abdominal wall.

4. Grasp the blue proximal ring in a manner that your hands are opposite to each other. Pull up on the blue proximal ring ensuring that the gray distal ring is fully seated against the peritoneal layer.

D) BLUE PROXIMAL RING

E) GRAY DISTAL RING

5. Roll the blue proximal ring inward by curling the edge over until the ring fully inverts.

6. Repeat step 5 until the incision is retracted and the SurgiSleeve™ wound protector is tight against the walls of the incision.

Removing the SurgiSleeve™ wound protector

7. If necessary, unroll the blue proximal ring of the SurgiSleeve™ wound protector.

8. Insert hand or finger through the opening, grasp gray distal ring and pull gray distal ring out of the incision.

3 SurgiSleeve™ wound protector with retraction ring

9. Prepare the abdominal cavity for the SurgiSleeve™ wound protector by making an incision adequate to accommodate the wound protector.

10. Remove the retraction ring from the SurgiSleeve™ wound protector.

A) RETRACTION RING

11. Remove the shipping wedge from the SurgiSleeve™ wound protector.

B) SHIPPING WEDGE

12. Insert the gray distal ring of the SurgiSleeve™ wound protector into the incision.

E) GRAY DISTAL RING

G) INCISION

13. Examine and sweep the incision area to ensure that tissue is not caught between the ring and the abdominal wall.

14. Grasp the blue proximal ring in a manner that your hands are opposite to each other. Pull up on the blue proximal ring ens[illegible] fully seated against the peritoneal layer.

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D) BLUE PROXIMAL RING

E) GRAY DISTAL RING

15. Roll the blue proximal ring inward by curling the edge over until the ring fully inverts.

16. Repeat step 15 until the incision is retracted and the SurgiSleeve™ wound protector is tight against the walls of the incision.

17. Insert the retraction ring into the blue proximal ring for additional support and to provide increased retraction of the incision. Ensure that the legs on the retraction ring is facing up.

A) RETRACTION RING

18. Ensure that the retraction ring is fully inserted 360 degrees around the proximal ring before proceeding.

Removing the SurgiSleeve™ wound protector with retraction ring

19. Remove the retraction ring from the blue proximal ring. Use the finger cutouts on the retraction ring to facilitate removal.

20. Unroll the blue proximal ring of the SurgiSleeve™ wound protector.

21. Insert hand or finger through the opening, grasp gray distal ring and pull gray distal ring out of the incision.

Incision Size Chart

Wound Protector Size	Incision Size
ExtraSmall	2 - 4 cm
Small	2.5 - 6 cm
Medium	5 - 9 cm
Large	9 - 14 cm
Extra Large	11 - 17 cm

STORE AT ROOM TEMPERATURE. AVOID PROLONGED EXPOSURE TO ELEVATED TEMPERATURES.

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