
Box size: 5.38 x 4.88
Overall size: 15 x 12
Folds to: 3 x 4

COVIDIEN™

Auto Suture™
Premium Extractor

1302756

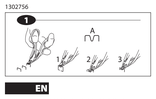EN

BEFORE USING PRODUCT, READ THE FOLLOWING INFORMATION THOROUGHLY.
IMPORTANT!
This handle is designed to assist in using the product. It is not a reference to surgical technique.
The device was designed, tested and manufactured for single patient use only. Reuse or reprocessing of this device may lead to its failure and subsequent patient injury. Reprocessing and/or resterilization of this device may create the risk of contamination and patient infection. Do not reuse, reprocess or resterilize this device.
MATERIALS AND PRECAUTIONS
This device was manufactured STERILE and is intended for use in a SINGLE procedure only. DO NOT ATTEMPT TO RE-STERILIZE.
1. Do not attempt to remove the staples.
INSTRUCTIONS FOR USE
NOTE: Do not pull up on staple before handle is fully squeezed.
1. Slip lower jaws of extractor under the backspan of the staple to be removed.
2. Close the instrument handles. You will unlatch the staple.
3. Lift the extracted staple away from the skin and discard. Proceed to remove the next staple from the skin.
4) STERILIZATION
AVOID PROLONGED EXPOSURE TO ELEVATED TEMPERATURES.

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
Covidien

Multifire Premium™

Pistol Grip Skin Stapler and Reload

1.822738

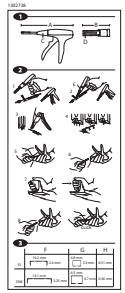EN

REVIEW STERILE PRODUCT IF SEAL OR PACKAGING INTEGRITY COMPROMISED.

IMPORTANT

Read the entire Instructions for Use carefully before using the product. Do not use if the package is opened or damaged.

The Multifire Premium™ Skin Stapler and reload are supplied sterile and are intended for single use only. Discard after use. Do not resterilize. Do not use after the expiration date. Store in a cool, dry place away from direct sunlight.

DESCRIPTION

The Multifire Premium™ Skin Stapler and reload are supplied sterile for external anastomotic purposes. The staples are manufactured from titanium alloy. Each reload contains 35 staples in 5 rows. The Multifire Premium™ Skin Stapler and reload are latex free.

EXCLUSIONS

Do not use the Multifire Premium™ Skin Stapler and reload for internal anastomotic purposes.

CONTRAINDICATIONS

Do not use if the package is damaged or has been opened.

WARNINGS AND PRECAUTIONS

- The Multifire Premium™ Skin Stapler and reload must be used in conjunction with the Multifire Premium™ Anvil. Do not use with any other anvil. Failure to use the Multifire Premium™ Anvil may result in staple malformation and/or failure to fire.
- Do not use the Multifire Premium™ Skin Stapler and reload if the package is damaged or has been opened.
- Do not use the Multifire Premium™ Skin Stapler and reload after the expiration date.
- Do not resterilize the Multifire Premium™ Skin Stapler and reload.
- Do not use the Multifire Premium™ Skin Stapler and reload for internal anastomotic purposes.
- Do not use the Multifire Premium™ Skin Stapler and reload if the sterile barrier is compromised.
- Do not use the Multifire Premium™ Skin Stapler and reload if the product is damaged.
- Do not use the Multifire Premium™ Skin Stapler and reload if the product is contaminated.
- Do not use the Multifire Premium™ Skin Stapler and reload if the product is expired.
- Do not use the Multifire Premium™ Skin Stapler and reload if the product is damaged.
- Do not use the Multifire Premium™ Skin Stapler and reload if the product is contaminated.
- Do not use the Multifire Premium™ Skin Stapler and reload if the product is expired.

DIRECTIONS FOR USE

1. Select the size of the Multifire Premium™ Skin Stapler and reload to be used.
2. Open the package and remove the Multifire Premium™ Skin Stapler and reload.
3. Remove the Multifire Premium™ Skin Stapler and reload from the package.
4. Remove the Multifire Premium™ Skin Stapler and reload from the package.
5. Remove the Multifire Premium™ Skin Stapler and reload from the package.
6. Remove the Multifire Premium™ Skin Stapler and reload from the package.
7. Remove the Multifire Premium™ Skin Stapler and reload from the package.
8. Remove the Multifire Premium™ Skin Stapler and reload from the package.
9. Remove the Multifire Premium™ Skin Stapler and reload from the package.
10. Remove the Multifire Premium™ Skin Stapler and reload from the package.

STORAGE

Store in a cool, dry place away from direct sunlight.

DISPOSAL

Dispose of in accordance with local regulations.

MANUFACTURER

Covidien

EC REPRESENTATIVE

Covidien

CONTACT INFORMATION

For more information, please contact your local Covidien representative.

LATEX STATEMENT

The Multifire Premium™ 35 reload is latex free.

SYMBOLS

	Sterile
	Do not resterilize
	Do not use if package is damaged
	Consult instructions for use
	Single use only
	Do not use if package is opened
	Do not use if sterile barrier is compromised
	Do not use if product is contaminated
	Do not use if product is expired
	Do not use if product is damaged

Autosuture™

APPOSE™ ULC

SINGLE USE SKIN STAPLER WITH AUTO RELEASE

10000-29718

295024

ENGLISH

BEFORE USING PRODUCT, READ THE FOLLOWING INFORMATION THOROUGHLY.

IMPORTANT!:
This booklet is designed to assist in using this product. It is not a reference to surgical techniques.

This device was designed, tested and manufactured for single patient use only. Reuse or reprocessing of this device may lead to its failure and subsequent patient injury. Reprocessing and/or resterilization of this device may create the risk of contamination and patient infection. Do not reuse, reprocess or resterilize this device.

DESCRIPTION
The APPOSE™ ULC single use skin stapler dispenses one staple each time the instrument trigger is activated. The stapler first penetrate the skin and are then formed, holding the tissue together.


INDICATIONS
The APPOSE™ ULC single use skin stapler can be used for the closure of skin in a wide variety of surgical procedures.

CONTRAINDICATIONS
1. Staples can be applied directly over bone or viscera; however, during application there must be a distance of not less than 6.5 mm from the surface of the skin to the underlying bone, vessel, or viscera. When using the regular skin staple, the distance required can be slightly less, but at least 4.5 mm. If these distances cannot be obtained by lifting the skin prior to application, use of the skin stapler is contraindicated.

WARNINGS AND PRECAUTIONS
1. The APPOSE™ ULC single use skin stapler is provided STERILE and is intended for use in a SINGLE procedure only. DISCARD AFTER USE. DO NOT RESTERILIZE.

1 SCHEMATIC VIEW
A) ANVIL
B) TRIGGER

2 INSTRUCTIONS FOR USE
1. For skin closure, approximate the skin edges in an everted fashion with the use of toothed forceps or a tissue approximating device. Skin may also be approximated by applying traction to both ends of incision. Make certain edges are everted.
2. The instrument nose is placed gently on the approximated skin at the desired site in the middle of the incision.
3. Squeeze the trigger completely and release. Failure to completely squeeze the trigger may result in incomplete staple formation and insufficient

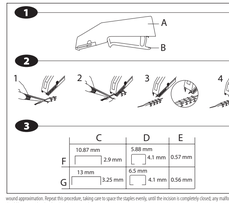wound approximation. Repeat this procedure, taking care to space the staples evenly, until the incision is completely closed; any malformed staples should be removed and replaced with a new staple.

3 STAPLE SPECIFICATIONS
D) STAPLE SIZE
D) APPROXIMATE CLOSURE
F) STAPLE DIAMETER
F) REGULAR
G) WIDE

STORE AT ROOM TEMPERATURE. AVOID PROLONGED EXPOSURE TO ELEVATED TEMPERATURES.

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