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Issue Date : ***24/03/2022***
Version : ***A/5***

Instruction for use

This product should be handled by professional dental technician, sit and adjust by dental specialist and follow doctor's instruction.

Product Name: Dental zirconia blank for aesthetic restoration

Intended use ：

The product is mainly used in prosthetic treatment, When the flexural strength of Dental zirconia blank for aesthetic restoration after sintering is between 300-500MPa, the product can be used for veneering, inlay, single crown, and substructure ceramic for three-unit prostheses not involving molar restoration, When the flexural strength of Dental zirconia blank for aesthetic restoration after sintering is over 500 MPa, it can be used for veneering, inlay, single crown, and substrate ceramic for three-unit prostheses. All blanks are processed through dental laboratories or by dental professionals.

Chemical compositions:

ZrO_2+HfO_2: 86.3%～94.2%
Y_2O_3: 5.8%～9.7%
Fe_2O_3: <0.5%
Al_2O_3: <0.5%
Er_2O_3: <2%
Other oxide: <0.5%

CTE（25-500℃）： （10.5±0.5）$\times 10^{-6}$ K^{-1}

Instruction for use:

1. Fix the zirconia blank in the specific clamping device.
2. Start the milling process.
3. Remove milled framework from blank and pre-trim with diamond disc.
4. Remove dust from the framework with air gun.
5. Color frame with coloring liquid, Final-dry frame under infrared lamp (see manual for coloring liquid).
6. Load into furnace.
7. Select sintering program and start sintering (see sintering schedule).
8. Remove sintered framework from furnace after the furnace temp drop to room temperature.
9. Check frame on working model for sitting and margin fit. Trim with grinder and polisher if required.
10. Blast using dia-50 micron Alumina at pressure 2-2.5 bar. Clean frame with brush using water. Do not steam-clean.
11. Frame is now ready for veneering, choose porcelain and glaze to finish the restoration.
12. Now the restoration is ready for the try-in.

Shelf period: 5-years is an estimated replacement time. The product is not expected to degrade during this time period.

Storage:

1. Keep container dry.
2. Store at room temperature.

NOTE: The product shall be used only by trained personnel. The manufacturer is not liable for any erroneous results, since we have no control if this instruction is properly followed or if the product is appropriately used. Any claims for damages therefore shall be limited exclusively to the commercial value of our products. Please make sure that the dentist has received all of the necessary information for the correct use of the product.

Medical Device Adverse Event: So far there is no medical device adverse event on this product.

Cautions*:* Federal Law restricts this device to sale by or on the order of a dentist.

Contraindications: Not found.

Safety instruction:

1. Avoid collision with hard material, avoid squeeze or shaking.
2. Don't apply to patient before the material is fully sintered.

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3. This product should only be handled by dental specialist.
4. Wear suitable face mask to prevent fine dust inhalation. Avoid fine dust air pollution by using the air blaster in short bursts only rather than a continuous flow.

Disposal: The material is not hazardous to the environment.

EC	REP

European representative:

Name: MedNet EC-REP GmbH

Address: Borkstrasse 10,48163 Münster, Germany.

DIMDI no.: *DE/0000048589*

Email: ecrep@medneteurope.com

Distributor: (please refer to labeling)

Manufacturer:

Name: Liaoning Upcera Co., Ltd.

Address: No.122 Xianghuai Road, Economic Development Zone 117004 Benxi,Liaoning China


Phone: +86 24-45565055

Fax: +86 24-45854457


Email: ln.upcera@upcera.com.cn

CE 0197			LOT	NON STERILE
	Date of manufacturer	Use-by date	Batch code	Non-Sterile
	#	UDI	MD	
Do not reuse	Model number	Unique device identifier	Medical device	

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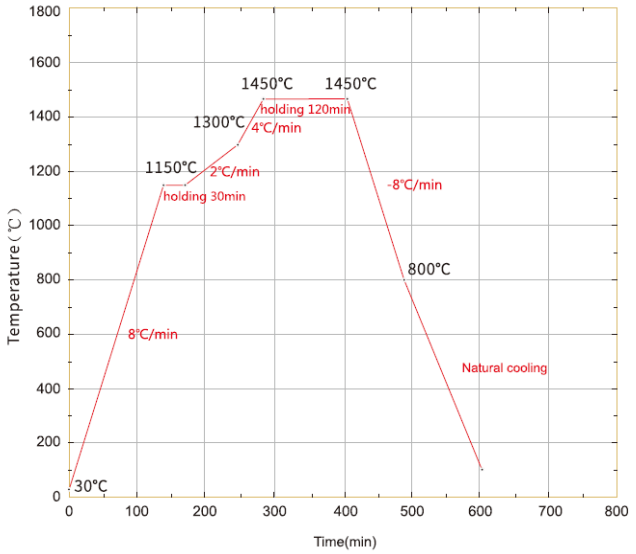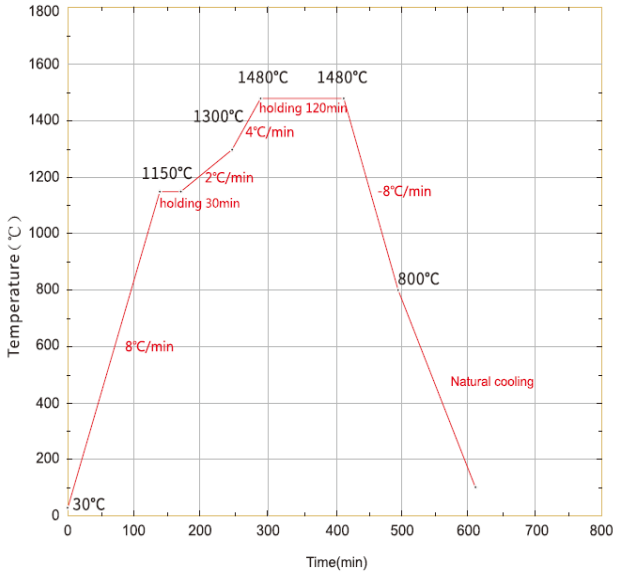Manufacturer:

Name: Liaoning Upcera Co., Ltd.

Address: No.122 Xianghuai Road, Economic Development Zone 117004 Benxi, Liaoning China

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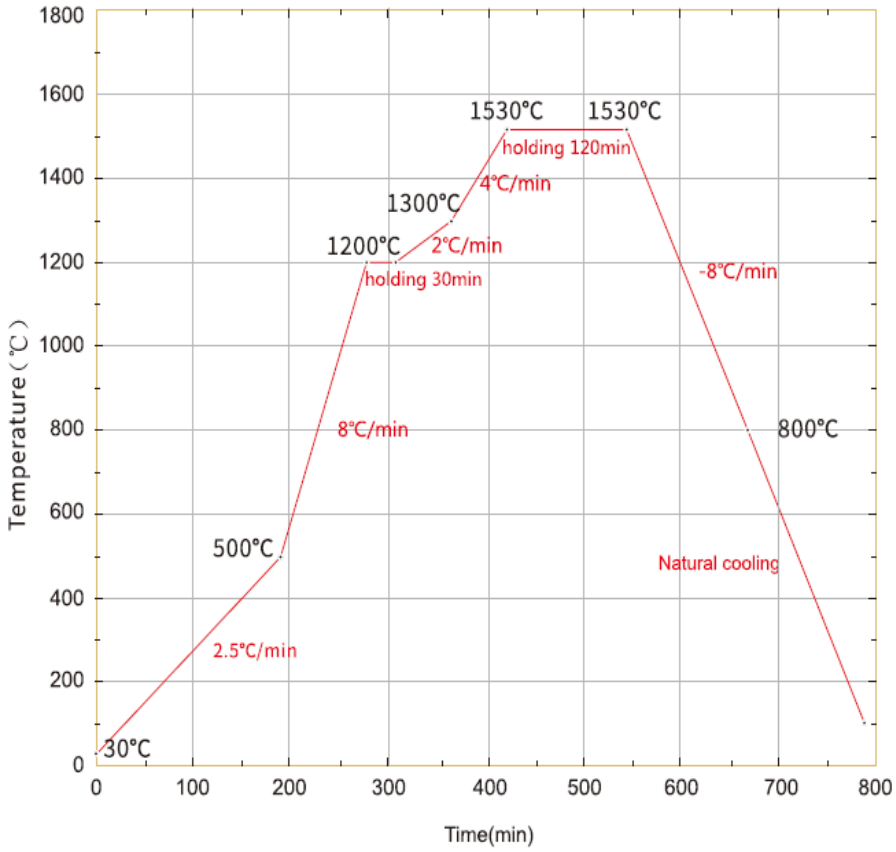
Product Name: Dental Zirconia Pre-shaded Blank

Intended use ：

Used for dental restoration using different CAD/CAM or manual milling machines. All Blanks are processed through dental laboratories or by dental professionals.

Chemical compositions:

Nanometer zirconia powder：＞98%

Fe_2O_3：＜0.3%

Pr_2O_3：＜0.2%

Er_2O_3：＜1%

Other oxide: ＜0.5%

CTE（25-500℃）： $10.5 \pm 0.5 \times 10^{-6}\ K^{-1}$

Instruction for use:

1. Fix the zirconia blank in the specific clamping device.
2. Start the milling process.
3. Remove milled crown or bridge from blank and pre-trim with diamond disc.
4. Remove the dust from the crown or bridge with air gun.
5. Color the crown or bridge with coloring liquid, Final-dry frame under infrared lamp (see manual for coloring liquid. Applicable to the zirconia which need to use coloring liquids).
6. Load into furnace.
7. Select sintering program and start sintering (see sintering schedule).
8. Remove sintered crown or bridge from furnace after the furnace temp drop to room temperature.
9. Check crown or bridge on working model for sitting and margin fit. Trim with grinder and polisher if required.
10. Clean with brush using water. Do not steam-clean.
11. Choose glaze to finish the restoration.
12. Now the restoration is ready for the try-in.

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NOTE: The product shall be used only by trained personnel. The manufacturer is not liable for any erroneous results, since we have no control if this instruction is properly followed or if the product is appropriately used. Any claims for damages therefore shall be limited exclusively to the commercial value of our products.

Please make sure that the dentist has received all of the necessary information for the correct use of the product.

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Cautions: Federal Law restricts this device to sale by or on the order of a dentist.

Contraindications: Not found.

Safety instruction:

1. Avoid collision with hard material, avoid squeeze or shaking.
2. Don't apply to patient before the material is fully sintered.
3. This product should only be handled by dental specialist.
4. Wear suitable face mask to prevent fine dust inhalation. Avoid fine dust air pollution by using the air blaster in short bursts only rather than a continuous flow.

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Email: ecrep@medneteurope.com
Distributor: (please refer to labeling)

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