AGU ADVANCED GROWING UP

Penguin

ASPIRATOR

NASAL

AGU NS19

en QUICK GUIDE

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en USER MANUAL

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NASAL ASPIRATOR **AGU NS19** EN

TABLE OF CONTENTS

1 Introduction 23
2 List of symbols 24
3 Important information before use 25
4 Product overview 27
5 Battery installation/replacement 28
6 Children's tunes 29
7 Instruction for use 29
8 Care and maintenance 31
9 Troubleshooting 32
10 Technical specification 33
11 Conditions of storage, transportation and operation 33
12 Technical data 34
13 Disposal 35
14 Certification 36
15 Manufacturer's warranty 36

1 INTRODUCTION

Nasal aspirator **AGU NS19** is specially designed to clear mucus in newborn babies and infants in home environment and/or while travelling.

Music function

12 tunes in the memory of this device may interest and comfort a baby and help to make the cleaning procedure more comfortable.

How to use

Insert aspiration Tip into a nostril, after that the device is ready to remove nasal secretions.

Safety

The soft aspiration Tip prevents the possibility for injury during operation.

Compact & easy design

This device is compact and easy to hold and operate.

Washable parts

Silicone Tips, Collection cup & head module can be removed from the device and are easily washed (Collection cup can also be disinfected in boiling water).

Light weight, low noise

The child calmly responds to the procedure of cleansing the nose.

Battery

This device uses two 1.5 V AA batteries; it is convenient to replace the batteries anytime, anywhere.

Please read this User's manual to the end before you start using the nasal aspirator **AGU NS19**.

2 LIST OF SYMBOLS

Symbols	Meaning
⚠	Warning/Note.
CE 0598	The product meets the basic requirements of: • The product meets the requirements of Directive 93/42/EEC concerning medical equipment; • Directive 2011/65/EU RoHS on the restriction of the use of certain hazardous substances in electrical and electronic equipment.
	WEEE (EC Directive on Waste Electrical and Electronic Equipment). This marking means that the device should not be disposed of together with household waste. In order to prevent a possible damage to the environment and harm to human health, such waste should be segregated from other types of waste and disposed of in accordance with the relevant accepted standards.
	Please read the instruction manual.
IP22	Device classification type BF. The degree of protection against insertion. The first figure (protection against insertion of foreign solid objects): 2 – protection against insertion of solid objects, the size of more than 12 mm, fingers or other objects with a length of not more than 80 mm, or solid objects. The second figure (protection against foreign liquids): 2 – protection against falling drops, falling objects from above at an angle to the vertical no more than 15° (the device in the normal position).
	Manufacturer.
EC REP	Manufacturer's authorized representative in the European Union.
CH REP	Authorized Representative in Switzerland.

3 IMPORTANT INFORMATION BEFORE USE


The device is designed to remove nasal secretion and mucus in children (age 0–12 years old) in at-home environment.

1. **WARNING!** Do not modify this device in any way without authorization of the manufacturer.
2. Do not expose the device to high / low temperature, very high humidity level. Avoid exposing the device to direct sunlight.
3. Avoid extreme stress, shock or accidental dropping of the device.
4. Please replace all batteries when the power and/or the sound of the pump weakens during the use of the device.
5. It is recommended to remove the batteries if the device is not used for a long period of time.
6. Please press and hold the Power button to start the device to generate vacuum for removing nasal secretion and mucus. If you want to stop the suction, simply release the Power button to turn off the device.
7. The nasal aspirator is for cleaning / removing the nasal secretion only. Do not use this device other than the intended purposes.
8. During the use of the device, the babies should not hold their breath, close their mouth or swallow their saliva.
9. If the baby does not feel well when using the device, stop aspiration and consult a doctor immediately.
10. Do not apply excessive force when placing the aspiration Tip into the nasal cavity.
11. This device is intended for use both at home and while travelling. Please consult your doctor before using the device in case of rhinitis and sinusitis and be aware of cleaning the aspiration Tip before each use.
12. The intended operators for this device shall be parents or adults who can follow the instructions and notices in User's manual to operate the device. The device is not allowed to be used by children themselves. Keep the nasal aspirator stored away from the children. Do not let them operate the device themselves.
13. **DO NOT** use the nasal aspirator when children have the following contraindications:
 - nasal bleeding;
 - acute head, facial, or neck injury;
 - coagulopathy or bleeding disorder.

14. Before using the device, please read this User's manual carefully. Incorrect use of the device may cause the following hazards/complications:
 - mechanical trauma;
 - nasal bleeding;
 - damage of mucous membrane of the nose.
15. Please keep the replaceable Tips out of the reach of children children and pets cannot access to avoid risk of them to be swallowed.
16. The Tip may cause potential allergic reactions.
17. Clean the device after use, avoid decontamination/ bacterial infection.
18. If any serious incident that has occurred in relation to the device should be reported to the manufacturer, local representative and the competent authority of the Member State in which the user is established.


⚠ This nasal aspirator must only be used with AGU NS series Tip. Please contact shops and pharmacies to buy those accessories.

4 PRODUCT OVERVIEW

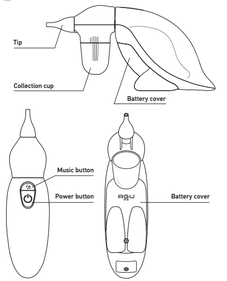5 BATTERY INSTALLATION/REPLACEMENT

1. Unscrew the screw.
2. Gently slide the Battery cover back.
3. Carefully remove the old batteries and properly dispose them.
4. Insert or replace batteries into the battery compartment, ensuring to match the indicated polarity symbols. Always use new batteries.
5. Replace the Battery cover.
6. Screw in the screw.

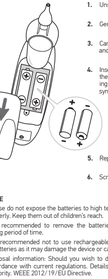NOTE

- Please do not expose the batteries to high temperatures and dispose them properly. Keep them out of children's reach.
- It is recommended to remove the batteries if the unit is not used for a long period of time.
- It is recommended not to use rechargeable, unqualified or different spec of batteries as it may damage the device or cause circuit shortcut.

Disposal information: Should you wish to dispose of the article, do so in accordance with current regulations. Details are available from your local authority. WEEE 2012/19/EU Directive.

6 CHILDREN'S TUNES

If you want to soothe your baby, press the Music button.

- Nasal aspirator **AGU NS19** offers 12 well-known melodies.
- Pleasant melodies will soothe your baby while the device is in use.
- To select another melody, press the Music button again.
- To turn off the music, press and hold the Music button.

7 INSTRUCTION FOR USE

7. Use the right Tip for your baby's nostrils.
 - The large Tip is recommended for use when the mucus is thick.
 - The small Tip is recommended for use when the mucus is loose and liquid.

- Dispose of the empty, fully discharged batteries at special collection points, recycling centers or consumer electronics retailers. You are legally responsible for battery disposal.

8. Make sure the Tip is properly attached to the suction unit. Before you start the device, check that the Tip is properly attached to the head module.

⚠ Check the Tip before each use. If you notice any signs of wear or damage, replace it immediately. The Tip is a consumable and its lifetime depends on the device's usage intensity. We recommend replacing the Tip every six months.

9. Press and hold the Start-Tast while the device is aspirating mucus.

10. Hold the baby upright and insert the Tip into one nostril. Be careful not to push against the nasal septum or insert the Tip too far into the nostril. Hold the nasal aspirator horizontally so that the Tip forms a 90° angle with the baby's nasal cavity. Mucus will be removed by suction without additional external force. Repeat the procedure in the baby's other nostril if necessary.

⚠ Mucus will not appear in the collection cup immediately. We recommend rinsing the Tip after each use by detaching it from the suction unit.

11. Release the Start-Tast to stop suction.

⚠ **CAUTION**

To prevent liquid from entering the device's main unit and causing it to malfunction:

- during the procedure, make sure the mucus in the collection cup does not exceed the maximum level (the maximum level is marked on the collection cup).
- after the procedure, detach the collection cup module from the main unit.

If used correctly, the device is protected against liquid entering the main unit. Please follow the instructions above. Keep in mind that liquid entering the device will void the warranty.

8 CARE AND MAINTENANCE

12. Separate the Tip and the head module from the suction unit by carefully pulling them apart.
13. Unscrew the collection cup by turning it counterclockwise.

CLEANING THE TIP, HEAD MODULE AND COLLECTION CUP

We recommend washing the collection cup, head module and Tips with mild soap in warm water for 2 minutes. Then rinse the accessories with warm tap water for 3 minutes. After cleaning, air-dry the accessories.

⚠ Do not wash the device in a dishwasher. Please clean the collection cup, head module and Tips with water and mild soap only.

CLEANING THE DEVICE

We recommend wiping the device's main unit with a clean, soft, dry towel and keeping it clean.

⚠ **Never immerse the main unit in water or other cleaning liquids.**

9 TROUBLESHOOTING

If you notice any problems while using the device, please check the following table first.

Problem	Possible cause	Action
Weak or no suction	Weak battery power or batteries are incorrectly installed	Check batteries or replace them
	No rubber ring (sealing) at the collection cup	Make sure there is a rubber ring (sealing) at the collection cup
	The collection cup is not properly attached	Make sure the vertical guides on the collection cup fit into the guides on the suction unit

10 TECHNICAL SPECIFICATION

Product name	Nasal aspirator **AGU NS19**
Type	**AGU NS19**
Effective vacuum pressure:	350–500 mmHg (±5%)
Noise level	55–60 dB (30 cm from the **AGU NS19**)
Tips	soft silicone
Motor	3 V/DC
Operating environment	16 °C–35 °C (60.8 °F–95 °F) relative humidity 15%–85% (non-condensing) atmospheric pressure 700–1060 hPa
Storage/transportation environment	-25 °C–+55 °C (-13 °F–+131 °F) relative humidity 15%–85% (non-condensing)
Power supply	2 x 1.5 V AA Alkaline Batteries
Weight	150 g (without Batteries)
Dimensions (L×W×H), mm	170×43×86

Battery life	approx. 90 minutes (depending on battery quality)
Expected lifetime	device – 3 years, accessories – at least 6 months
Package contents*	nasal aspirator, 1.5 V AA batteries (2 pcs), silicone tips (2 pcs) small tip and large tip, transparent collection cup, storage bag, user manual, warranty card
Warranty	2 years

NOTE

- If the device does not work properly even after the above-mentioned measures, please contact your dealer or the nearest service center.
- Do not use the device if any of its parts is damaged or malfunctioning.

11 CONDITIONS OF STORAGE, TRANSPORTATION AND OPERATION

Nasal aspirator **AGU NS19** is intended for use in the electromagnetic environment specified below.

The customer or the user of **AGU NS19** should assure that it is used in such environment.

EMC table

NS19 is intended for use in the electromagnetic environment specified below. The customer or the user of **NS19** should assure that it is used in such environment

Guidance and manufacturer's declaration – electromagnetic emissions		
Phenomenon	Professional healthcare environment a)	HOME HEALTHCARE ENVIRONMENT a)
RF emissions CISPR 11	CISPR 11 Group 1 Class A (non-BLE function) Group 2 Class A (BLE function)	CISPR 11 Group 1 Class B (non-BLE function) Group 2 Class B (BLE function)
Harmonic distortions	Not applicable (Note: Power supply via battery or DC input)	
Voltage fluctuations/flicker emissions	Only AC input must be tested Not applicable (Note: Power supply via battery or DC input) Only AC input must be tested	

a) The device is intended for use in home healthcare environments and professional healthcare environments, limited to patient room and breathing apparatus in hospitals or clinics. The remaining requirements of CISPR 11 Group 1 Class B were observed and applied. The device is intended for use in the mentioned environments if connected directly to the public power supply network.


b) The test is not applicable in this environment, unless the used ME-EQUIPMENT and ME-SYSTEM is connected to the PUBLIC POWER NETWORK and the power consumption is within the EMC basic standard requirements.

Guidance and manufacturer's declaration – electromagnetic immunity – enclosure port			
Phenomenon	EMC basic standard or test method	Immunity test levels	
		Professional healthcare environment	HOME HEALTHCARE ENVIRONMENT
ELECTROSTATIC DISCHARGE	IEC 61000-4-2	± 8 kV Contact ± 2 kV, ±4 kV c, ±8 kV, ±15 kV Air	
Radiated RF EM fields	IEC 61000-4-3	ai	

The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.).


6 CHILDREN'S TUNES

If you want to comfort your baby, just press the Music button.

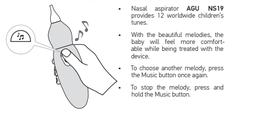- Nasal aspirator **AGU NS19** provides 12 worldwide children's tunes.
- With the beautiful melodies, the baby will feel more comfortable while being treated with the device.
- To choose another melody, press the Music button once again.
- To stop the melody, press and hold the Music button.

7 INSTRUCTION FOR USE

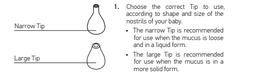1. Choose the correct Tip to use, according to shape and size of the nostrils of your baby.
 - The narrow Tip is recommended for use when the mucus is loose and in a liquid form.
 - The large Tip is recommended for use when the mucus is in a more solid form.
2. Be sure to place the Tip securely on the aspirator. Make sure the Tip is completely covering the rib before using.

Check the Tip before each use. If you notice any signs of wear or damage, change it immediately. The Tips belong to consumables, their service life depends on the frequency of use of the device. We recommend replacing Tips every six months of use.

4. Press and hold the Power button while the device is clearing the mucus.

5. Hold the baby upright and insert the Tip into one nostril of the baby. Be careful not to press against the inside wall of the nostril or plug the Tip too far into the nostril. Be sure to hold the aspirator horizontally keeping the Tip of the aspirator at a 90° angle to the nasal passage. Mucus is cleared by the aspirator without any additional external force. Repeat the same procedure in the other nostril of the baby. If necessary, the procedure may be repeated.

The mucus may not immediately appear in the Collection cup. After each use it is recommended to do a quick cleaning of the Tip by removing the Tip from the aspirator.

6. Release Power button to stop the suction.

ATTENTION


In order to prevent getting liquid inside the main unit of the device and its failure:

- during the procedure, make sure that the amount of mucus in the storage tank does not exceed the maximum level (the maximum level is marked by a special inset in the tank and the device clean.
- after completion of the aspiration procedure, immediately disconnect the head unit with the storage tank from the main unit.

If the device is used correctly, liquid will not get into the main unit. Please follow the recommendations above. Remember, getting liquid inside the device will lead to loss of warranty.

8 CARE AND MAINTENANCE

1. Separate the head module with the Tip from the base of the aspirator by gently pulling it apart.
2. Unscrew Collection cup by turning it clockwise.

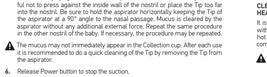CLEANING PROCESS OF COLLECTION CUP, HEAD MODULE AND ASPIRATION TIP

It is recommended to wash the Collection cup, head module and aspiration Tip with mild soap in warm water for 3 minutes. Then rinse each part with running hot tap water for 3 minutes. After the cleaning process, these parts must be completely dried after.

DO NOT wash this device with any kinds of of strong soap or cleaner. Please just wash the Collection cup, head module, and Tip with water and detergent.

DO NOT expose the device details to high temperatures (more than 60°C).

UNIT CLEANING

It is strongly recommended to wipe the unit with a clean soft dry cloth and keep the device clean.

DO NOT wash the housing body of this device with water or any kind of liquids.

9 TROUBLESHOOTING

If any trouble occurs while using the device, please check the following list first.

Problem	Possible Cause	Solution
The suction is weak or there is no suction at all	Batteries are low or improperly installed	Check or replace batteries
	Lack of gray rubber ring in the place where the Collection cup is screwed in	Make sure there is a gray rubber ring in the place where the Collection cup is screwed in
	The Collection cup is not properly screwed in	Make sure the three guides on each side of the Collection cup line up with the guides on the housing body of the aspirator
The music will not play when the Music button is pressed	Batteries are low or improperly installed	Check or replace batteries
The music will not stop when the Music button is pressed	Pressing of the button is not long enough	To stop the music, keep pressing the button until the music stops playing

10 TECHNICAL SPECIFICATION

Product	Nasal aspirator **AGU NS19**
Model	**AGU NS19**
Effective Vacuum Pressure	350–500 mmHg (±5%)
Sound pressure level	55–60 dB (30 cm away from the **AGU NS19**)
Aspiration Tip	soft silicone
Motor	3V/DC
Operating environment	16 °C–35 °C (60.8 °F–95 °F) relative humidity 15%–85% (non-condensing) atmospheric pressure 700–1060 hPa
Storage / Transportation environment	-25 °C–+55 °C (-13 °F–+131 °F) relative humidity 15%–85% (non-condensing)%
Power source	2x1.5 V AA alkaline batteries
Weight, g	150 (without batteries)
Dimensions (LxWxH), mm	170x43x86
Battery operation time	about 90 minutes (depends on the quality of batteries)
Typical service life	3 years
Completeness*	nasal aspirator, 1.5 V AA batteries (2 pieces), silicone tips (2 pieces), transparent collection cup, storage bag, user's manual, warranty coupon
Warranty period	2 years

Notice: the manufacturer has the right to introduce changes into the device design.

11 CONDITIONS OF STORAGE, TRANSPORTATION AND OPERATION

- The device should be transported at -25 °C–+55 °C (-13 °F)–+131 °F), relative humidity 15%–85% (non-condensing).
- The device can be operated at a temperature of 16 °C–35 °C (60.8 °F–95 °F), relative humidity 15%–85% (non-condensing), atmospheric pressure 700–1060 hPa.

WARNING

- If the device does not work even after taking the above-mentioned actions, please contact the store where it was purchased or the nearest dealer.
- Never turn the device on if any part of it is broken or does no longer work properly.

- Nasal Aspirator AGU NS19 should be stored at a temperature between 25 °C–+55 °C (-13 °F–+131 °F), relative humidity 15%–85% (non-condensing).

ATTENTION!

After transport or storage at low temperatures, it is necessary to keep the device at room temperature for at least 2 hours before switching on.

12 TECHNICAL DATA

Guidance and manufacturer's declaration electromagnetic emissions.

The **AGU NS19** nasal aspirator is intended for use in the electromagnetic environment specified below. The customer or the user of the **AGU NS19** nasal aspirator should assure that it is used in such an environment.

EMC Table

NS19 is intended for use in the electromagnetic environment specified below. The customer or the user of **NS19** must make sure that it is used in such an environment.

Guidance and manufacturer's declaration – Electromagnetic emissions		
Phenomenon	Professional health-care facility environment a)	HOME HEALTHCARE ENVIRONMENT a)
Conducted and radiated RF EMISSIONS	CISPR 11 Group 1 Class A (Net BLE Function) Group 2 Class A (With BLE Function)	CISPR 11 Group 1 Class B (Not BLE Function) Group 2 Class B (With BLE Function)
Harmonic distortion	Not applicable (Note: Power by Battery or DC Input)	
	Only the AC input needs to be tested	
Voltage fluctuations and flickering	Not applicable (Note: Power by Battery or DC Input)	
	Only the AC input needs to be tested	

a) The equipment is suitable for use in Home Health Environments and Professional Health Care Environments limited to patient rooms and respiratory treatment facilities in hospital or clinics. The more restrictive acceptance limits of Group 1 Class B (CISPR 11) have been applied to the Home Health Environment. The use in the mentioned environments when directly connected to the Public Mains Network.

b) The test is not applicable in this environment unless the ME EQUIPMENT and ME SYSTEM used will be connected to the PUBLIC MAINS NETWORK and the power input is otherwise within the scope of the Basic EMC standard.

Guidance and manufacturer's declaration – Electromagnetic immunity – Enclosure port		
Phenomenon	Basic EMC standard or test method	Immunity test levels
		Professional health-care facility environment
		HOME HEALTHCARE ENVIRONMENT
ELECTROSTATIC DISCHARGE	IEC 61000-4-2	± 8 kV contact ± 2 kV, ±4 kV, ± ±8 kV, ±15 kV air
Radiated RF EM fields	IEC 61000-4-3	10 V/m 80 MHz – 2.7 GHz 80% AM at 1 kHz
Proximity fields from RF wireless communications equipment	IEC 61000-4-3	COMPLIANT NOTE: Further information about distances to be maintained between portable and mobile RF communications equipment (transmitters) and the **NS19** can be requested from supplier of the device. The device is designed for use in a home environment. However, it is advisable to keep the equipment at an adequate distance of, at least, 0.5 m from mobile phones or other RF communications transmitters to minimise possible interference
RATED power frequency magnetic fields	IEC 61000-4-8	30 A/m c) 50 Hz or 60 Hz

a) The equipment is suitable for use in Home Health Environments and Professional Health Care Environments limited to patient rooms and respiratory treatment facilities in hospital or clinics. The more restrictive IMMUNITY acceptance limits have been considered and applied. b) Before modulation is applied.

c) This test level assumes a minimum distance of at least 15 cm between the ME EQUIPMENT or ME SYSTEM and sources of power frequency magnetic fields.

13 DISPOSAL

The symbol on the product or on its packaging indicates that this product is not classified as household waste.

- Disposing of this device in a proper way you will help to prevent possible adverse consequences for the environment and people's health.
- For environmental protection purposes, this device should not be disposed of together with domestic (household) waste. The disposal should be performed in accordance with local regulations and laws.

- The device should be disposed of in accordance with the EU Directive 2012/19/EU WEEE (Waste Electrical and Electronic Equipment) for waste electrical and electronic equipment.

If you have any questions, please contact the local utility responsible for waste disposal.

14 CERTIFICATION

Nasal aspirator **AGU NS19** complies with the main requirements of the EU Directive 93/42/EEC concerning medical equipment, and Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment.

15 MANUFACTURER'S WARRANTY

The warranty period for the nasal aspirator **AGU NS19** is indicated in the warranty card.

The warranty is valid subject to availability of a warranty card filled out by a duly authorized representative confirming the date of sale and also of the relevant receipt. The warranty does not cover the consumable parts, components and constituent elements, namely: batteries, silicone narrow and large tips , transparent collection cup, storage bag.

The warranty and free maintenance service are not provided in case:

- the device is used with violation of the requirements set by the operation manual;
- damage is caused by consumer's willful or wrong actions as a result of improper or negligent use;
- traces of mechanical effect (dents, cracks, chips, etc.), tampering, dismantling, traces of repairs outside the authorized technical service center, traces of ingress of moisture into the case or effect of corrosives, or of any other outside interference with the device design are observed on the device case, and also in other cases of violation by the consumer of the rules of storage, cleaning, transportation and maintenance of the device stipulated by the operation manual;
- penetration of oil, dust, insects, liquids (not intended for the use with devices) and other foreign matters inside the device.

The warranty does not cover defects or malfunctions caused by:

- natural wear and tear of life-limited component parts;
- damage to accessories and tips that are not an integral part of the product (narrow silicone tip, large silicone tip, collection cup, storage case, package);
- using low-quality, worn out accessories and tips;
- force majeure (accident, fire, flood, electrical failure, etc.).

WARNING

Strictly follow the instructions in order to ensure reliable and long-term operation of the device.

Apply to the sales assistance in case the device does not work properly.

For the repair and maintenance apply to a specialized service department.

The manufacturer reserves a right to change the design of the device.

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